USE & CARE MANUAL

ZENET

Foot Massage Device

Model: ZET-762

Thank you for purchasing our foot massager. Before use, please read this Use and Care Manual carefully and pay special attention to the safety precautions in order to operate and use this product to properly.

SAFETY REGULATIONS

a) Description of important safety issues

- In order to ensure safe and proper use of the product, be sure to have read this Manual carefully before operation.
- Please, unplug from the socket after use or before cleaning the product.
- Use the heating function with caution to avoid possible serious burns.

Do not use on insensitive or overly sensitive skin or areas with poor blood circulation. Never allow children or people with physical or mental disabilities to use the heating function themselves. Improper use may be dangerous.

- In order to prevent fire, electric shock, or injury due to overheating, do not use the product under pillows or blankets.
- Do not use it closed in order to avoid motor overheating, which may cause malfunction or damage.
- If the product is used or placed near children, people with physical or mental disabilities or disabilities, they need someone else to guide them and be responsible for their safety.
- If you feel discomfort during the massage, stop the massage immediately.
- Do not use this product if the following conditions occur: the power cord or plug is damaged, operation is abnormal, the device has been dropped on the ground or damaged, or it has been dropped into water. Please contact your technical service center for more details, inspection and repair.
- Do not use replacement parts not recommended by the manufacturer.
- Do not insert or drop objects into openings of the product.
- If you want to turn off the power, turn off all controls, and then unplug the power cord from the socket.
- This product should not be used for self-diagnosis or self-treatment. If diagnosis or treatment is required, it should be performed by professional medical personnel.
- Do not use the product in a humid or dusty environment in order to avoid product malfunction.
- Do not place the product near an electric heater, stove, or in direct sunlight.

b) The use of this product is not suitable for the following groups of people


If you have the following conditions, consult your doctor before use.

- People suffering from malignant tumors
- People using a pacemaker or an implanted medical device
- People suffering from acute illnesses

c) Precautions when switching

- Do not operate the control switch with wet hands.
- Do not allow the control switch to come into contact with liquids such as water.
- Do not fall asleep while the device is operating.
- Do not use this product if you are intoxicated or feel unwell.
- Do not press the massage head directly against your arms, hand joints, feet, knees, chest or abdomen.
- Do not use this product continuously for more than thirty minutes.
- Do not use the product if its fabric cover is damaged.
- Do not use the product an hour before or an hour after meals in order to avoid side effects.
- When using the product at low room temperature, do not suddenly turn on the heating function. Before use it is recommended to gradually bring the temperature back to the normal range. If the product is moved from a cold storage location to a warm location for use, water vapor may accumulate in components and this will affect the normal operation of the product or even cause a component failure. After moving from a cold location to a warm location the product is recommended to be kept at a normal room temperature for 1 hour before using it.

COMPOSITION OF THE PRODUCT

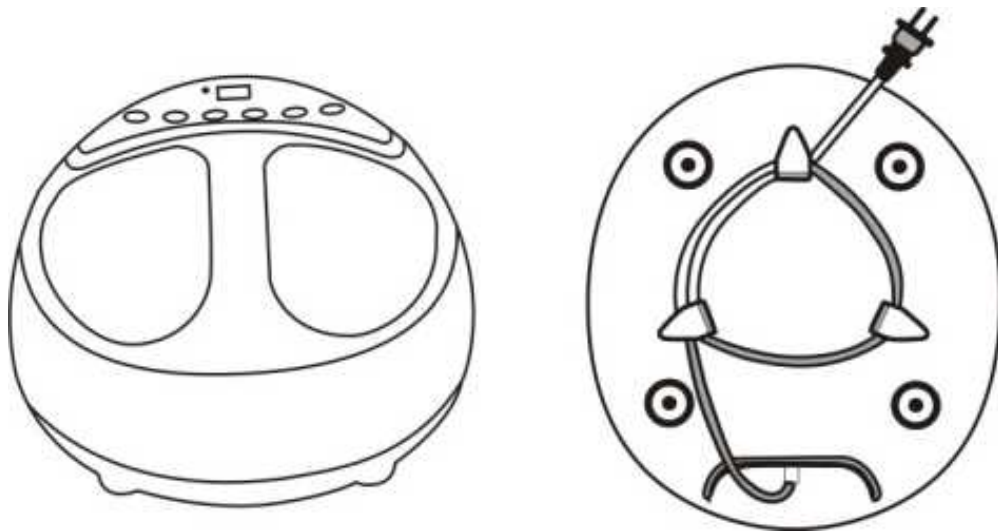- (1) Central surface: 6 independent control buttons to control the massage function of the product one by one.
- (2) Front and bottom center: Rotating massage head - for massaging acupuncture points in the foot arch area, which can be heated.
- (3) Front and lower front part: rotating massage head for acupuncture massage of feet one by one.
- (4) Front part, bottom part and rear part: rub the top massage head - in order to achieve proper massage of acupuncture points in the heel area.

- (5) Surface and interior: air cushions alternately press on the instep and heel for air pressure massage.
- (6) Front bottom surface: removable washable fabric cover which can be cleaned one at a time.
- (7) In the bottom: The power cord storage device allows you to store power cords one by one.

METHODS OF OPERATION AND DESCRIPTION OF FUNCTIONS

TIME – time setting button:

Adjust the massage time.

ROLL – button to turn compression massage on and off.

MODE – setting massage modes. The three massage modes of this massager are indicated by means of red lights.

STRENGTH – adjust the massage intensity, intensity levels can be switched by pressing the button.

Control the massage intensity of the whole machine and mark it with red light (middle light is on for medium range, left light is on for weak level, and right light is on for strong level)

HEAT - Heating keys:

When the massage function is turned on, the heating function can be controlled independently, and the indicator light is red during heating (default heating when turned on)

POWER – On/Off key:

Control the entire machine function to turn it on or off.

a) Preparation before use

- Check the product for damage.
- Make sure the placement is stable and comfortable.

b) How to use

1, Unplug the power cord of the wiring harness and connect it to a 220V/50Hz outlet.

2, Method of operation

- After connecting to the power supply, the whole machine is in standby mode and the key indicator on the control panel is flashing.
- Press the **POWER** button, the on/off button indicator light will turn red, the **MODE** button will enter mode 1 (the mode indicator light on the left above the mode button will turn on)

PRODUCT SPECIFICATIONS

Name of the product	Foot massager with heating function and an airbag
Product number	ZET-762
Rated voltage	220 V
Rated power	25 W
Timing	15 minutes

FAULT TRACING

If you have problems when operating this product, please use the following tips and possible solutions. If the product still has malfunctions, send it to a service center or to an authorized technician for inspection and repair.

Emergency situation, possible causes, remedy	Emergency situation, possible causes, remedy	Emergency situation, possible causes, remedy
The product does not work	1. The product is not switched on 2. Functional buttons are not switched on	1. Connect to the power supply or check if the power supply is well connected. 2. On/off button.
Sudden stoppage	1. Power outage 2. The massage time exceeds 15 minutes. 3. If the continuous use time is too long, the product will overheat. 4. The force is too strong and the massage ball is blocked.	1. Check whether the power interface is connected correctly. 2. Restart the on/off button. 3. Turn on the on/off button again after the product has cooled down. 4. Reboot on/off button

STORAGE OF THE PRODUCT

a) Storage

- When the product is not in use, please place it in a color box and place it in a safe, dry, ventilated and cool place.
- Do not twist the switching power supply or wires.
- Do not allow sharp objects to scratch the surface of the product.

- Do not place the product in direct sunlight or high temperatures.

b) Cleaning

- Before cleaning, make sure the main power switch is turned off and unplug the power cord.
- Use a dry cloth to clean the surface of the product.

Notification

- Never use harsh detergents such as benzene or thinner to clean the product.
- Keep the product clean and away from places with high humidity.

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