

142mm

90mm

Hot &Cool Skin Care Device Manual▶



Relieve fatigue
Lifting & tightening
Apply to Face,eye,neck;easy to carry,suitable for home and travel use

Diastolic & shrink pores
Anti wrinkle & whitening
D818

INTRODUCTION

1.Acoustic vibrations beauty principle

Each person's face has at least three million pores,while cleansing liquid and hands only clean skin surface,it is difficult to remove the oil pores,only clean out the pores of the oil residue.can the skin recover from the deep clean.Lessons learned 46 years specializes in skin care experience,the whole development process guided by skin specialists,combined with our advanced engineers to sonic micro-vibration technology,fitting face of intelligent design,no dead comprehensive clean your skin.8000 times per minute,high-frequency vibration,vibrate out of deep pores oil and deeply clean every pores.Pores lightweight transparent,delicate skin smooth,once again change the skin more through cleanliness.

2.Hot and cool beauty principle

A Japanese professor at the university shows that repeatedly reaching the skin at a temperature of 42 °C helps to produce HSP protein, repair skin cells and strengthen skin elasticity to resist the formation of wrinkles. The warm 42°C treatment helps to dilute the pores for deep cleansing, while the cold treatment with a cold effect of up to 6°C instantly tightens pores and locks in nutrition. Alternating warm and frozen skin care, it stimulates the liquid circulation at the bottom of the skin and speeds up the metabolism of the cells. Besides effectively reducing the elimination of eye circles and simple eye lines, it also brighten skin. Even without makeup, the skin is also revealed. Combined with micro-vibration, it is beneficial to enhance the skin's firmness and relieve fatigue.

3. EMS beauty principle (Electrical Muscle Stimulation)

EMS is used to activate the skin cell, stimulating rejuvenation of collagen, to improve skin elasticity, soften wrinkles and leave skin looking smooth and healthy.

Ice sense effect:

- 1 -

- 1.Help skin firmness,maintain flexibility,prevent allergies and other problems.
- 2.Shrink pores,lock moisture and reduce wrinkles.
- 3.Can soothe muscles,relieve skin fatigue.

Warm effect:

- 1.Warming,makes people feel warm and ,helps to relax muscles,relieve pain;
- 2.Stimulate the blood circulation acceleration,so the pores open,can help the skin to quickly import all kinds of skin care products absorb double skin,make the skin more elastic and moist.
- 3.Putting a sponge to deep cleansing,better able to export the facial toxins,blackheads, dirt.


4.LED photon beauty principle

- 1.Red light:wavelength 620nm(±10nm)18 LEDs deeply reaches into the 630nm skin to promote blood circulation and the collagen regeneration enhance the skin elasticity,to make the skin to be smooth and moisture.Be suitable for the slack & wrinkly skin.
- 2.Blue light:wavelength 470(±10nm) 18LEDs deeply reaches into the 480nm skin,to kill the bacteria,promote the blood circulation,can penetrate into 2.5mm under the cortex,destroy the bacteria living environment,anti- inflammatory,bactericidal, remove acne,acne scars and so on effectively.

NAME OF PARTS

- | | |
|---------------------------------|----------------------------|
| ①. Hot and cool panel | ②. LCD display |
| ③. Charging jack | ④. USB charging cable |
| ⑤. Power & Mode switch | ⑥. LED photon button |
| ⑦. Level of EMS and Temperature | ⑧. Circle of makeup cotton |
| ⑨. Base | |

- 2 -

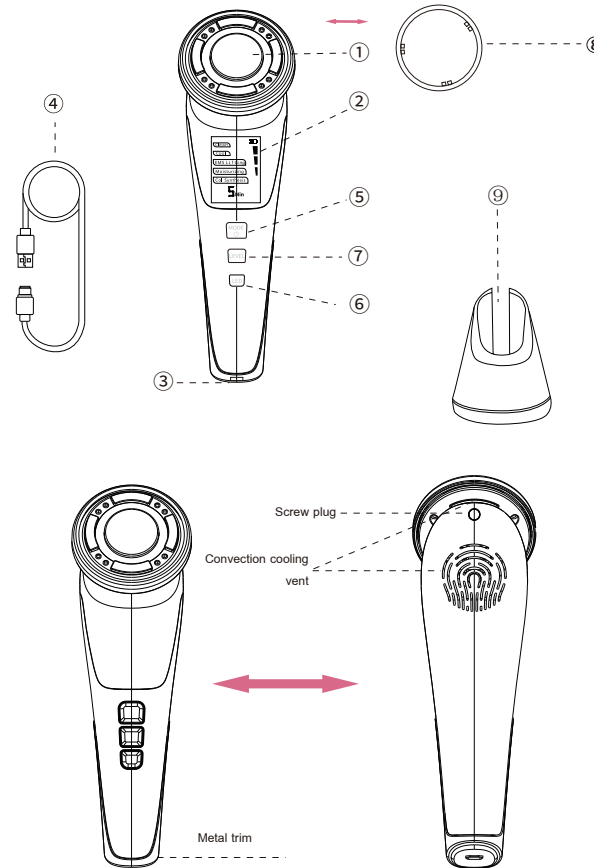

- 3 -

SPECIFICATION

Name:EMS Hot &Cool Skin Care Device	Model:D818
Input : DC 5V	Temperature : 6-42°C(±5°C)
Vibration:8000 (±10) rpm	Charging time :About 3 hours
Charging current:480mA	Power:MAX 8.5W
Battery capacity:1000mAh	Standby time:45 days
Net weight:145g	Use time:5min/time or 2 min/time
Machine size:175*55*67mm	Box size:215*150*84mm


OPERATION INSTRUCTIONS IN DETAIL

1. Please take full charge before the first using. The charging time is about 3 hours. Long press with 3 seconds on power key 5 to start machine, defaulting clean mode in 1 level, then the beauty instrument work with vibrate and heat. According to personal skin needs, temperature can be adjusted by press key 7. You can adjust LED Photon light by press key 6 in each mode (5min smart timing for mode 1).
- 2.Press the power key to the Cool mode, the beauty instrument work with vibration and cools. According to personal needs, you can press 7 to adjust the temperature (2min smart timing, after power off the fan still works for 30S to cool down)
3. Press the power key 5 to the EMS lifting mode.The beauty instrument vibrates and production high EMS micro current. According to individual needs, you can adjust the EMS level by press key7 (5min smart timing for mode 3).
4. Press key 5 to switch the moisturizing mode. The beauty instrument will vibrate & heat and sent EMS micro current to skin. According to individual requirements, you can adjust the EMS intensity by press key 7(5min smart timing for mode 4).
5. Press the power key 5 to switch the Collagen Synthesis mode. The beauty instrument vibrates and cools while emitting EMS micro current to skin. (Same 2 min smart timing as model 2)


- 4 -

According to individual requirements, you can adjust the EMS intensity by press key 7 (2min smart timing for mode 5).


6. Using Method as following pics and texts:

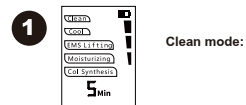

Need to use with remove makeup cotton.Take off the Circle, attaching the cotton, pour with makeup remover water. Then turn on the clean mode, warm and open the pores, shatter the dirt and export the deep layer of rubbish, clear pores, say goodbye to blackheads!

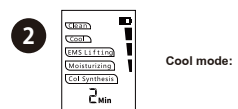

Put the mask on the face and turn on the Cool mode. The cool function can effectively shrink the pores and lock the skin moisture.

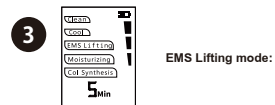


After cleansing, apply skin care products on the face and turn on the EMS mode. EMS stimulates the muscles and makes the skin more energetic and elastic. It improves fine lines and remove crow's feet and wrinkles on face, and restores the skin's youthfulness, smoothness and tenderness.



- 5 -

After cleansing, apply skin care product such as lotion and essence to your face, and activate moisturizing mode to promote double absorption of skin care products! (On each mode, you can select red light / blue light according to your skin needs)

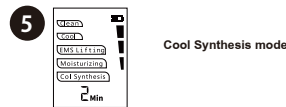

After the skin care products were absorbed, turn on the collagen synthesis mode, cool to shrink pores, lock the nutrition. The EMS micro-current stimulates the skin cells to synthesize collagen and rejuvenate the muscles anti-wrinkle. (On each mode you can Select red/blue light according to your skin needs)

SAFETY INSTRUCTIONS ⚠

- 1.Read all safety instruction before using.
 - 2.Follow the operating instructions to use this device.
 - 3.Do not pull down or remake the device,please consult the distributor for repairing.
 - 4.Only use the device for its intended purpose as described in this manual.
 - 5.Do not store or use the device in the high temperature and humidity place such as under sunshine or near fire or in the bathroom,etc.To avoid fire or electric shock.Do not use or place it in or near the high temperature appliance.
 - 6.Do not use and place the device in the humid condition, in case of any short circuit, electrical shock or fire.
 - 7.Keep away from children and physical disabled people allow to use, close supervision is necessary when this device is operated by user under 18 years of age.
 - 8.When operate cool function, the power consumption is relatively large, the fan will have a delay in heat dissipation. It needs to be placed for about 10 minutes for next time use after use cool function.
- WARNING:**
1. Please stop using when product abnormal heating to avoid burns.
 2. Please do not drop or clash the device.
 3. Please apply serum constantly during use.
 4. Please do not use it if the power cord is damaged.
 5. Below conditions are not allowed to use.

- 6 -

Heart disease,fever,communicable disease,cancer,open cuts,acute disease, hypotension,hypertension,install pacemaker,artificial heart or lung,during medical

WARRANTY CARD

1. Free maintain for the breakdown caused by product quality problems.
2. One year guarantee from the date of purchase.
3. Free to change broken parts, freight is collected by customers.
4. If breakdown caused by the human factor, free maintain won't be given and freight shall be collected by customers.
5. Please fill the guarantee clearly, without any revision. Please take good care of guarantee. If lost, won't be re-provided.

Product Name		Item No.	
Customer Name		Phone	
Purchase Date		Means Of Purchase	
Address			
Records			
Date	Content	Results	Service Provider

Implementation of standards: GB 4706.1-2005
GB 4706.15-2008

