

zenwire

MagniSonic
Oral Irrigator WF238

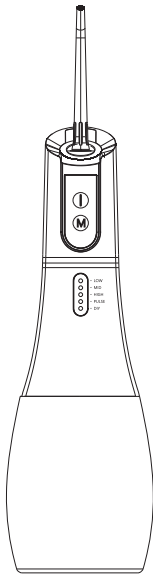User Manual

Product descriptions

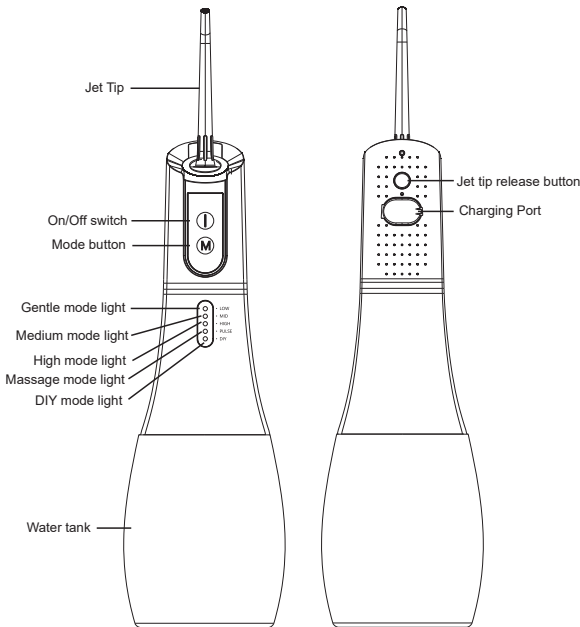
Product Accessories


Carry pouch x1

Charging Cable USB-C x1

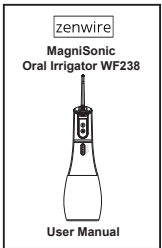
User Manual x1

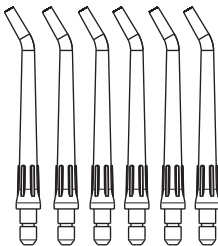
Jet Tip x6

How to Use

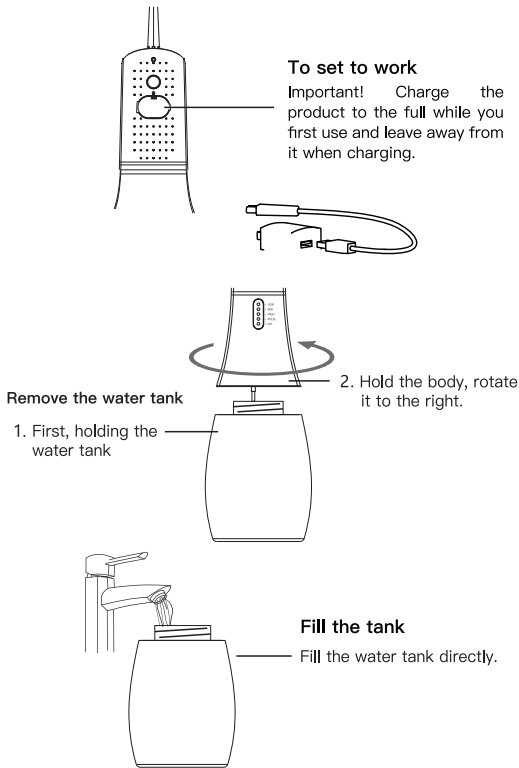Operating Procedures

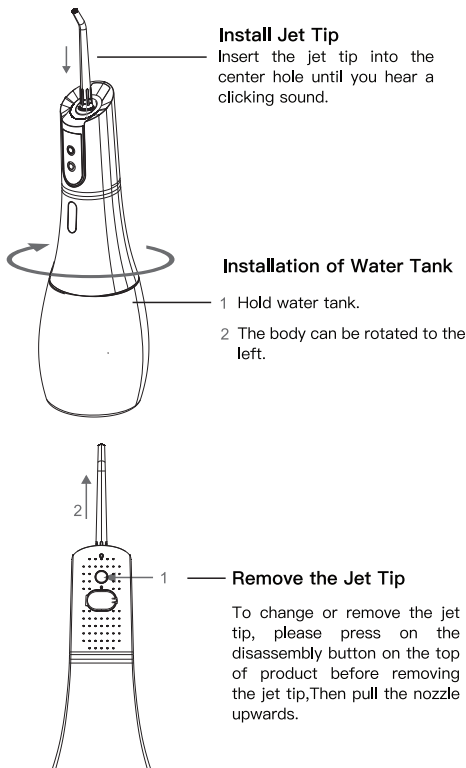
Instructions for Use

- ◆ The device has a memory function. When turned back on, it will resume the same water pressure mode used before it was switched off.
- ◆ When powered off, press the mode button to switch between the 5 available modes: Gentle, Medium, High, Massage, and DIY.
- ◆ You can change the mode while the device is off. Select your desired mode before turning it on.
- ◆ The device will automatically power off after 120 seconds of operation, regardless of the selected mode.
- ◆ If the Massage Mode is used for more than 3 seconds, the remaining battery level will be displayed.
- ◆ During charging: All five indicator lights will flash. When fully charged: All five lights will remain on steadily.
- ◆ If all five lights flash continuously during use, the battery is low and needs recharging.
- ◆ The device cannot be used while charging.
- ◆ Rotate the jet tip dial to adjust the direction of the water flow.
- ◆ DIY Custom Mode Setting:
 1. While the device is off, press the mode button to switch to DIY mode. The indicator light will turn on.
 2. Press and hold the power button. The device will start operating and will increase speed every 4 seconds. When the desired speed is reached, release the power button to save your custom setting.Note: The device has a memory function and saves the custom setting once.

Recommendation skills

1. Hold the device at a downward angle. Place the jet tip in your mouth and aim it at your teeth before switching on the device.
2. Make sure the water jet and gums form a 90° angle. Slightly close your mouth to prevent water from splashing out and to ensure effective cleaning.
3. Start by cleaning the back teeth (molars) and move toward the front teeth (incisors).
4. Guide the water flow along the gum line. Pause briefly between teeth to ensure the area around the gums and between the teeth is thoroughly cleaned.

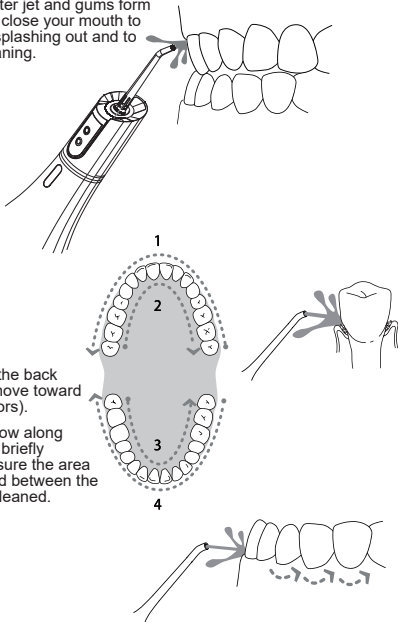Common Problems and Solutions


Problems	Reason	Method of handling
Leakage of water between jet tip and screw part	The seal ring on the fixed socket for the jet tip is damaged, and the jet tip cannot be installed properly	Contact the retailer for repair, and install the jet tip properly
Water tank leaks water	Seal ring is loose	Re-install correctly
Inadequate pressure	Low voltage usage	Use after recharging
Device cannot be switched on	Low voltage usage	Use after recharging
Product is not discharging water	Filtration net is clogged	Clean the filtration net
Can't draw water	The inlet filter screen is blocked	Remove the washer

If you still can't solve the problem according to the above methods, please consult the sales store you purchased or go to the designated after-sales maintenance point for repair.

Warnings

1. Keep out of reach of children and animals. This product is not a toy and should be operated by adults only.
2. Do not use while charging.
3. Do not disassemble or attempt to modify the device.
4. Avoid direct contact with eyes or ears during use.
5. Use clean water only. Do not use mouthwash, saltwater, or other additives unless specified.
6. Stop using immediately if the device malfunctions or shows visible damage.
7. Do not expose to extreme temperatures or fire.
8. IPX8 Waterproof: Suitable for use in the shower or wet environments. Do not submerge the device for extended periods.
9. Fully charge before first use (approximately 4 hours).
10. If storing for an extended period, charge fully in advance.
11. Recharge promptly when power is low to preserve battery life.
12. Contains a built-in non-replaceable lithium battery.
13. At the end of battery life, dispose at an approved recycling center. Do not throw in regular trash.

E-mail: sklep@zenwire.eu
Manufacturer: Zenwire Sp.z o.o. Kołobrzaska 4, Poznań 60-418, Made in PRC.

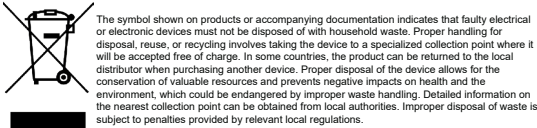Business users in the European Union

If you need to dispose of electrical or electronic devices, please contact your nearest point of sale or supplier for additional information.

Disposal of waste outside the European Union

This symbol is only valid in the European Union. If you need to dispose of this product, please contact local authorities or the seller for information on the correct method of disposal.

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