ITCH HEALER

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Výrobca:

Zhejiang Sorbo Technology co., Ltd.
No. 2, Daqing Road, Xinqiao Street,
Ouhai, Wenzhou, Zhejiang,
China 325016
sb@wzshengbo.com

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Weight: 40g

Principle:

In the event of a sting, insect saliva activates a defense reaction.

Immune cells are activated, histamine is released, and sent to the nerve cell receptors.

This creates an immune response, where the area may become irritated and itch.

Itch Healer can sooth this skin reaction, using "local thermotherapy" which is based on two operational principles:

1. Heat acts a counter-stimulus. Applying heat inhibits the itch signaling pathway to the nerve cells of the skin and relieves itching.
2. Heat slows the histamine release. Studies show heat activates various positive immune reactions including a decrease of histamine.

Functions and Features:

- 1.CHEMICAL-FREE TREATMENT: Itch Healer is an electronic device that provides symptomatic relief for bites and stings caused by insects such as mosquitoes, wasps and bees.
- 2.FAST AND EASY TO USE: Choose seconds setting for fast relief from the itching, pain, and swelling caused by insect bites and stings.
- 3.Safe Bite Relief for Child: Temperature Control Exactly, Sting healer has 4 different temperatures to choose from, which are suitable for everyone.

Operation Method:

4 Different Temperatures Choose :

Position1 40℃

Position2 43℃

Position3 46℃

Position4 51℃

Keep the Temperature Switch switch for 2 seconds into choose Temperature Control(automatic shutdown after 60 seconds) ,then press Healer Switch start working(automatic shutdown after 30 seconds).

Input Charge

Type C cable connect with AC Adaptor, PC etc., Temperature Indicator go on sparking when input charging and go on lighting after fully charge.

Fast,chemical-free treatment for insect bites and stings

mosquito

honeybee

wasp

Ants


Lice

spider

Bedbug

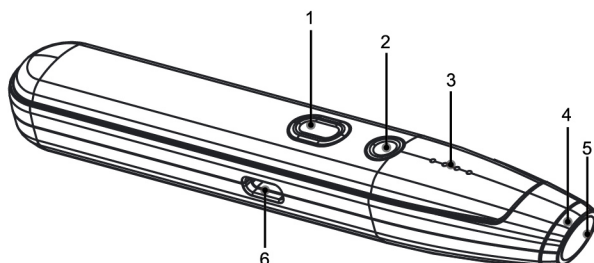1. Temperature Switch
2. Healer Switch
3. Temperature Indicator
4. Working Indicator
5. Ceramic Healer
6. Type C Input Charging

MADE IN CHINA