TWO610

PL / INTLIGENTNA WAGA OSOBOWA **EN /** SMART PERSONAL SCALE
DE / INTELLIGENTE PERSONENWAAGE **RU /** ИНТЕЛЛЕКТУАЛЬНЫЕ ЛИЧНЫЕ ВЕСЫ
ES / BALANZA PERSONAL INTELIGENTE **FR /** BALANCE PERSONNELLE INTELLIGENTE
IT / BILANCIA PERSONALE INTELLIGENTE

Eldom Sp. z o.o. • ul. Pawła Chromika 5a • 40-238 Katowice, POLAND
tel: +48 32 2553340 • fax: +48 32 2530412 • www.eldom.eu

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/EN/ The disposal of used electrical and electronic equipment (applies to European Union countries and other European countries with separate waste-collection systems).

This symbol on the product or its packaging indicates that it should not be classified as household waste. It should be handed over to an appropriate company dealing with the collection and recycling of electrical and electronic equipment. The correct disposal of the product will prevent potential negative consequences for the environment and human health resulting from hazardous substances present in the product. Electrical devices must be handed over to restrict their re-use and further treatment. If the device contains batteries, remove them, and hand them over to a storage point separately. DO NOT THROW EQUIPMENT INTO THE MUNICIPAL WASTE BIN. Material recycling helps to preserve natural resources. For detailed information on how to recycle this product, please contact your local authority, the recycling company, or the shop where you bought it.

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SAFETY INSTRUCTIONS

Carefully read the entire contents of this manual before first use.

Keep all warnings and safety instructions for future reference.

1. The equipment is intended for domestic use only.
2. This equipment may be used by children of at least 8 years of age and by persons with reduced physical and mental capabilities and persons with inexperience and lack of familiarity with the equipment, if supervision or instruction is provided on how to use the equipment in a safe manner so that the associated risks are understood. Children should not play with the equipment. Unsupervised children should not carry out cleaning or maintenance of the equipment.
3. Do not exceed the maximum allowable weight (180kg) - this may damage the scale or cause it to calibrate.
4. Handle the scale with care. Do not step on the scale when it is wet. Do not stand on the scale with wet feet. Do not jump on the scale or drop the scale.
5. Use a damp, soft cloth without chemicals to clean, then wipe dry.
6. Do not immerse the equipment in water.
7. Repairs to the equipment may only be carried out by an authorised service centre. Any upgrades or the use of other other than original spare parts or components of the device is prohibited and jeopardises the safety of use.
8. Store the scale in a shaded and dry place.
9. Do not press on the scale with other objects during storage.
10. The accuracy of the scale indication depends on the way it is used.
11. Eldom sp. z o.o. shall not be held liable for any damage resulting from improper use of the device.

WARNING: Improper use of the appliance may cause injury.

CAUTION: The plastic bag can be a hazard - to avoid suffocation with the bag, keep it away from babies and small children.

ADDITIONAL INFORMATION

CAUTION

1. When in use, the scales generates a low-intensity current, which passes through the body. This can cause electro-medical devices (e.g. pacemakers) to malfunction, which poses a risk to the health of the user.
2. The use of the scales by pregnant women or patients on dialysis is not recommended.

CAUTION

- A. The indicators estimated by the device are only given as approximated values for the needs of the user. They are not intended for medical purposes. If these percentage values are lower or higher than the provided reference values, contact your doctor for consultation.
- B. For accurate measurement results, do not place the scales within an electromagnetic field (e.g. near a microwave oven).
- C. The surface of the scales can be slippery if wet.

Hereby, Eldom Sp. z o.o., 40-238 Katowice, ul. Pawła Chromika 5a declares that the radio equipment type smart personal scales TWO500, TWO550, TWO600 is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: **www.eldom.eu/wsparcie/deklaracje**

OPERATING INSTRUCTION

SMART PERSONAL SCALE TWO610

EN

• GENERAL DESCRIPTION

1. Display
2. Electrodes


• TECHNICAL DATA

- maximum weight: 180kg
- accuracy of graduation: 5g
- accuracy of measurement: 50g
- power supply: 2xAAA battery (included)

• INTENDED USE

The scale is designed for weight measurement, as well as for determining the level of fat, water, muscle, bone mass, visceral fat, BMR, protein, obesity, biological age, and LBM.

Additionally, the device displays the ambient temperature.

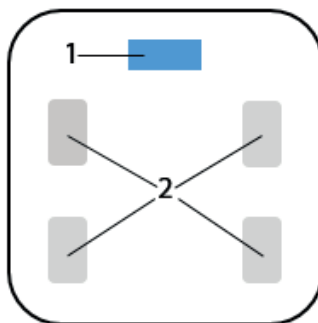• TIPS

- The indicators are determined based on such parameters as weight, height, age, etc., and they facilitate the determining of the overall condition of the body.
- It is recommended that you take measurements at the same time daily (e.g. between 7:00-9:00 p.m.).
- Measurement results can be erroneous after exhaustive physical exercises, when on a fortified diet, or in a state of serious dehydration.
- Place the scales on a firm and flat surface.
- Stand still while taking the measurement.
- The measurements of body-fat and water levels can give invalid results, if taken by:
 - adults over 70 years of age (the normal weighing mode can be used),
 - bodybuilders and professional athletes.

• USE

- Before use, check that the unit is complete and not damaged.
- Remove all transport protection.
- Insert the batteries.


• BATTERIES

The scale is powered by two AAA batteries.

Battery replacement

- Open the battery cover located on the bottom of the unit.
- Gently lift up and remove the batteries.
- Insert new batteries. **Pay attention to the polarity.**
- Close the cover.

Used batteries should be disposed of in special containers.

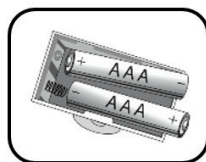• OPERATION

Taking measurements without using the app – standard weighing mode

- Place the scale on a stable, dry surface.
- Stand on the scale, evenly distributing your feet (the scale will turn on automatically).
- Wait a few seconds for the weight indication to stabilize.
- Read the display.

Taking measurements without using the app

NOTICE Due to the continuous development of our products, the app and its appearance may change over time.

The scale connects with devices that support wireless communication with at least Bluetooth 4.0 and operating systems iOS 12.0/Android 6.0 or newer.


You can create up to 8 user accounts.

1. Download the app using the code located on the cover of the manual or search for the keyword "OKOK-International" in the "App Store" or "Google Play".
2. For the first use, create a new user account.
3. Ensure that Bluetooth is enabled on your smartphone.
4. Launch the app.
5. Stand barefoot on the scale and evenly distribute your feet on the electrodes (2).
6. Measurements will appear on the scale's display and in the app.
7. All measurements will be automatically saved in the app's data history.

The app allows tracking changes over time for 18 parameters. Detailed explanations of the indicators are available in the app.

Automatic switch-off

A few seconds after stepping off the scale, the device will automatically switch off. If you do not step on the scale after switching it on, the unit will also switch off after several seconds.

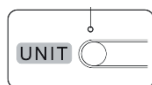Changing the weighing units

You can change the weighing unit (kg/lb) using the button on the bottom of the scale.

Weight overload message

The maximum permissible weight (180 kg) must not be exceeded. This can cause damage to the device.

Low battery indicator

Indicates that the batteries need to be replaced.

Error message

When the error message "LCER" appears on the display or the scale does not switch off, the battery can be removed for a few seconds. Incorrect scale operation and the above-mentioned message can also be caused by the proximity of devices such as mobile phones, microwave ovens, CB radios, etc. If neither recharging the battery nor moving the scale away from these devices has solved the problem, please contact service.

• CLEANING AND MAINTENANCE

- Do not immerse in water.
- Use a damp, soft cloth without chemicals to clean, then wipe dry.
- Wipe the housing with a moistened cloth.

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