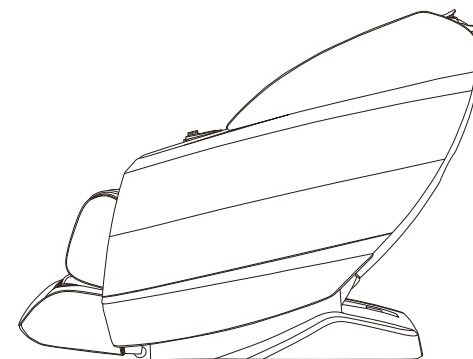
**iRest**

**SAVE THESE INSTRUCTIONS**

iRest Queentri(A357)  
**MASSAGE CHAIR**  
— Operation Instructions —

Thanks for purchasing this product. Please read this manual instruction carefully in order to know how to use this product correctly before operation. Please pay more attention to some important safety warning and keep this manual book properly for future use.

Note: Our company reserves the right to change and explain the design of the product. If there is some difference between the product function in market circulation and the description in the manual, it is a normal update of the product, and our company will not inform you otherwise. The specific product function and color shall prevail in kind.

## CONTENTS

### Before use

---

Safety precautions	02-05
Installation parts	06
Place massage chair	07-08
Before use	09
Product structure	10-15

### Usage method

---

Operation Instructions	16-31
Massage close	31

### Maintain

---

Cleaning and Maintenance	32
--------------------------	----

### Important

---

Troubleshooting	33
Product Specifications	34

**Safety Precautions** (must be observed)

- Please read all relevant instructions carefully before use.
- The notes marked here will help you to use the product safely and correctly.
- In order to clearly indicate the extent of harm and damage, precautions are divided into **"Warning"** and **"Caution"** that may be caused by incorrect use. Both of them are important for safety and must be followed.

 <b>Warning</b>	Indicates that serious injury or death may result from improper use.
 <b>Caution</b>	Indicates that misuse may cause minor injury or property damage.

- Please pay attention to the symbol before each safety item.

 <b>Forbidden</b>	Absolute prohibition
--	----------------------

**The following populations should not use this product**

 <b>Forbidden</b>	<ul style="list-style-type: none"> <li>■ The device is not intended for use by persons (including children) with physical disabilities, sensory or neurological impairments or who lack experience and common sense, unless they are supervised and instructed in the use of the device by those responsible for their safety. Children must be supervised to ensure they do not use tools for play.</li> <li>■ Patients with osteoporosis.</li> <li>■ Patients with severe heart disease and wearing implantable medical devices such as pacemakers and defibrillators that are susceptible to electromagnetic waves.</li> <li>■ People who have a fever (above 38 degrees ).</li> <li>■ Pregnant women or shortly after childbirth and during their menstruation period.</li> <li>■ People who follow medical advice and are prohibited from receiving massage due to thrombosis, aneurysm, severe varicose veins, or other circulatory diseases. Within two weeks after the appearance of blood clots, it is strictly prohibited to massage any form of body parts that have blood clots.</li> <li>■ Patients who are injured or have a skin disease.</li> <li>■ When unsupervised, any children are under the age of 14 and those unconscious are not allowed to use this product.</li> <li>■ People who require rest or feel unwell as instructed by a doctor.</li> <li>■ People with spinal abnormalities or scoliosis are prohibited from receiving this product due to medical advice due to other diseases.</li> <li>■ People whose bodies are wet are forbidden to use this product.</li> <li>■ This product has a hot surface. Persons insensitive to heat must pay attention when using this product.</li> </ul>
--	--

**Safety Precautions** (must be observed)

**Important safety warnings**

 <b>Warnings</b>	<ul style="list-style-type: none"> <li>■ Keep children away from the movable parts of this product.</li> <li>■ Use the well-grounded power supply unit suitable for this product.</li> <li>■ Pull out the plug after use or before cleaning to avoid any damage to physical injury or the product.</li> <li>■ Operate this product as instructed in this manual.</li> <li>■ Do not use any parts or accessories which are not recommended.</li> <li>■ Do not use this product outdoors.</li> <li>■ Do not use near open flames such as stoves, and do not smoke while using.</li> <li>■ Do not use heating equipment such as electric blankets on this product.</li> <li>■ Do not sit on the back of the chair, leg massagers, shoulders, armrests, or stand on the seat.</li> <li>■ Please read the Operation Instructions carefully before operation.</li> <li>■ Never use this product for any other purpose not listed in the manual .</li> <li>■ 20 minutes of usage each time is recommended.</li> <li>■ Do not use this product if the cover or leather is broken or damaged.</li> <li>■ Do not drop anything into this product.</li> <li>■ Never fall asleep while using this product.</li> <li>■ Do not use this product if you are drunk or feel unwell.</li> <li>■ Do not use this product within one hour after a meal.</li> <li>■ Do not use this product with excessive force to avoid any injury.</li> <li>■ Do not use for purposes other than massage.</li> <li>■ If the supply cord is damaged, to avoid danger, it must be replaced by a professional from the manufacturer, its maintenance department or similar department</li> </ul>
---	---

**Operating environment**

 <b>Caution</b>	<ul style="list-style-type: none"> <li>■ Do not use this product in a very moist environment such as the bathroom.</li> <li>■ Stop using this product immediately when the ambient temperature changes sharply.</li> <li>■ Do not use this product in any heavily-dusty or caustic environment.</li> <li>■ Do not use this product in a small space or a place without good ventilation.</li> </ul>
--	---

**Safety Precautions** (must be observed)

Storage conditions	
 <b>Caution</b>	<ul style="list-style-type: none"> <li>■ Storage room temperature: 5 °C to 40 °C.</li> <li>■ Relative humidity (RH) between 20% and 80%; free of corrosive gases, and a well-ventilated room.</li> </ul>

Product service and maintenance	
 <b>Caution</b>	<ul style="list-style-type: none"> <li>■ This product should be maintained by the service center designated by the manufacturer only. No user should disassemble or maintain it without permission.</li> <li>■ Do not forget to turn off the main power switch after usage.</li> <li>■ Do not use this product if the power socket becomes loose.</li> <li>■ If this product will be left unused for a long period of time, it is necessary to roll up the power cord of this product and store this product in a dry and dust-free environment.</li> <li>■ Do not store this product at a high temperature or near open flame. Avoid prolonged exposure to direct sunshine.</li> <li>■ Please clean this product with dry cloth. Never use thinner, benzene or alcohol.</li> <li>■ The mechanical components of this product are specially designed and manufactured, so no special maintenance is needed.</li> <li>■ Do not use a sharp object to stab at this product.</li> <li>■ Do not roll or pull this product on any uneven ground. Instead, it should be lifted before moving.</li> <li>■ Please use this product intermittently. Do not keep it working continuously for a long period of time.</li> </ul>

Solutions to ordinary malfunctions	
 <b>Caution</b>	<ul style="list-style-type: none"> <li>■ If the motor sounds during operation of this product, this is the normal operation sound.</li> <li>■ If the controller cannot work normally, check whether the plug and the power socket are connected reliably, and whether the power switch is turned on.</li> <li>■ If the rated working time is over, this product will automatically activate the power-off button; If this product keeps operating for a long period of time, the temperature protector will stop this product automatically. The product can be used again after half an hour.</li> </ul>

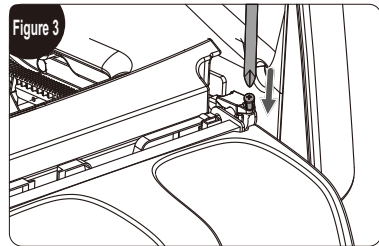
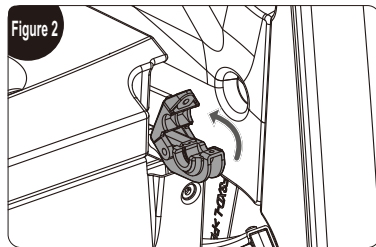
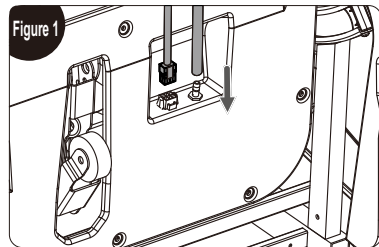
**Safety Precautions** (must be observed)

Safety precautions	
 <b>Caution</b>	<ul style="list-style-type: none"> <li>■ Check whether voltage is consistent with the specification of this product.</li> <li>■ Never insert or pull out the plug with a wet hand.</li> <li>■ Do not let water flow into this product to avoid electric shock or cause damage to this product.</li> <li>■ Do not pull the power cord when plugging or unplugging it. Rough handling is forbidden.</li> <li>■ Do not damage the wire or modify the circuit of this product.</li> <li>■ Do not clean live parts of this product with a wet cloth such as a power switch and plug.</li> <li>■ Be away from this product in case of power failure to avoid any injury if power is restored suddenly.</li> <li>■ Stop using this product at once while it works abnormal, and consult your local supplier.</li> <li>■ Stop using this product if you feel unwell, and consult your health care practitioner.</li> </ul>

## Installation parts

### Calf installation instructions

1. Insert the wires of calf and the air tube into the front socket of the steel frame assembly respectively (Figure 1).
2. After the wire docking is completed, open the two calf fixing covers on the steel frame assembly (Figure 2).
3. Lift the calf, put it into the notch, cover the fixing cover, and finally fix it in place with 2M4\*10 screws (Figure 3).


## Place massage chair

### Installation placement

**Make sure there is enough space for the massage chair can be tilted at least 12 cm.**

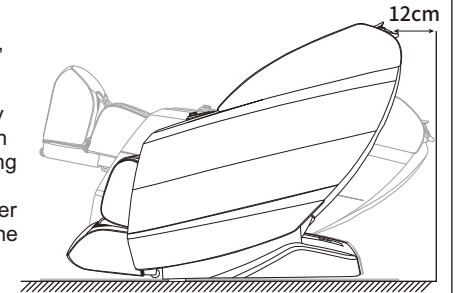
Do not expose the massage chair directly to high temperature environment due to sunlight (such as: in front of heating equipment), and place it to where can avoid high temperature and direct sunlight. It is recommended to lay a soft cushion for long-term use in a designated location. When laying the cushion, the size of the cushion should be enough to cover Keep the range where the machine touches the floor and the range where your feet touch the floor.

### Product lift instruction

**Require team work to lift and move the chair to avoid injury.**

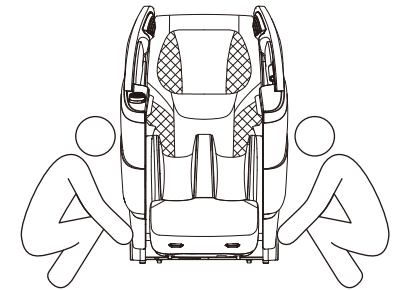1. Please lift the massage chair to move it on vulnerable ground such as wooden floor.
2. When moving on a similar wooden floor, it is recommended to lift it. Avoid damage to the casters. When the floor is lifted and moved to the designated position, be careful to keep your hands out of the way to avoid injury, and be careful of your feet being crushed. Do not let go of the unit until the unit is completely flat, the floor may be damaged, therefore, it is recommended to place mats, etc. on the floor.


## Place massage chair

### Move method

Use the rear casters to move the massage chair. Move the chair by tilting it down and backward about 45° as shown in the diagram.


**NOTE: Excessive force to title will cause the chair to tip over completely and product damage.**

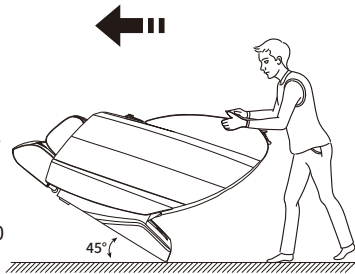1. Before moving, please shut down the machine properly to reset the whole machine, and then cut off the power.
2. Please keep the chair vacant when moving. (People, pets or clutter will not be able to stay on the chair.)
3. When moving, do not cross a drop barrier higher than 2mm or a gap greater than 5mm.
4. When using casters to move this product, single continuous movement is not supported; if the distance exceeds 50 meters, please make this product move slowly and evenly.

### How to disassemble and use headrest cushion and backrest cushion

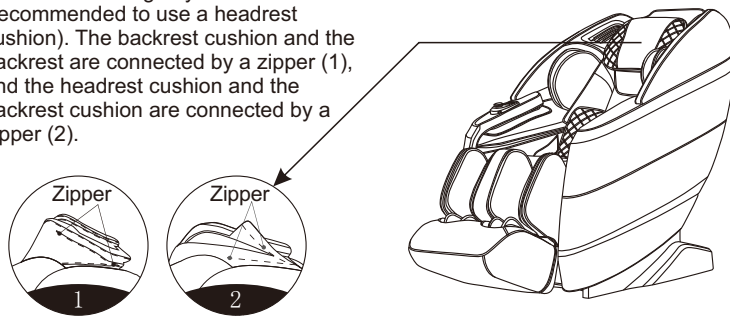Using a headrest cushion can reduce the intensity of kneading massage on the neck and shoulders, and you can decide whether to use a headrest cushion according to your needs (recommended to use a headrest cushion). The backrest cushion and the backrest are connected by a zipper (1), and the headrest cushion and the backrest cushion are connected by a zipper (2).


## Before use

### Check around

Caution

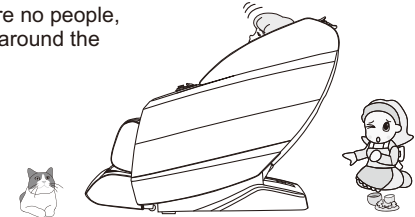Make sure that there are no people, pets and other objects around the unit.


### Make sure there are no foreign objects in the gaps of the machine

Caution

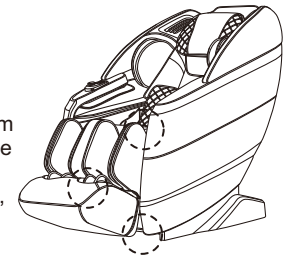
Before massaging, make sure that no foreign objects are caught in the legs, feet, etc.  
Under no circumstances allow hands, feet, and head to penetrate into the area between the lower leg mechanism and main body, or the area between the armrest and the capsule.  
Under no circumstances should hands, feet, or head sink under the fabric or leather covering of the massage area.


### Turn on the power

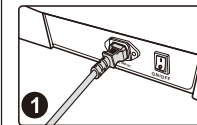
Caution

1 Plug in power cord

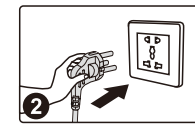2 Schematic diagram of power on

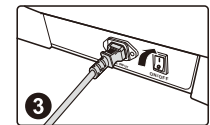3 Turn on the power switch ("I" position)


Warnings

- Please check the power cord and plug for damage before use.
- Please supply power to this product in accordance with the working power requirements of this product's calibration.
- It is strictly forbidden to use a power supply other than the nominal power supply of this product.
- It is strictly forbidden to connect various power conversion devices without authorization to supply power to this product.

## Product structure
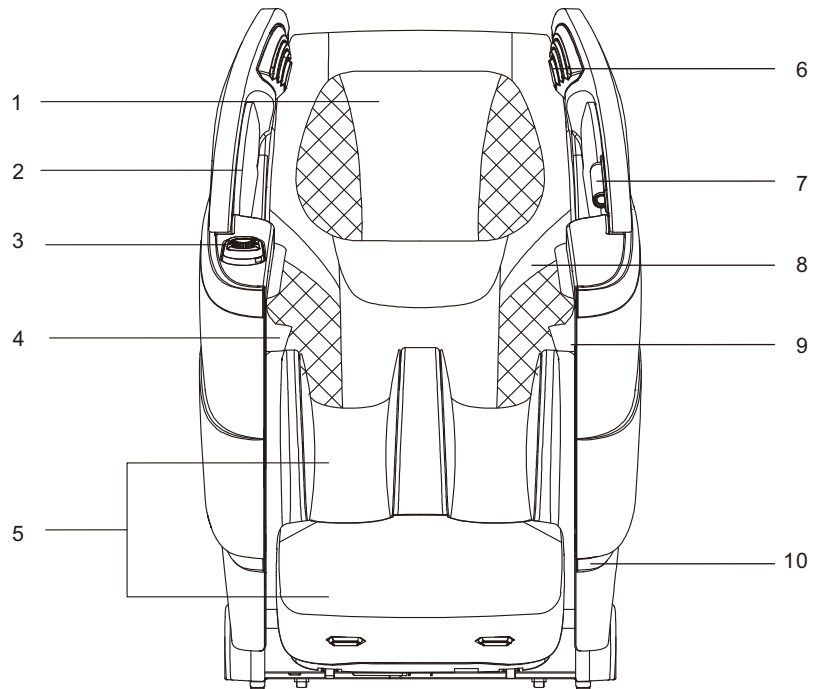
### Part description


- 1.Head cushions
- 2.Shoulder air pressure assembly
- 3.Handrail shortcuts
- 4.Armrest air pressure assembly
- 5.Leg massage mechanism

- 6.Flared hood
- 7.Reading light
- 8.Back cushions
- 9.Remote control
- 10.Mood lights

## Product structure
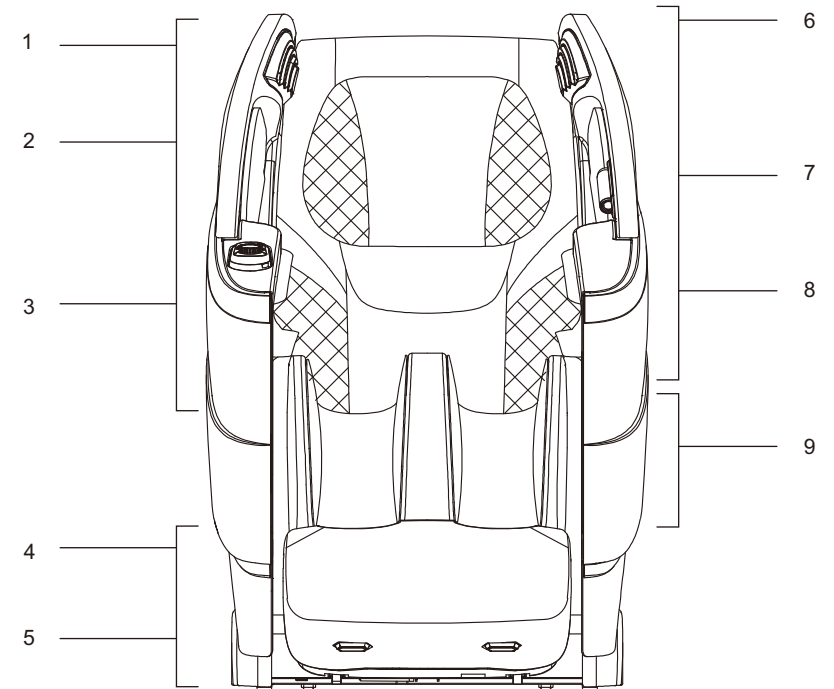
### Schematic view of functional distribution


- 1.Built-in voice Control microphone
- 2.Shoulder Air bag massage
- 3.Palm Air bag massage
- 4.Foot Air bag massage
- 5.Sole Air bag massage
- Roller massage

- 6.3D digital audio music playback on both sides of head
- 7.Mechanism massage from head and neck to buttock
- 8.Back heating Heat function
- Waist air bag
- 9.Foot Air bag massage

## Product structure
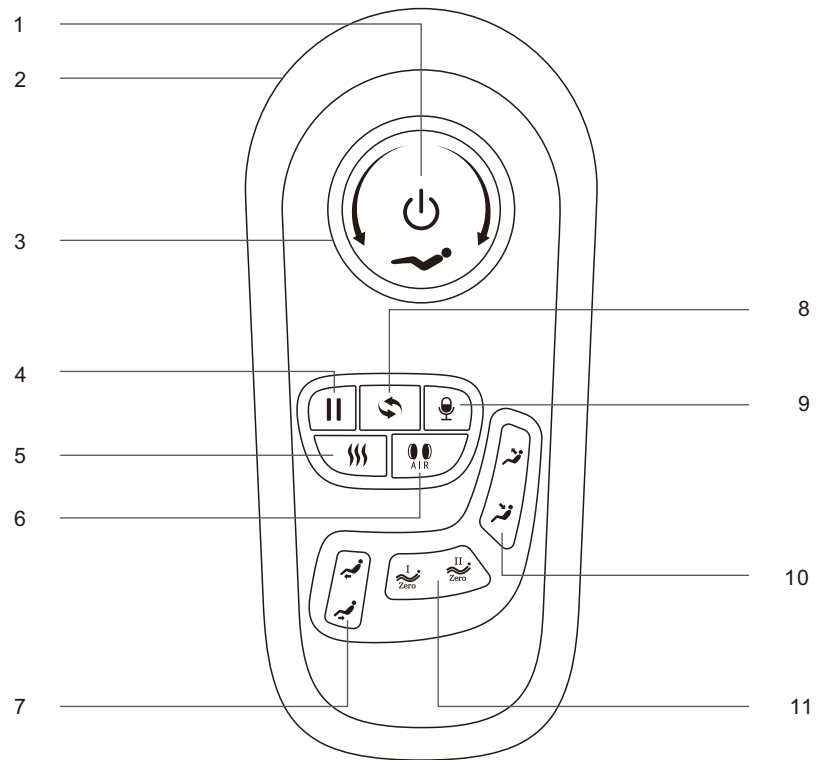
### Armrest shortcut key description

- 1.ON/OFF key: Turn the massage on or off.
- 2.USB charging port:18W fast charge.
- 3.Linkage adjustment key: Linkage adjustment of backrest frame and calf frame lifting.
- 4.Pause key: Pause all massage functions.
- 5.Heat therapy key: Turn the heating function on or off.
- 6.Air pressure key: Turn the air pressure massage function on or off.
- 7.Leg mechanism lift key: Adjust the calf frame to rise or fall.
- 8.AUTO key: Switches the automatic program.
- 9.Voice key: Short press to wake up the voice assistant, long press to turn on or off the voice function.
- 10.Seated lift and lower key: Adjust the backrest frame to rise or fall.
- 11.Zero gravity key: Adjust to the massage chair zero gravity angle.


## Product structure
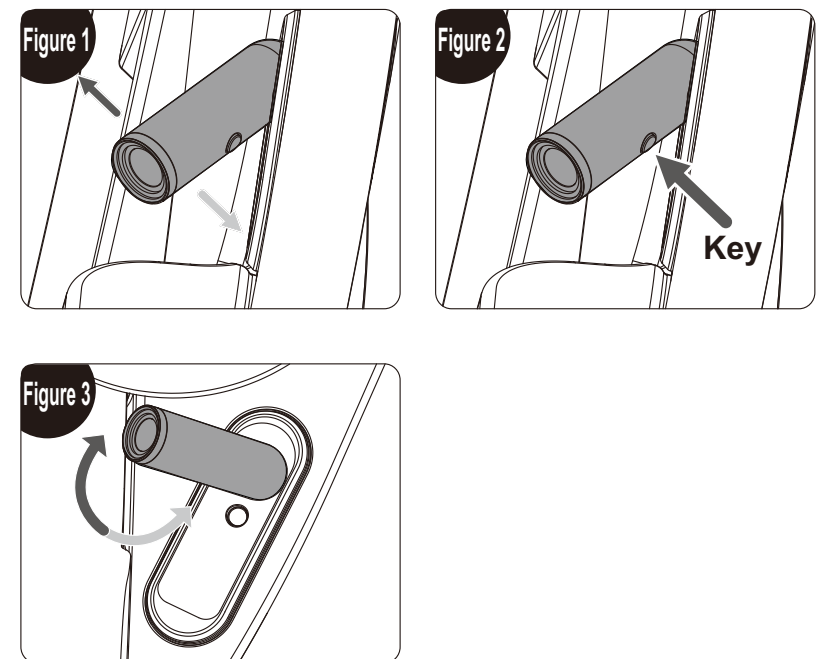
### Reading light instructions (can be used without power on after charging to electricity)

Note: After the user is sits, rotate the reading light out (Figure 1), and rotate the reading light back after use.

**Operation instructions:** Rotate the reading light outward, press the key (Figure 2) to turn the light on, and then rotate the reading light up and down (Figure 3) to adjust to the desired position. After use, rotate the reading light to its original position.


**Key description (Figure 2):** The first time the key is pressed, the light is on; Press the key for the second time, and the brightness of the lamp is increased to the second level; Press the button for the third time, and the brightness of the lamp is increased to third level; The fourth time you press the button, the light turns off; And so on and so forth. After setting the brightness each time, the light will be turned off automatically after 1 hour.

## Product structure
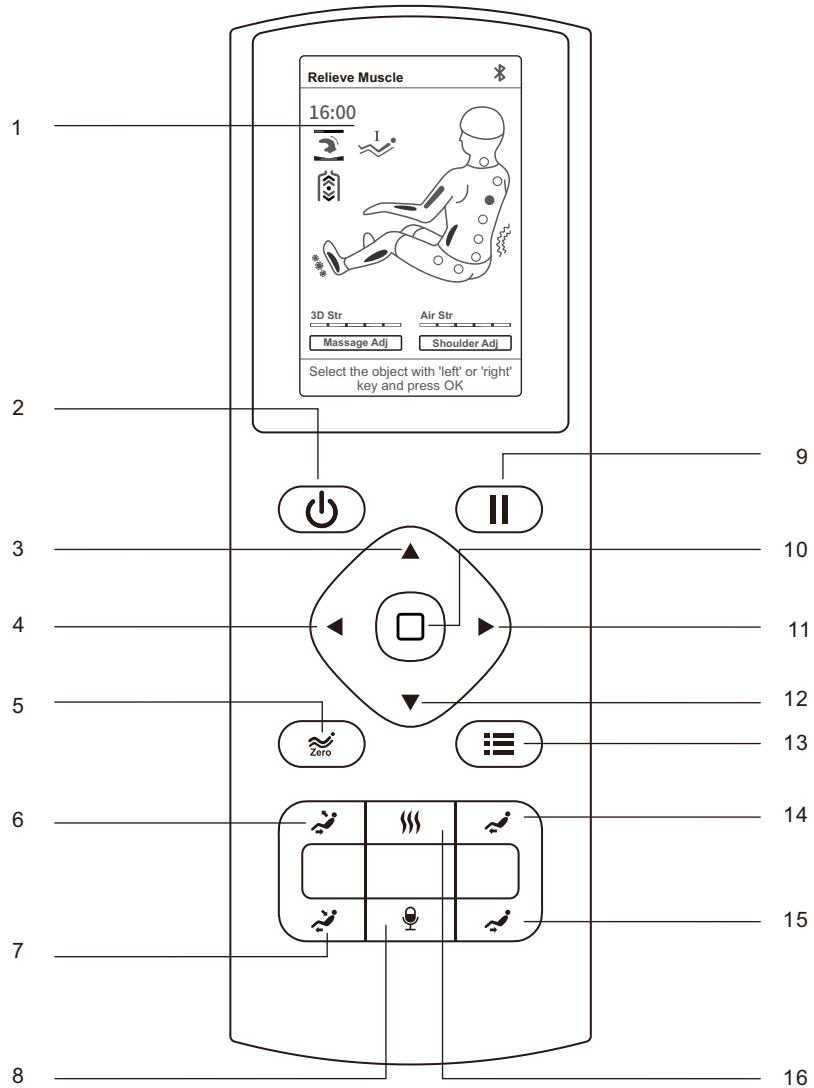
### Remote control key description


## Product structure

1. Display screen
2. ON/OFF key: Turning the massage chair on or off
3. Up key: Navigation up key, select the object up in the menu, or adjust it up.
4. Left key: Navigation left key, select an object to the left, or return to the previous menu.
5. Zero gravity key: Zero gravity angle setting, zero gravity I zero gravity I zero gravity III and zero gravity off, cycle adjustment.
6. Sitting up key: Adjust the massage chair angle upwards.
7. Sitting down key: Adjust the massage chair angle downwards.
8. Voice key: Wake up the voice assistant.
9. Pause key: Pause massage or continue massage.
10. Confirm key: To confirm selected massage function or cancel selected massage function.
11. Right key: Navigation right key, select an object to the right or go to the next menu.
12. Down key: Navigation down key, select the object down within the menu, or adjust it down.
13. Menu key: Into main menu selection screen.
14. Leg mechanism rise key: The leg mechanism adjust upwards.
15. Leg mechanism down key: The leg mechanism adjust downwards.
16. Heat therapy key: Turn the heating function on or off.

## Operation Instructions

### VI. Personalization program (please use this function after the chair is turned on)

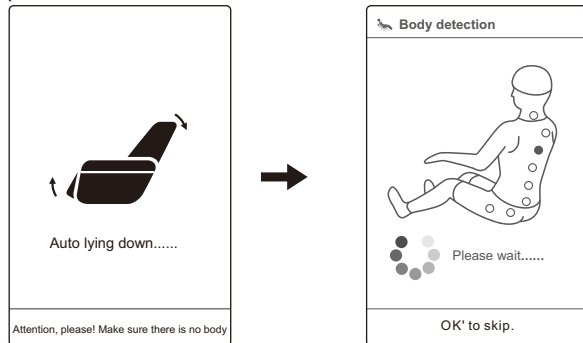
1. Press the remote control ON/OFF key or long press the ON/OFF key on right armrest to turn on massage.
2. Wait for the automatic program to be selected, select the automatic program according to the remote control or voice prompts, and automatically shut down if it has not been operated for more than 20 minutes.
3. Automatic angle adjustment, after selecting the automatic program, the massage chair will automatically adjust to the zero gravity angle, and wait for its adjustment to be completed. You can press the confirm button to skip and go directly to the next step.

Automatic angle adjustment diagram

Body shape detection diagram

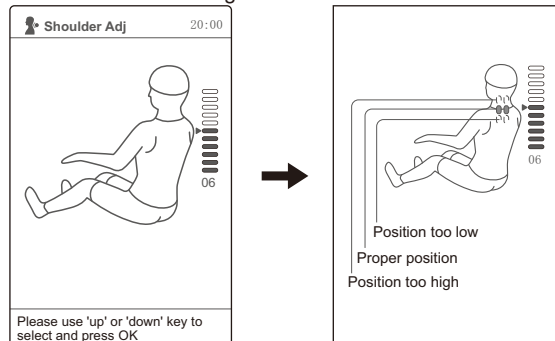
4. Wait for the body shape test to be completed, you can press the confirm key to skip the test and start the massage directly.
5. Shoulder position fine-tuning, after body detection is completed, the remote controller or voice assistant will prompt the user to fine-tune the shoulder position, if the massage wheel stops at your shoulder position, it does not need to be adjusted directly and press the confirmation key to start massage, if not, please use the upper key , or the lower key to adjust the massage wheel to your shoulder position and then press the confirm key to start massage. If more than 20 seconds without operation directly default the current position of the massage wheel for the shoulder position and start massage.

## Operation Instructions

### V. Pause function (please use this function after power on)

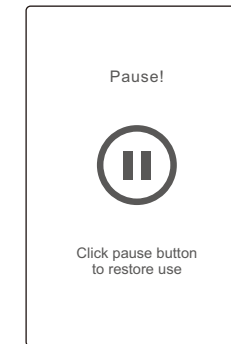
- After the massage chair is turned on, press the pause key or armrest pause key and will enter the pause state, all massage functions stop, press the pause key or armrest pause key again to resume the previous massage state. Automatic shutdown after pausing for more than 20 minutes. Note: In the suspended state, all other functions except shutdown cannot be used.

Message pause state diagram


### VI. Automatic program (please use this function after power on)

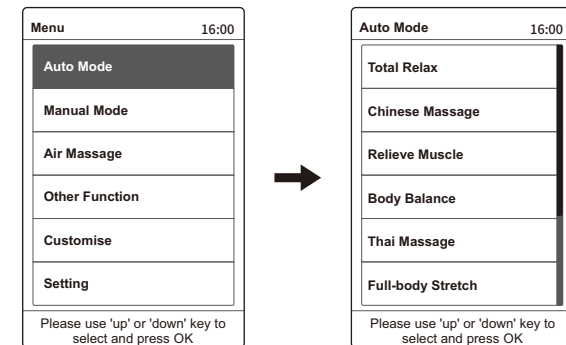1. Select automatic program
  - (1) Press the remote control menu key to enter the main menu.
  - (2) Select the Automatic Program option and press Confirm to enter the Automatic Program menu.

## Operation Instructions

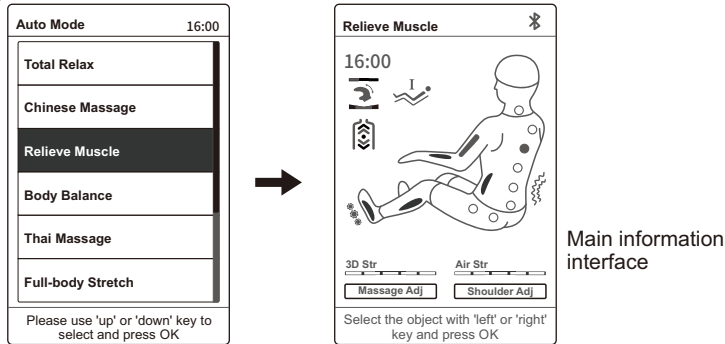(3) Select the desired automatic program and press the confirm key  to start the massage.

### 2. Automatic program introduction

Function	Description
Total Relax	With the aid of flexible guide rails, let the body posture angle closer to lying flat; with the effective anti-arch stretching of 3D movement to the back, achieve zero pressure stretching, effectively relieve lumbar muscle pain, and eliminate fatigue, thus to get a better massage experience, relieve the stress on the body and mind.
Chinese Massage	By means of Traditional Chinese medicine massage technique mainly, apply 3D massage on the back and waist so as to promote the blood circulation of the back and waist, soothe the meridians and relieve the blockage of the back and waist through massage. It is suitable for people with back and waist pain, and lumbar muscle strain. Massage intensity is gentle.
Relieve Muscle	Apply deep pressure on the neck, back and waist, buttocks, and legs by means of kneading, patting, finger pressing techniques mainly to relieve the fatigue on body, allowing the body to better relax.
Body Balance	Apply full air massage mainly, with the use of airbags for cross circulation of kneading and pressure, so as to promote blood circulation, alleviate pain in various body parts, and deeply relieve muscle fatigue. (Remarks: full air massage mainly, with frequent air massage action)
Thai Massage	Mainly use calf and backrest electric push rod to clamp legs, hands and shoulders with the aid of air massage, and then apply stretching massage, while the movement holds against the back and waist for 3D medical massage, so that the body is in the Thai stretching state. It is suitable for people under high stress to improve their head-down on chest and humpback. Massage intensity is strong.

## Operation Instructions

Full-body Stretch	By means of kneading, patting, finger pressing and medical massage techniques mainly, clamp the legs, hands and shoulders with the aid of air massage, and then apply stretching massage to soothe the meridians, promote blood circulation, relieve blockages in the body.
Neck & Shoulder Care	Massage the neck and shoulders by means of kneading and stretching techniques mainly, and also the upper back by medical massage, and then press the Jianjing acupoint to relieve cervical spondylosis and achieve the effect of relaxing the muscles of the neck and shoulders. Massage intensity is gentle, suitable for people working in office and suffering from cervical fatigue.
Spine Traction	Massage the waist by means of kneading, patting and medical massage mainly, while the movement holds against the back and waist for 3D medical massage, so as to relieve lumbar spine pain; combined with the thermostatic heat configuration around the waist, maintain a comfortable temperature to bring warm massage experience; with gentle and moderate massage intensity and comfortable temperature, bring warm massage experience.
Leg&foot Relieve	Massage the shoulders and waist by means of kneading, patting, kneading and patting mainly, apply air massage on the legs and feet for squeezing massage, allowing more soothing leg and foot massage, suitable for sedentary people with low stress on the legs and feet, not only relieving the cervical spine, lumbar spine pain of sedentary people, but also soothing the legs and feet. Massage intensity is gentle.
Superior Enjoyable	Apply air massage on the whole body, with the foot roller and calf kneading massage. Calf stretching function is canceled, suitable for the elderly and people suffering stress on back. The massage intensity is gentle.
President Cosy	Massage the shoulders, neck, back and waist by means of kneading, patting, kneading and patting mainly, to relieve fatigue, promote blood circulation, and release stress by such deep massage. It is suitable for men, and the massage intensity is moderate.
Keep Fit	Massage the shoulders, neck, waist and hips by means of kneading and finger pressing techniques, through the massage on the whole body, to relax the body and mind, and enjoy the queen-level treatment. It is suitable for women, and the massage intensity is gentle.

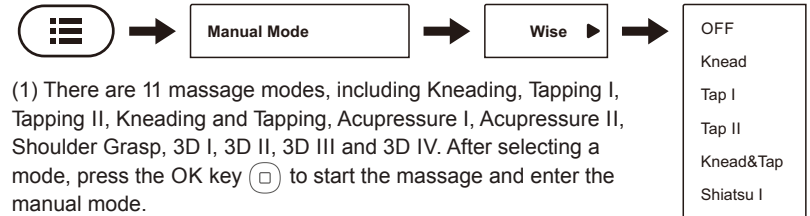
## Operation Instructions

Office Regimen	Massage the shoulders, neck and waist by means of kneading, patting, kneading and patting techniques mainly, to relieve shoulder and neck pain, improve lumbar muscle soreness and maintain lumbar spine through such deep massage. It is suitable for office white-collar and sedentary people, and the massage intensity is moderate.
Health Recharge	Deeply massage the whole body by means of kneading, patting and finger pressing techniques, to promote blood circulation throughout the body and prevent blood clots. It is suitable for people who live at home for a long time and do not like to go out, and the massage intensity is moderately gentle.
Energy Boost	Massage the waist, shoulders and neck by means of kneading, patting, kneading and patting techniques, to awaken the body's vital energy, stimulate physical energy through such deep massage of the whole body with the massage intensity being heavy, which is suitable for people with strong stress, such as athletes, for stimulating the potential before a race.
Lunch Break	Massage the back and waist by means of kneading and patting techniques, and press the Shenyu acupoint, to improve the quality of lunch break, so that the user is more energetic in afternoon work. It is suitable for office workers, and the massage intensity is gentle.
Back Spa	Deeply massage the whole body by pressing related acupoints by means of kneading, patting, finger pressing and other massage techniques, to let the stimulation fully reach the deep layer of muscle tissue, feel warm locally combined with the function of hot compress, and regulate the Qi and blood in a bid to achieve the effect of warmth and damp expelling. (Note: the whole body massage together with hot compress)
Sweet Dream	Massage the shoulders, neck, back and waist by means of kneading and medical massage techniques, and press the acupoints of Fengchi and Shenyu, to improve the sleep quality and help sleep. It is suitable for people suffering from poor sleep, and the massage intensity is gentle. Note: There is background music, which can be turned on or off through the Voice Prompt option in the settings.

## Operation Instructions

### VII. Manual program (please use this function after the power is turned on)

1. Massage modes: as shown in the figure below, open the mode selection menu to select a massage mode.

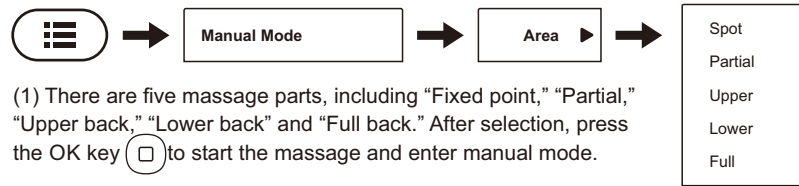
(1) There are 11 massage modes, including Kneading, Tapping I, Tapping II, Kneading and Tapping, Acupressure I, Acupressure II, Shoulder Grasp, 3D I, 3D II, 3D III and 3D IV. After selecting a mode, press the OK key  to start the massage and enter the manual mode.

(2) Massage mode display on the information interface.

Function	Massage roller adjustment	
OFF	Width.	Adjustable speed in five gears Adjustable width in five gears
Knead 	Speed.	
Tap I 	Speed and width.	
Tap II 	Speed and width.	
Knead&Tap 	Speed.	
Shiatsu I 	Speed and width.	
Shiatsu II 	Speed and width.	
Grasp 	Speed.	
3D I <b>3D</b>	Speed.	
3D II <b>3D</b>	Speed.	
3D III <b>3D</b>	Speed and width.	
3D IV <b>3D</b>	Speed and width.	

## Operation Instructions

2. Massaged parts: As shown in the figure below, open the massage part selection menu to select massage parts.

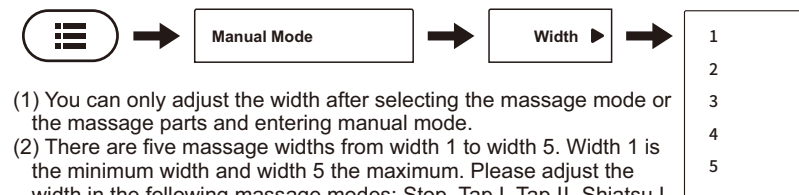
(1) There are five massage parts, including “Fixed point,” “Partial,” “Upper back,” “Lower back” and “Full back.” After selection, press the OK key  to start the massage and enter manual mode.

(2) Massage part display on the information interface.

Function	Description	
 Spot	Fixed-point massage at fixed-point position.	In the selection of spot and local massage techniques, the massage area can be adjusted by pressing “up key” and “down key” in the information interface status.
 Partial	Small-scope local massage back and forth.	
 Upper	Upper back massage back and forth.	
 Lower	Lower back massage back and forth.	
 Full	Full back massage back and forth.	


3. Massage width: as shown in the figure below, open the width selection menu to select a massage width.

(1) You can only adjust the width after selecting the massage mode or the massage parts and entering manual mode.


(2) There are five massage widths from width 1 to width 5. Width 1 is the minimum width and width 5 the maximum. Please adjust the width in the following massage modes: Stop, Tap I, Tap II, Shiatsu I, Shiatsu II, 3D III and 3D IV.

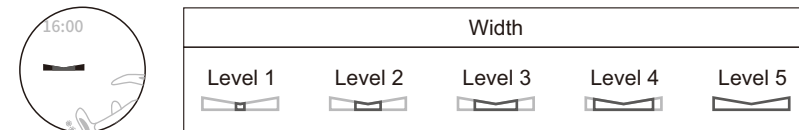(3) Massage part display on the information interface.

## Operation Instructions

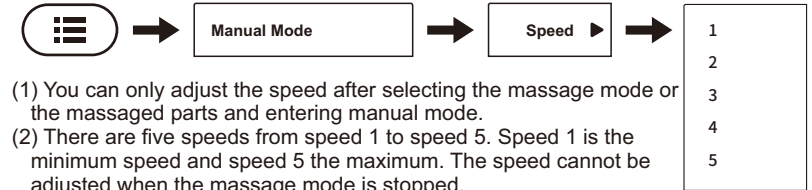
4. Massage speed: as shown in the figure below, open the speed selection menu to select a massage speed.

(1) You can only adjust the speed after selecting the massage mode or the massaged parts and entering manual mode.

(2) There are five speeds from speed 1 to speed 5. Speed 1 is the minimum speed and speed 5 the maximum. The speed cannot be adjusted when the massage mode is stopped.


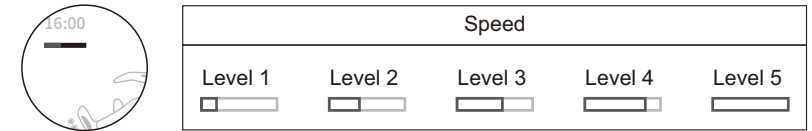
(3) Speed display on the information interface.

## VIII. Massage adjustment (please use this function after the chair is turned on)

### 1. 3D intensity adjustment

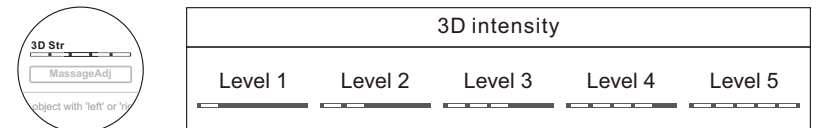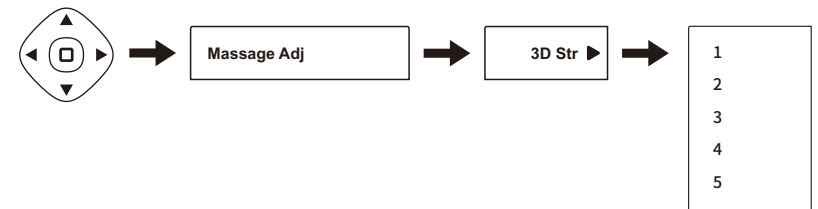
As shown in the figure below, open the 3D intensity adjustment menu to adjust the 3D intensity. There are five levels from Level 1 to Level 5. Level 1 represents the minimum intensity and Level 5 the maximum. The stronger the 3D intensity, the higher the massage roller is lifted and the greater the massage intensity. (In the main information interface)

## Operation Instructions


### 3. Time adjustment

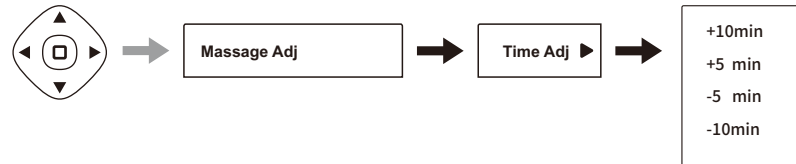As shown in the figure below, open the time adjustment menu to increase or decrease the massage time and the adjustment range is 0–40 minutes. When the massage time is counted down to 0, the massage chair will be turned off automatically, the massage rollers will be reset and the posture of the massage chair will remain unchanged. (In the main information interface)

### 4. Shoulder position adjustment

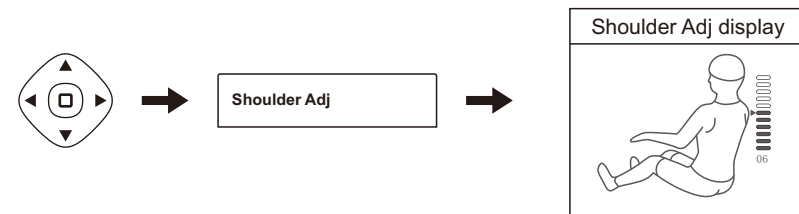
As shown in the figure below, open the shoulder position menu to finely adjust the shoulder positions upwards or downwards. When the mechanical arms press your shoulders, if you feel the positions are inaccurate, you can adjust the shoulder positions. (In the main information interface)

## IX. Air pressure massage (please use the function after the power is turned on)


### 1. Air pressure mode

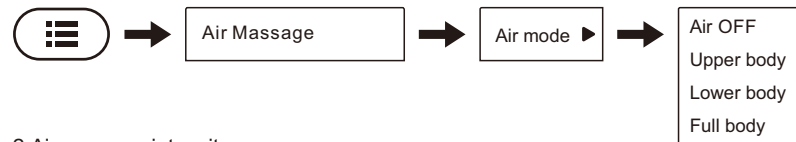Air pressure mode selection can be chosen by entering into air pressure menu or armrest air pressure key  as shown below. Air pressure mode include close, upper body air pressure, lower body air pressure, full body air pressure.


### 2. Air pressure intensity

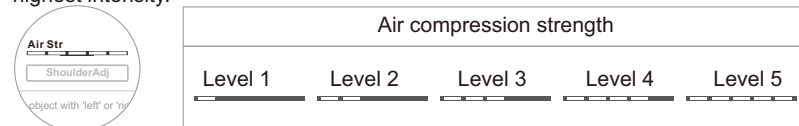Air pressure intensity can be adjusted after opening the air pressure function. There are 5 adjustable levels from 1 to 5, 1 indicates the lowest intensity, 5 indicates the highest intensity.

## Operation Instructions

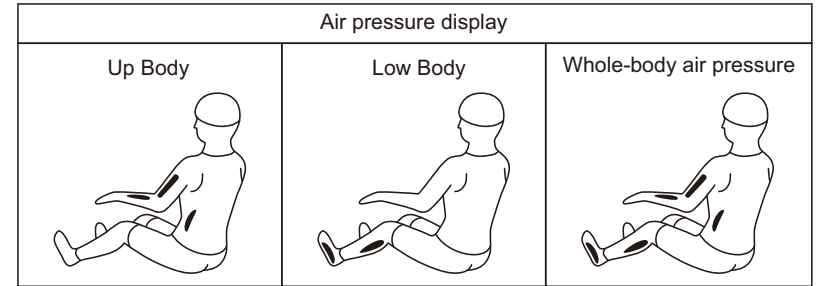
Air pressure display on the information interface.


## X. Auxiliary function (please use this function after powering on)

### 1. Foot roller

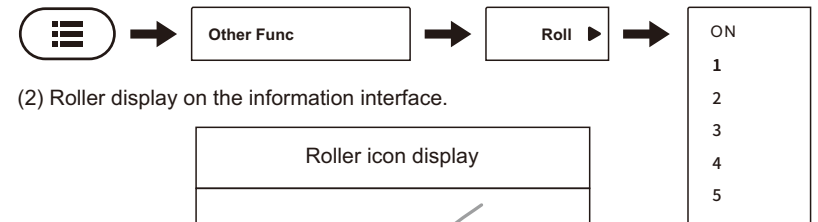As shown in the figure below, enter the foot roller adjustment menu can adjust the strength of the foot roller. There are 5 adjustable levels from 1 to 5, 1 indicates the lowest intensity, 5 indicates the highest intensity. The stronger the strength of the foot roller massage wheel top out the higher the massage intensity. (Under the main information interface)

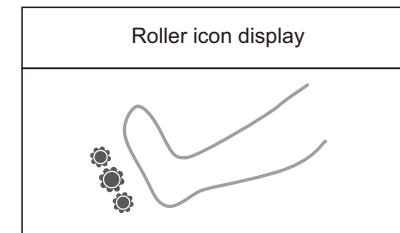(2) Roller display on the information interface.


### 2. Back heating

Back heating intensity can be chosen by entering the auxiliary function menu as showed below. There are 2 adjustable levels. Level 1 means the weakest intensity, level 2 means the strongest intensity.

## Operation Instructions

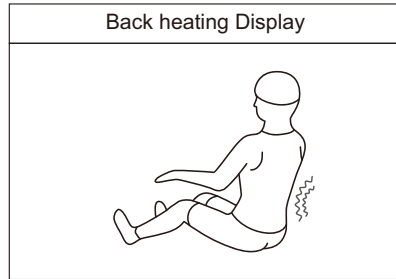Message interface Heating display


### 3. Toe heating

Toe heating intensity can be chosen by entering the auxiliary function menu as showed below. There are 2 adjustable levels. Level 1 means the weakest intensity, level 2 means the strongest intensity.

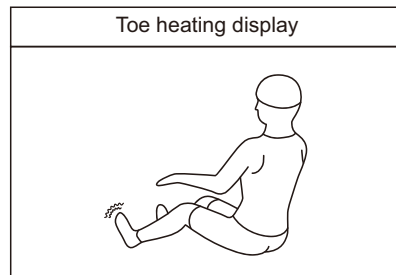
Heating display on information interface

### 4. Led light

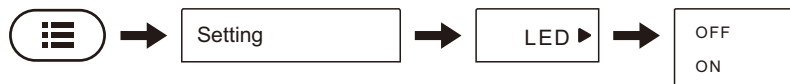Led light can be turned on/ off by entering into settings function menu as shown in the following figure.


## Operation Instructions

### XI. Angle adjustment (please use this function after the chair is turned on)

#### 1. Zero-gravity

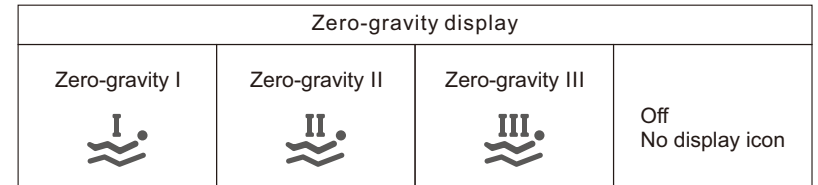
There are three levels corresponding to the three set angles of the massage chair. You can press the Zero-gravity key to switch repeatedly between the three zero-gravity levels or turn off the zero-gravity.

#### 2. Whole body angle

Press and hold the down key on hand controller to adjust the massage chair angle, then the backrest go down and footrest up, release the key will stop adjusting. When the adjustment reaches the limit position, the adjustment will stop and a "tick" tone will be emitted.

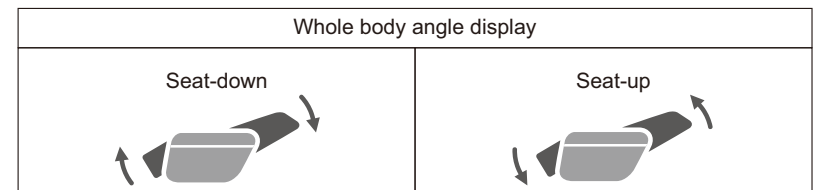
Press and hold the up key on hand controller to adjust the massage chair angle, then the backrest go up and footrest down, release the key will stop adjusting. When the adjustment reaches the limit position, the adjustment will stop and a "tick" tone will be emitted.

#### 3. Footrest angle





Press and hold the down key on hand controller or shortcut Keys on the armrest to adjust the footrest angle down, release the key will stop adjusting. When the adjustment reaches the limit position, the adjustment will stop and a "tick" tone will be emitted.

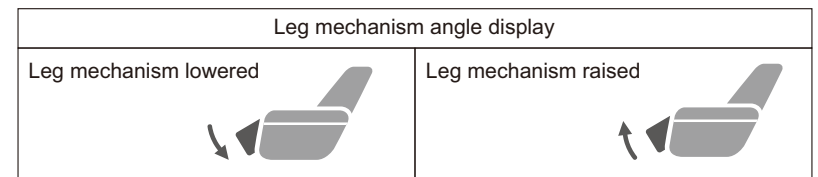
Press and hold the up key on hand controller or shortcut Keys on the armrest to adjust the footrest angle up, release the key will stop adjusting. When the adjustment reaches the limit position, the adjustment will stop and a "tick" tone will be emitted.

## Operation Instructions
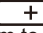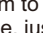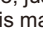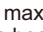
### XII. Memory storage (please use this function after the power on)


As shown below, enter the memory storage menu interface, select the  option and press the OK key  can add the current running massage program to the memory storage. When you use the memory storage interface next time, just select the previous collection of programs and press the OK key  to call this massage program. The same massage program can not be repeatedly added, a maximum of 3 kinds of massage programs can be collected. Select the program has been collected and long press the confirmation button  to delete this collection.

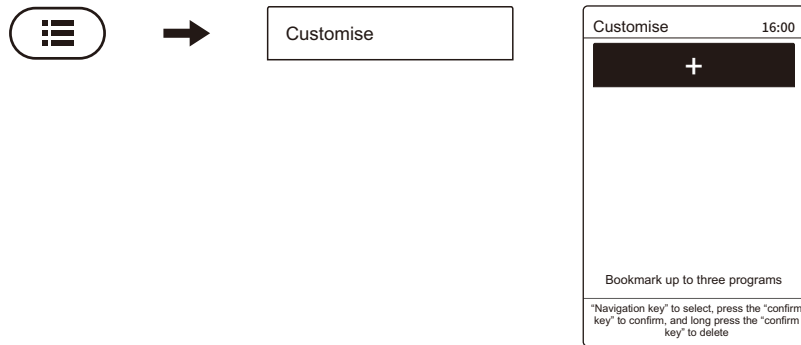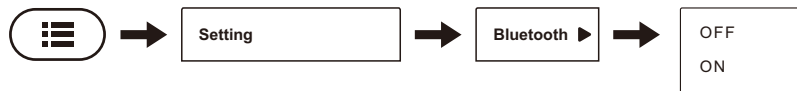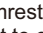
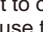
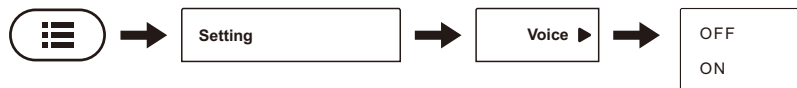
### XIII. System Settings (Please use this function after powering on)

1. Bluetooth: as shown in the figure below, open the Bluetooth Settings menu to enable or disable the Bluetooth speaker function. After the Bluetooth speaker function is enabled, you can play music by connecting to a smart device (mobile phone, tablet, etc.). The connection method is as follows:

- (1) Turn on the Bluetooth function of the smart device, search for the Bluetooth name with the words "IMCM-XXXXXX" and connect.
- (2) Open the music player on the smart device to play music.

2. Voice control function can be turned on or off by entering into settings function menu. Short press voice key  on remote control or the voice key  on armrest to wake up the voice assistant function. Long press the voice key  on armrest to open or close the voice control function. After opening the voice control, you can use the voice assistant function to control the massage chair. Using method is as follows.

## Operation Instructions

You can say "Hi Alice" or "Hey Alice" near the speaker position to wake up the voice assistant, or press the armrest voice key  to wake up the voice assistant, and the voice assistant will reply "I am here" to indicate that it has been woken up. After waking up the voice assistant, please say the command within 6 seconds to achieve the function of the massage chair control. You can continue to say the commands after the voice assistant responds to them. The interval between two adjacent commands can not exceed 6 seconds or you need to re-awaken the voice assistant. Detailed voice command table is as follows.

VOICE CONTROL ENTRY		
Entry: "Hi Alice" or "Hey Alice". Answer: I am here.		
No.	Voice control entry	Answer entry
1	Massage on	Answer: Ok, Massage on
2	Massage close	Answer: Ok, Massage close
3	Total Relax	Answer: Ok, Total Relax
4	Energy Boost	Answer: Ok, Energy Boost
5	Sweet Dream	Answer: Ok, Sweet Dream
6	Neck & Shoulder Care	Answer: Ok, Neck & Shoulder Care
7	Lunch Break	Answer: Ok, Lunch Break
8	Thai Massage	Answer: Ok, Thai Massage
9	Open the air pressure	Answer: Ok, Open the air pressure
10	Close the air pressure	Answer: Ok, Close the air pressure
11	Up the seat position	Answer: Ok, Up the seat position
12	Down the seat position	Answer: Ok, Down the seat position
13	Change the other mode	Answer: Ok, Change the other mode
14	Go little down	Answer: Ok, Go little down
15	Go little up	Answer: Ok, Go little up

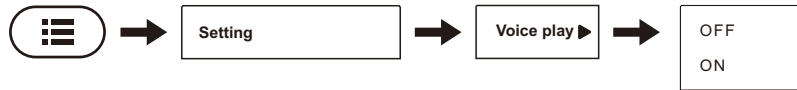
Note: commands 14 and 15 can only be used after command 13.


## Operation Instructions

3.Voice broadcast can be turned on or off by entering into he voice announcement settings menu as shown in below.


Voice prompts: After opening the voice broadcast function, in the process of using the massage chair, the user will be prompted how to use the massage chair correctly, helping new users to become familiar with the massage function faster.

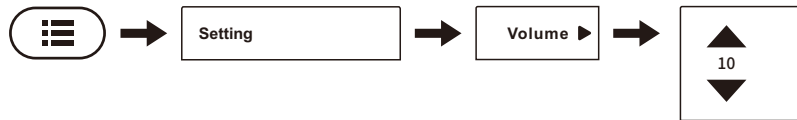4.The system volume can be set by entering into volume settings menu as shown in the following figure. Volume level can be adjusted through pressing the upper or lower key. There are total adjustable levels from 0 to 15.

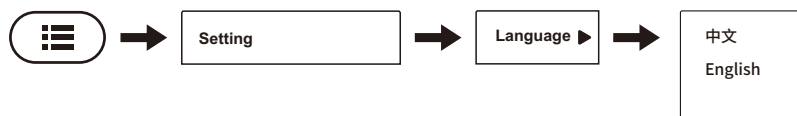
5.Language can be set by entering the language settings menu as shown below. There are simplified Chinese and English to choose from.

6.Serial number

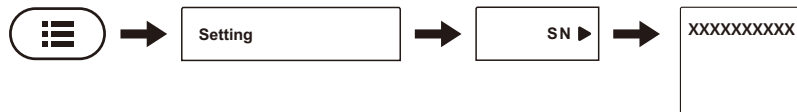You can view the serial number by entering to the Setup function menu as shown below.


## Operation Instructions

### XIV. End massage


During the massage, press on/off key  on hand control or press the on/off key  on armrest to turn off the machine.

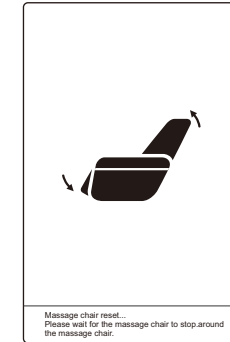Immediately close all massage functions, backrest and footrest will be reset, or massage timing time to, all massage functions will be closed automatically, backrest and footrest is not reset.


Cut off the whole machine power supply, as shown in the figure (Schematic diagram of power switch position of the whole machine).

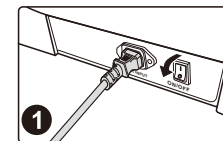
Turn off the power switch ("0" position)

Schematic diagram of pulling out the power plug

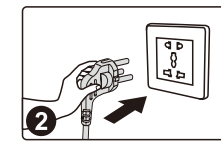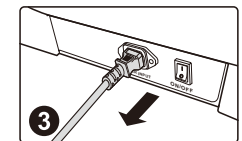Remove the product connection power cord


## Cleaning and Maintenance

### Cleaning synthetic leather and plastic parts

Wipe the leather with a soft, dry cloth.

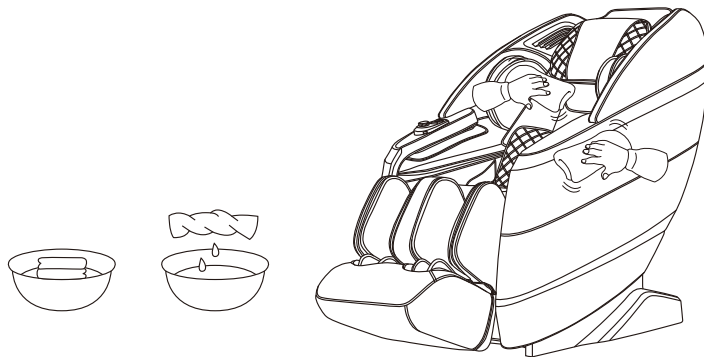- Do not use chemicals such as thinner, gasoline, alcohol, etc. to clean the unit.
- When using commercially available leather care products (rags), please follow the instructions for use.
- If synthetic leather is particularly dirty, wipe it in the following way.
  - ① Immerse a soft cloth in water or hot water containing 3% to 5% neutral detergent, then wring it out thoroughly.
  - ② Use the above soft cloth to wipe the leather surface.
  - ③ Wipe off the cleaning solution with a cloth soaked in clean water and thoroughly wrung out.
  - ④ Wipe with a wrung soft cloth.
  - ⑤ Let it air dry naturally.
- When it is difficult to remove stains, please immerse commercially available melamine foam in neutral detergent, and then wipe the machine with it.
- Do not use a hair dryer to dry the surface.
- Be careful when wiping the machine with denim and colored fabrics, as the color of the fabric will stain the surface of the synthetic leather.
- Do not keep these parts in contact with plastic for a long time, as this will cause discoloration.
- Synthetic leather may be discolored, so when using hair dye, please cover the touched part with a towel.

## Troubleshooting

<b>Fault</b>	When the calf mechanism or the backrest reaches a certain position, press the corresponding adjustment button and it will no longer change, and a continuous "beep, beep" sound will be emitted.
<b>Solution</b>	Due to product structure design and ergonomic requirements, this product has a limited protection design for the limit state of attitude adjustment. At the same time, a reminder sound will sound when the limit is reached.
<b>Fault</b>	The function does not work or some functions do not work after the machine is turned on.
<b>Solution</b>	After the machine is turned on by pressing the switch key, it is necessary to manually select an automatic massage program or other working state. If it is not operated, the machine will automatically shut down after 20 minutes. If the manual operation mode is selected, the functions that are not selected may be disabled, and the working mode needs to be manually set in the manual mode.
<b>Fault</b>	The machine makes a leather rubbing sound when adjusting its posture.
<b>Solution</b>	This product adopts high-grade anti-corrosion leather that meets environmental protection requirements. Due to the inherent characteristics of leather materials, there will be a normal slight noise when rubbing at the relative movement, which is a normal phenomenon. Note: If you hear fabric or leather tearing sound when the machine is running, please stop using it immediately and seek professional maintenance.
<b>Fault</b>	The massage wheel did not reach the shoulders or neck.
<b>Solution</b>	If the head does not touch the head cushion or the back does not touch the backrest, the shoulder position may be lower than the actual position during the body shape detection process. Sit in the deepest position of the seat, rest your head on the headrest, and start over from the beginning.
<b>Fault</b>	The unit is damaged. The power cord or power plug is abnormally hot.
<b>Solution</b>	To prevent accidents, be sure to contact an authorized service center.

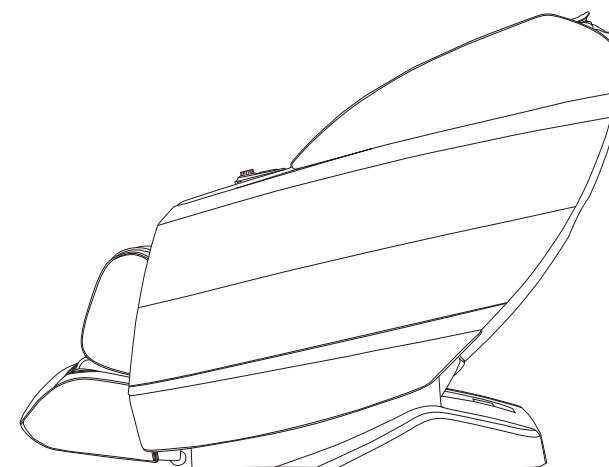
Please feel free to contact the designated service personnel for any questions or doubts about this product.

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