**BLOW**

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**Wave 4000**

## Masażer

Massager

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## DISCLAIMER

**PLEASE DO NOT USE THE MASSAGE GUN OR ANY VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY:**

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device but we advise you to consult a doctor first.


Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above. Practical experience has shown that integration of vibration exercise into a treatment plan is advisory in a number of cases. This must be done on the advice of and in accompaniment with a doctor, specialist or physiotherapist.

## USER MANUAL

For first time using , please must fully charge the battery. Please note that when screen capacity LED show number 99, you still need to keep charging until the number 99 stop flashing.

1. Switch the handle botton Power button ⑧ to ON, then touch the screen Power button ④, the massage gun will start working.
2. After turn on the massager, short press button ⑤ to increase massage speed, press button ③ to reduce massage speed. The LED ① will show the working speed, highest speed is 30.
3. To charge this massage gun, we suggest using a 5V,2A(or higher power) mobile charger and TYPE-C cable. The LED ② will show the left battery capacity, when battery is fully charged, the LED ② will show 99.

## MESSAGE GUN DESCRIPTION
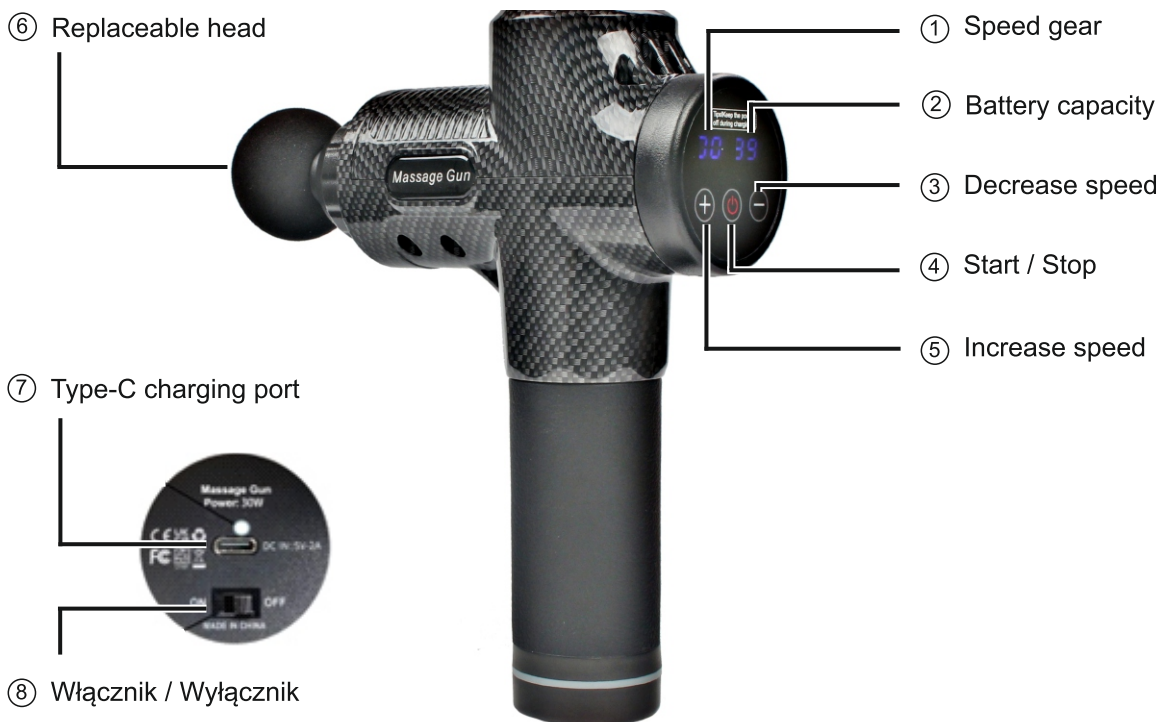
## SPECIFICATION

- 30 speeds - from 1200 to 3800 rpm
- Type-C port DC IN 5V/2A
- power: 100W
- battery: lithium-ion
- stroke amplitude 12mm
- charging time: 60-90 min
- working time: 1-5 hours
- noise level: 45~65dB
- touch LCD display
- charger 5V/2A

## INTERCHANGEABLE HEADS

1. flat head - reaches deep muscles, with a strongly compact structure, such as the back or quadriceps muscles, and other deep tissues.

2. pneumatic (air) head - for a very gentle massage of the most tense muscle groups. Perfect for massaging the bone area. Perfect for a relaxing massage.

3. large arc head - Tip for massaging large muscle groups. Tiny protrusions on the head stimulate circulation, which allows for faster regeneration and better results. This tip allows you to prepare sensitive areas.

4. paddle head - intended for larger muscle groups, isolates pressure on a smaller surface.

5. arc-type head - for massaging tightly packed muscles, thanks to its structure and the protrusions used, it stimulates the muscles very well and prepares them for training. It is also great for massaging deep-seated trigger points or deep muscles.

6. bullet head - the best solution for a strong spot massage. Muscle stimulation in this way effectively eliminates radiating muscle pain.

7. U-shaped head - for massaging the muscles around the spine, the head is designed to massage only the muscles, such as the back extensors, calves or forearms.

8. point head - ideal for gentle massage of very tense muscle groups. The best choice if you need gentle stimulation of fascia and muscles. It provides a very pleasant and relaxing massage.

9. ball head - standard tip, with strong pressure, intended for larger muscle and body parts, e.g. abdomen or arms.

10. cork head - a head with the most flexible structure, intended for massaging muscles around bones,

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The product meets the requirements of the directives:  
2014/30/EU, 2014/35/EU, 2011/65/EU, (EU)2015/863, (EU)2017/2102  
The declaration of conformity is available for download from the website [www.prolech.com.pl](http://www.prolech.com.pl)

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