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GENERAL INFORMATION

Scope of delivery

- 1 x Massage gun
- 11x Massage heads
- 1 x Power supply unit including charging cable
- 1 x Bag
- 1 x Instructions for use

Technical data

Power:	20-30W
Voltage :	24V DC
Frequency:	50-60Hz
Charging time:	ca. 2 hours
Operating time:	ca. 4-6 hours
Speed:	1 100-3500rpm
Battery:	2600mAh
Volume:	45-55Db

Product data

Dimensions:	25.12 x 24.4 x 5.1 cm
Weight:	1.1 kg
Material:	Aluminium + PP + ABS + EVA
Colour:	Black

11 of the attachments shown here are included with the massage gun.

- 1. Flat head**
suitable for large areas of the body, such as hips, legs etc.

- 2. U-shaped head**
suitable for both sides of the spine.

- 3. Lotus head**
suitable for the soles of the feet.

- 4. Round head**
can be used for the shoulder and neck.


- 5. Shovel head**
suitable for the spine.

- 6. Club head**
suitable for large muscle groups, such as hips, legs etc.

- 7. T-head**
suitable for back muscles.

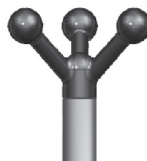- 8. Three-jaw head**
suitable for arms and hands.

- 9. Cylindrical head**
can be used for local pain points or small muscle groups.

- 10. 5-arm attachment**
can be used for muscle protrusions, e.g. calf muscles.

- 11. Mushroom head**
suitable for pectoral muscles, shoulder and neck.

When to use

- Before or after sport
- If experiencing muscle tension
- To relax

How to use the massage gun

1. Ensure that the massage gun battery is charged.
2. Insert the desired attachment.
3. Turn on the device using the LCD touch display.
4. Set the speed and time using the LCD touch display by pressing '+' or '-'.
5. Hold the massage gun by the silicone handle.
6. Massage the area of the body to be treated.

Instructions for battery removal

Only remove the battery when disposing of the device. Before removing the battery, ensure that the device is not plugged in and that the battery is completely empty. Take the necessary precautions when using tools to open the device and when disposing of the battery.

1. Check whether any screws are inserted to the front or rear of the device. If so, unscrew them.
2. Detach the back/front device cover using a screwdriver. Dismantle any other covers and parts until the circuit board containing the battery is visible.
3. Remove the battery.

Warning information

To charge the battery, only use the power supply unit provided with the delivery.

Never cut open the power supply unit in order to attach a different socket as this can be dangerous.

Do not open the device to replace the battery. Do not use the power supply unit in or near wall sockets containing electric air fresheners as these can cause irreparable damage to the power supply unit.

Safety instructions:

Carefully follow these instructions to ensure the correct and safe installation, use and function of the product. Store them in a safe place for future reference.

Always switch off the power supply before carrying out any maintenance or repair work.

Keep all electronic and electrical products away from water and other liquids unless the item is specifically intended for such use.

Keep packaging out of the reach of children.

Warning

Attention, risk of suffocation!
Loose cables present a suffocation hazard for babies and toddlers. Therefore, ensure that all cables are firmly and securely attached. Keep babies, toddlers and pets away from cables!

Disposal guidelines for electronics

 If this symbol is found on the device, please observe your local regulations for the disposal of electrical and electronic products. Please observe the applicable regulations in your country, and do not dispose of disused equipment in domestic waste.

Guarantee

We hereby guarantee that this product is free from material and processing defects. This guarantee is valid for two years if the defective product is returned to Amazon within 30 days.

Intended use

This product is intended for household use only and is not designed for commercial or industrial use.

Disclaimer

We accept no liability for damages which may arise as a result of incorrect or improper use of this product.

Care instructions

Clean the product using a dry cloth. Do not use any abrasives or detergents and avoid contact with liquids.

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Lajd Concept GmbH