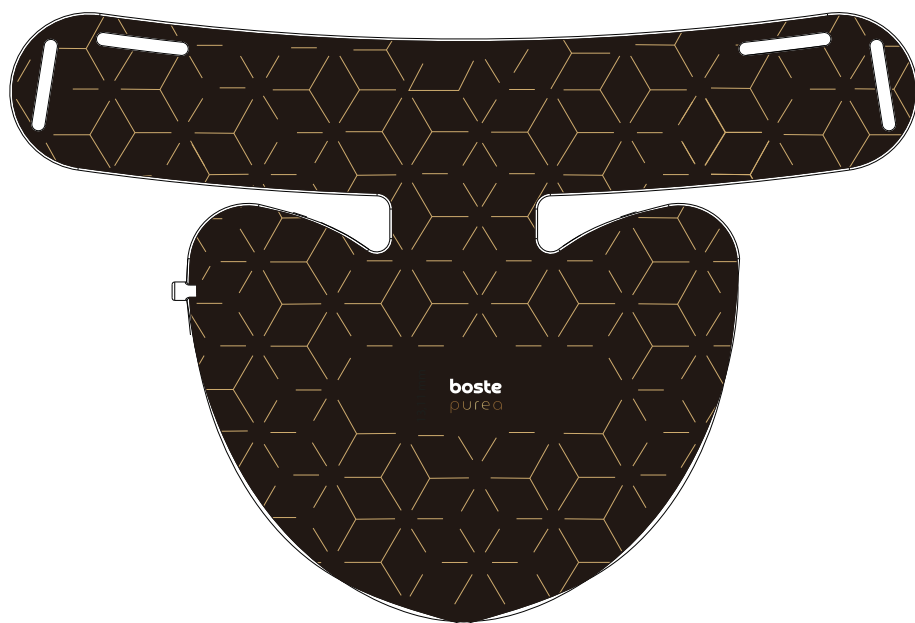
SPECTRA LED Neck Mask

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English

Product Specifications

Product name: Spectra LED Photon Neck
Tender Black
Model: SM2327

Dimensions: 338 × 220 × 47.7 mm

Net weigh: 180 g

Included in the set::

- LED neck mask
- Remote control
- DC Charging Cable
- Adjustable nylon strap

Charging voltage: 5V/1A

Operating current: ≤ 3.5A

Rated voltage: ≤ 4.7V

Battery capacity: 2200 mAh (lithium)

Number of LEDs: 62 groups

Water resistance rating: IPX7

Country of origin: Chiny


Product Description

This is a modern LED neck mask designed for home-use phototherapy skincare. It is lightweight, foldable, and easy to use. The housing is made of flexible medical-grade silicone that adapts to the shape of the neck, ensuring wearing comfort. The adjustable nylon strap allows easy customization for different users.


The neck care model is equipped with 62 LED groups, totaling 248 high-quality LED diodes, each with four chips built into one unit. The LEDs emit light in four different wavelengths:

- Red light (630 nm) – improves skin elasticity and firmness
- Blue light (460 nm) – anti-inflammatory effect, reduces acne and breakouts
- Yellow light (590 nm) – brightens and evens out skin tone
- Infrared (850 nm) – penetrates deeply to support regeneration


Using the included controller, the user can select from three intensity levels and adjust session duration based on personal needs. Light therapy based on these wavelengths supports the skin's natural regenerative processes, stimulates collagen and elastin production, and helps improve the overall appearance and condition of the neck.

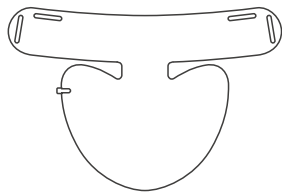
LED neck mask

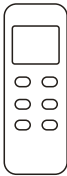Remote control

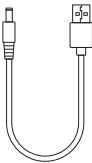DC Charging Cable

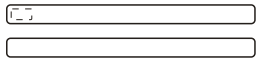Adjustable nylon buckle

Light Modes

-  Red Light (630 nm): Level 1, 2, 3
-  Blue Light (460 nm): Level 1, 2, 3
-  Yellow Light (590 nm): Level 1, 2, 3
-  Mixed Light (Combination): Level 1, 2, 3

Infrared light can be controlled independently or used in combination with other light modes.

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Red Light – 630 nm

Anti-aging effect, supports collagen production, firms the skin, and improves elasticity.

Blue Light – 460 nm

Sterilizes the skin, eliminates acne-causing bacteria, reduces inflammation and scars, and increases skin oxygenation.

Yellow Light – 590 nm

Evens skin tone, boosts blood circulation, reduces redness, puffiness, and fine lines.

Infrared – 850 nm

Regenerates the skin, improves blood circulation, enhances hydration, and supports deep skin repair.

Important Information

- Before first use, make sure to fully charge the mask battery.
- During charging, the LED indicator on the controller will start blinking. Once charging is complete, the display will stay lit.
Estimated charging time is approximately 4 hours.
- The adjustable Velcro strap allows for a comfortable and secure fit.
To properly wear the mask, thread the strap through the openings on both sides, wrap it around your neck, and fasten it securely.
- The operating modes and treatment effects of the neck device are identical to those of the face mask.

Device Controller

Device Operation

- During charging:** The power icon signals this with a "breathing light" (blinking). After full charging, the indicator will stay steadily lit.
When the battery voltage drops to 3.0V, the indicator will blink and the device will shut down automatically after approx. 5 seconds, with a beep sound reminding the user to recharge.
- Press and hold the power button for 2.5 seconds:** to turn the device on or off. When turned on, the red indicator (R) will appear. By default, the device starts at power level 1.

Press the E button to change the light intensity level: level 2 and level 3 are brighter.
Short press the M button to change the light color mode – e.g., yellow light (Y).
After selecting yellow, the device returns to level 1.
Press the E button again to switch to level 2 or level 3.
Then you can press M again to activate blue light (B), or press it again to activate a combination of three light colors (G).
Press M once more to switch to independent infrared mode (NIR) – all other lights will be turned off and only the infrared function will remain.
- When infrared (NIR) light is off press the NIR button to activate it.**
Press the button again to turn it off. The infrared function can be used alone or combined with other light modes.
- To turn off the device, press and hold the power button for 2.5 seconds.**
- When connecting or disconnecting the charging cable, a buzzer will emit a prompt sound (beep).
- The device can remain on while charging. After use, press and hold the power button – the device will shut down and automatically switch to charging mode.

Precautions

- 1. Use the device strictly according to the manual.**
Improper use may lead to equipment damage or unexpected reactions.
- 2. Do not use the device for purposes other than those intended.**
This product is designed exclusively for cosmetic skincare treatments.
- 3. Not suitable for pregnant or breastfeeding women, or children.**
- 4. Avoid using if you have skin conditions or other serious health issues.**
If you suffer from autoimmune diseases, light sensitivity (photosensitivity), lupus, eczema, albinism, or neurological/mental disorders — the use of this device may be contraindicated. Always consult a doctor (e.g. dermatologist) before use.
- 5. Do not use if you have light sensitivity (photosensitivity).**
It may cause severe skin reactions in individuals with this condition.
- 6. Always protect your eyes during use.**
Use protective goggles and keep your eyes closed. If you have a genetic eye disease — do not use LED masks.
- 7. If you feel any discomfort — stop using the device immediately.**
- 8. The device is not waterproof.**
Do not immerse in water or rinse. If water contact is suspected — do not use.
- 9. Keep out of reach of children and infants.**
Store safely.
- 10. Do not use if the product is damaged.**
Damaged components can cause injury or malfunction.
- 11. For home use only.**
Do not disassemble, modify, or repair the device yourself — it may lead to malfunction or damage.
- 12. Recommended usage frequency:**
3–4 times per week, for 15–25 minutes per session.

Cleaning and Maintenance

- After each use, wipe the surface of the mask with a dry, soft cloth (e.g. microfiber) to remove sweat, dust, and residue.
- Do not use water, detergents, alcohol, or chemicals — they may damage the silicone and LED surface.
- Do not immerse the controller or mask in water. The product is not waterproof.
- Store the product in a cool, dry place, away from heat sources, sunlight, or moisture.
- Do not bend the mask sharply — it may damage internal LED circuits.
- Do not pull on cables or press hard on buttons or display components.
- If the device is not in use for a longer time — recharge the battery at least once every 1–2 months to avoid deep discharge.

Frequently Asked Questions

Can I use the mask every day?

– It is recommended to use it 3–4 times a week, for 15–25 minutes per session. Using it too frequently may irritate the skin.

Can I use the mask after applying cosmetics?

– Do not use the mask directly after applying products containing acids, retinol, or other photosensitizing ingredients. Wait a few hours or perform the session on clean skin.

Can I use the mask around the eyes?

– It is recommended to close your eyes during the session and use the included protective goggles. Do not shine the light directly into your eyes.

Is the device waterproof?

– No. The mask and controller are not water-resistant. Do not immerse in water or rinse under running water.

How long does charging take?

– A full charge takes approximately 4 hours. The LED indicator will stop blinking when the battery is fully charged.

Can pregnant or breastfeeding women use the mask?

– Use is not recommended for pregnant or breastfeeding women without prior consultation with a doctor.

Can children use the mask?

– No. The device is intended for adult use only.

What do the light colors mean?

– Each color has a different effect: Red – anti-aging, Blue – anti-acne, Yellow – brightening, Infrared – regenerating

This product complies with the requirements of the following EU directives:

Directive 2014/30/EU – Electromagnetic compatibility

RoHS Directive 2011/65/EU – Restriction of hazardous substances

The crossed-out wheeled bin symbol means

that within the European Union, used electrical and electronic equipment should not be disposed of with household waste.

It should be taken to appropriate collection points for recycling.

THE BOSS IS HERE.

Magnum RF Sp. z o.o.
ul. Wojska Polskiego 71
18-300 Zambrów Poland
NIP:7231638410
REGON:386610073
KRS:0000894731

www.boste.eu

Deutsch

Produktparameter

Produktname: Spectra LED Photon Neck –

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