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

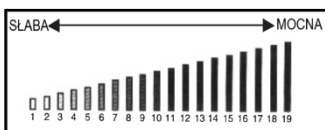
Product features: Foot electrostimulator. The set includes: mat, USB cable, LCD display. Number of massage modes: 8. Number of massage intensity levels: 19. Mat dimensions: 29 x 27 cm. Voltage: 3.7V. Frequency: 1~100Hz. Working time: 15 min. Output power: max. 130 mA. Battery capacity: 120 mAh.

How to use: Electromuscular stimulation (EMS), also called neuromuscular electrical stimulation, is the use of electrical impulses to induce repeated muscle contractions.

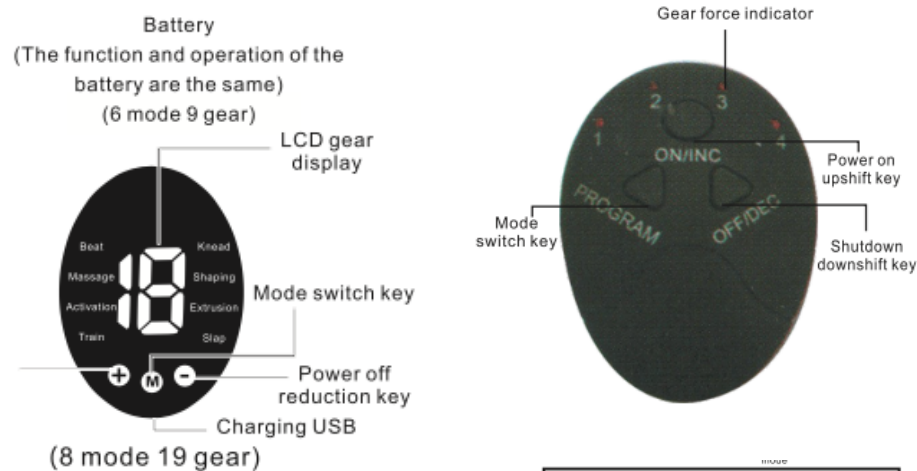EMS has proven effectiveness as a tool for active muscle strengthening. Uses: muscle strengthening, pain relief, preventing atrophy of unused muscles, muscle relaxation, muscle conditioning for beauty. First use and charging: Before use, the device must be charged and disconnected from the power source. The red light is on while charging, and when the device is fully charged, the light will turn off. Charging time is 1-2 hours and after that the device will work for approximately 160 minutes. The higher the power setting of the device, the higher the energy consumption. Electrostimulation: Before using the electrostimulator, remove your socks. Place the electrostimulation mat on the floor. Place both feet on the device and press the power button. Taking your foot off the device will stop the operation. The electrostimulator has 19 intensity levels and 8 modes. It is recommended to start the stimulation at the lowest intensity level and, after getting used to it, slowly move to a higher intensity level, then to the next one, and so on. The device will automatically turn off 15 minutes after switching to intensity level 1-19. If the machine is set to level 1-19 and is not in use, it will automatically switch to gear 0 within 5 seconds. While the battery level gradually decreases, the operating intensity of the device will also gradually decrease. Before further use of the device, please charge the battery.

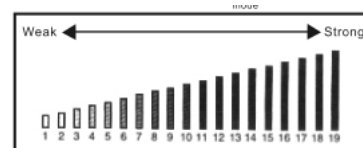Problem solving		
Problem	Reason	Solution
device cannot be turned on	Discharged battery.	Charge your device.
	Ambient temperature too high or too low.	recommended ambient temperature (15~35°C).
Electrical stimulation can not be felt.	Power is not turned on	Switch button until the LED light comes on.
stimulation is too weak.	intensity is set to low.	Increase the intensity of the device's operation.
stimulation is too high.	intensity is set to high.	Reduce the intensity of operation of the device.
Sudden interruption of electrical stimulation	remaining battery power has been exhausted.	Charge your device..

Information for users on the disposal of electrical and electronic equipment (for households)

The symbol shown on the products, or the accompanying documentation indicates that faulty electrical or electronic equipment cannot be disposed of together with household waste. The correct procedure, if necessary, to dispose of electrical or electronic equipment, to dispose of, reuse or recover components, is to transfer the equipment to a specialized collection point, where it will be accepted free of charge. In some countries, you can give the product back to your local distributor when you purchase another device. Proper disposal of the equipment preserves valuable resources and avoids negative health and environmental impacts that may be compromised by inadequate waste management. Please contact your local authority for details on the nearest collection point. Improper disposal of waste is punishable by penalties provided for in the relevant local legislation.

Users in European Union countries: If you need to dispose of electrical or electronic equipment, please contact your nearest point of sale or your supplier for further information. Disposal of waste in countries outside the European Union: This symbol is only valid in the European Union. If you need to dispose of this product, please contact your local authorities or dealer for information on how to proceed correctly.

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