RY100B/C

/ PL / BEZPRZEWODOWY IRYGATOR DO ZĘBÓW **/ EN /** WIRELESS ORAL IRRIGATOR **/ CZ /** AKUMULÁTOROVÁ ÚSTNÍ SPRCHA **/ DE /** DRAHTLOSE MUNDSDUSCHE **/ RU /** БЕСПРОВОДНОЙ ИРРИГАТОР ДЛЯ ПОЛОСТИ РТА **/ SK /** BEZDRÔTOVÁ ÚSTNA SPRCHA **/ HU /** VEZETÉK NÉLKÜLI SZÁJZUHANY **/ ES /** IRRIGADOR DE DIENTES INALÁMBRICO

Eldom Sp. z o.o. • ul. Pawła Chromika 5a • 40-238 Katowice, POLAND
tel: +48 32 2553340 • fax: +48 32 2530412 • www.eldom.eu

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/EN/ Disposal of used electrical and electronic equipment (applies to the European Union countries and other European countries with separate waste collection systems).

The symbol on the product or its packaging indicates that the product should not be classified as household waste. It should be handed over to a collection point that deals with the collection and recycling of electrical and electronic devices. Correct disposal of the product prevents potential negative consequences for the environment and human health that could be caused by improper disposal of the product. Material recycling helps to preserve natural resources. For detailed information on recycling of this product, please contact your local authority, a cleaning service provider or a shop where the product has been bought.

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SAFETY INSTRUCTIONS

- Before using the appliance for the first time, please read carefully these Instructions for Use.
- The device can be used by children aged 8 and older and by persons with limited physical, sensory or mental capabilities, or persons without sufficient experience and knowledge only under supervision or if previously instructed on the safe use of the appliance and the possible risks. Children shall not play with the appliance. Cleaning and maintenance that are the responsibility of the user shall not be carried out by children without supervision.
- Do not place the appliance in the vicinity of heat sources.
- The appliance is intended for domestic use only.
- Do not immerse the cord or the appliance in water.
- The appliance must always be cleaned after use. The appliance must always be turned off before cleaning and maintenance as well as in case of malfunction and after operation.
- Do not use with a broken or otherwise damaged hairdressing head as it may lead to injury.
- Do not use in the open air.
- Use only with the originally included accessories and feeder.
- When in operation, do not leave the appliance unattended.
- when charging, do not place the appliance where it can fall into a bath, a washbasin or a shower tray.
- do not touch the appliance if it has fallen into one of the above mentioned fixtures containing liquid while charging or with the plug inserted. If it were to happen, immediately disconnect the appliance from the socket.
- The appliance can be repaired only in an authorised Service Centre. Any modernisation or application of spare parts or elements other than the original ones is forbidden and dangerous for the user.
- The company Eldom Sp. z o.o. shall not be liable for any damage caused by inappropriate use of the appliance.

SPECIFIC RECOMMENDATIONS

- Do not fill the appliance with water warmer than 50°C.
- Use the appliance as instructed in the manual or as directed by a dentist.
- Do not point the water jet directly under the tongue, inside the ear, or to other sensitive areas. The pressure produced inside the appliance may cause damage to such areas.
- People with pacemakers should consult with a specialised physician before using the appliance.
- Do not use mouthwashes/preparations containing iodine, bleaches, or tea tree oil.
- Remove any jewellery from the mouth before use.
- If you have received dental treatment within the

last 2 months, consult your dentist before using the irrigator.


- If you experience bleeding after using the irrigator, consult your dentist. Do not use if there are any wounds found in the mouth.
- Do not use the appliance each time for more than 5 minutes and more often than every 2 hours.
- Charge the appliance by only using the charger with the parameters indicated in the manual.
- Do not turn on the appliance without the attached tip (1).
- Do not press the tip (1) too hard against the gums, teeth, or tongue.

ENVIRONMENTAL PROTECTION

- the appliance is made of materials which can be recovered or recycled
- the appliance should be handed over to the appropriate centre of collection and recycling of electrical and electronic equipment.

WARRANTY

- the appliance is intended for personal use at a household,
- it is not intended for professional purposes,
- in case of misuse, the warranty is rendered null and void.

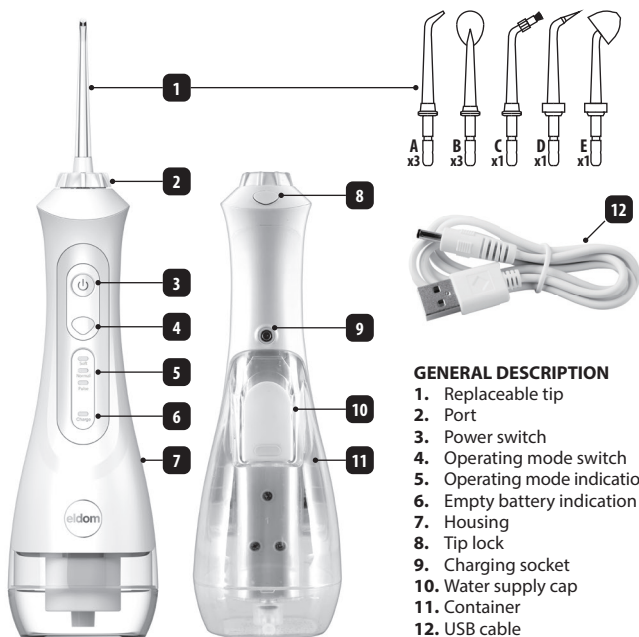GENERAL DESCRIPTION

1. Replaceable tip
2. Port
3. Power switch
4. Operating mode switch
5. Operating mode indication
6. Empty battery indication
7. Housing
8. Tip lock
9. Charging socket
10. Water supply cap
11. Container
12. USB cable

TECHNICAL DATA

Power: 7.5W

Pressure: 80 - 160 PSI (5.5 - 11 Bar)

Pulsation: 1500-2000 pulsations per minute

Volume: 200 ml

Battery: DC 3.7 V, 1500 mAh, Li-Ion

Charging time: approx. 6 hours.

Operating time: up to 14 days
(twice a day)

Irrigator water resistance: IPX7

Charging method: cable

Connect the cable (6) to a charger with the following parameters: DC 5 V, 0.6 mA
(not included)

INTENDED USE

The irrigator is intended for daily oral hygiene. The different tips enable you to clean both the interdental spaces, hard-to-reach places, gaps in between the braces, and the tongue surface. The appliance is also used for massaging the gums and for clearing the nose.

Please note that the irrigator cannot replace a toothbrush and toothpaste when brushing your teeth. It is meant to compliment them.

It is worth remembering that using the appliance regularly can reduce the development of gum diseases and significantly improve oral hygiene.

APPLIANCE OPERATION

Read the entire instructions manual carefully before using the vacuum cleaner for the first time. After unpacking the appliance, make sure that it has not been damaged during transport. In case of any doubts, do not use the appliance until after contacting the service station. **Before first use, rinse the appliance with warm water.**

CHARGING THE DEVICE

- The device may only be used after it has been charged.
- Switch off the device before charging.
- Plug one end of the cable (12) into the socket (9) and the other end into the power supply unit (not included) with the parameters provided above, and then into the mains.
- When charging, the light (6) will be on and will go off when the device is charged.
- The device cannot be used while charging.
- When charging, the light (6) will light up red, and when charged, green.
- Do not use the appliance when charging.

NOTE! In case of a fully discharged battery, the light may light up (6) after approx. 5 minutes. This is normal.

TYPES OF TIPS AND THEIR USE

Replace the tips (1) every 3 months for hygienic reasons, and whenever they get damaged.

- A** standard tip for irrigation.
- B** tongue cleaning tip. Place the tip in the middle of the tongue; press gently and then switch on the irrigator. Start the cleaning with the lowest level and gradually increase it, adjusting it to the individual needs.
- C** orthodontic tip. Designed for cleaning hard-to-reach places between braces.
- D** periodontic tip. Designed for cleaning periodontal pockets. It is also used to apply medicinal fluids between the gum line and the periodontal pocket. Use the tip only with the lowest pressure. To use the tip features, place it on the tooth at an angle of about 45 degrees and move to the tooth root (the place where the gums and the tooth meets).
- E** nose clearing tip. Place the tip onto one nostril and bend over the washbasin. Choose "**SOFT** or **PULSE**" mode. Then, switch on the irrigator pumping the irrigation liquid inside the nostril. Let the fluid out with the second nostril. Repeat for the other side. Clean the nose between the two cleaning procedures. Remove any retaining liquid from the nose after finishing the irrigation.

Tips (1) are consumables and are not subject to a warranty repair.

Replace the tips (1) with only the original tips available in the set or purchased separately in **RYZ101** sets (standard tips - 2 pcs.) or **RYZ102** (9 tips as in this appliance).

USE

- open the water supply cap (10),
- fill the tank (11) with water or a special liquid - use dental fluid in 1:1 ratio; use other preparations according to the manufacturer's recommendations,
- close the water supply cap (10),

NOTE! Fill the appliance with water (liquids) not exceeding 50°C.

- attach the appropriate tip (1) into the port (2),
- select the operating mode using switch (4),
SOFT - cleans sensitive teeth and gums,
NORMAL - for everyday use,
PULSE - for gum massage,
- place the tip (1) in the mouth and turn on the appliance with the switch (3),
- to replace the tip, press and hold the button (8), then gently pull out the tip (1),

NOTE! The tip can only be replaced when the irrigator is switched off.

- to turn the appliance off, press the power switch (3) again,
- after the irrigation is complete, fill the tank (11) with warm water and rinse the appliance - the tip (1) should be pointed at either the washbasin, bath, or shower tray. Then, turn on the irrigator and wait until the tank is completely empty.

NOTE! Rinsing the irrigator with warm water is particularly important when using a special liquid for irrigation.

TIPS

- hold the appliance so that the tip is facing upwards,
- the mouth should be open so that water can flow freely into the washbasin,
- place the tip **(1)** on the teeth at a 90 degree angle and move it clockwise,
- wash sensitive areas between the teeth carefully,
- move the water jet along the gum line - this will effectively rinse any food residue
- from both the surface and between the teeth, while providing a gentle massage of the gums.

Problem	Cause	Possible solution
Leakage from the container	Filler is not closed (10)	Close the filler (10) correctly
	Internal gasket is defective	Contact the manufacturer's maintenance service
The appliance does not work	Low battery	Charge the battery
Abnormal pressure	The tip (1) is clogged by limescale	Replace the tip (1)
Leakage around the port (2)	Incorrectly mounted tip (1)	Remove and properly attach the tip (1) into the port (2)

CLEANING AND MAINTENANCE

Clean the appliance after each use.

- Before cleaning, make sure the device is turned off.
- Do not immerse the cable **(12)** in water.
- Do not use strong agents, cleaning agents or cleaning products which may damage the surface.
- Wipe the cable **(12)** and housing **(7)** with a damp cloth.
- Rinse the tip **(1)** thoroughly.
- Remember to descale the appliance at least once every 3 months. To do this, use only descaling agents for teeth irrigators.
- Do not clean the unit in water, e.g. by submerging it or placing it under a stream of water, as this may cause irreparable damage to the appliance.

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