VEVOR®
TOUGH TOOLS, HALF PRICE
Technical Support and E-Warranty Certificate
www.vevor.com/support

Red Light Face Massager Model: RL05-B

We continue to be committed to provide you tools with competitive price. "Save Half", "Half Price" or any other similar expressions used by us only represents an estimate of savings you might benefit from buying certain tools with us compared to the major top brands and does not necessarily mean to cover all categories of tools offered by us. You are kindly reminded to verify carefully when you are placing an order with us if you are actually saving half in comparison with the top major brands.

VEVOR®
TOUGH TOOLS, HALF PRICE

Red Light Face Massager

Model: RL05-B

NEED HELP? CONTACT US!

Have product questions? Need technical support?
Please feel free to contact us:

Technical Support and E-Warranty Certificate
www.vevor.com/support

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

Product Parameters

Material: Aluminum alloy+PC

Input Voltage: 3.7V

Power: 1.67W

Net Weight: 60g

Wavelength: 630nm+ Other 6 kinds

Battery Capacity: 350 mAh


Charging Time: 90~120mins

Power Supply: USB charging

Product size: 146x43x12 mm

Product Features

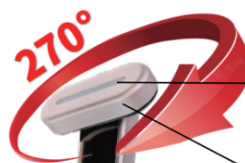Red Light
620nm-630nm


Inlet Surface
270°rotating head
touch to start,Smart
power saving

Body
Aluminum alloy

Power Supply
USB charging port

How It Works ?

The 630nm wavelength is ideal for targeting all manner of skin concerns. In one study, the 630nm wavelength was combined with other beauty treatments such as microneedling to create a statistically significant improvement in photodamage, fine lines and sagging.

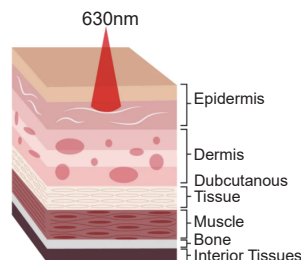Product Functions

The steps on using beauty device

Step 1

Before use, wipe with a clean cloth or damp tissue (30%humidity), please don't put it water clean!

Step 2

Firstly, clean your face, then apply water lotion to keep your face or neck moist.

Step 3

Apply an appropriate amount of eye cream or serum around the eyes to keep the skin moist.

Step 4

Gently slide the massage head, lift and massage the eye area, repeat 8-10 times. (Note: When holding device the massage head touches the skin, it starts to work, and stops when it not touch.)

Step 5

After use, wipe the massage head with a clean cloth or damp tissue (30% humidity). Don't put it water clean.

CORRECT DISPOSAL

This product is subject to the provision of european Directive 2012/19/EU. The symbol showing a wheelie bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to acollection point for recycling electrical and electronic devices.

BATTERY DISPOSAL

To minimize hazards to health and the environment at the end of this product's life, laws dealing with Waste Electrical and Electronic Equipment (WEEE) and The Waste Battery Directive require you to dispose of this product at a suitable collection facility where it will be sent in order to remove the batteries and for appropriate recycling. Please contact your local authorities for more details on recycling and safe disposal of these in your area.

FCC Information:

CAUTION: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment!

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:


- 1) This product may cause harmful interference.
- 2) This product must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

Note: This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules, These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the distance between the product and receiver.
- Connect the product to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for assistance.

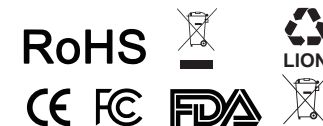Charging Method: Insert one end of the USB charging cable into the charging port at the bottom of the eye beauty instrument, and connect the other end to the charger or computer USB port. The indicator light is in the state of breathing and flashing. It is recommended that the charging time should not exceed 24 hours.

Charging Time: 100mins fully charge, can be use for 60mins if continuous use.

About Charging Power Reminder

- When the battery is low, the indicator light flashes quickly
- When charging, the indicator light is breathing and flashing
- When fully charged, the indicator light is always on

Warm Notices

- Do not disassemble and modify this product without permission.
- Do not select the guide head in the opposite direction.
- When using it around the eyes, there is a light flickering sensation, which is a normal biological micro-electricity effect, please do not worry.
- Due to the deviation of the resistance value of the human body, individual people have a strong electric shock when using it. Please choose whether to continue using it according to the individual's tolerance level. If there is any discomfort, please stop using it immediately.
- If you find any abnormality such as vibration during use, please check whether the battery is low or shut down.
- When not in use for a long time, please turn off the machine and store it in a dry place to avoid sunlight exposure.
- Please place this product out of the reach of children. This product cannot be used as a toy.

VEVOR®
TOUGH TOOLS, HALF PRICE
Technical Support and E-Warranty Certificate
www.vevor.com/support

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