
Little Doctor

LD Special
LD SteTime


ENG STETHOSCOPE LD. *Instruction Manual.*
POL STETOSKOP LD. *Instrukcja obsługi.*
HUN LD SZTETOSZKÓP. *Használati útmutató.*
ROU STETOSCOPI LD. *Manual de instrucțiuni.*
BGR СТЕТОСКОП LD. *Ръководство за експлоатация.*

FIGURES / ILUSTRACJE / RAJZOK / FIGURI / ФИГУРИ


LD SPECIAL

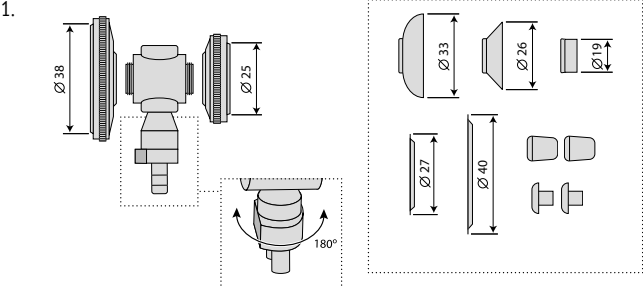
LD STETIME

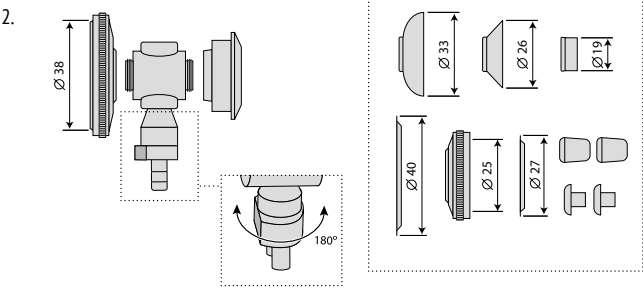Change of power cell / Wymiana baterii / Az elemek cseréje /
Înlocuirea bateriei / Смяна на батерията

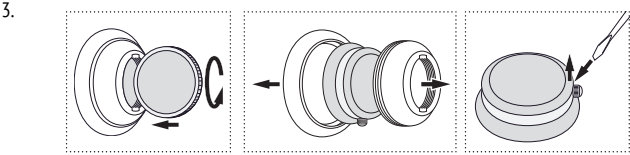Setting the clock / Ustawianie zegara / Az idő beállításá /
/ Setarea ceasului / Настройка на часовника

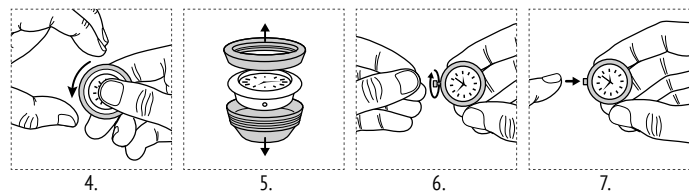ENG

LD SPECIAL

Multifunctional high-quality stethoscope (type SPRAGUE RAPPAPORT) is designed for auscultation of heart, lungs, vessels, gastrointestinal tract of adults and children.

Special completing set provides for 4 working combinations of stethoscope from ordinary to neonatal.

Complete set:

- Twin metal head (LD-S018), made by special design, with big and small (children) diaphragms – 1 pc. Working diameter of big diaphragm – 38 mm, working diameter of small diaphragm – 25 mm. Each diaphragm is provided with transparent plastic membrane. Design feature of head makes it possible to change any diaphragm for bell of needed size. Diaphragms and bells are attached to the head with the help of thread. Changing over between working sides is carried out by turn of head round rod valve by 180°. Flat side of rod valve points out to working side of head.
- Two thickened sound-conducting tubes of 56 or 72 cm long (indicated in box) – 1 set.
- Binaural metal tubes with spring and ear tips – 1 set.

Additional accessories entering the set:

- Plastic big bell (diameter 33 mm) – 1 pc.
- Plastic small bell (diameter 26 mm) – 1 pc.
- Plastic neonatal bell (diameter 19 mm) – 1 pc.
- Big membrane (diameter 40 mm) – 1 pc.
- Small membrane (diameter 27 mm) – 1 pc.
- Plastic ear tips – 2 sets.
- Plastic bag for additional accessories – 1 pc.

Size with package, mm: 300 x 130 x 40.

Weight (with package), g: 380.

LD STETIME

Multifunctional high-quality stethoscope (typ SPRAGUE RAPPAPORT) with quartz clock built-in into head is designed for: for auscultation of heart, lungs, vessels, gastrointestinal tract of adults and children.

Special completing set provides for 4 working combinations of stethoscope from ordinary to neonatal).

Complete set:

- Twin metal head (LD-S041), made by special design, with big diaphragm and extension piece with quartz clock – 1 pc. Working diameter of diaphragm – 38 mm. Diaphragm is provided with transparent plastic membrane. Design feature of head enables to change both diaphragm and clock to any extension piece from the set of additional accessories. Diaphragms, clock and bells are attached to the head by means of thread. Change-over between working sides is carried out by head turn round rod valve by 180° Flat side of rod valve points out to working side of head.
- Two thickened sound-conducting tubes of 56 or 72 cm long (indicated in box) – 1 set.
- Binaural metal tubes with spring and ear tips – 1 set.

Additional accessories entering the set:

- Plastic big bell (diameter 33 mm) – 1 pc.
- Plastic small bell (diameter 26 mm) – 1 pc.
- Plastic neonatal bell (diameter 19 mm) – 1 pc.
- Big membrane (diameter 40 mm) – 1 pc.
- Small diaphragm (working diameter 25 mm) – 1 pc.
- Small membrane (diameter 27 mm) – 1 pc.
- Plastic ear tips – 2 sets.
- Plastic bag for additional accessories – 1 pc.

Size with package, mm: 300 x 130 x 40.

Weight (with package), g: 398.

Change of power cell

Power cell, type LR 626/377, 1.5 V, is used in clock of stethoscope LD SteTime.

Change of power cell is necessary on the average once a year.

Sequence of power cell change:

1. Unscrew attachment piece with clock from stethoscope head.
 2. Unscrew with the help of coin rear part of attachment piece.
 3. Remove clock.
 4. Using thin screwdriver, pry off and remove clock lid.
 5. Change power cell (type LR 626/377, 1.5 V).
 6. Make assembly in reverse order.
- Guarantee for power cell is not applied.

MAINTENANCE, STORAGE, SERVICE AND DISPOSAL

1. THE CLOCK IS TURNED OFF DURING TRANSPORTATION AND STORAGE TO PREVENT BATTERY DISCHARGE. PERFORM THE SEQUENCE OF ACTIONS TO TURN THE CLOCK ON (for LD SteTime).

1.1. Hold the clock by the central part with two fingers and turn its case counterclockwise (Fig. 4). Release and remove the clock mechanism (Fig. 5).

- 1.2. Set the time by rotating the crown (Fig. 6).
- 1.3. Turn on the clock by pressing the crown until it clicks (Fig. 7).
- 1.4. Put the clock mechanism into the case and repeat the procedure the other way round.
2. Do not keep or use the stethoscope close to heaters and open fire.
3. Keep the stethoscope out of contamination.
4. Make sure that the stethoscope does not contact with any aggressive solutions.
5. When required, repair this device at authorized service centers only.
6. For disposal of the device follow the local regulation applicable at your region. No special disposal provisions are stipulated by the producer.
7. Operating conditions: +10°C to +40°C; relative humidity: 85% or below. Storage and transportation conditions: -34°C to +50°C; relative humidity: 85% or below.
8. Change clock power cells (for LD SteTime) in specialized clock workshops.

WARRANTY

The purchased device is warranted for 2 years from the date of sale. Faults detected during the warranty period will be repaired free of charge by a specialized service center within 21 business days from the date of device delivery to the service center. The costs for the defective device transportation are borne by Little Doctor Europe Sp. z o. o. (door-to-door service), except for faults that are not covered by the warranty. For the addresses of the warranty service providers, see MANUFACTURER AND AUTHORIZED REPRESENTATIVES or go to www.littledoctor.sg.

The party submitting the claim must properly prepare the device for shipment, protecting it from damage during transportation and enclose a completed warranty card (with the seller's stamp, date of sale and seller's signature) and proof of purchase (invoice or receipt). The warranty does not cover damages resulting from inadequate protection of the device during transportation.

The warranty period is extended for the time the device is under repair. Post-warranty repair of the device is performed at the user's expense. The warranty does not cover the following: consumables and accessories for individual use (tubes, nebulizers, masks, mouthpieces, etc.); products with mechanical, thermal or chemical damage, or damage caused by violations of the rules specified in the operating instructions; products with signs of modification, opening and/or repair by an unauthorized service center (by an individual); products that have been connected to unrecommended power supplies or used with unrecommended consumables; products that have been operated with non-remedied defects; products that are defective through no fault of the manufacturer, such as: actions of third parties, natural phenomena and natural disasters, animals, insects, foreign objects or liquids entering the device, etc.; internal and external contamination, scratches, cracks, scrapes and other mechanical damage to the equipment occurring during its operation; products with illegible or altered serial numbers or a damaged warranty seal; batteries, glass parts, light bulbs, as well as operations such as adjusting, cleaning and other maintenance of the product.

The buyer has the right to replace the device with a serviceable one if an authorized service center finds a manufacturing defect that cannot be eliminated or after 3 unsuccessful repairs of the same malfunction.

MANUFACTURER AND AUTHORIZED REPRESENTATIVES

SYMBOLS:

CE Compliance with Directive MDR (EU) 2017/745

Important: Read the manual

Authorized Representative in the EU

Manufacturer

UDI Unique Device Identification

Protect from moisture

Serial number

Medical Device

Importer

The revision date of this Instruction Manual is shown on the last page as XX-XXXX-YYMM-NN, where YY is the year, MM is the month, and NN is the revision number.

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

This unit manufacturing is certified according to international standard ISO 13485. Stethoscope comply with the requirements of European Directive MDR (EU) 2017/745.

Manufactured under control and for Little Doctor International (S) Pte. Ltd. (7500A BEACH ROAD #11-313 THE PLAZA SINGAPORE 199591. Postal Address: Sengkang Central PO Box 960 Singapore 915411).

Manufacturer: Little Doctor Electronic (Nantong) Co. Ltd. (No. 8, Tongxing Road Economic & Technical Development Area, 226010 Nantong, Jiangsu, People's Republic of China).

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