LUMINEO®

SONIC | EMS | ION+ | ION- | LED

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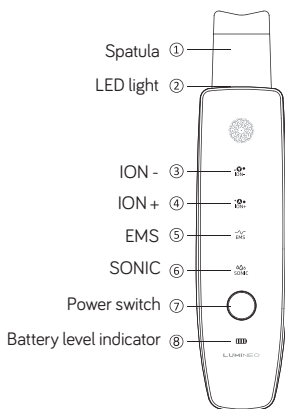
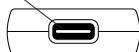
FUNCTIONS OF THE DEVICE

SONIC (Cavitation Peeling)

It is a modern, effective and gentle way of face deep-cleansing that you can use at home. It gently exfoliates the skin without disrupting it, which is why it can also be used on sensitive or couperose skin. The method uses ultrasound that breaks down and removes dead skin cells. The treatment can be conducted in any season. It perfectly prepares the skin for absorbing nutrients from cosmetics.

EMS (Electrical Muscle Stimulation)

The device is intended for conducting treatments using the EMS technology at home. EMS is an electrotherapy facial using low-frequency electric impulses. Electrotherapy stimulates muscle contractions, thus relaxing muscles, boosting blood circulation and removing redundant metabolic products. In the beauty industry, electrostimulation is most commonly used to generate the lifting

effect on the face skin. The effects of EMS are visible as soon as after several treatments. The skin is firmer, less rough and its structure is less granular.

ION + ION - (Iontophoresis)

This treatment uses direct (galvanic) current. This technology is considered to be one of the most effective anti-aging methods as it improves the condition of the skin, including smoothing out and reducing the appearance of wrinkles. The method boosts the deep layers of the skin by generating gentle microcurrents (both negative and positive), which initiates the cell regeneration process. Thanks to the microcurrents, the skin is visibly smoother, softer and firmer.

TURNING THE DEVICE ON AND OFF

Press the power switch ⑦ to turn the device on. Press and hold the power switch ⑦ for 2 seconds to turn the device off. The device will turn off automatically after 5 minutes.

SETTING OPERATING MODES

After switching on, the device will start operating in the SONIC mode (cavitation peeling) and the red LED light will turn on automatically.

In order to switch to the EMS mode, press the power switch ⑦ again. Blue LED light will turn on automatically.

In order to switch to the ION+ or ION- mode, press the power switch ⑦ again. Red/blue LED light will turn on automatically.

HOW TO USE THE DEVICE?

PREPARING THE SKIN BEFORE TREATMENT

1. Remove make-up and cleanse your face before using the device. For cleansing, we recommend using the LUMINEO facial cleansing brush and applying the LUMINEO SUPREME CLEANSER.
2. Before using the device, test it out on a small patch of the skin and check both the device and the spatula. Make sure that there are no rough or defective surfaces and no bubbles. If you notice any irregularities, do not use the device and contact the manufacturer.


SONIC (Cavitation Peeling + Red LED Light)

Apply the LUMINEO cavitation peeling toner provided with the device to your face.

Turn on the SONIC mode (cavitation peeling).

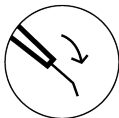Place the curved part of the spatula against your skin at a 45° angle and start sliding the spatula along your face without keeping it too long in one place. If your face starts to get dry, reapply the toner.

Once the treatment is over, remove the remaining toner gently or wait until it is absorbed.

Cavitation peeling can be conducted not only on the face but also on the neck, hands or back. Recommended time of the treatment is 5 minutes, 1-2 times per week.

EMS (+ Blue LED Light)

Apply a cream of your choice, e.g. LUMINEO BOOSTER*, that is recommended for EMS treatments on your face.

* available separately, not included in the set.

Turn on the EMS mode.

Place the flat part of the spatula against your skin at a 45° angle and start sliding the spatula along your face without keeping it too long in one place. If your face starts to get dry, reapply the cream during the treatment.

Once the treatment is over, gently rub the remaining cream into your face.

EMS treatment can be conducted not only on the face but also on the neckline. Recommended time of the treatment is 2 minutes. You can conduct treatments everyday during the morning or evening routine using LUMINEO products (recommended for EMS treatments).

ION+ ION- (Iontophoresis)

Apply a cosmetic of your choice, e.g. LUMINEO Me-so-Serum*, that is recommended for iontophoresis in a correct order.

* available separately, not included in the set.

• V •
ION-

Turn on the ION- mode and apply Meso-Serum no. 1.

Place the flat part of the spatula against your skin at a 45° angle and start sliding the spatula along your face without keeping it too long in one place. If your face starts to get dry, reapply the serum during the treatment. Continue the treatment for 2 minutes.


Once the treatment is over, rinse the remaining Meso-Serum no. 1 off your face.

• A •
ION+

Turn on the ION+ mode and apply Meso-Serum no. 2.

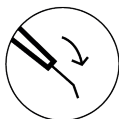Place the curved part of the spatula against your skin at a 45° angle and start sliding the spatula along your face without keeping it too long in one place. If your face starts to get dry, reapply the serum during the treatment. Continue the treatment for 3 minutes.

Once the treatment is over, rinse the remaining Meso-Serum no. 2 off your face.

Iontophoresis treatment can be conducted not only on the face but also on the neckline. Recommended time of the treatment is 5 minutes. The treatment can be conducted 2–3 times per week.

COUNTERINDICATIONS

1. Metal objects in the body (heart pacemaker, implants, braces, metal plates, splints).
2. Nervous system disorders (epilepsy).
3. Thyroid diseases.
4. Thrombophlebitis.
5. Multiple sclerosis.
6. Postoperative period.
7. Fever.
8. Inflamed rosacea, cystic acne.
9. Any skin conditions: caused by bacteria, fungi and viruses (including herpes).
10. Numerous telangiectasias.
11. Cancer.
12. Pregnancy, lactation, period.
13. Never use the device on damages skin, beauty marks, raised moles or if you have a skin condition.
14. Do not use the device if you are allergic to anything, have skin inflammation or if your skin is prone to irritation or allergies.
15. The device may only be used on healthy skin. It should not be used if you have open wounds, excessively sensitive skin, skin prone to rosacea, cystic acne or if you are allergic to metal.
16. If you suffer from another disease, you should consult your physician before using the device.

CHARGING

1. Before you start charging the device, make sure that the charging cable and the device are dry.
2. Do not use the charging cable in a humid environment, by the bathtub, shower, swimming pool or other places with water.
3. While charging, the battery level indicator flashes red. Once the battery is fully charged, the indicator is red.
4. Charge the device for 2 hours before the first use.
5. The device cannot be turned on while it is charging.
6. For charging, you should only use the cable provided with the device.
7. Use a charger with rated voltage of 5V and rated output power of 1A or a USB port. Chargers with different specifications may slow down the charging

process or result in overheating the device.

8. The device can be charged many times by connecting the USB cable to the wall mount charger or computer.
9. Disconnect the cable after the device has been charged.
- 10. If you are not planning to use the device for more than 3 months, charge it before putting it away. Otherwise, the battery may get damaged and it may be impossible to charge it again.**

CLEANING AND MAINTENANCE

1. Clean the spatula and the device on a regular basis.
2. Turn off the device while cleaning it.
3. Use a damp towel first and then a dry towel to clean the device carefully and remove all the remaining skincare product from the spatula.
4. Do not use chemicals to clean the device or the spatula.

NOTES ON THE STORAGE

1. Switch off the device when it is not used.
2. Store the device and the charging cable in a clean, dry place, ensuring they do not fall into water.
3. Do not disassemble the device on your own.
4. Keep the device out of the reach of children.

NOTES ON THE USE

1. Read the user's manual and the notes before your start using the device.
2. The device should be used in line with the user's manual.
3. Never use the device on dry skin face. Always remember to wet the skin using skincare products.
4. Do not use the device too intensively on the same skin area.
5. Remove any metal jewelry before using the device.
6. Slight redness may occur after the treatment, it is a normal symptom.
7. If you feel any pain or discomfort during the treatment, stop it immediately.
8. Never place the device into water, use it in the bathtub or shower and rinse it under running water. The device is not waterproof and may be damaged.
9. Do not share the device with other persons for hygienic reasons.
10. Do not use the device if you suspect it may be damaged – in this case, contact the distributor or manufacturer.

11. Do not expose the device to direct sunlight, high temperature or boiling water as it may lead to device failure.
12. Avoid using the device at high temperatures, in places with large amounts of dust and in places where flammable substances are kept, such as kitchen or a car.

WARRANTY

1. The device is covered by the limited warranty for one year from the date of purchase.
2. In the case of filing a complaint, it may be necessary to show an original receipt of purchase indicating a date of the device delivery.
3. The manufacturer reserves the right to review the condition of the device.
4. The warranty does not cover the damage caused by:
 - improper use or storage of the device by the customer
 - using unauthorized spare parts or disassembling the device by the customer
 - using the device in a way that is contrary to the user's manual
 - mechanical damage to the device
 - natural wear of the device
 - having the device fixed in an unauthorized service
 - warranty expiry

SERVICING

Do not try to fix the device on your own, otherwise the warranty will be void.

DISPOSAL OF THE DEVICE

This device is marked with a symbol of a crossed-out wheeled bin in accordance with the Polish Waste Electrical and Electronic Equipment Act and the WEEE Directive 2012/19/EU. This mark means that the device, after its life period, may not be disposed of together with other household waste. The user is obligated to give it to a facility that collects waste electrical and electronic equipment. The collecting entity, including local collecting points, shops and communal units, set up an appropriate system enabling the user to dispose of this equipment. Proper management of the waste electrical and electronic equipment helps to avoid consequences that are harmful for people and natural environment and that stem from the presence of hazardous

components and inappropriate storage and recycling of such equipment.

TECHNICAL SPECIFICATION

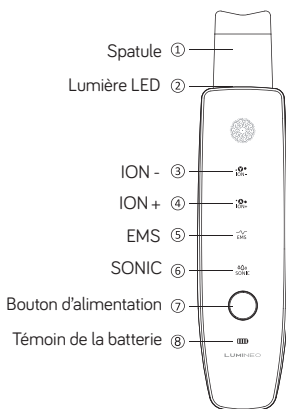
Multifunction device: SONIC | EMS | ION+ | ION- | LED

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