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
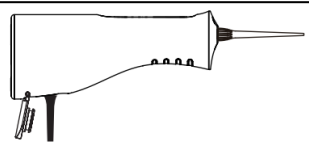
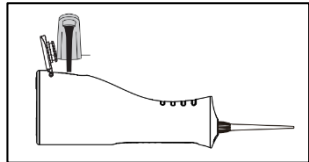
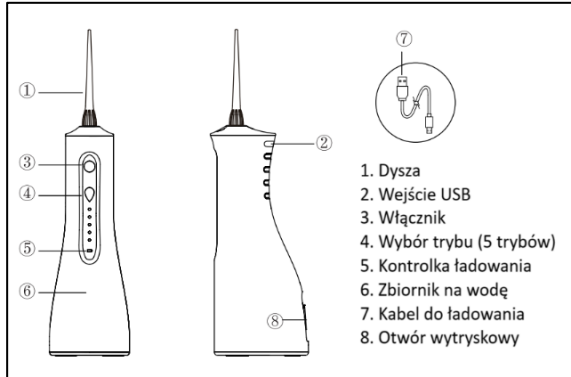
USER MANUAL U264

Importer: VITOLOG.PL Sp. z o.o, ul. Miedzyleska 6, 50-514 Wrocław, Poland
Product features: Cordless dental irrigator. Number of tips in the set: 4. Number of operating modes: 5. Charging station: none. Tank capacity: 300 ml. Power: 5W, 5V/1A. Waterproof: IPX7. Power supply: USB cable (1800 mAh/18650 battery). Operating time: 4 hours (on one charge). Material: ABS + PC. White color.

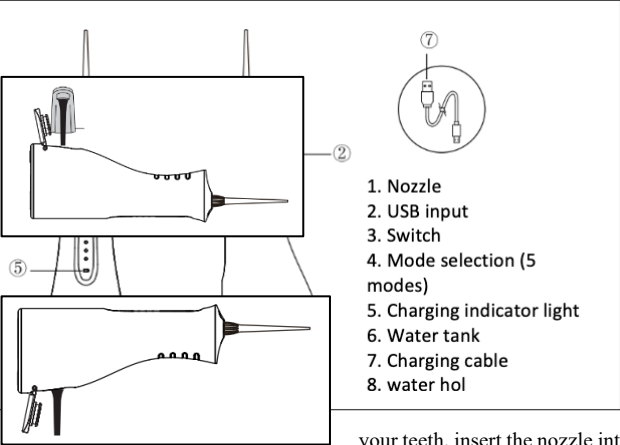How to use: A strong stream of water can thoroughly cleanse the oral cavity of bacteria, discoloration, tartar and food remains. A very gentle water pulse can deeply clean the gap between your teeth, improving periodontal health. The device deeply cleans the tooth surface and works better than traditional dental floss. Filling the tank with water: open the water tank cover, place the device horizontally and pour water. Method of correct tooth cleaning: select the appropriate operating mode of the device. It is recommended to set the power to the lowest setting at first and gradually increase it until the desired power is achieved. When brushing

your teeth, insert the nozzle into your mouth, aiming at the teeth or gums, opening your mouth slightly to ensure a smooth flow of water while rinsing your teeth. After selecting the appropriate operating mode of the device, start rinsing your teeth from the outside in one direction and when finished, rinse your teeth from the inside in the same direction (as shown in the image below). Note: Do not use strong mode to wash the damaged gum area or bleeding area continuously to prevent inflammation. The intensity of water pulsation should be adjusted to the condition of the oral cavity.

Removing water from the tank: if there is water in the tank after finishing work, pour it out. To do this, open the water tank cover and place the product horizontally with the tank facing down to pour out the water.

Cleaning the device: after use, rinse the device under the tap or wipe it with a wet cloth, and then wipe it with a dry cloth or dry paper towel.

Battery disassembly: disassembly of the device should only be done before disposing of the device. Before removing the battery, discharge it completely and turn off the device. Open the side cover compartment. Cut off the positive and negative electrodes, then pull out the battery.

Warnings: Please read the instruction manual before use and keep it for future reference. Use the product as intended. The product is intended for home use. Do not cover the head during use. Do not use the device for more

than 10 minutes at a time. Do not open or disassemble the product yourself, as this may damage it or result in electric shock. Do not fill the tank with water at a temperature higher than or equal to 40°C, as this may cause burns. Do not immerse the device in water while charging or while it is connected to the power supply. If the device falls into water while charging, disconnect the plug from the power source before removing the device from the water. Do not immerse the product in water as this may damage it. Regularly check the technical condition of the


!

ATTENTION

Read this manual before using the device for the first time or if it has not been used for a long time.

If your device is not responding, see below:
1. Lift them for a few seconds;
2. If this doesn't work, blow into the nozzle.

device and charger. Do not use a charger with a damaged cable or plug. A damaged charger must be replaced. Do not use a damaged product. Do not expose the product to extremely low or high temperatures. Protect the product from direct sunlight. Protect the product against strong vibrations, falls and mechanical damage. Regularly check the technical condition of the device. Any repairs may only be performed by an authorized person. Do not open the device. Do not modify the device or perform any repairs. If any problems occur with your device, stop using it. Do not use the product for any purpose other than cleaning your mouth. Do not press the nozzle against the gums or teeth to avoid damaging the gums. Do not use the product in damp rooms. Do not use the product at high temperature. Do not use the irrigator without the nozzle in place. Some mouthwashes may contain corrosive agents, so do not

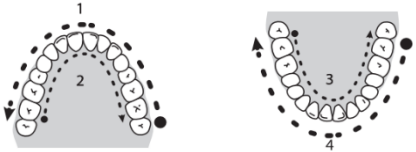pour them into the irrigator tank. Cleaning with running water or a wet cloth. Do not soak the products in water. Do not wash the product with hot water above 40 degrees. Do not bend or stretch the nozzle. It is recommended to replace the nozzle every 3 months. Empty the water tank and place it in a cool and dry place

Troubleshooting		
Problem	reason	Solution
The device does not work or will not turn on	The device is new or has not been used for more than 3 months.	The device should be charged for more than 4 hours.
	Low water pressure.	Replace nozzles.
After charging, the device only works for a moment.	The nozzle is blocked by foreign objects or deformed.	Clean the filtering net.
	The filter net is blocked.	Skontaktuj się z serwisem.
The charging light does not light or flashes randomly.	The battery life is ending.	Charge the device for at least 4 hours.
	Charging time was too short	Connect the charger correctly
The device stopped dispensing the water	The device is not properly connected to the charger.	Change the adapter or connect directly to the USB port.
	The adapter (charging head) is incorrect.	Contact service.
	The product circuit is malfunctioning.	
	The new device had not been used for some time.	Remove the nozzle and clean the body of the device, especially the nozzle opening.

Information for users on the disposal of electrical and electronic equipment (for households) The symbol shown on the products or the accompanying documentation indicates that faulty electrical or electronic equipment cannot be disposed of together with household waste. The correct procedure, if necessary, to dispose of electrical

or electronic equipment, to dispose of, reuse or recover components, is to transfer the equipment to a specialized collection point, where it will be accepted free of charge. In some countries, you can give the product back to your local distributor when you purchase another device. Proper disposal of the equipment preserves valuable resources and avoids negative health and environmental impacts that may be compromised by inadequate waste management. Please contact your local authority for details on the nearest collection point. Improper disposal of waste is punishable by penalties provided for in the relevant local legislation. **Users in European**

Union countries: If you need to dispose of electrical or electronic equipment, please contact your nearest point of sale or your supplier for further information. **Disposal of waste in countries outside the European Union:** This symbol is only valid in the European Union. If you need to dispose of this product, please contact your local authorities or dealer for information on how to proceed correctly.

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