
lovio

BEAUTY

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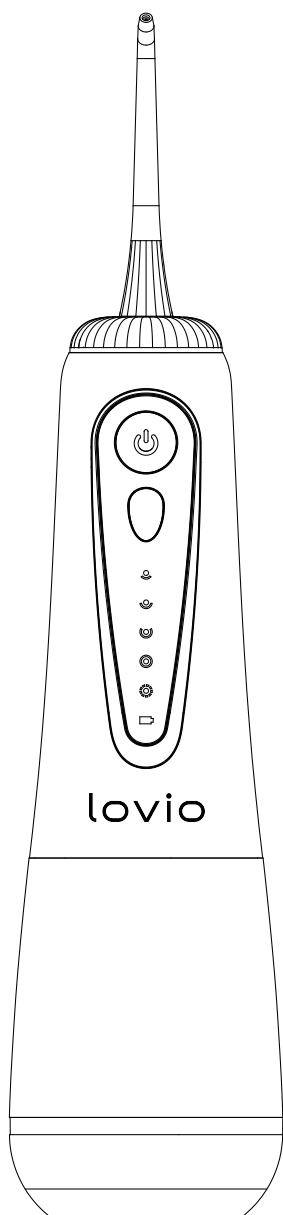
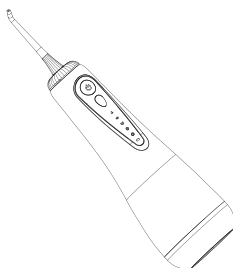
SPECIFICATION

Product Model:	Lovio Water Flosser LVWF124WH		
Material:	ABS	Battery Capacity:	2000mAh
Power:	2W	Water Pressure:	80--160psi
Charging Time:	5 hours	Waterproof:	IPX8
Water Tank Capacity: 340mL			

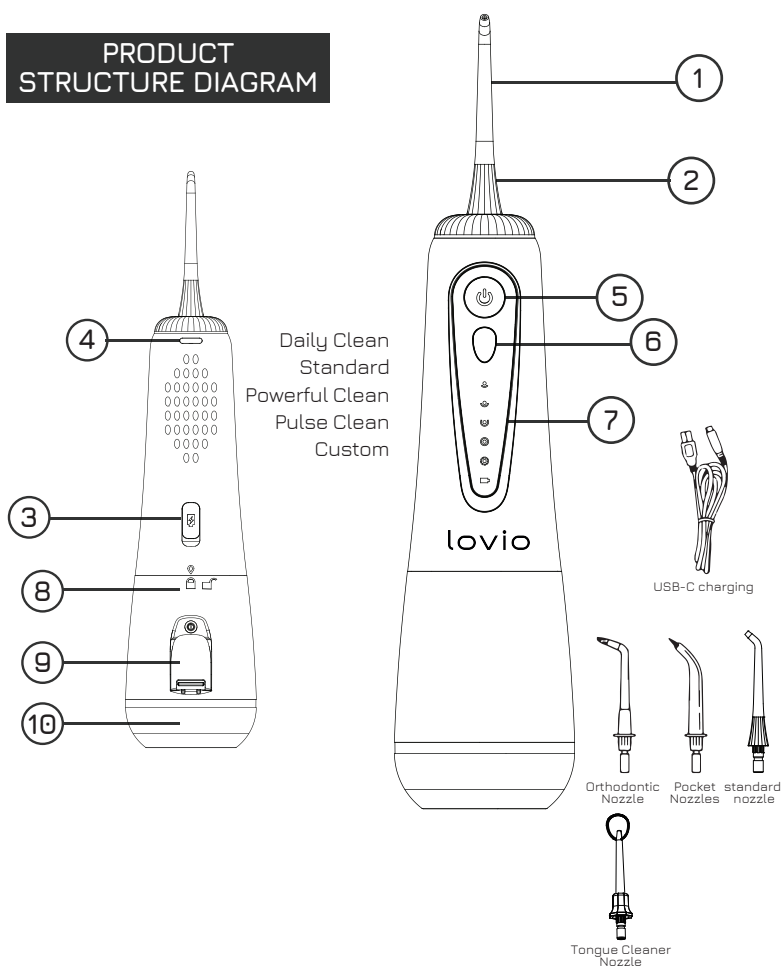PRODUCT
STRUCTURE DIAGRAM

- | | | |
|------------------------|-------------------|------------------------|
| 1. Nozzle | 2. Nozzle Rotator | 3. USB-C Charging Port |
| 4. Nozzle Eject Button | 5. Power Switch | 6. Mode Button |
| 7. LED Indicator Light | 8. Water Tank | 9. Water Inlet |
| 10. Water Tank Cap | | |

IMPORTANT SAFETY TIPS

Please read all instructions before use.

DANGER:

To reduce the risk of electric shock, please note the following:


- Do not throw or insert any foreign objects into any openings of the product.
- Keep power cord away from heated surfaces.
- Do not operate where aerosol (spray) products are being used or where oxygen is being used.
- Remove all oral jewelry before using this product.
- Do not use if there are open wounds in the mouth or tongue.
- Do not open your mouth or rinse your teeth from a distance, and avoid splashing water on your body or face.

WARNING:

- Keep away from water while charging.
- Always make sure the handle is completely dry before connecting the USB 3.5mm cable. Never charge this product outdoors or near hot surfaces.
- Please use the official charging cable and use 5V safe voltage to charge when charging.
- Keep away from fire sources when charging, using, and storing. Do not expose it to direct sunlight. The ambient temperature should be between 0°C and 40°C.
- Do not allow children to use this product without supervision.
- Do not let children play with this product.
- Please only use LOVIO original accessories or consumables.
- Never fill the sink with hot water. Hot water may cause burns.
- Never apply excessive pressure to the tip of the nozzle.
- To avoid injury, never fill the sink with hot water. Never put foreign objects into the sink.
- Please stop using damaged nozzles. Please replace the nozzle every 6 months or earlier if there are signs of wear.
- This product does not contain any user-serviceable parts. If this product is damaged, please stop using it and contact the service department.
- Do not put any part of this product in the dishwasher.

- This product is designed to clean teeth and gums only. Do not use it for any other purpose. (e.g. eyes, nose and ears).
- If your gums bleed profusely after using this product and continue to bleed, feel unwell or painful after 1 week of use, stop using this product and consult your dentist/doctor.
- If you have had oral or gum surgery in the previous 2 months, please consult your dentist before using this product.
- If you have a pacemaker or other implanted device, please consult your physician or the manufacturer of the relevant implanted device before use.
- If you have a medical problem, please consult a physician before using this product.
- This product is for personal use and should not be used by dental clinics or institutions for multiple patient use.
- Use this product only for its intended use as described in this user manual. Misuse may result in hazard or serious injury.
- Before charging, be sure to empty the remaining water in the water tank, slide the water tank to the bottom of the main body to put it away, and then place the machine stably.
- When not in use for a long time, please be sure to empty the residual water in the water tank.
- If the product is abnormally hot or smelly, discolored, or takes longer than usual to charge, please stop using it and stop charging it, and contact after-sales customer service.
- Never place the product in a microwave oven or on an induction hob.
- To prevent the battery from heating up or releasing toxic and harmful substances, never open, modify, puncture, damage or disassemble the product or battery. Do not short circuit, overcharge or reverse charge the battery.

Charge the unit before first use

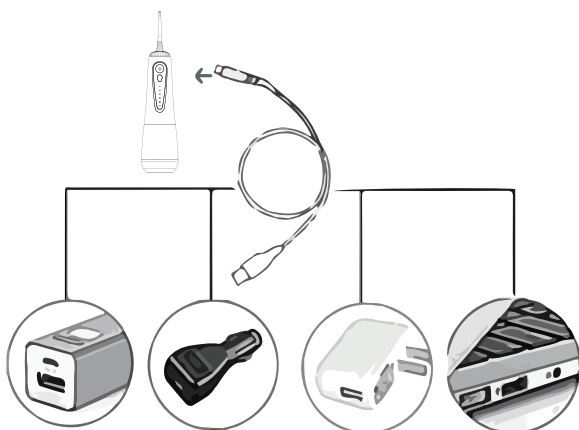
Tips: When charging this product, the user can use the USB 3.5mm charging cable provided with the product to connect any 5V 1A mobile phone charging adapter, computer and charging treasure, etc.

- **Main body charging:** When the product charging, it automatically stops working; The charging light flashes when charging, and stays on when fully charged.
- **Low battery reminder:** When the charging light at the bottom flashes during working, it indicates low-voltage and reminds the user that charging is required.
- If you have not used the oral irrigator for 6 months, please charge it according to the charging requirements of the purchased oral irrigator model until it is fully charged.

Fill the tank with water

Open the lid of the water tank, place the oral irrigator horizontally and fill the water tank with water

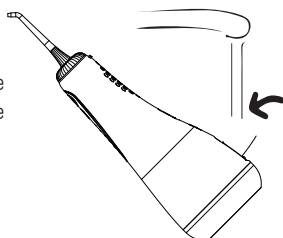
Inserting and removing nozzles

Insert the nozzle into the top of the oral irrigator. Press down firmly until the nozzle "clicks" into place. To remove the nozzle from the body, press the nozzle release button and pull the nozzle out of the body.

Note: Do not press the nozzle release button when the machine is running

When it's over

Turn off the appliance and empty any remaining liquid from the water tank.

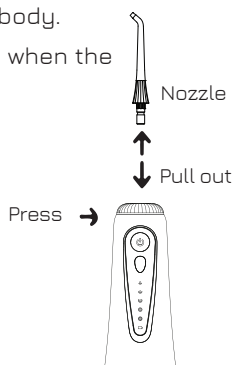Model Introduction


Daily Mode: Gentle water pressure, suitable for first-time users.

Standard mode: efficient cleaning, satisfying most people's daily life clean.

Powerful Mode: Strong water pressure to remove stubborn stains.

Pulse Mode: Alternating strong and weak, deep cleans and supports gums.


Custom mode: Long press to adjust the flow rate of water without any adjustment.

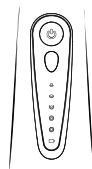ACCESSORY USAGE


Orthodontic Nozzle

Rinse hard-to-reach areas around braces, braces, etc. use When cleaning the braces, the braces and the gap between the teeth with the orthodontic nozzle, supplemented by jet of water.

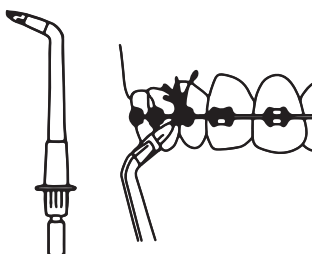
Standard nozzle

When using, lower the head slightly, put the nozzle in the mouth, open or close the mouth slightly, align the nozzle with the gap between the teeth, and then press the switch button to rinse, rinse each tooth along the gum line, pits and grooves, wisdom teeth and other hidden parts You can focus on rinsing, let the water flow out or spit out from the corner of the mouth naturally, press the switch button after use then remove the nozzle from your mouth.

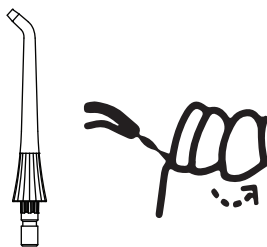Deep/Pocket Nozzles

Specifically designed to deliver water or antiseptic solutions into deep periodontal pockets. Always set the unit to minimum pressure when using deep/periodontal pocket nozzles. Aim the flexible nozzle at a tooth at a 45-degree angle, and move the nozzle gently until it is placed under the gum line and into the periodontal pocket.

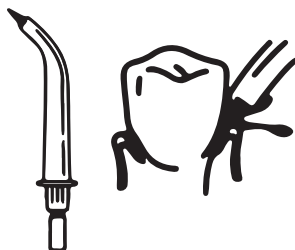Tongue cleaner nozzle

This tip is specifically designed to clean the tongue by combining a gentle scraping and rinsing action. While rinsing, it removes bacteria and sulphur compounds that are found on the tongue and cause bad breath.

Warm tips for product use

Dental experts suggest: the first time users use the soft mode for a week, and then use the normal mode after getting used to it. Press the switch button to use it. When using it, you need to pay attention to the following points:

- (a) The direction of water flow is perpendicular to the teeth or gums;
- (b) Move slowly along the teeth and clean the teeth sequentially;
- (c) The nozzle is aligned with the gum line and perpendicular to the gum;
- (d) Thoroughly clean teeth, trays, crowns and dentures.

For the first time use, do not spray water directly along the periodontal pocket for a long time, so as not to cause tooth discomfort.

Please use an adapter conforming to national standards with an output of 5V 1A and an output power not exceeding 15W to charge the product.

Removal of hard water deposits/ general cleaning

Hard water deposits may build up in the unit depending on the mineral content of the water. If you ignore this, it will affect the performance. Clean the internal parts: fill the water tank with warm water, add 1 tablespoon of white vinegar anchor point. Point the body spray head into the washbasin. Turn on the appliance and let it run until the water tank is empty. Then fill the water tank with clean warm water and repeat the above operation.

Battery Care and Disposal

This device contains a lithium battery. When the battery life has expired, please send it to your local recycling center. The battery is not replaceable.

Note: If the machine stops during use, it is because of excessive heat. Please wait two hours before using it!


Common faults and solutions

Phenomenon	Possible reason	Measure
Water tank leaking	Abnormal sealing ring of water tank cover	Disassemble the sealing ring and reinstall
The water pressure is too low	Not enough power	Fully charged
	Nozzle clogged	Replace nozzle
It can only work for a few minutes after being fully charged	The battery life is exhausted (approximately 500 charge and discharge times)	The battery life is exhausted (approximately 500 charge and discharge times)
	Charging cable failure	Replace the charging cable
	Charging head failure	Replace the charging head
There is sound but no water	Nozzle clogged	Replace nozzle
	Coming out Internal disc tilted	When the product is turned on, pour water from the nozzle
No Sound No Water	Low battery (not used for 3 months)	Continuous charging for more than 16 hours to activate the battery
Suddenly Stopped Working	Product 2-minute antimistouch design	Reboot to use
	Motor temperature control protection at 75 degrees (maximum gear for continuous use of 5 boxes of water)	Stop for 2 hours and turn on again

If you still can't solve the problem according to the above methods, please consult the dealer you purchased and go to the designated repair center for repair.

When oral irrigator is used for the first time or has not been used for a long time, this phenomenon may occur; after pressing the power switch, it needs to wait for a while before the nozzle can spray water. In this case, please proceed as follows.

After adding an appropriate amount of water to the installation port of the nozzle, install the nozzle again, and then turn on the power switch.

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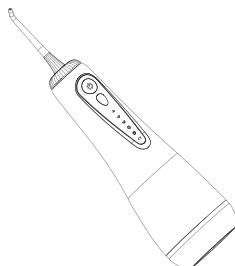tel: +48 22 123 96 60 lub e-mail: serwis@manta.com.pl

Monday to Friday, 9 a.m. to 5 p.m.

The manufacturer reserves the right to make changes to the product specifications without prior notice.

Manta S.A. hereby declares that this product complies with the essential requirements and other provisions of the so-called “New Approach” Directives of the European Union.

The device is marked with the symbol of a crossed-out waste container, in accordance with European Directive 2012/19/EU and the Polish Act of 11 September 2015. (Journal of Laws 2015 item 1688) on Waste Electrical and Electronic Equipment (WEEE). Products marked with this symbol should not be disposed of or thrown away with other household waste at the end of their useful life. The user is obliged to dispose of waste electrical and electronic equipment by delivering it to a designated point where such hazardous waste is recycled. The collection of this type of waste in separate locations and the proper recovery process contribute to the protection of natural resources. Proper recycling of waste electrical and electronic equipment has a positive impact on human health and the environment. For information on where and how to dispose of waste electrical and electronic equipment in an environmentally safe manner, the user should contact the relevant local authority, the waste collection point or the point of sale where they purchased the equipment.

Attention!

Product design and technical specifications are subject to change without prior notice. This applies in particular to technical specifications, software, pictures and the user manual. This user manual serves as a general orientation for the operation of the product. The manufacturer and distributor assume no liability for any inaccuracies resulting from errors in the descriptions contained in this user manual.

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