
Purea S6 Pro Hand Massager

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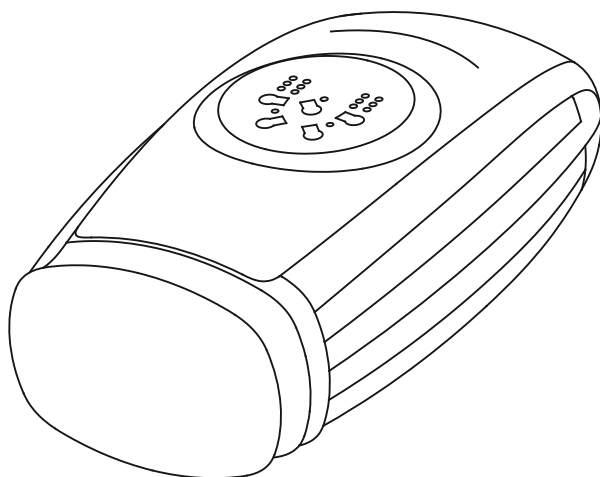Thank you for choosing Boste.

User Manual

Read carefully and make sure you understand the instructions before use.
Pay particular attention to the safety precautions.

After purchasing the device, charge it first — this will extend the battery life.
Use a 5V, 1A (or 2A) charger, like a phone charger.

While charging, the power indicator will flash — this means the device is charging.

When charging is complete, the indicator will turn off.


Target Group

- People who work on computers for long periods.
- People with poor blood circulation in their hands.
- Women who want to maintain the appearance of their hands.

Product Features

- Ergonomic shape that matches hand acupressure points and natural curves.
- Pneumatic, segmented compression massage.
- Full-palm acupuncture-style massage.
- Heating function for cold environments (optional feature).
- Micro-vibration increasing massage intensity (optional feature).


Control Panel

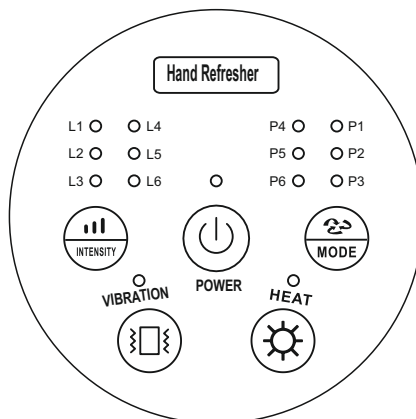
| Functions |


Air Pressure

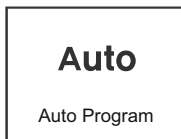Auto Program

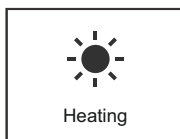
Heating

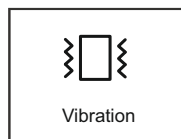Vibration

| Device Components |
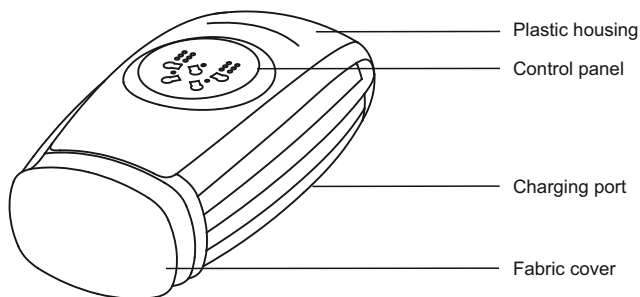
| How to Use the Control Panel |

Power Button

Turns the device on and off.
Hold for 2 seconds to activate the power supply, then press the mode button — the device will start operating.
The default massage time is 15 minutes, after which the device switches off automatically.

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The default intensity is L3. There are 6 intensity levels available.
You can switch through them in the sequence: L3 - L4 - L5 - L6 - L1 - L2.
L1 is the weakest level, L6 is the strongest.

| How to Use the Control Panel |

Mode Button

The default mode is P1. There are 6 massage modes in total.
You can choose your preferred mode by switching in sequence:
P1 - P2 - P3 - P4 - P5 - P6 - P1.

Vibration Button

Press to turn the vibration on or off.
There are 2 vibration modes — choose according to your needs.

Heating Button

Press to turn the heating function on or off. After about 5 minutes, the internal temperature will noticeably increase.

| Troubleshooting |

Behaviour	Possible Causes	Solution
The product will not start	Battery is discharged	Charge the device
	The device was not started correctly	Hold the power button for 2 seconds
Faint sounds during operation	Sound of the air pump (buzzing)	Normal operation, not a fault
	Sound of the solenoid valve (clicking)	

| Specifications |

Model	Purea S6 Pro Hand Massager	Power	5W
Model number	SM-001	Working time	15 min
Voltage	DC 3.7V	Weight	1.1 kg

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