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Product features: Buttock muscle electrostimulator. Dimensions: 20 x 26 cm. Packed in a box. The set includes a control panel and a band. Panel for 2x AAA batteries (not included). High stimulation power. Fast and visible results. High strength.

Before use: Insert the battery before use: open the battery compartment and insert 2 x AAA battery (not included) with correct polarity, then close the battery compartment. Put the device on your buttocks and turn it on using the control panel. The device can be used during exercise or walking, but also while cleaning, watching movies, reading books, etc.

How to use:

1. Disassemble the packing bag and install the main motor on the massage patch.
2. Remove the transparent protective film from the massage patch.
3. Attach the massage patch to the part of the body that needs to be massaged.
4. Startup: Press the "ON/INC" button to start and keep pressing

Press the "ON/INC" button to increase the strength.

5. Press the "OFF/DEC" button to decrease the strength. The default strength is level

11 after changing the mode

6. Mode: Press the "PROGRAM" button to select the mode and then press "ON/INC" to start working, the mode can also be changed during operation.

7. Power off: Keep pressing the "OFF/DEC" button to weaken the force until the red light is off.

8. After turning off, take off the massage patch and cover it with a transparent film cover and place it in a cool place.

Good practices:

1. To protect the life of the hydrogel, please clean your skin before using this product.
2. It is recommended to use 12-15 minutes, 2-3 times a day. If you continue

Using it for a long time may cause fatigue or sickness.

3. Do not exercise within half an hour after eating, which may cause discomfort.


Warnings: Please read this manual in its entirety before using the product and keep it for the lifetime of the product. The owner is responsible for ensuring that all users of this product are properly informed on how to use it safely. Use the product as intended. This product is for home use only. Protect against moisture and dirt. Do not use or store the product in wet conditions. Protect the product from direct sunlight. Avoid contact with water and other liquids. Do not allow the product to heat up or stay near sources of fire. Do not disassemble or modify the product. Any repairs may only be carried out by an authorized person. The product may only be maintained by an authorized person. Never use a damaged product. If the product does not operate as intended, turn it off immediately. If you are not sure whether the stimulation will not harm your health, consult your doctor before using it. Places with neoplastic lesions or their vicinity must not be stimulated. Do not use the device in the bath or shower. Do not use the device while sleeping. Do not use stimulation while driving, operating machinery, or engaging in other activities that present a risk of injury. Only attach the device to normal, undamaged, clean and healthy skin. Stimulation should be stopped or the intensity level should be lowered when we begin to feel discomfort during it. If this feeling doesn't go away, see your doctor. Do not use the stimulator in the vicinity of shortwave or microwave therapy devices as this may reduce its output power. The stimulator must not be used together with electronic health monitoring devices (e.g. pulse oximeter, ECG alarm unit), as in this case their indications may be incorrect (noisy). Skin irritation or burns may occur where the electrodes come into contact. Do not use this device if you have a pacemaker, defibrillator or other metal or electronic equipment. If this happens, there is a risk of electric shock, burns, electrical interference, or death. Do not mix old and new batteries. Batteries that have not been used for a long time should be removed. Batteries must be inserted with the correct polarity. The item is not a toy. Keep away from children.

Information for users on the disposal of electrical and electronic equipment (for households) The symbol shown on the products or the accompanying documentation indicates that faulty electrical or electronic equipment cannot be disposed of together with household waste. The correct procedure, if necessary, to dispose of electrical or electronic equipment, to dispose of, reuse or recover components, is

to transfer the equipment to a specialized collection point, where it will be accepted free of charge. In some countries, you can give the product back to your local distributor when you purchase another device. Proper disposal of the equipment preserves valuable resources and avoids negative health and environmental impacts that may be compromised by inadequate waste management. Please contact your local authority for details on the nearest collection point. Improper disposal of waste is punishable by penalties provided for in the relevant local legislation. **Users in European Union countries: If you need to dispose of electrical or electronic equipment, please contact your nearest point of sale or your supplier for further information. Disposal of waste in countries outside the European Union: This symbol is only**

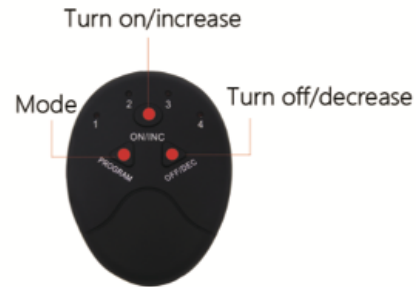valid in the European Union. If you need to dispose of this product, please contact your local authorities or dealer for information on how to proceed correctly.

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