



SNAILAX SL-52A2 Deep Kneading Foot Massager with Heat User Manual

[Home](#) » [SNAILAX](#) » SNAILAX SL-52A2 Deep Kneading Foot Massager with Heat User Manual 

Contents

- [1 SNAILAX SL-52A2 Deep Kneading Foot Massager with Heat](#)
- [2 Product Information](#)
- [3 Product Usage Instructions](#)
- [4 Warranty Information](#)
- [5 IMPORTANT INFORMATION](#)
- [6 Safety Instructions](#)
- [7 Before using the unit](#)
- [8 Maintenance and cleaning](#)
- [9 Useful Information And Contents](#)
- [10 Technical Data](#)
- [11 Setup and Operation](#)
- [12 Cleaning and maintenance](#)
- [13 Warranty](#)
- [14 Extend Warranty for Free](#)
- [15 Documents / Resources](#)
 - [15.1 References](#)
- [16 Related Posts](#)



SNAILAX SL-52A2 Deep Kneading Foot Massager with Heat



Product Information

- **Product Model:** SL-52A2
- Detachable foot sleeve
- Machine washable

Product Usage Instructions

1. Detach the foot sleeves from the product.
2. Wash the foot sleeves in a washing machine or by hand.
3. Allow the foot sleeves to air dry.

Warranty Information

This warranty does not include the following:

- Refund procedures (refer to <https://www.snailax.com/pages/refund-policy> for detailed instructions)
- Any warranty or affirmation of fact, express or implied, other than as outlined in this warranty statement
- Implied warranties of merchantability and fitness for a particular purpose
- Any warranty or affirmation of fact related to misuse, improper selection, recommendation, or misapplication of any product
- Any warranty or affirmation of the fact that the catalogs, literature, and websites accurately illustrate and describe products

Extend Warranty for Free

1. Visit our Snailax online service center at <http://support.snailax.com> or scan the QR code below to extend your warranty to 3 years.
2. Alternatively, like our Snailax Facebook page at <https://facebook.com/snailax> or scan the QR code below, and message us with your purchase proof.

3. Complete your application for a free warranty extension within 60 days of purchase.

IMPORTANT INFORMATION

RETAIN FOR FUTURE USE

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.



- **Explanation of symbols** This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.



- **WARNING** These warning notes must be observed to prevent any injury to the user.



- **CAUTION** These notes must be observed to prevent any damage to the device.



- **NOTE** These notes give you useful additional information on the installation or operation.



- **Protection category II**



- **Lot number**



- **Manufacturer**

Safety Instructions

about the power supply

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your main supply.
- Only connect the mains plug to the wall socket when the device is switched off.
- Only use the mains adapter that is provided with the foot massager.
- Keep the power supply, the cable, the control unit, and the foot massager away from heat, hot surfaces, moisture, and liquids. Never get hold of the power supply with wet hands.
- Do not take hold of a device if it falls into the water. Pull out the power plug and the power supply unit from the socket immediately.
- Always unplug the equipment immediately after use.
- To disconnect the unit from the power supply, always pull the plug out of the socket. Do not pull on the cable!
- Do not carry, pull, or twist the foot massager on the cable, or the power supply unit.
- You must not continue to use the unit if the lead or the mains adaptor is damaged. For safety reasons, these parts may only be replaced by authorized service personnel. To prevent hazards, always send the unit to the service center.

- Ensure that the cable will not cause tripping. The cable must not be kinked, trapped, or twisted.
- Do not use this product for medical treatment.

Specific groups of people

- The surface of the appliance is hot. People who are insensitive to heat must be careful when using the appliance.
- Children must be supervised to ensure that they do not play with the device.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse
- You should refrain from using the FOOT MASSAGE or should consult your doctor before using it, if
- you are pregnant,
- you have a cardiac pacemaker, artificial joints, or electronic implants,
- you suffer from one or more of the following illnesses and complaints:

Circulation disorders, varicose veins, open wounds, bruises, broken skin, and inflammation of the veins.

- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Never leave the unit unattended if it is connected to the power supply.
- Take care when using the heater. Do not use it if you have circulatory problems or on areas not sensitive to heat. Risk of burning!
- Consult your doctor if you have any questions concerning the therapeutic use of a foot massager.
- If you have unexplained pain, are receiving treatment from a doctor, and/or use medical equipment, consult your doctor before using the foot massager.

Before using the unit

- Always carefully inspect the power supply, the cable, the control unit, and the foot massager for damage before use. Do not attempt to use a defective unit.
- Do not use the unit if there are visible signs of damage on the unit or the cable components, the unit does not work properly or the massager or the control unit has been dropped or become damp. To prevent hazards, always send the unit to the service center.

operating the unit

- The foot massager must only be used for its intended purpose as described in the instruction manual.
- Using the device for any other purpose invalidates the warranty.
- Never leave the unit unattended if it is connected to the power supply.
- The maximum operating time for one session is 15 minutes.
- Do not allow the foot massager to come into contact with pointed or sharp objects.
- Never put the unit down next to an electric stove or other sources of heat.
- Do not stand on the unit.
- Do not touch the power supply unit when you are standing in water and always take hold of the plug with dry hands.

WARNING: The device is not to be used for commercial use or medical purposes.

- If you have health concerns, consult your doctor before using the foot massager.
- The driver may not use the foot massager when driving.
- Only use the foot massager in a parked vehicle.
- Only use the foot massager indoors.
- Do not use the foot massager in damp rooms (e.g. when bathing or showering).

Maintenance and cleaning

- The unit is maintenance-free. If a fault does occur, just check whether the mains power adapter is properly connected.
- You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all guarantee rights. Ask your specialist dealer and have repairs carried out by the service.
- Do not use the unit if there are visible signs of damage on the unit or the cable components, the unit does not work properly or foot massager or the control unit has been dropped or become damp. To prevent hazards, always send the unit to the service center.
- Do not immerse the device in water or other fluids.
- If, nonetheless, the liquid should penetrate the device, unplug the main cable immediately.
- Cleaning and user maintenance shall not be done by children without supervision.



Useful Information And Contents

Thank you very much

- Thank you for purchasing the Snailax foot massager.
- With normal care and proper treatment, it will provide years of reliable service.

Please read all instructions carefully before using this product.

Contents



Technical Data

- **Name and model:** Foot massager SL-52A2
- **Voltage supply:** 100-240 V~50/60 Hz
- **Power supply unit:** 12V = 2000 mA
- **Power consumption:** 24W
- **Timer:** 15 minutes
- **Operating conditions:** only use in dry rooms
- **Storage conditions:** in a cool, dry place
- **Dimensions:** approx. 34.8×32.5×17.3 cm
- **Weight:** approx. 2.4 kg
- **Adapter cable length:** approx. 1.8 m

WARNING: Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

Setup and Operation

Operating

1. Place the massage machine on the flat floor.
2. Connect the adapter cable to the corresponding port on the massager
- 3.



Plug the home adapter into an electrical outlet.

4. Turn on the device through the operating control panel.


WARNING

- Before plug please check and make sure the power cord is undamaged.
- Do not use the wet hand to plug in the adapter.
- Please don't use the product by standing posture or put any heavy object on the product to avoid any broken.

- Please don't stand or sit on the product.
- When moving the product, please don't pull the power cord or drag the product directly to avoid any damage caused.
- Please don't let children and pets play with the products to prevent any accidents from happening.
- Do not use the wet hand to pull out the power cord to avoid any electric shock.
- Do not use this product in a damp environment to avoid electric shock.

Setup and Operation

- Please choose a chair with the desired height and be seated, then put your feet in the massager.
- Operate below the control panel to get foot massages.

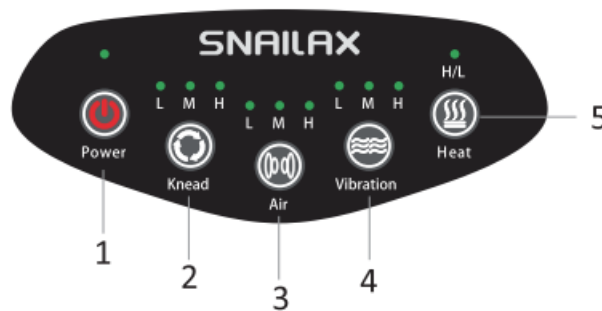

1. Power button:

Press this button to turn on or off the massager;

Turn on the massager and it automatically will work in the default mode: i.e. The rolling +compression +vibration+heating functions will work at Low speed

2. Rolling/Kneading button:

After powering on, press this button to turn on the kneading massage function, continuously press to cycle through Low/Medium/High Speed or Off setting, and the corresponding indicator L/M/H will light up.

3. Compression button:

After powering on, press this button to turn on the compression function, continuously press to cycle through the Low/Medium/High level of intensity of Off setting, and the corresponding indicator L/M/H will light up.

4. Vibration Button:

After powering on, press this button to turn on the Vibration function, continuously press to cycle through the Low/Medium/High level of intensity of Off setting, and the corresponding indicator L/M/H will light up.

5. Heat button:

After powering on, press this button to turn on the Heat function, and continuously press to choose Low/High heat level or Off setting.

- **Low level:** the temperature is 45-50 C (113-122 F), the indicator light is blue;
- **High Level:** the temperature is 50-55 C (122-131°F), the indicator light is red.

WARNING:

- Please don't put water or other liquid on the control panel.
- Please don't press the button too long time at a time to avoid the program being confused.
- When power failure happens during a massage, please turn off the power and pull out the plug from the socket.

- Please do not let your hand or fingers close to and insert into the gap of products to prevent any injuries caused.
- Please do not let children and pets play with the products to prevent any accidents from happening.

Cleaning and maintenance

- Before cleaning, make sure the power supply is turned off and the plug has been pulled out from the socket.
- It is forbidden to dip the whole product in water or other liquid when cleaning.
- Use a soft cloth to wipe the surface of the product and control panel.
- Use a piece of cloth that has been dipped with cleaning agents to wipe the main unit of the foot Massager, and then use a piece of dry and soft cloth to wipe it.
- Do not pull the power plug with damp hands.
- Do not use water washing products, and do not spray the water or other liquids to the interior of the product.
- Do not use corrosive cleaning agents like benzene or thinner to clean the Foot Massager, otherwise, the color of the upholstery will fade or the Foot Massager will crack or be eroded.
- Do not iron the upholstery.
- When a product is wet, let the product dry naturally and thoroughly.

Removable & Washable Cover

- Detach the foot sleeves and wash them in the machine or by hand. Allow to air dry.

Maintenance

- Make sure the power supply is turned off and the plug has been pulled out from the socket.
- Sweep off the dust. Keep the foot massager clean and away from damp places.
- Store the massager in a cool, dry place out of reach of children.
- Avoid contact with sharp edges or pointed objects that might cut or puncture the cover of the massager.

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

By our policy of continual product improvement, we reserve the right to make technical and optical changes without notice.

Warranty

Warranty Limitations(this warranty does not include):

Normal wear and tear; Any condition resulting from other than ordinary residential wear or any use for which the product was not intended, such as use in rental or contract trade or commercial use; Any condition resulting from incorrect or inadequate maintenance or care; Damage resulting from misuse, abuse, negligence, accidents or shipping damage, exposed to environmental conditions temperatures, water, improper operation, improper installation, improper use of electrical supply; transportation damage; theft, vandalism or; Damage resulting from the attachment or modification of any unauthorized accessories; Fail to show the valid proof of purchase: Lost or stolen products: Items that have expired their warranty period; No quality-related issues (after 30 days of purchase); Free product/ Gift offer by Snailax; unauthorized repairs or modification; Replaced or refurbished

product. Dissatisfaction due to buyer's remorse (after 30 days of purchase).

Refund Procedures:

(Refer to <https://www.snailax.com/pages/refund-policy> for detailed instructions)

Warranty Disclaimer

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. SNAILAX SHALL HAVE NO LIABILITY FOR, AND EXPRESSLY DISCLAIMS ANY WARRANTY OR AFFIRMATION OF FACT, EXPRESS OR IMPLIED, OTHER THAN AS SET FORTH IN THIS WARRANTY STATEMENT, INCLUDING, WITHOUT LIMITATION

1. THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE;
2. ANY WARRANTY OR AFFIRMATION OF FACT RELATED TO MISUSE, IMPROPER SELECTION, RECOMMENDATION, OR MISAPPLICATION OF ANY PRODUCT; AND
3. ANY WARRANTY OR AFFIRMATION OF FACT THAT THE CATALOGS, LITERATURE AND WEBSITES IT PROVIDES ACCURATELY ILLUSTRATE AND DESCRIBE PRODUCTS.

Snailax reserves all the right to final interpretation of these Terms of Warranty and Return

Policy above and the right to change, modify, add, or remove portions of these terms at any time without prior notification.

Extend Warranty for Free

1. Extend your warranty to 3 years by visiting our Snailax online service center at <http://support.snailax.com> or by scanning the QR code below:



2. Alternatively, like our Snailax Facebook page at <https://facebook.com/snailax> or scan the QR code below, to message us with your purchase proof



Your application for a free warranty extension must be completed within 60 days of purchase.

SNAILAX CORPORATION

- Email: supporteu@snailax.com
- Internet: www.snailax.com

TANMET INT'L BUSINESS LTD 9 Pantygraigwen Road, Pontypridd, Mid Glamorgan, CF37 2RR, UK
tanmetbiz@outlook.com

Documents / Resources

 SNAILAX Foot Massager SL-52A2 CE	SNAILAX SL-52A2 Deep Kneading Foot Massager with Heat [pdf] User Manual SL-52A2 Deep Kneading Foot Massager with Heat, SL-52A2, Deep Kneading Foot Massager with Heat, Kneading Foot Massager with Heat, Foot Massager with Heat, Massager with Heat, Foot Massager, Massager, Massager SL-52A2
--	---

References

-  [OnlineServiceCenter](#)
-  [Snailax Massager | Buy Snailax Shiatsu Neck & Back Massager | Shop Shiatsu Massage Mat & Cushion at Snailax.com](#)
-  [Refund Policy - Snailax](#)
- [User Manual](#)