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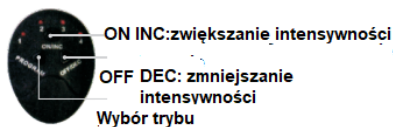
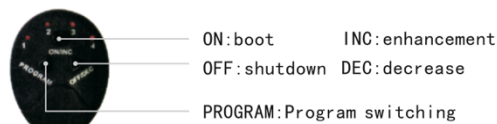
Mode indicator light

Installation guide

Product features: muscle electrostimulator; material: ABS, PU; set includes: 1 x belly pad, 2 x shoulder pad, 3 x remote control; 6 operating modes (A-F); working temperature: 0-45°C; power supply: 2 x 1.5V AAA batteries; frequency: 1-100Hz; output: max 9.8mA; dimensions: 21 x 10 x 2.3 cm; weight: 57.4g; automatic shutdown; single workout time: 12 minutes

How to use: Remote Control Buttons: ON/INC: Power On/Intensity Increase Button. OFF/DEC: button to turn off/decrease the intensity.

PROGRAM: switching between modes (choose between 6 modes (A-F)). Mode Indicators: When you press the ON/Intensity Up button, the light is intense. When the OFF/intensity down button is pressed, the light is a shade darker. There are 10 intensity levels

Before use, install the battery according to the polarity in the remote control. Select the overlay and place the remote control in the control unit. Clean and dry the area of your body where you will use the stimulator. Remove the vinyl cover and place the pad on the body where the muscles are stimulated (abdomen, arms). Turn on the device by pressing the ON button. Select the operating mode by pressing the letter button. Press the +/- button to increase/decrease the frequency. To adjust the strength of the pulses, press the +/- button.

To turn off the device, press the OFF button or the device will automatically turn off after 12 minutes.

Warnings: It is recommended to use the lowest output power at the beginning. Do not use the device immediately after a meal, shower or bath. The product is intended for home use. The device is not a medical device. Before using the product, consult your doctor and make sure that there are no health contraindications. Do not use for any injuries or injuries to the limbs. Do not use in case of problems with the cardiovascular, respiratory and circulatory systems, diseases of the spine, problems with pressure and other health problems. Do not use by people with a pacemaker, pregnant women. Do not use the product while driving a car or operating machinery. Do not stick electrodes on diseased skin. Do not use the product in a place with high humidity. Do not use the product while sleeping. Do not use the product together with electronic health monitoring devices. Do not place the electrodes on the neck or in the chest area. Exercise extreme caution when using the product. If you feel any discomfort while using the product, stop using the product immediately. Failure to do so may result in serious injury. Protect against falls, shocks, mechanical damage. Do not immerse in water. Do not handle the product with wet hands. Keep the product clean. Do not stick overlays on wet, oily skin. Do not touch the adhesive part. Place the protective caps on the electrodes during storage. Clean with a soft, damp cloth to remove any remaining sweat or dirt from the surface of the device. Do not use bleach, corrosive cleaners, sharp objects to clean.

Do not expose the device to extremely high or low temperatures. Protect the product

from direct sunlight. Store the product in a dry, shaded place away from heat sources such as radiators. If you do not use the product for a long time, remove the batteries. Use the product as intended. Do not use if you notice that the product is damaged. The product is not a toy. Keep the product away from children.

Information for users on the disposal of electrical and electronic equipment (for households) The symbol shown on the products, or the accompanying documentation indicates that faulty electrical or electronic equipment cannot be disposed of together with household waste. The correct procedure, if necessary, to dispose of electrical or electronic equipment, to dispose of, reuse or recover components, is to transfer the

equipment to a specialized collection point, where it will be accepted free of charge. In some countries, you can give the product back to your local distributor when you purchase another device. Proper disposal of the equipment preserves valuable resources and avoids negative health and environmental impacts that may be compromised by inadequate waste management. Please contact your local authority for details on the nearest collection point. Improper disposal of waste is punishable by penalties provided for in the relevant local legislation. **Users in European Union countries: If you need to dispose of electrical or electronic equipment, please contact your nearest point of sale or your supplier for further**

information. **Disposal of waste in countries outside the European Union:** This symbol is only valid in the European Union. If you need to dispose of this product, please contact your local authorities or dealer for information on how to proceed correctly.

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