E-SHAWL 60W

PL

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

Mission Air

ENG

ELECTRIC HEATING PILLOW FOR SHOULDER AND NECK

You can find instructions in other languages on our website.

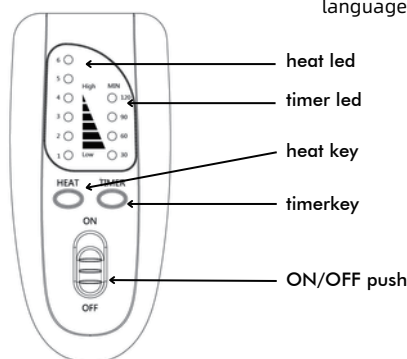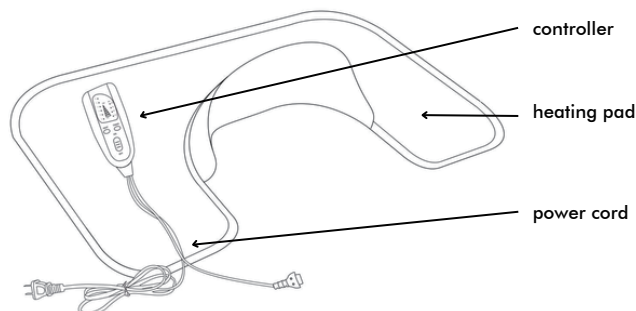CONTROL

1. Make sure the plug is connected to the power supply.
2. Slide the ON/OFF switch to the ON position. The heating pad will start working in the default mode (48°C, 60min timer).
3. Press the heating button to adjust the temperature to your needs. The LED indicator will show the current heating level. 6 temperature settings: 40 / 45 / 48 / 52 / 56 / 60°C.
4. Press the timer button to set the working time. The LED indicator will show the currently selected time program. 4 TIMER function settings: 30 / 60 / 90 / 120min.
5. Slide the ON/OFF switch to the OFF position to turn off the device.

CLEANING AND STORAGE

Turn off the device and unplug it from the power source. Unplug the controller from the pillow. The pillow can be washed by hand or by machine in a protective bag at low temperature. It is recommended to clean the heating pad only when necessary. If it has been washed frequently, the product may wear out quickly. We recommend using a cloth or damp sponge to remove small stains and add a little liquid detergent to clean the fabric. Do not use cleaning agents or chemicals. Only natural drying is allowed, do not tumble dry, do not iron. The switched off device should be stored in a dry, clean, cool and safe place.

CLEANING AND STORAGE

This is an electrical product. Before using the neck and shoulder heating pad, please make sure you read and understand the following instructions. Failure to follow them may result in electric shock or fire causing serious injury.

1. Read all instructions before using this product.
2. For indoor use only.
3. DO NOT use if any damage is visible (cracks, tears or exposed wires and heating element).
4. Immediately turn off the heating pad and unplug it from the power supply if it overheats or becomes charred during use.
5. Serious heart patients, pacemaker users and pregnant women should consult a doctor before use.
6. DO NOT use on unconscious and heat-sensitive individuals.
7. DO NOT make any modifications to the switch button, wires and heating element.
8. DO NOT use while sleeping; do not place the heating pad under a mattress or quilt during use.
9. DO NOT use the wire as a handle.
10. DO NOT pierce or use any sharp objects to secure the heating pad.
11. DO NOT use another electronic heating device in conjunction with the heating pad.
12. DO NOT roll up the heating pad during use.
13. Cover the heating pad only with your body, DO NOT sit or lie down on the heating pad.
14. DO NOT place or cover anything on the heating pad during use.
15. NOT for use by children.
16. Unplug when not in use.
17. DO NOT use if wet.
18. DO NOT dry clean this heating pad. Dry cleaning solvents may damage the wires.
19. DO NOT use bleach or cleaning fluids.
20. DO NOT use a wringer.
21. DO NOT iron the heating pad.
22. DO NOT use mothballs or any type of spray.
23. DO NOT reconnect the heated product or use it until it is COMPLETELY dry.
24. DO NOT use pain relieving cream or other body cosmetics in combination with the heating pad.