
zenwire


## MagniSonic Portable Oral Irrigator WF309
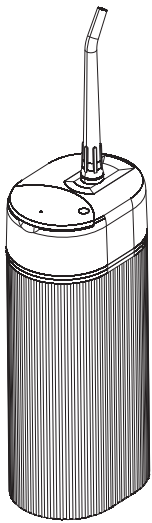
### User Manual

01


### Product Descriptions

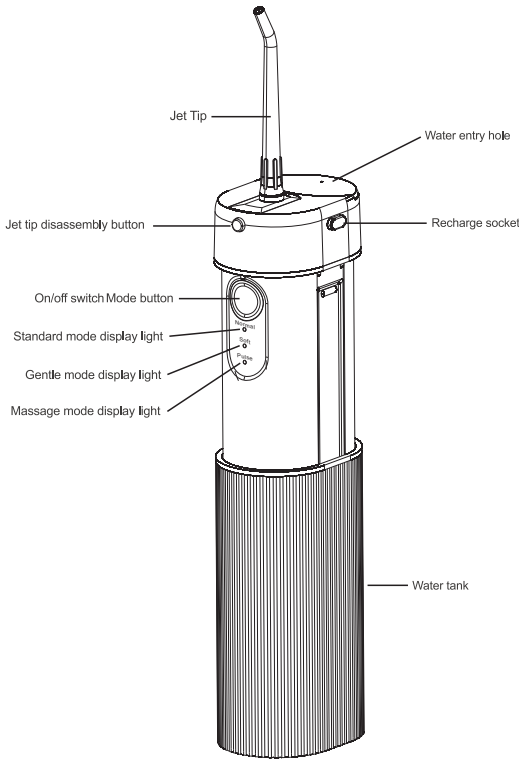02

### Product Accessories

#### Carry pouch x1

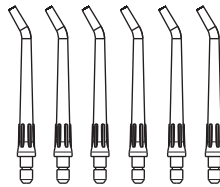Jet tip x6


Cable USB-C

### Operating Procedures

#### Preparation


Important: Before using the device for the first time, please fully charge it. Insert the USB plug into a power socket and ensure it is connected to a power source. Note: The device cannot be used while it is charging

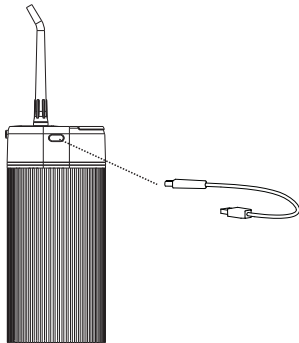03

### Installation and Removal of Jet Tip

Insert the jet tip into the spray head hole by pushing or twisting it in until you hear a click indicating it's securely in place.

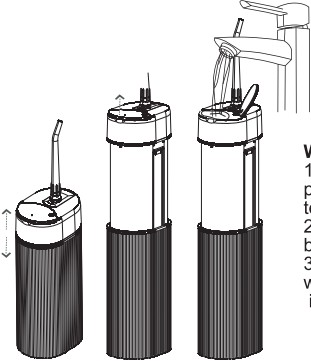To remove or replace the jet tip, press the release button on the handle, then pull the tip out.

04

### Usage Instructions


1. Press the button repeatedly to cycle through the following modes: Power On → Standard Mode → Soft Mode → Massage Mode → Power Off
2. Memory Function: When powered on, the device will start in the last mode used before automatic shutdown.
3. Battery Level Display: When turning off, the mode lights will briefly show the approximate remaining battery level before fading out.
4. Auto Shutdown: Any mode will automatically turn off after 60 seconds of operation.
5. Mode Switching: Press the button again within 15 seconds to switch modes. If no button is pressed for more than 15 seconds, the device will shut down.
6. Low Battery Indicator: When the massage mode light flashes, the battery is low and needs charging.
7. Charging Indicator: When charging is complete, all three lights will stay on continuously.

**Important:**  
The device cannot be used while charging.


05

### Recommended Techniques

Hold the device at a downward angle and place the jet tip inside your mouth, aiming it toward your teeth before turning the device on.

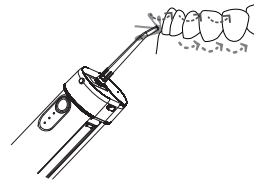Position the water flow so it forms a 90-degree angle with your gums. Slightly close your mouth to prevent water from spraying out and to achieve the best cleaning effect. Start by gently clenching your rear teeth, then move the jet tip to clean your front teeth (incisors). Slide the water flow along your gums, pausing briefly between each tooth, until you have thoroughly cleaned the gumline and the spaces between your teeth.

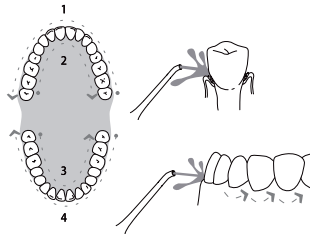
06

### Disassembly Example


#### After Use

1. Turn off the device.
2. Open the water tank door and tilt the tank so that the water injection hole faces downward.
3. Place the water tank back into the device.
4. Turn the device on to flush out any residual water inside.
5. Keep the device running until no more water comes out.
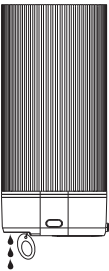
#### Removing the Water Tank for Cleaning

1. To separate the water tank from the main body, hold the water tank with one hand and the device body with the other.
2. Press upward gently and pull the water tank out.
3. The removable water tank can now be cleaned inside.

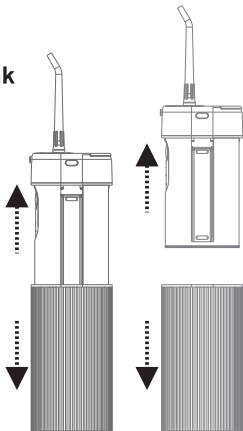07

### Warnings

1. Keep out of reach of children and animals. This product is not a toy and should be operated by adults only.
2. Do not use while charging.
3. Do not disassemble or attempt to modify the device.
4. Avoid direct contact with eyes or ears during use.
5. Use clean water only. Do not use mouthwash, saltwater, or other additives unless specified.
6. Stop using immediately if the device malfunctions or shows visible damage.
7. Do not expose to extreme temperatures or fire.
8. The device is not intended for use in water.
9. Fully charge before first use (approximately 4 hours).
10. If storing for an extended period, charge fully in advance.
11. Recharge promptly when power is low to preserve battery life.
12. Contains a built-in non-replaceable lithium battery.
13. At the end of battery life, dispose at an approved recycling center. Do not throw in regular trash.
14. Avoid using extremely cold or hot water. It is recommended to use water at a temperature of 35–40°C.

08

### Common Problems And Solutions

| Problems                                     | Reason                                                                         | Solution                                                          |
|----------------------------------------------|--------------------------------------------------------------------------------|-------------------------------------------------------------------|
| Water leakage between jet tip and screw part | Seal ring on the fixed socket is damaged, or jet tip is not installed properly | Contact retailer for repair, and reinstall the jet tip correctly. |
| Water tank leaks water                       | Seal ring is loose or not installed properly                                   | Reinstall the seal ring correctly.                                |
| Device cannot be switched on                 | Low battery voltage                                                            | Recharge the device before use.                                   |
| Product is not discharging water             | Filtration net is clogged                                                      | Clean the filtration net.                                         |
| Water leaks from sides of water tank         | High air pressure causing water leakage (normal phenomenon)                    | This is normal due to air pressure. No action needed.             |

09

E-mail: sklep@zenwire.eu  
Manufacturer: Zenwire Sp.z.o.o. Kołobrzaska 4,  
Poznań 60-418, Made in PRC.

The symbol shown on products or accompanying documentation indicates that faulty electrical or electronic devices must not be disposed of with household waste. Proper handling for disposal, reuse, or recycling involves taking the device to a specialized collection point where it will be accepted free of charge. In some countries, the product can be returned to the local distributor when purchasing another device. Proper disposal of the device allows for the conservation of valuable resources and prevents negative impacts on health and the environment, which could be endangered by improper waste handling. Detailed information on the nearest collection point can be obtained from local authorities. Improper disposal of waste is subject to penalties provided by relevant local regulations.

#### Business users in the European Union

If you need to dispose of electrical or electronic devices, please contact your nearest point of sale or supplier for additional information.

#### Disposal of waste outside the European Union

This symbol is only valid in the European Union. If you need to dispose of this product, please contact local authorities or the seller for information on the correct method of disposal.

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