VEVOR[®]

TOUGH TOOLS, HALF PRICE

Technical Support and E-Warranty Certificate www.vevor.com/support

VEVOR Derma Pen Model:M8

We continue to be committed to provide you tools with competitive price. "Save Half", "Half Price" or any other similar expressions used by us only represents an estimate of savings you might benefit from buying certain tools with us compared to the major top brands and does not necessarily mean to cover all categories of tools offered by us. You are kindly reminded to verify carefully when you are placing an order with us if you are actually saving half in comparison with the top major brands.

VEVOR[®]

TOUGH TOOLS, HALF PRICE

Derma Pen

Model:M8

NEED HELP? CONTACT US!

Have product questions? Need technical support? Please feel free to contact us:

Technical Support and E-Warranty Certificate
www.vevor.com/support

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

Parts Name

How to operate

Wireless Model:

The instrument used with both built-in battery or AC power

Wired Model:

The instrument used with with AC power only

- 1.Plug in DC adaptor to the jack
- 2.Plug in DC adaptor to AC outlet
- 3.To charge the battery:Plug the adapter into an outlet and connect it to the power interface

Charging Status


- The power indicator stays red while charging
- The power indicator turns green after charging completed

How to operate

Preparation

- 1.Take out needle cartridge from the sterilized packaging
 - 2.Fitting slot of cartridge into the pen correctly, push needle cartridge horizontally up to the end.
 - 3.Lock the needle cartridge by turning it to the counterclockwise
 - 4.Turn the needle cartridge clockwise and then pull out horizontally if you want to remove it. [Fig1]
 - 5.Start the device
- Press the power button for 2-3 seconds to start or shut down the device.

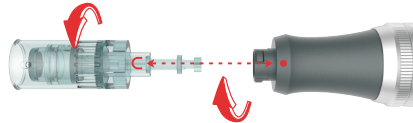[Fig1]

Before assembling or disassembling the needle cartridge, make sure the U-shaped point of the needle cartridge is aligned with the point of the handpiece

Attention

In order to protect the pen well,do not invert the instrument to prevent the liquid backflow

VEVOR[®]

TOUGH TOOLS, HALF PRICE

Technical Support and E-Warranty Certificate
www.vevor.com/support

How to operate

Adjustment of needle length

Adjust the depth of needle by turning the scale adjustment ring. The depth of needle should be adjusted after start the pen.

Turn it clockwise,needle depth will be shorter
Turn it counterclockwise,needle depth will be longer

Speed Display

6 Speeds Adjustment

The pen starts from low to high speed

The speed changes once each time you press the button

Correct Disposal:

This product is subject to the provision of european Directive 2012/19/EU. The symbol showing a wheelie bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to acollection point for recycling electrical and electronic devices.

FCC Information:

CAUTION: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment!

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:


- 1) This product may cause harmful interference.
 - 2)This product must accept any interference received, including interference that may cause undesired operation.
- WARNING: Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.
- Note: This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules, These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception,which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the distance between the product and receiver.
- Connect the product to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for assistance.

Needle cartridge introduction

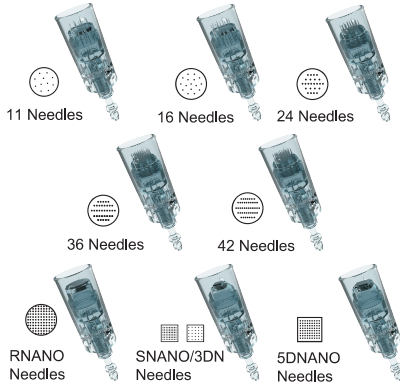Attention

Different models of the product vary in needle cartridge configuration. Please refer to the actual product for details.

The needle cartridge is disposable.

The needle cartridge is designed for M8 only.

Please sterilize the needle cartridge with 75% alcohol for a few minutes before use.

Please dispose the used needle cartridge to the trash container.

Reference chart

REGION	NEEDLE DEPTH (mm)
Forehead	0.25 - 0.5
Between the Eyebrows	0.25 - 0.5
Nose	0.25
Around the orbital rim	0.25
Cheek and Chin	0.5 - 1.0
Surgical and Scar Maximum	2.5

Specification

Product name : VEVOR Derma Pen

Model name : M8

Power adapter input: AC100-240V 50/60Hz

Power adapter Output: DC5V 1A

Battery Capacity: 380mAh (lithium-ion battery)

Packing List

Item	Quantity
Hand Piece	1 pc
Adapter	1 pc
Manual	1 pc
Power Cable	1 pc

Precaution

	Please be sure to read instructions before use and follow as described.
	Do not apply it to the man as follow: 1.Who has acute disease or bloody disease over 38℃ 2.Who has allergy and hemophilia 3.Who has melanosis of the face 4.Who has high internal heat 5.Who regularly visits hospital 6.Who takes medicines 7.Who is under 12 years old
	Do not use it on the power other than AC 220V. It may cause fires or malfunction
	Do not use it near the place with high temperature, dust or humidity (Near the humidifier).It may cause malfunction.
	Do not repair, disassemble or alter the controller and the body. It may cause malfunction.
	Do not touch the device or detach the power cords with wet hands. It may cause electric shock.
	Do not use the stripped cord or propped cord. It may cause fires(Be sure to replace it with another).
	Do not use the propped cord which has crack or scares. It may cause scares on the skin.
	When you don't use it for a long time, unplug it from the outlet.
	Be sure to have an inspection for test. If any troubles are found, get it repaired.
	No one other than the operator can touch it. It may cause accidents or malfunction.

	Alternating current		Serial No
	Direct current		EU Representative
	Class I equipment		CE mark
	Attention, consult accompanying documents		Manufacturer
	Date of manufacture		Floating applied part

Warranty

The warranty period of device is 12 months.
We offer repair to our products from the day they are purchased. If any of our products has damages or natural trouble. We offer repair or replacements

General maintenance

Make sure to inspect the device and adapter before use

Keep it with its power off.Wipe its surface out with a dry cloth

Keep it away from direct sunlight or high temperature and humidity

Do not leave the measurement, repair or alter of the device to non-professional staff

When you use the device which has not been used for a long time, make sure to check if it works properly with safety

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