BERDSEN

GALAXY
SONIC
TOOTHBRUSH

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DE SCHALLZAHNBÜRSTE

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







322000 / 322001 / 322002 / 322003 / 322020 / 322021 / 322022 / V. 1.1

www.berdson.pl

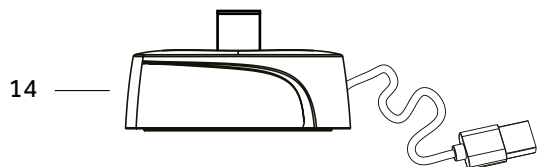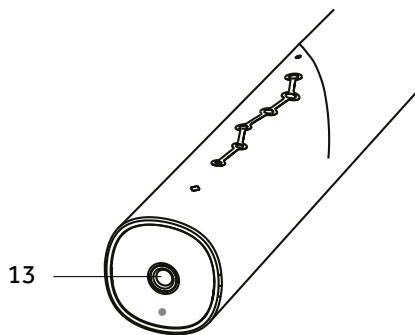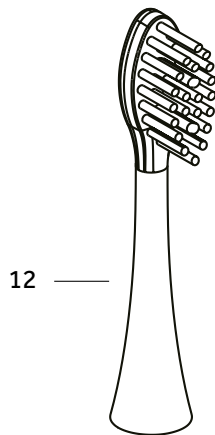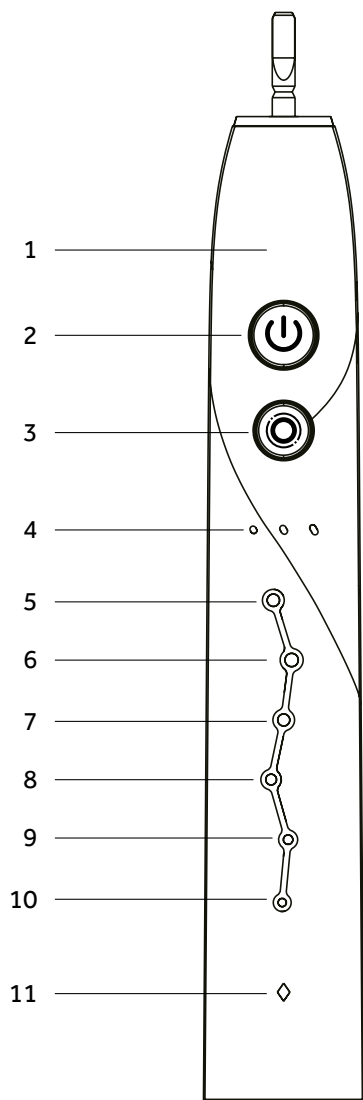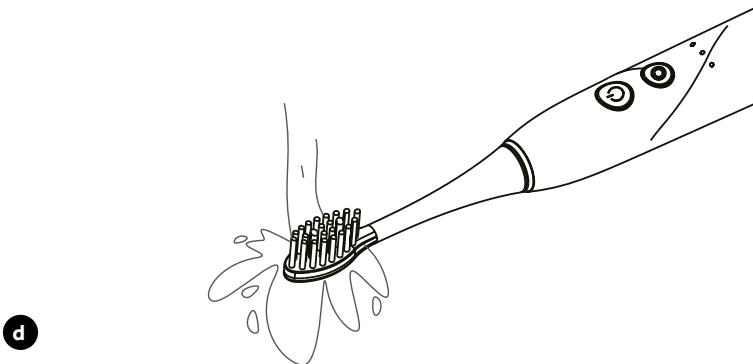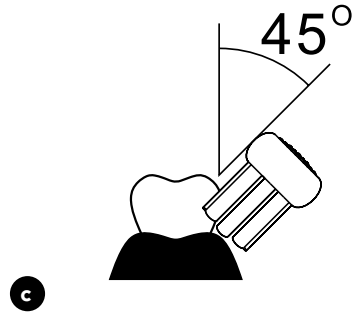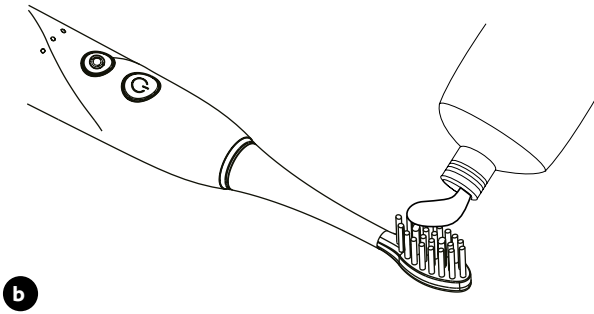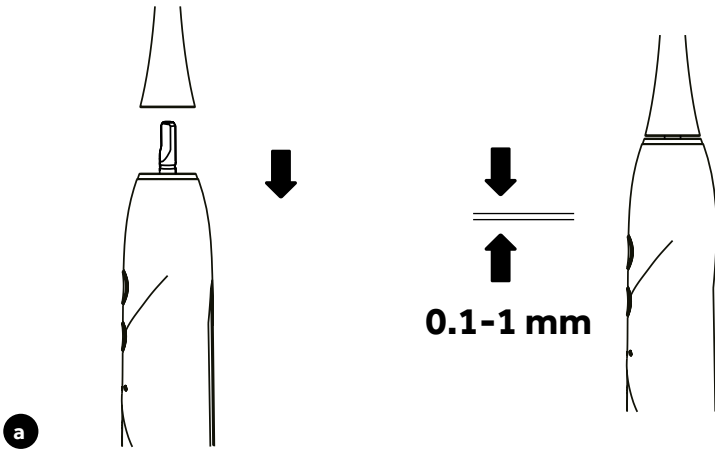Thank you very much for choosing a Berdsen product. This manual contains important safety, usage, and disposal guidelines for the device.

WARNINGS!

- Use the product according to its intended purpose and the recommendations in this manual. The device is designed for cleaning teeth, gums, and the tongue. Do not use any accessories or attachments other than those provided and recommended by the manufacturer.
- The toothbrush is for personal home use only and is not intended for use in a dental office.
- Never use the device if the cord, plug, or charger is damaged.
- Avoid contact of the charger with water or other liquids. Risk of electric shock!
- Do not place the charger in water, near a bathtub, sink, or in their vicinity.
- If the charger falls into water or another liquid, do not attempt to retrieve it! Immediately disconnect the power cord.
- Disconnect the charger before cleaning. Before plugging the charger into an outlet, ensure it is dry.
- Unplug the charger by pulling the plug, not the cord. Avoid placing cords under furniture, on sharp edges, or near heat sources. If the cord is damaged, replace it at an authorized service center. Connect the charger only to an outlet with the appropriate voltage.
- Do not open the device casing or replace the battery yourself—only authorized personnel may perform replacements.
- The device may be operated by children over 8 years of age and by persons with reduced physical, sensory, or mental capabilities only under supervision or after proper training.
- At the beginning of product use, you may experience irritation or tingling of the teeth due to the high vibration speed, and even gums may begin to bleed. This is a normal phenomenon. If this condition persists after two weeks of continuous use, stop using the product and consult a dentist.
- If you have undergone oral or gum surgery, had dental issues, or received extensive oral treatment within the last two months, consult a dentist before using the product.
- If you experience pain, stop using the product immediately and consult a dentist.
- If the toothbrush has not been used for an extended period, the indicator may not light up immediately during charging. Wait patiently for a few minutes until the device returns to normal operation.
- Heating during charging is normal, but if the toothbrush feels too hot, disconnect it from the power supply immediately and take it to a service center for inspection.
- Use the product only indoors. Do not expose the toothbrush to direct sunlight. Store it at room temperature or in a dry, ventilated place.
- The battery should be disposed of safely.

DEVICE STRUCTURE (FIG. I)

- | | |
|----------------------|-----------------------------|
| 1. Handle | 8. Gum care mode |
| 2. Power button | 9. Sensitive mode |
| 3. Mode/Power button | 10. Deep cleaning mode |
| 4. Power indicator | 11. Charging indicator |
| 5. Cleaning mode | 12. Brush head (3 pcs.) |
| 6. Whitening mode | 13. Charging port |
| 7. Polishing mode | 14. Inductive charging base |

OPERATION (FIG. II)

- Make sure the brush head is correctly attached to the handle (Fig. a).
- Wet the brush head and apply toothpaste to the bristles (Fig. b).

STARTING/CHANGING MODE/TURNING OFF

- When off, press the power button to turn on the toothbrush.
- **Changing mode:** After turning on, briefly press the mode/power button to switch between different brushing modes.
- **Changing power:** After turning on, hold the mode/power button for 1.5 seconds to change the brushing power to low, medium, or high.

BRUSHING MODES

You can choose the most comfortable brushing mode. This toothbrush has 6 modes (each mode additionally has three power levels: high, medium, low):

- **Cleaning mode:** Suitable for most people.
- **Whitening mode:** Suitable for people with teeth discoloration caused by stains from smoke, tea, and other dark or yellow pigments.
- **Polishing mode:** Suitable for people with tartar and plaque.
- **Gum care mode:** For those with a tendency to bleed gums.
- **Sensitive mode:** Suitable for people with sensitive gums and teeth.
- **Deep cleaning mode:** Suitable for people with tartar and plaque.

TURNING OFF

- When the device is operating, touch the power button to turn it off.
- After 3 seconds of operation, if you turn off the toothbrush, it will remember the previous mode and power (strong, medium, or low), setting them as default for the next use.
- After each shutdown, the power indicator will light up for 2 seconds in each operating mode, displaying the remaining battery level.

PROPER BRUSHING TECHNIQUE

- Hold the bristles at a 45-degree angle from the gum line (Fig. c). Start with circular motions, slowly and gently rubbing the space between the gums and teeth.
- Continue moving the brush head to clean all teeth.
- During brushing, the toothbrush's vibration effect will briefly pause every 30 seconds to remind you to change brushing position.
- The best time to brush teeth is 2-3 minutes.
- The toothbrush will automatically turn off after 2 minutes. If you want to continue brushing, turn it on again.

CLEANING THE BRUSH HEAD

- After brushing your teeth, clean the brush head bristles under running tap water (Fig. d).
- After cleaning the brush head, place the electric toothbrush in a dry and cool place.
- It is recommended to replace the brush head every three months.

CHARGING

- When the battery level is low, the charging indicator flashes red, or the toothbrush cannot be turned on.
- Connect the device for charging.
- During charging, the charging indicator will pulse red.
- Full charging takes about 8 hours.
- When the battery is fully charged, the charging indicator will light up green for 30 seconds, then turn off.

Note: The original charging cable is dedicated to charging this sonic toothbrush model; do not use it for other purposes.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Cannot turn on the product	Low battery level	Charge the toothbrush
Cannot attach the brush head to the handle	Contamination in the brush head opening	Remove the contamination
	Contamination on the handle's drive shaft	Clean the contamination
	Inappropriate brush head	Select the correct brush head
Short operating time	Excessive use	Charge the toothbrush
Weaker brush head vibrations	Improper assembly of the brush head and handle	Remove any contamination and reattach the brush head
	Low battery level	Charge the toothbrush

ENVIRONMENTALLY FRIENDLY USAGE PERIOD

The device has a ten-year environmentally friendly usage period, during which there is no risk of leakage or mutation of toxic substances, making it safe for health and the environment. This period also applies to all included accessories. After this period, the product should be returned to a recycling collection point.

SPECIFICATIONS

Power	2 W
Input	5 V $\pm$ 0.4 A
Battery	Lithium 3.7 V 800 mAh
Connection	USB
Charging time	about 8 hours
Water resistance rating	IPX7

INFORMATION FOR USERS ON HOW TO DEAL WITH WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT.

If the device, packaging, instruction manual and the like are marked with the symbol of a crossed-out wheeled waste container, it means that the device is subject to separate waste collection in accordance with Directive 2012/19/EU of the European Parliament and of the Council. The used device should not be disposed of with other household waste and should not be treated as municipal waste. By disposing of electro-waste in the container you pose a threat to the environment. Take the used device to a separate collection point for electrical and electronic equipment organized by the public administration. By segregating and submitting used electrical and electronic equipment for processing, recovery, recycling and disposal, you protect the environment from pollution and contamination, contribute to reducing the use of natural resources and lowering the cost of producing new equipment. Proper disposal and scrapping helps eliminate the adverse effects of scrapped products on the environment and human health. For details on the recyclability of this product, please contact your local municipality, city cleaning services or the store where you purchased the product.

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