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
Eye Beauty Massager

ZH-HK

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EN

Instruction Manual

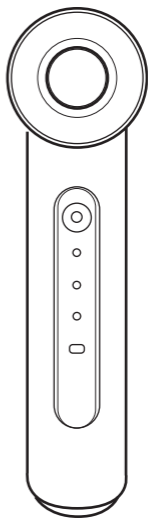**Ragic** 朗見

<b>ZH-HK</b>	<b>說明書</b> .....	<b>2</b>
<b>EN</b>	<b>Instruction Manual</b> .....	<b>8</b>

---

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---

Congratulations on the purchase of your Rogic Eye Beauty Massager. We hope it brings you many years of satisfaction.

Should you encounter situations that have not been properly addressed in this Instruction Manual, please contact our Customer Service Center where a technician can address your specific problem or question.

For more information, please visit the official website:  
[www.rogic.com.hk](http://www.rogic.com.hk)

The Company reserves the right to make technological and/or design changes to this product for continuous improvement.

Thank your for your choosing!

# Contents

---

1. Important safety Instructions . . . . .	8
2. Package Contents and Technical Specifications . . . . .	9
2.1 Technical Specifications . . . . .	9
2.2 Product Design . . . . .	9
3. Usage . . . . .	10
4. Problem Shooting . . . . .	12
5. Customer Service . . . . .	13

# 1. Important Safety Instructions

**READ ALL INSTRUCTIONS BEFORE USING THIS APPLIANCE.  
PLEASE KEEP THE INSTRUCTION MANUAL**

When using our product, basic precautions should always be followed, including the following:


## **i** Attention

01. Do not store it within the reach of children and infants, and do not allow children or infants to play with the appliance.
02. Please unplug the USB and main body while not using it for a long time.
03. Do not use alkaline detergent, alcohol, thinner, gasoline, nail polish remover, etc.
04. Please clean the massage head after each use.
05. Keep it out of reach of children and in a well-ventilated place.
06. When not in use for a long time, remove the power adapter and USB cable from the device and socket for storage.
07. Do not place in hot and humid places (near fire, under hot sun, bathroom, etc.)
08. If the product has not been used for a long time, please make sure that the product is in normal and safe condition before use.
09. In case of abnormal operation of the product, please contact the dealer for maintenance.
10. Please avoid smearing skin care products with your hands to keep the skin care products clean or avoid slipping of the device due to hand slip.
11. If you feel not smooth in the process of using, please add appropriate skin care product.
12. Please completely touch the electrode tip to the skin when using.
13. Be sure to move it when use on skin.
14. When you feel hot, move a little faster or lower the intensity.
15. Please move a little faster when you feel hot.
16. Rechargeable batteries are consumables.
17. After charging, the use time is very short, the battery life can be considered. (Normally about 300 charges)
18. If you need to replace the battery, contact customer service hotline.
19. If the product is not used for a long time, please charge it at least once every 3 months.
20. It is normal that product heats up during use or charging.
21. Be sure to take out built-in rechargeable battery before waste.
22. When removing the battery, turn off the power and pull out the plug.
23. When removing the battery, use a screwdriver to avoid injury.
24. When removing the battery, do not scratch it with hard objects. May cause batteries to heat up and catch fire.

## **i** Remove the battery

Remove the battery by referring to steps ①~② in the following figure.

- ① Use a tapered tool to pry open the back shell.

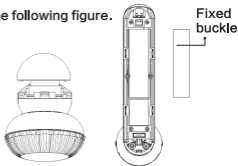
## 2. Package Contents and Technical Specifications

### 2.1 Technical Specifications

Model	RE100
Voltage	5V == 1A
Battery	Li-ion battery 750mAh
Charging Time	~ 60 minutes (Due to the charging environment, the charging time is slightly different.)
Lasting Time	Mode 1: about 100 minutes ; Mode 2: about 360 minutes ; Mode 3: about 70 minutes (full battery state begins, depending on the environment)

\* Package Contents and Technical and design specifications may be changed in the course of continuous product improvement.
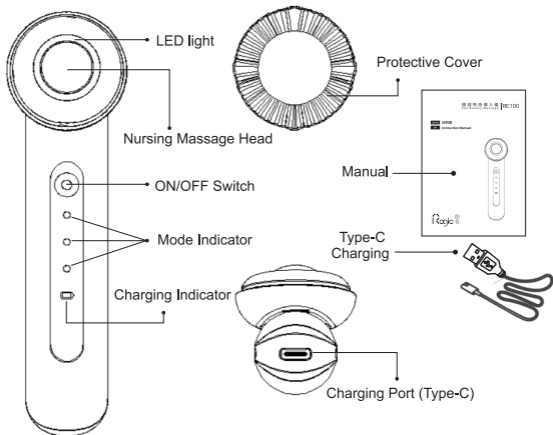
## 3. Usage

### Before use

1. Plug the Type C connector of the charging cable into the charging port of the device.
  2. Insert the USB plug of the charging cable into the USB charger (mobile phone charger can be used).
    - For new products purchased for the first time, please charge them before use.
    - During use, the charging indicator **blinks**, indicating that the power is low and the machine is about to be powered off. Please charge it.
    - During charging, the charging indicator **blinks**.
    - After charging, the charging indicator **will be on for 3 minutes and then off**. Please remove the USB cable.
    - At room temperature around 25°C, it takes about 1 hour to **fully charge and lasts about 70 minutes**. The charging time and usable time are **slightly different** due to the remaining power and charging environment.
- \* If it still cannot charge normally after 10 seconds, please confirm whether the socket of the USB cable is loose

### Skin fitness test

Be sure to test the product's compatibility with the skin in advance when using for the first time.

- Be sure to test beforehand to prevent skin problems.
  - Stop the test immediately if any abnormal conditions occur at the test site, such as diagnosis, redness, itching or irritation.
  - Do not use when skin feels abnormal.
- ① Clean the inner side of the arm with soap and other products, and then wipe it with a towel.
  - ② Directly squeeze the appropriate amount of skin care products into the central part of the nursing massage head, paste it on the washed arm, and use it for about 3 minutes.
  - ③ Massage the excess skin care products until absorption, and then confirm whether the skin is abnormal.
  - ④ If the skin is not flushed or hot within 24 hours, you can use this product.
- Make sure the electrode tip is completely attached to the skin of the hand and moved while using.

### Safety features

- Temperature sensor  
When the temperature of the skin or massage head is detected to be high, the main machine will automatically reduce or turn off the output for safety.
- Intelligent function  
Skin contact sensor intelligent sensing, item start to work soon as skin contact.
- Automatic power off function  
The machine will automatically shutdown (with sound "beep beep") if no operation after 3 minutes. More energy saving and environment protection.
- Charging tips  
When the power is on, the power indicator indicates the battery level. When the white light is blinking, please charge it before using it. Battery power varies slightly depending on the state of the battery and the environment in which it is used.

### Basis function

- **Constant Thermal Care**  
42°C constant warm temperature promote blood circulation, relax and brighten skin around eyes.
- **High-frequency Micro Vibration**  
High frequency micro vibration about 210 times per second can activate dermal fibroblasts, promote collagen regeneration and improve elastic fibers.
- **Iontophoresis**  
High-tech intelligent ion introduction technology, promote skin care essence penetrate into the deep muscle bottom, revitalize skin cells, smooth away age traces.
- **630 (±10) nm red light**  
Red upgrade energy improves dullness and brightens skin tone.

### Using Method

Best used to use it after removing make up and cleansing. It is recommended to use with eye cream and the skin care products.

- ① **ON/OFF** Long touch the power button for about 2 seconds, the default mode 1-indicator lights up, constant temperature thermal care start running.
- ② **Mode Selection**

Modes	Mode Instruction	Mode Duration
	Mode 1, indicator lights up thermal care red light therapy, ion induction	5 minutes
	Mode 2 indicator lights up high frequency micro vibration ion induction	5 minutes
	Mode 3 indicator lights up thermal care high frequency micro vibration red light therapy, ion induction	5 minutes

- ③ Squeeze skin care products directly in the center of massage head, and move on the skin.
- ④ skin care after use

General use recommendations, cleaning face with water after skin caring, and then do daily skin care and maintenance.

### Maintenance

- Keep it out of reach of children and in a well-ventilated place.
- When not in use for a long time, remove the power adapter and USB cable from the device and socket for storage.
- Do not place in hot and humid places (near fire, under hot sun, bathroom, etc.)

## 4. Problem Shooting

If you feel abnormal during using this product, please confirm the following items before submitting it for maintenance.

The following state	Possible causes and handling methods
Unable to boot	Touch the power button for about 2 seconds
	The battery is empty, please charge it
Even after charged, it only lasts a few minutes	The machine has the function of automatic power off, it will automatically turn off the power after 3 minutes
	The battery is low, please charge it
	Battery life issues. Please contact the dealer
Can't be charged	Insert the original USB cable into the device. Check whether it is loose
Constant Temperature overheating	Use on the same area for a long time. Be sure to use it while moving
	Please move a little faster
Constant temperature is not felt	The machine shuts down after 3 minutes. Please press the power supply to start the machine again
The white power indicator blinks	Indicates that the remaining battery is low. Please charge it before using it

### FAQ

Q: What is the standard of usage time and frequency?

A: Use 5 days a week for better results (no more than 6 minutes once a day)

※ According to your skin condition, increase or decrease the frequency of use.

※ When skin feels abnormal, please stop using immediately.

Q: Can it be used continuously?

A: It is recommended to use this product once a day for 6 minutes.

Q: How many skin care products should be used at one time?

A: Please refer to the dosage diagram of P.14. If you feel less skin care products, you can add appropriate amount at any time.

Q: Can you use it in the bathroom?

A: This product is not waterproof. Please avoid using in the bathroom.

Q: Skin hot after use.

A: This is caused by the constant temperature warming effect.

## 5. Customer Service

---

### WARRANTY CONDITIONS

Thank you for purchasing ROGIC products, for the sake of excellent quality of our products, we certify that ROGIC products on this warranty is free of defect in materials and workmanship the end user is entitled to this warranty subject to the terms and conditions listed below:

- (1) For any defect, in the judgement of our technicians, caused under normal usage, is responsible for repairing and replacing parts free of charge within the first year free warranty period (except accessories and some parts), any replaced defective parts will become the company's property.
- (2) For the information of designated product, please visit [www.rogic.hk](http://www.rogic.hk) for details.
- (3) This warranty is invalid in the following cases:
  - (a) If the product has been damaged out of misuse, negligence, accident and natural disaster.
  - (b) If the product has been repaired, modified by persons other than PRO LOGIC INTERNATIONAL LIMITED.
  - (c) If the serial number has been altered, effaced and/or removed.
  - (d) Fail to follow the installation instructions provided.
  - (e) Machines for commercial use.
- (4) The company will not be liable for any loss directly or indirectly caused from breakdown of the equipment.
- (5) Customer will be requested to bring the product to and from our service centre whenever service is required.
- (6) This warranty is applicable in HongKong and Macau areas only.
- (7) Please keep and present this warranty card together with official purchase invoice to our technicians for inspection during each service.

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Unit 106-107, 1/F, Heng Nagai Jewelry Centrem 4 Hok Yeun Street East,  
Hung Hom, Kowloon, HK

Tel: (852) 2365 0655

分銷商 Distributor:

愛家名品有限公司 iPlusHome Limited

[www.iplushome.hk](http://www.iplushome.hk)