BLACKROLL®

THERMO PAD

Instruction Manual

EN | DE | FR | NL | IT | PL | CN

Hi!

Congratulations on purchasing your new **BLACKROLL® THERMO PAD!**

This innovative heat technology supports your recovery, relaxation, and overall well-being—whenever and wherever you want.

Please read these instructions carefully before using the product for the first time to ensure safe and optimal use.

Keep the user manual in an easily accessible place, ideally directly with the product, so that you can refer to it at any time if necessary.

The **BLACKROLL® THERMO PAD** is a rechargeable heat pad that generates pleasant warmth of around 43 °C for up to 8 hours. At the touch of a button, it provides targeted, soothing warmth – for pure relaxation, comfort, and recovery.

How does the THERMO PAD work?

The device emits even heat, ensuring a pleasant, relaxed feeling in the desired target area. It is easy to use and reusable – ideal for use at home, in the office, or on the go.

Thanks to its durability and rechargeability, it is also a sustainable alternative to disposable heat products.

Note: The BLACKROLL® THERMO PAD is not a medical device.

The BLACKROLL® THERMO PAD Set includes:

1 x THERMO PAD (integrated heating unit with rechargeable battery)

3 x adhesive films (up to 15 uses possible – may vary depending on skin type)

1 x USB-C charging cable

1 x Quick Start Guide

Safety Instructions

- + Use the product only on clean, dry, and intact skin.
- + Do not use while sleeping or during charging.
- + Do not apply on open wounds, irritated, or inflamed skin.
- + Protect from water, moisture, and heat sources. Do not use while showering or bathing.
- + Use only the included 5V USB charging cable.
- + Do not open, puncture, or throw the device into fire.
- + Stop use immediately if redness, irritation, or discomfort occurs.
- + Not suitable for children or individuals unable to remove the device independently.
- + Clean only with a slightly damp cloth — do not use solvents, cleaners, or alcohol.
- + Improper use voids the warranty.

Contraindications

There are absolute contraindications that exclude the use of the device entirely, and relative contraindications where individual assessment is required.

“If one or more of these points apply to you, consult your doctor before using the device.”

Absolute Contraindications

- + Diseases such as diabetes mellitus, multiple sclerosis, peripheral vascular disease, spinal cord injuries, or severely impaired circulation.
- + Open wounds, heavily irritated, or inflamed skin.
- + Electronic or metallic implants.
- + Limited ability to perceive heat or temperature changes (e.g., diabetes, neurological conditions, advanced age).

Relative Contraindications

- + Current injuries, inflammation, or tendency for bruising.
- + Conditions affecting heat perception, sensitivity, or circulation.
- + Very sensitive or damaged skin (e.g., sunburn, skin conditions).
- + Pregnancy (do not use on the abdomen, lumbar, or pelvic region).

Implants

Wearing metallic or electronic implants is generally not a contraindication because no electrical energy or electromagnetic fields are used.

Caution is advised when:

- + the implant is near the surface
- + there are sensitivity issues in the application area
- + local inflammation is present

In these cases, prior medical clarification is recommended.

Product overview

The BLACKROLL® THERMO PAD generates even, pleasant heat of approx. 43 °C.

This heat application promotes relaxation, regeneration, and general well-being through localized heat in the desired target area.

Features


- + Up to 8 hours of operation
- + Constant deep heat ≈ 43 °C
- + Easy to use
- + Rechargeable battery via USB-C
- + Soft, skin-friendly, water-repellent silicone surface
- + Reusable adhesive pads for up to 3–5 applications per pad (depending on skin type)
- + Automatic temperature safety shut-off

Focused heat can help with:

- + recovery after physical exertion
- + Relaxation of muscle and tissue tension
- + a pleasant feeling of well-being in the muscles and joints
- + stimulation of local blood circulation and metabolism
- + general mobility and flexibility of tissue

Components

- A** Power Button (On/Off)
- B** USB-C Charging Port
- C** Status LED
- D** Controller
- E** Soft-touch silicone surface
- F** Heating layer

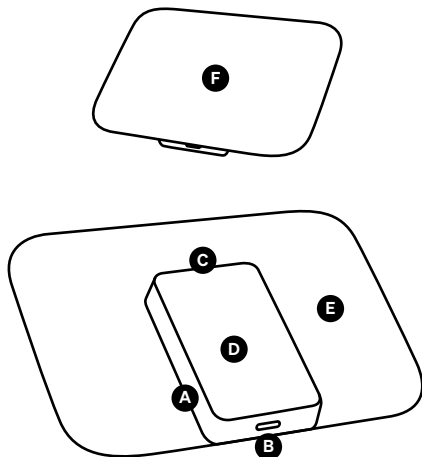
Application Notes

- + Ensure the pad sits firmly but comfortably on the skin.
- + Before first-time use, test for a few minutes on a non-sensitive body area to check personal heat tolerance.
- + Do not bend or fold.
- + Use only the included adhesive film.
- + Do not use simultaneously with other heat sources (heating pad, hot water bottle, sauna).
- + Do not apply creams, oils, or ointments directly before use.
- + Recommended duration: as needed, max. 8 hours (battery life).
- + Safety feature: device shuts off after ~8 hours or when temperature limits are exceeded.

Quick Start

1. Ensure the controller is fully charged.
2. Attach the adhesive film to the heating layer:
 - 2.1 Peel off the first half of the film.
 - 2.2 Stick on the first half of the THERMO PAD.
 - 2.3 Peel off the second half of the film.
 - 2.4 Stick on the second half of the THERMO PAD.

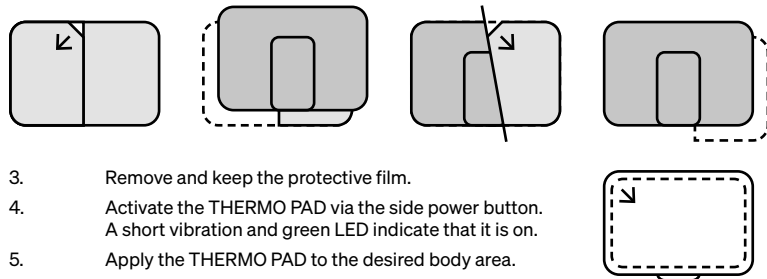3. Remove and keep the protective film.
4. Activate the THERMO PAD via the side power button. A short vibration and green LED indicate that it is on.
5. Apply the THERMO PAD to the desired body area.
6. After use: remove the pad from the skin and place it back on the protective film to prolong adhesive lifespan.
7. Turn the THERMO PAD off. A short vibration confirms shutdown. Recharge fully.

Troubleshooting

If problems occur during use, the following tips can help you quickly troubleshoot:

No heat generation

Cause: Battery empty
Solution: Fully charge the device. The status LED flashes red; once fully charged, the LED goes out.

Status LED does not light up

Cause: Device switched off
Solution: Press the power button to switch on the device. The status LED lights up green. If the status LED lights up red, please charge the battery.

Heat generation stops

Cause: Automatic safety shutdown active
Solution: Allow the device to cool down and then restart it. If the status LED is red, please charge the battery.

Battery does not charge

Cause: Incorrect cable or incorrect power supply
Solution: Use the supplied USB-C cable and a 5 V power adapter in accordance with IEC 62368; otherwise, contact support.

Surface damaged

Cause: Wear or improper cleaning
Solution: Stop using immediately and contact support.

Device not responding

Cause: Software error
Solution: Press and hold the power button for 10 seconds to reset the device.

Cleaning

- + Turn off the device and allow it to cool completely.
- + Clean the surface with a soft, slightly damp, non-abrasive disinfectant wipe.
- + Do not immerse in water or clean under running water.
- + Allow to dry completely before reuse.

Storage

- + Store in a cool, dry place away from direct sunlight.
- + Keep out of reach of children.
- + Do not place sharp or heavy objects on the product.

Technical specifications

Maximum temperature:	Approx. 43 °C
Operating time/battery life:	Up to 8 hours
Ambient temperature:	5 °C – 35 °C (optimal room temperature below 25 °C/relative humidity 10% to 90%)
Charging time:	Approx. 90 minutes
Soft-touch surface material:	Silicone
Battery type:	Lithium polymer (Li-Po), rechargeable
Nominal voltage in (V DC):	3.7
Maximum charging current in (A):	2
Maximum charging voltage in (V DC):	5
Energy content (Wh):	11.1
Battery capacity (mAh):	3000
Number of battery cells:	1
Battery weight (g):	Approx. 60g
Battery life:	Min. 1100 applications
Type of charging port:	USB-C
Is a charging cable included?	Yes
Power adapter included?	No, use a 5V DC USB power adapter that complies with the IEC 62368-1 standard for charging.
Film skin compatibility:	In vitro cytotoxicity (ISO 10993-5)
Film material:	Polyethylene terephthalate (PET)

Notes on the film

To date, there are no known allergies to the materials used in the adhesive tape or adhesive film. However, if skin reactions occur after contact with the heating film, discontinue use and consult a doctor.

Do not share the adhesive films with other people, as this may lead to skin irritation or infection. The films are intended for personal use only.

Disposal

This product must not be disposed of with household waste. Return it to a certified collection point for electronic devices and batteries to protect the environment. Note: The adhesive tape can be disposed of in a normal waste bin and does not require special recycling.

Disclaimer

Use of the device is at your own risk. It is recommended that you consult a specialist if problems arise before or after use, or if recurring symptoms occur.

Warranty

The seller grants the buyer, provided that the buyer is a consumer, a statutory warranty for a period of 2 years from the date of purchase of the item in accordance with the statutory provisions. The seller's warranty obligations are initially limited to subsequent performance, whereby the seller is free to choose whether to deliver a defect-free item or repair the defective item. If subsequent performance is unsuccessful, the buyer is entitled to further legal claims.

The seller rejects warranty claims for damage caused by accident, obvious misuse, abuse, or improper transport; damage caused by repairs carried out by unauthorized personnel; damage caused by use and operation not in accordance with the operating instructions; damage attributable to improper handling of the purchased item.

Furthermore, claims for damages by the buyer are excluded. This does not apply to claims for damages by the buyer arising from injury to life, limb, or health, as well as liability for other damages based on an intentional or grossly negligent breach of duty by the seller, its legal representatives, or vicarious agents. Liability for breaches of duty, the fulfillment of which is essential for the proper execution of the contract and on the observance of which the buyer may regularly rely, also remains unaffected.

In the event of a slightly negligent breach of such contractual obligations, the seller shall only be liable for foreseeable damage typical for this type of contract, unless the customer's claims for damages arise from injury to life, limb, or health.

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The device contains electronic components and batteries and must not be disposed of with household waste. Dispose of it at appropriate collection points according to local regulations.

The product complies with applicable EU directives for consumer products, including basic safety and health requirements.

Li-ion

Indicates that a rechargeable lithium-ion battery is included. Recycle the battery properly and do not dispose of it in household waste.

LOT

Indicates the production batch or serial identification to ensure product traceability.

The packaging and printed insert are recyclable and must be collected separately.

This product complies with the requirements of Council Directive 2014/30/EU on electromagnetic compatibility.

This device, its cables, and its battery are recyclable. Take the product to an approved recycling center.

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Indicates that the packaging is made of paper and can be recycled according to local regulations.

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BLACKROLL AG
WSVK Oederan GmbH
Produced by CALOPAD AG