


Available in 5  
different colors:  
Orange, Clear, Blue,  
Yellow, and Green

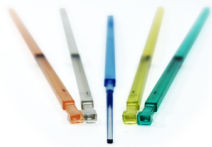

## PRODUCT SPECIFICATIONS

Product information is identical unless otherwise noted.

### Description:

- The body is a square shape stick made of medical grade resin, has a fine concave tip where the embryos are placed.
- The cap is made of the same resin, provides an airtight seal by the coupling of two tapered surfaces in a 0.250" of sealing surface.

### Dimensions:

- Cryolock®  
Body 4.56"L x 0.118"W x 0.118"H  
Tip width 0.050"  
Tip thickness 0.01"  
Cap 1.78"L x 0.118"W x 0.118"H
- S-Cryolock®  
Body 4.56"L x 0.094"W x 0.094"H  
Tip width 0.037"  
Tip thickness 0.01"  
Cap 1.78"L x 0.094"W x 0.094"H

### Performance:

- Cryolock®  
Cooling rate  $\cong$  -19,800°C/min  
Warming rate  $\cong$  27,000°C/min
- S-Cryolock®  
Cooling rate  $\cong$  -29,500°C/min  
Warming rate  $\cong$  40,500°C/min

### Certificate of Analysis:

- Available upon request
- 1 cell MEA  $\geq$  80% expand blastocysts within 96 h.
- Endotoxin LAL  $\leq$  2 EU/device.
- Sterility: 25-40 kGy (SAL10<sup>-6</sup>).

### Long-term storage:

- Biotech supported up to 3.5 years under liquid nitrogen, beyond that time is unknown.

### Notice to the user:

Any serious incident that has occurred in relation to this device should be reported to the manufacturer and the competent authority of the Member State in which the user is established.

For more information go to:

[www.cryolock.info](http://www.cryolock.info)


















Manufacturer by

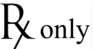

BioTech


5975 Shiloh Rd, Suite 101  
Alpharetta, GA 30005 USA  
1-800-313-7793

### Glossary of Symbols:

	Catalog Number: CL-R-CT for Cryolock and S-CL for S-Cryolock
	Lot Number
	Expiration date: Year-Month-Day
	Sterilized using radiation
	Single sterile barrier
	Do Not use of package is damaged
	Do Not- Re-Sterilize
	Do not Re-use
	US Caution: Federal law restricts this device to sale by or on the order of a licensed healthcare practitioner
	Medical Device
	Consult instructions for use
	Unique Device Identifier
	Country of Manufacturer
	Manufacture By
	CE Mark
	UKCA Mark
	EU Authorized Representative
<p>Atlantico Systems Ltd 34 Oldfield, Kingston Galway, Ireland +35391443609</p>	
<p>UK Responsible Person: Compliant Medical Devices Ltd 58 London Street, Whitchurch, Hampshire, RG28 7LN England +44(0)7581 317 331</p>	

EN ISO 15223-1:2021 – Symbols to be used on the medical device labels.  
(FDA only symbol: Rx only).

# CRYOLOCK® S-CRYOLOCK® CRYOLOCK FAMILY DEVICES

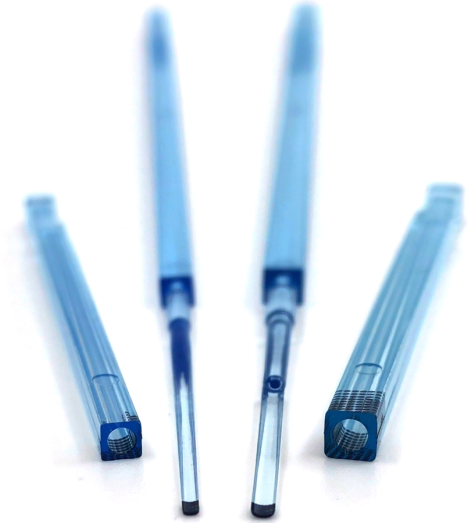

Open System

Vitrification of  
Oocyte and/or  
Embryo

### Intended For Use:

Cryolock® Family Devices are cryopreservation storage devices that are intended for use in vitrification procedures to contain and maintain human 1-Cell stage embryos.

For non US-countries: For Oocytes and/or Embryos.

## INSTRUCTIONS FOR USE

### Warnings


- All procedures must be performed under aseptic laboratory conditions.
- To avoid injuries with LN<sub>2</sub>, wear protective gloves and glasses.
- **Do not use device if:** (a) Pouch or package is open or damaged, (b) Gamma indicator is yellow or missing, or (c) Expiration Date has expired.
- Before loading MII oocytes or embryos, verify integrity of device under microscope view, discarding any device with cracked tips, scratched brittle, with flash, bubbles, presence of foreign material or abnormal shape.
- For better survival rates, use MII Oocytes or good quality embryos. Use only with licensed vitrification media for the embryo stage being vitrified.
- Avoid direct contact of the tip of device at any time; with any surface or material different to vitrification/warming media or pipettes holding the specimens.
- Always use the device with its corresponding cap as it was originally packaged.
- To prevent accidental loss of embryos, perform loading and unloading of 1-cell embryos under microscope view, avoiding contact of the tip against other surfaces. (i.e. edge of petri dishes, or liquid nitrogen containers)
- Load specimens with a maximum of 1 µL of vitrification media, excessive media may cause low survival rates as well as attachment of the tip to the inner cavity of the device cap and possible breakage of tip or cap when warming.
- Immediately and after immersing the device into LN<sub>2</sub>, carefully insert the device Tip into the Cap twisting tightly until secure
- To avoid accidental rushing, or inappropriate time of exposure of specimens to vitrification solutions during loading and plunging into LN<sub>2</sub> perform ONLY 1 or 2 specimens at a time.
- When plunging device into LN<sub>2</sub> always use a separate fresh aliquot LN<sub>2</sub> per patient. Be careful when releasing the device under LN<sub>2</sub>, don't throw devices into LN<sub>2</sub>, place them gently into the corresponding goblets previously equilibrated under LN<sub>2</sub>.
- It is important that the container holding LN<sub>2</sub> be filled no less than 20cm. Not doing so could cause the user to add unnecessary stress to the device and potentially causing the device to break.
- Never leave several devices with oocytes or embryos loaded uncapped under LN<sub>2</sub>.
- Do not re-sterilize or re-use Cryolock® or S-Cryolock® devices. Device properties may change decreasing device performance. Possible contamination, low survival rates, lysis and/or embryo degeneration may occur.
- If device is dirty, discard it, DO NOT clean or wipe device tips with alcohol or equivalents, material properties may change.
- The long-term vitrification safety of oocytes and embryos on children born following this procedure is unknown.

### Precautions


- The correct use of the device is responsibility of the user. For exclusive use of embryologists, biologists or laboratory technicians duly trained on cryopreservation techniques and vitrification protocols.
- For vitrification and warming purposes, have all necessary materials, tools and equipment ready and handy before starting procedures.
- For Laboratory use only. Not for diagnostic use.
- **Storage Instructions:** Store at room temperature
- **Disposable:** After each package containing 5 devices is opened, all devices need to be used or discarded. Cryolock® and S-Cryolock® is for single use only.
- For infectious patients: a) Never reuse liquid nitrogen when plunging specimens from infectious patients. b) Always keep infectious specimens on a separate "dewar" exclusively assigned for "infectious patients", follow laboratory procedures for infectious patients.

## LOADING AND CLOSING


1. Use a liquid nitrogen-resistant label to identify oocytes or embryos of the patients, using the label on the same surface where Cryolock® or S-Cryolock® is engraved.

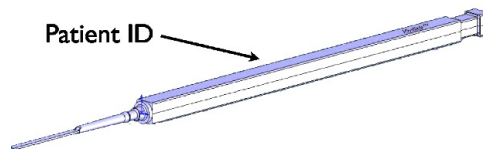


2. Prepare the sample for vitrification according to laboratory vitrification protocol.
3. Using a micropipette, carefully load a maximum of 2 specimens on the concave surface of the tip (same side of Cryolock® or S-Cryolock® logo) and about 3mm (1/8") from the edge of tip (use black mark as a reference) removing any excess of cryoprotectant solution leaving as minimum volume of vitrification media as possible ( $\leq 1 \mu\text{L}$ ).

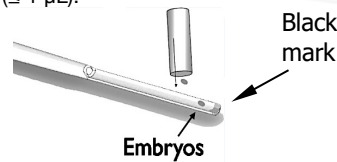

4. Immediately immerse Tip and use the black mark of the Cap as a guidance for capping under LN<sub>2</sub>. Allow equilibration until stop bubbling. Carefully insert the tip into the cap twisting tightly enough until secure. Never throw uncapped devices containing specimens and leave unattended

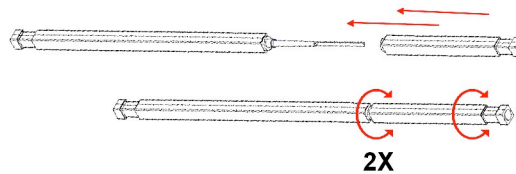



5. Storage specimens on dewars following the laboratory vitrification protocol. Always store the device with the cap facing down.

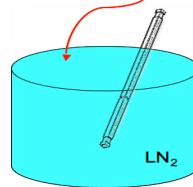


**Note:** Perform steps 3 to 5 in less than 1 minute. After vitrification, device must be kept under liquid nitrogen at all times.

## WARMING

1. Prepare the warming solutions according to laboratory vitrification protocol.
2. Identify the sample to be thawed.
3. Place the warming solution under microscopic view.
4. Using forceps hold the upper extreme of the device body and then quickly remove the cap with a gentle twist pulling the cap straight and away from the device body until release, never bend device body or cap.

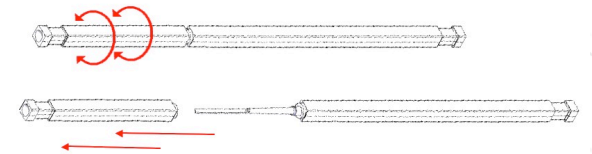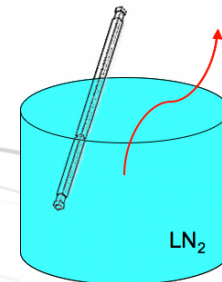

5. Immediately plunge the tip of Cryolock® or S-Cryolock® with specimens facing up into the warming solution at 37°C.

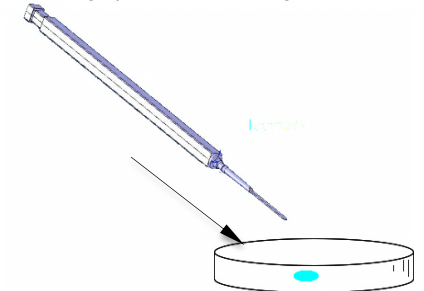

6. Under microscopic observation, gently move the device until embryos are released from the tip.
7. Continue the warming according to laboratory vitrification / warming protocol.
8. Discard device after completion of procedure, in accordance with the applicable medical/hazardous waste regulations.

**Note:** Transition between steps 4 to 5 should be no longer than 5 seconds.