Affordable. Reliable. Home Improvement.

BED HANDLE

MODEL:LQX-110055

VEVOR Support Center

Affordable. Reliable. Home Improvement.

BED HANDLE

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

Warning-To reduce the risk of injury, user must read instructions manual carefully.

SAFETY INSTRUCTIONS

WARNING:

Read this material before using this product. Failure to do so can result in serious injury.

Use precautions

TO PREVENT SERIOUS INJURY AND DEATH FROM TIPPING:

1. DO NOT SIT OR STAND ON THIS ITEM.
2. This product is not a toy. Do not allow children to play with or near this item. Please keep small spare parts out of children's reach.
3. Use as intended only.
4. Inspect before every use; do not use if parts are loose or damaged.
5. Keep away from sharp points, blades and other items, so as not to scratch the product.

Technical Parameters

Model	maximum load capacity
LQX-110055	300lb
FDA Registration No.: 3016734231	

Instructions for Use

APPLICATION OF BED SAFETY RAIL


This Bed Safety Rail is an assistive balance device. It is designed to assist people who are mobile, but require some assistance to get in and get out of bed, or reposition in bed safely. The bed rail is NOT intended for people with very poor mobility, or who are permanently bedridden. It is NOT designed to prevent people from falling out of bed.

CAUTIONS and WARNINGS

1. Please make sure that you read and follow these instructions carefully. Failure to comply may result in injury.
2. Ensure all parts are assembled securely and all spring buttons are properly engaged in position prior to use.
3. The area between an assistive bed rail and mattress creates a potential for serious, life-threatening entrapment. It is possible for a person's head, neck, arm (s), leg (s) or body to become wedged between the side of the mattress and any of the bed rail.
4. The mattress must be tight against the assistive handrail, with no gaps at all times. This should be checked frequently-at least daily.
5. The mattress must be tight against the assistive handrail, with no gaps at all times. This should be checked frequently-at least daily.
6. NEVER use unless the product is tight against the mattress without gaps, and at least 12.5 in. from headboard and foot board.
7. NEVER use with children.
8. NEVER use on toddler, bunk, water, or inflatable beds, beds with soft compressible pads, articulated beds, or on beds with a mattress height taller than 13in (33cm) or shorter than 7 in (17.8cm). ONLY use with twins, queen and king size beds.
9. The product is an auxiliary product. Patients without autonomous ability must be under supervision.

Instructions for Use

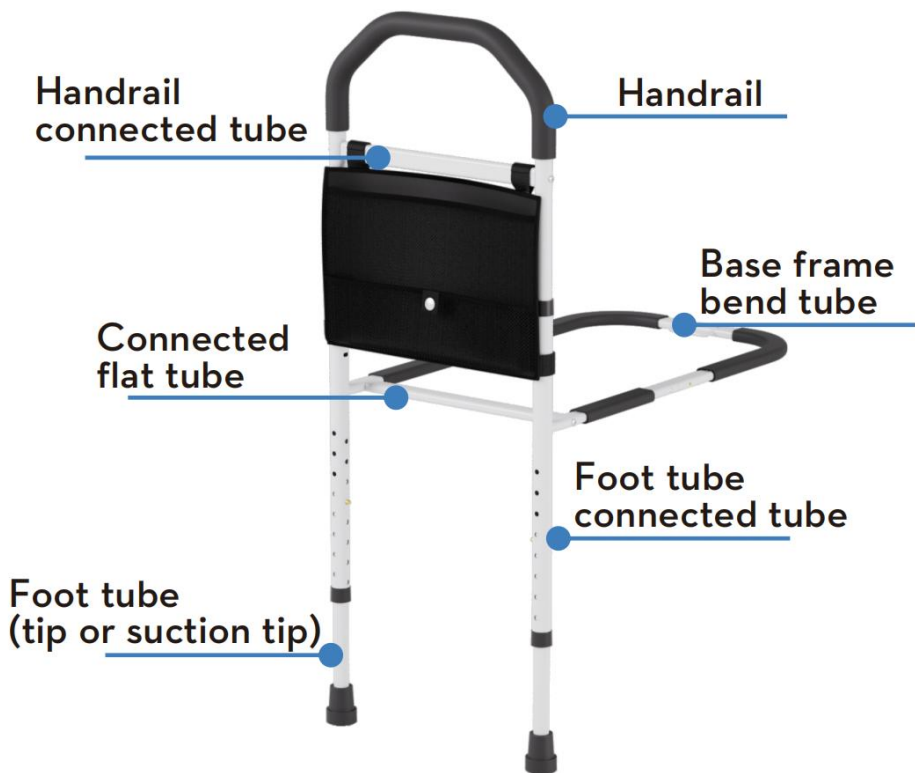[Application] Designed for the elderly and those with leg injuries, providing the necessary support and stability for a safe and sound sleep.

[Structure] It is mainly composed of Handrail tube, handrail lower tube, reinforced tube welding body, base frame bend welding body, base frame welding body.

[Feature] T-style bedside rail is made of high-strength carbon steel tube with white powder-coated surface. The height of the handrail and foot tube can be adjustable. The height of handrail has 4 levels to set. the height of foot tube has 16 levels to set. Very convenient to adjust and use.

[Specification]

1. The appearance is bright and clean, the color is uniform, no obvious scratches, cracks, edges, burrs.

2. Handrail height adjustable range: 390-465mm.

Foot tube height adjustable range: 330-517mm.

[Instruction]

1. Open the package and check whether all the parts are included, as below picture shows:

①.Foot tube(tip or suction tip)

②.Handrail tube

③.Connected flat tube

④.Foot tube connected tube-Left

⑤.Foot tube connected tube-Right

⑥.Base frame bend tube

⑦. Storage bag

⑧.Straps

⑨.Opp bag:

Opp bagx1

Tapping screwsx4


Screwsx2

Hexagonal wrenchx1

2. Take out the foot tube connected tube-Left and Right, connected flat tube, tightening them by using the hexagonal screws, as shown in below picture ①:


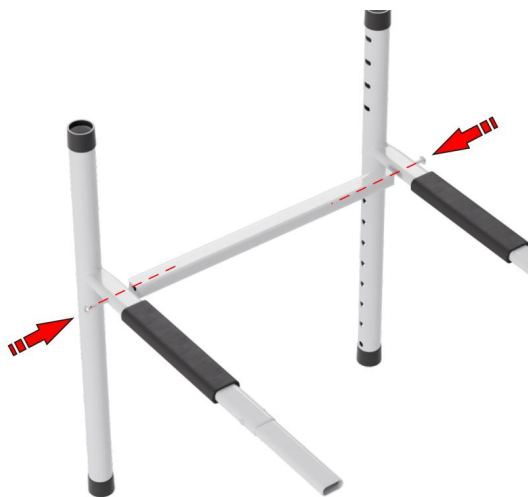
Picture ①

3. After the handrail tube is inserted into the handrail connection tube, the pin pops out, indicating a secure fixation. as shown in below picture ②-1:

4. After the foot tube is inserted into the foot tube connection tube, the pin pops out, indicating a secure fixation, as shown in below picture ②-2:


5. After the Foot tube connected tube is inserted into the foot tube connection tube welding body, the pin pops out, indicating a secure fixation, as shown in below picture:

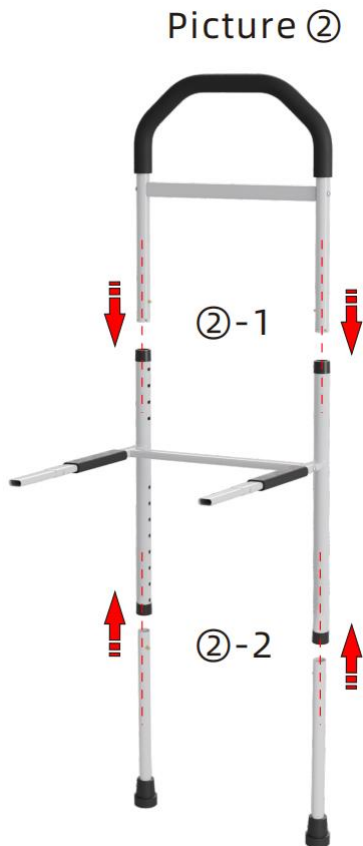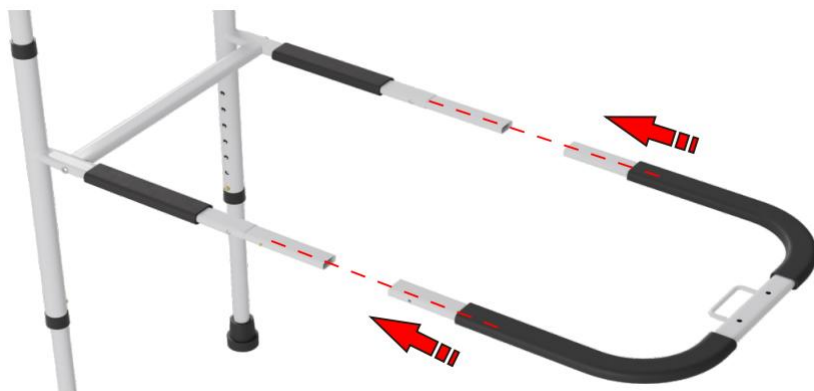6. Install the storage bag.

[Height adjustment] Adjust the position of the pins and adjust the handrail to the desired height, as shown in below picture:

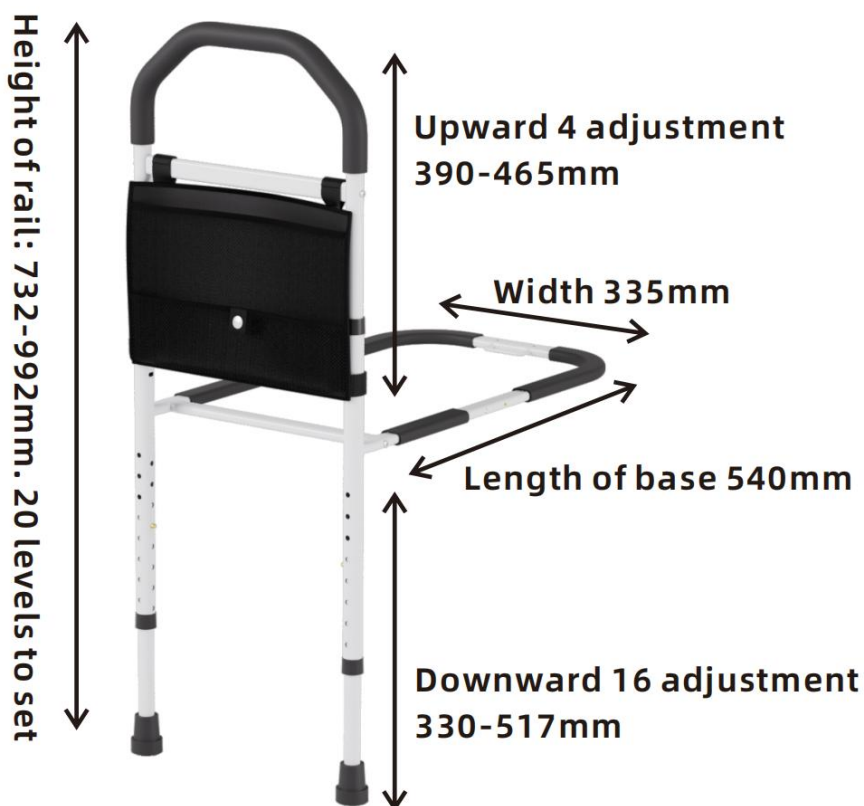[Storage] Forbidden to be pressed horizontally, handle with care.

[Scraped] Disposal in accordance with local laws.

Option -Fix the Bed Rail on Bed Frame (it' s non-mandatory)

The Bed Safety Rail could be fixed on base of the bed with screws provided and provide increased stability for users, but it is not necessary to do.

CLEANING THE BED SAFETY RAIL

Note: Please check the suitability of the cleaning products before use. Do not use abrasive cleaners or cloths as these could permanently damage the surface and invalidate any warranty.

1. Clean regularly using soapy water or warm water. Then wipe it with dry towel.
2. Allow product to fully dry prior to use/storage.

Our After-Sales Service:

1. If you discover any product damage or missing parts, please contact our customer service team for assistance.
2. Our dedicated customer service team is available to address any questions or concerns you may have regarding the product.
3. We strive to provide prompt and efficient support to ensure the satisfaction of our customers.
4. For any inquiries, warranty claims or assistance, please contact our customer service.

Usage Guidelines for Bed Rail:

1. Ensure that the bed rail is properly installed and securely attached to the bed frame before use.
2. When adjusting the height or position of the bed rail, make sure it is securely locked in place before use.
3. Do not allow children to play or hang on the bed rail. It is intended for adult use only.

4. Use caution when getting in or out of bed, ensuring that the bed rail is properly positioned to provide support.

Manufacturer: MOBILITY SOURCE MEDICAL TECHNOLOGY CO., LTD.
Address: 160 HEHENG ROAD, YANGHE TOWN, GAOMING DISTRICT, Foshan, Guangdong, 528513, CHINA.
Imported to AUS: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australia
Imported to USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730

UK	REP	YH CONSULTING LIMITED. C/O YH Consulting Limited Office 147, Centurion House, London Road, Staines-upon-Thames, Surrey, TW18 4AX
EC	REP	E-CrossStu GmbH Mainzer Landstr.69, 60329 Frankfurt am Main.

Affordable. Reliable. Home Improvement.

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Importato in AUS: SIHAO PTY LTD. 1 ROKEVA STREET EASTWOOD NSW
Importato negli USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place,