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PLEASE READ WHOLE INSTRUCTION
CAREFULLY BEFORE USE

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Thank you for choosing our product. We are deeply convinced that it will meet your expectations. Please read this document before first use. This manual is passed in order to provide the necessary information for the safe use of the device. The correct application of the recommendations will contribute significantly to extend the durability and aesthetics of the product. We want to emphasise that constantly we are trying to improve the quality of our products, which may result in minor changes that are not included in the manual

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The shower stool is a light construction made of aluminum. Stool legs with regulated height assure stable support. The seat has a function of rotation which allows the patient to rotate in any direction without stand up.

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STOOL'S EQUIPMENT:

1. Polyurethane seat,
2. Aluminum legs,
3. Telescopic regulation of height,
4. Rubber straps on the top of stool's legs

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For patients with significant movement limitations, generally poor physical status, aging persons; using a bath or shower for hygienic practices. For use by patients who cannot themselves perform actions and hygiene procedures. Stool secure the stable position of the patient, protects against the risk of falling. Protection against wound infection or dipping during hygiene procedures. Stool can be used during hospitalization and home care.

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
Contraindications to a sitting position, the lack of trunk stabilization. Patients with major coordination disorders, balance and intellectual disabilities should not be left without supervision.

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- After unpacking the package, check if it contains all parts in accordance with the list of elements set out above.
- Ensure that all components have no visible damage. If any doubts, please contact the manufacturer.
- Do not use the stool until all parts are properly connected.
- Stool assembly, should start with the connection two frames of stool.
- In order to do that should prepare pad. Insert screw into mounting hole and then tighten it. (Illustration 5; 6.)

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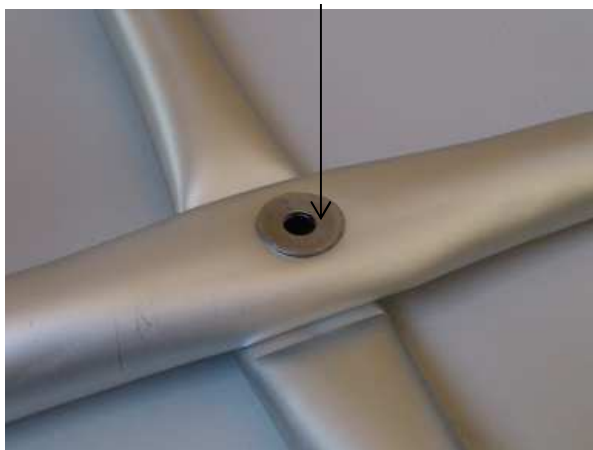- Then should connect stool frame with seat by using four screws. Between the stool frame and the screw should be placed a pad. (Illustration 7.)
- At the ends of the stool frame, should insert the legs, in order to do this, press latches and then set the desired height of legs.
- Make sure all legs are set to the same height and that they are properly fixed.

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- a) Before using the shower stool, refer to this manual and keep it for future reference.
- b) The shower stool is designed to be used on flat, straight surfaces.
- c) The use of shower stool on the rough surface (of the slope more than 5%) may result in falling down and injury.

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d) It is strictly forbidden to stand on the stool.

e) Maximum load of the stool – 100 kg.

f) Use the stool in case of observing any visual or detectable damages (such as frame split) it may result in falling down.

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Keeping all the elements clean is essential for preserving the product working parameters. Keeping the shower stool for a long time in a very moist room or outdoor may result in loss of esthetic and functional parameters of the product.

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
Weight: 2 kg

Frame material: Aluminum,

Serial number/ production date: on the product

Product durability: 5 years from production

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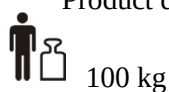Repair

- In case of reparation, please contact the producer.

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Customer service

- If you have any questions or need help, contact the manufacturer.

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