

Public Lab's

Papercraft Spectrometer v2.0.8

Open Source CERN OHL v1.2

This open source hardware design was developed by contributors like you. Have questions? Ideas for improvement? Want to collaborate with others?

Visit: PublicLab.org/paper

This template can be photocopied for reuse! For best results, use cardstock or a similar heavy paper. Make sure to do double-sided copies so the black box and slit cutout will print on the opposite side.

What you'll need:

- scissors
- tape or glue
- a butter knife or ballpoint pen
- a DVD-R

Instructions


more info at PublicLab.org/paper

1. Cut out the spectrometer along the red lines (these may be grey if photocopied in black & white)


2. Use a butter knife or ballpoint pen to score the dotted lines for folding



3. Fold over and cut out the "windows" at A1, A2, & A3

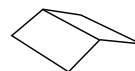

Mountain fold

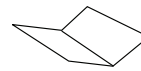

Valley fold

4. "Mountain" fold along the dotted lines: -----

5. "Valley" fold along dashed lines: - - - - -

6. Cut, peel, and attach DVD-R fragment (see A3 for directions)



7. Assemble the slit cutout on the opposite side. Tape the slit over window A2 and trim as needed.

8. Align tabs B, C, D & E and tape or glue them down.

9. Visit PublicLab.org/makespec to learn how to calibrate and use your new spectrometer!


DIY Spectrometry Paper Slit

Instructions: PublicLab.org/makespec
Open Hardware - CERN OHL 1.2

1. cut out bar and U-shape below, keeping lines as straight as possible



2. move upper piece down to form slit; glue or tape in place

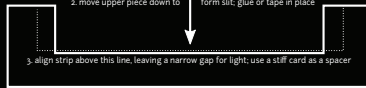

3. align strip above this line, leaving a narrow gap for light; use a stiff card as a spacer

aim for this narrow - 0.4mm - or less; hold up to light to compare