



# PrismView™ Blades

## English & American Profile

### Laryngoscopes

### With Conventional Illumination

#### PRISMVIEW™

The NEW English PrismView™ and American PrismView™ blade may be used:


- when the view of the vocal cords is impaired, especially the most anterior larynx.
- during nasal intubation
- for post operative examination of the larynx.

Both contain an optically polished prism which provides a refraction of 30° without image inversion. To prevent condensation, warm the prism before use by immersing the complete unit in warm water. Blades made of surgical stainless steel.



#### REPLACEMENT LAMPS

5-0234-14  
NEW! Parabolic Reflector Lamp  
Provides 40% More Light!

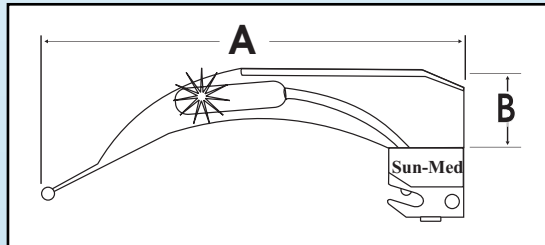

CAT#	SIZE	DESCRIPTION	A	B
------	------	-------------	---	---

#### ENGLISH PRISMVIEW™

5-5028-02	2	Child	100mm	17mm
5-5028-03	3	Medium Adult	130mm	24mm
5-5028-35	3.5	Ext. Med. Adult	144mm	25mm
5-5028-04	4	Large Adult	155mm	26mm

#### AMERICAN PRISMVIEW™

5-5029-02	2	Child	100mm	17mm
5-5029-03	3	Medium Adult	130mm	24mm
5-5029-35	3.5	Ext. Med. Adult	144mm	25mm
5-5029-04	4	Large Adult	155mm	26mm