



Product Catalogue

Aesthetic, Revision and Total Breast Reconstruction™ Surgery



| A Confident Choice for Life™

Contents

GC Aesthetics®	3
A Confident Choice for Life™	3
Pre-Surgery	13
Eve™ 4.0	17
GCA®Consultation Kit™	19
Surgery	35
Implant Insertion Device	
HydroCone™	43
Round Breast Implants	
The Round Collection™ Smooth	45
PERLE™	51
The Round Collection™	55
Anatomical Breast Implants	
Luna XT™	65
Nipple-Areola Implant	
FixNip™ NRI	71
Breast Tissue Expanders	
Blossom™ IVTE	73
Post-Surgery	77
Noa™	81
Healena™	83
GCA Comfort Plus™ Warranty	85
Face & Body Implants	87
Chin Implants	93
Gluteal Implants	95
Testicular Implants	97
Calf Implants	99
Body Tissue Expanders	
Renuvé™ RVTE	101



GC Aesthetics®

A Confident Choice for Life™

GC Aesthetics® is an established global medical device company focused on women and their journey through Aesthetic, Revision and Reconstruction procedures.

We provide the most comprehensive implant portfolio in the industry with breast, nipple-areola, face, body implants and tissue expanders¹ to give surgeons quality surgical options.

Our breast implants are manufactured in accordance with the internationally recognized Quality Management System Standard ISO 13485:2016 and are continuously monitored and verified, meaning our products fulfil all relevant standards required for medical devices in Europe and by local notified bodies worldwide².

A Confident Choice for Life™.

+40

Years of expertise
in Aesthetic and
Reconstructive surgery²

+70

Countries world-wide
with GCA® Products³

MDR

First breast implant
on the market under
the latest regulations
(Luna XT™)⁴

10Y

Prospective and
multicentric long-term
proven safety data⁵

17Y

Extended long-term
high patient
satisfaction data⁶

100%

Designed and
manufactured
in Europe²



The Total Breast Care™ Company

GCA® has collaborated with the medical community to develop new and exciting high-quality products that complement each other from the initial consultation all the way through to patient aftercare.



PRE-SURGERY



SURGERY



POST-SURGERY






Long-term Data Safety and performance

Results Kaplan-Meier 10-year cumulative risk*


~ 1,000 Breast Implants Cristalline™ ParaGel™

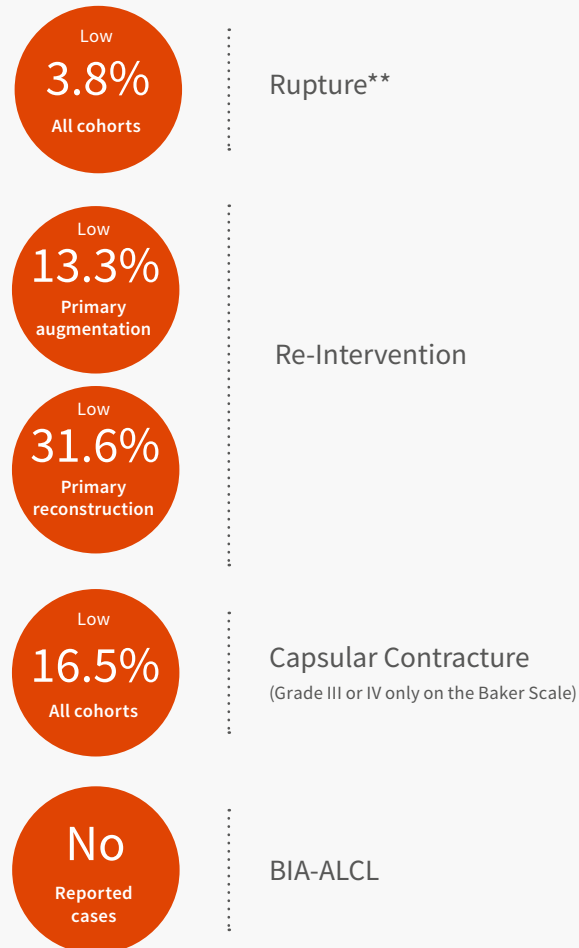

10 Year Prospective & Multicentric Clinical Study

Surgeons and patients across the globe, have trusted in GCA® products. Over 500 patients and almost 1,000 breast implants were monitored for 10 years in a prospective and multicentric clinical study⁵.

The outcome from our 10 year data study confirms that GCA® Products*** **have low rates of rupture, re-intervention and capsular contracture⁵**.

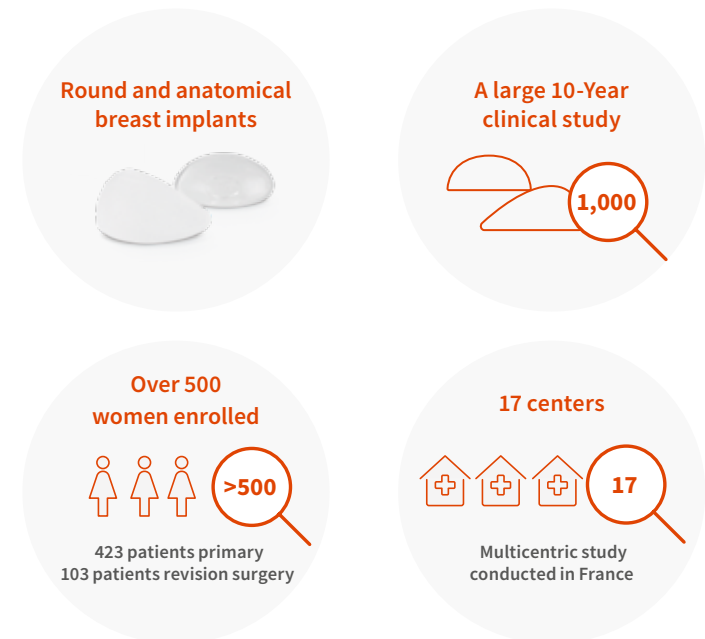

* For more information please consult the GCA® 10 years safety and performance brochure.


** GC Aesthetics®: Implant rupture was diagnosed by surgeon examination only (Not MRI cohort).

*** Study conducted with GCA® breast implants: The Round Collection™ and The Matrix™.



Extended long-term follow-up patient satisfaction study⁶

The longest follow-up study in a single centre

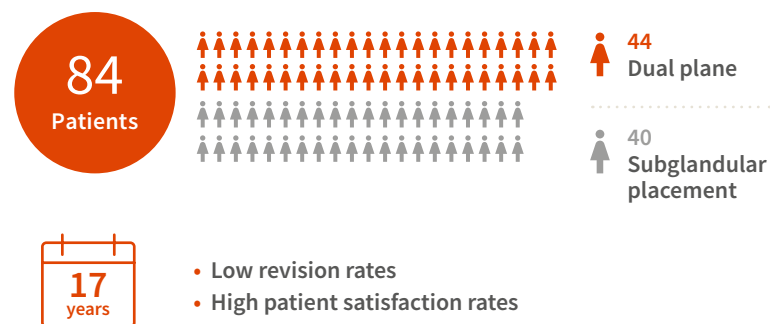

Results:

Quality of Life Results

Physical well-being score of
100/100

Psychosocial well-being score of
74/100

Sexual well-being score of
67/100

Satisfaction Results

Satisfaction with breasts score of
71/100

Satisfaction with outcome score of
81/100

Satisfaction with care score of
100/100



GCA Academy™

The GCA Academy™ is dedicated to providing a progressive learning forum for surgeons, with a focus on advancing scientific knowledge, surgical approaches and clinical management in order to aid surgeons, drive clinical excellence and improve patient outcomes.

Register for free at: gcaacademy.com

- Join our digital 24/7/365 platform to access exclusive content.
- From breast augmentation to breast reconstruction and other procedures.
- Access expert talks, live surgeries, and business boost webinars.

Our long-term learning partnership

Learning Path



Find out more about the GCA Academy™ at: gcaacademy.com



GCA
ACADEMY™



gcaacademy.com

More than
>2,500
Surgeons
registered




PRE-SURGERY

Make your consultation process more inclusive and reassuring by encouraging your patients to make an informed decision about their breast surgery.

Using a digital approach gives an added value to your practice that makes your patients' experience comfortable and unique.

Many benefits are potentially associated when using a pre-surgery tool during the first visit:

- Patients may be more satisfied with the results after the surgery and this can lead to less reoperations.
- Makes the consultation process more inclusive and reassuring, as patients are involved in the decision-making process.
- Enables patients to experience the look of breast augmentation prior to surgery.
- Helps surgeons build the required trust in the initial consultation.
- Easier to meet patients' expectations.
- Improves conversion rates and may reduce time from consultation to surgery.
- Allows surgeons to proceed with confidence.

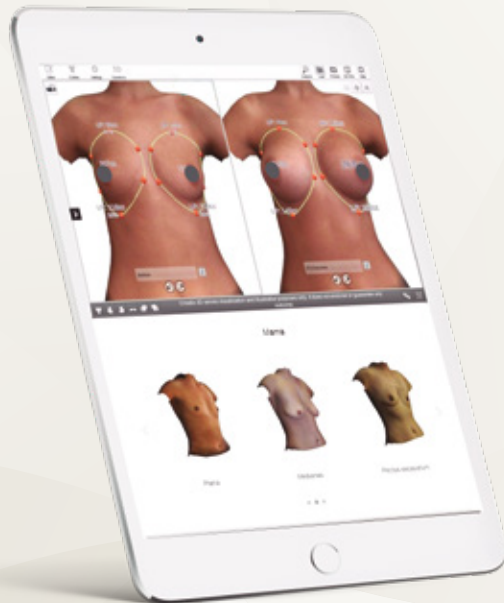

Eve™ 4.0

The New Digital Experience

A full business platform to improve the dialogue between surgeons and patients. The software includes a 3D/Augmented Reality simulator and video consultation integrated in the same platform. With this technology, patients will see the possible achievable results before the surgery.

Main benefits of 3D:

- Standardise & document your before and after pictures.
- Upgrade your existing picture database in 3D.
- Share your results in animated 3D on social media channels.

Eve™ 4.0 Purchase options*

- **Gold:** 3D Simulation limited to 60 patients.
- **Platinum:** 3D Simulation for unlimited patients.
- **Diamond:** 3D and Augmented Reality Simulation for unlimited patients.

eve 4.0
By GC Aesthetics®

*All options include video consultation. Hardware required: iPad Pro.



GCA® Consultation Kit

For pre-surgical planning

This kit focuses on achieving a higher level of patient satisfaction after breast implant surgery, a better and more realistic preview of the desired results, enabling a shared responsibility in the final decision between the surgeon and the patient.

Content GCA® Consultation Kit:

- 1 x Outer case/Inner Case
- 12 x External sizers
- 2 x Stackable sizer
- 2 x Noa bras (S&M sizes)
- 1 x Tool case



EXTERNAL SIZERS

Product Name	Width (cm)	Height (cm)	Projection (cm)	Volume (g/cc)
SIZE S	10.9	14.8	3.4	180
SIZE 0	11.7	16.1	4.3	240
SIZE 1	12.3	16.5	4.5	270
SIZE 2	12.8	17.0	4.8	320
SIZE 3	13.0	17.3	4.9	340
SIZE 4	13.4	17.9	5.3	400
Stackable	9.90	13.2	1.0	110

*These are not medical devices and are intended for an approximate estimation of sizing only. These are not intended to replace a review with the physician on the most appropriate size implant to be used for each individual patient as this can vary due to many factors such as patient anatomy and implant location. They are not intended to replace the use of internal sizers which provide a more accurate representation of the size of the implant.

Size S

Width (cm)	Height (cm)	Projection (cm)	Volume (g/cc)
10.9	14.8	3.4	180



ROUND

Corresponding Breast Implant*

The Round Collection™ Smooth					ES 261
					Medium Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
261180	10.6	3.2	180	6.8	35R180
261200	11.1	3.3	200	7.1	35R200

PERLE™					SOR-MR
					Moderate Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	
SOR-MR 190	9.9	3.6	190	6.8	

The Round Collection™					ES 802N
					Medium Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
802180N	10.6	3.2	180	6.8	35R180
802200N	11.1	3.3	200	7.1	35R200

The Round Collection™					ES 805N
					Medium - Extra High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	
805195N	10.5	3.5	195	6.9	

The Round Collection™					ES 810N
					High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	
810180N	9.5	3.8	180	6.8	

The Round Collection™					ES 801
					Medium Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
801180	10.6	3.2	180	6.8	35R180
801200	11.1	3.3	200	7.1	35R200

The Round Collection™					ES 81
					High Profile 85% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
81200	9.9	3.7	200	6.8	37M200



SMOOTH



SMOOTH
OPAQUE



MICRO
TEXTURED



ANATOMICAL

Corresponding Breast Implant*

Luna XT™ 					MIA-HR1 High Height, Low Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-HR1/185	10.5	11.0	3.4	5.1	185	7.7

Luna XT™ 				MIA-MR1 Moderate Height, Low Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-MR1/195	11.0	10.5	3.6	5.2	195	7.7

Luna XT™ 					MIA-LR1 Low Height, Low Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-LR1/175	11.0	9.5	3.6	5.0	175	7.3

Luna XT™					MIA-MR2		
					Moderate Height, Moderate Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-MR2/175	10.0	9.5	3.9	5.2	175	7.6	GS-XM2 165

*See footnote on page 19.

Size 0

Width (cm)	Height (cm)	Projection (cm)	Volume (g/cc)
11.7	16.1	4.3	240



ROUND

Corresponding Breast Implant*

The Round Collection™ Smooth					ES 28
					High Profile 85% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
28280	10.9	4.3	280	7.7	37M280

PERLE™					SOR-MR
					Moderate Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
SOR-MR 280	11.0	4.2	280	7.7	GS-MR 270
SOR-MR 300	11.4	4.3	300	7.9	GS-MR 300

PERLE™					SOR-HR
					High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
SOR-HR 235	9.7	4.6	235	7.4	GS-HR 240

The Round Collection™					ES 805N
					Medium – Extra High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	
805240N	11.4	3.7	240	7.5	
805335N	12.1	4.3	335	8.2	

The Round Collection™					ES 810N
					High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	
810240N	10.3	4.2	240	7.4	
810265N	10.6	4.4	265	7.7	

The Round Collection™					ES 812N
					Extra High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
812240N	10.2	4.4	240	7.5	37R240

The Round Collection™					ES 81
					High Profile 85% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
81280	10.9	4.3	280	7.7	37M280



SMOOTH



SMOOTH
OPAQUE



MICRO
TEXTURED



ANATOMICAL

Corresponding Breast Implant*

Luna XT™ 					MIA-HR1 High Height, Low Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-HR1/275	12.0	12.5	4.0	5.9	275	8.9

Luna XT™ 				MIA-MR1 Moderate Height, Low Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-MR1/245	12.0	11.5	4.0	5.7	245	8.5

Luna XT™ 					MIA-LR1 Low Height, Low Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-LR1/230	12.0	10.5	4.0	5.5	230	8.1

Luna XT™					MIA-HR2		
					High Height, Moderate Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-HR2/245	11.0	11.5	4.3	6.0	245	8.8	GS-XF2 230

Luna XT™					MIA-MR2		
					Medium Height, Moderate Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-MR2/230	11.0	10.5	4.3	5.8	230	8.4	GS-XM2 215

Luna XT™ 					MIA-LR2 Low Height, Moderate Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-LR2/245	11.5	10.0	4.5	5.9	245	8.5

Size 1

Width (cm)	Height (cm)	Projection (cm)	Volume (g/cc)
12.3	16.5	4.5	270



ROUND

Corresponding Breast Implant*

The Round Collection™ Smooth					ES 28
					High Profile 85% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
28280	10.9	4.3	280	7.7	37M280
28325	11.4	4.5	325	8.0	37M325

PERLE™					SOR-MR
					Moderate Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
SOR-MR 300	11.4	4.3	300	7.9	GS-MR 300
SOR-MR 340	11.9	4.4	340	8.1	GS-MR 330

PERLE™					SOR-HR
					High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
SOR-HR 270	10.0	4.9	270	7.8	GS-HR 270

The Round Collection™					ES 805N
					Medium - Extra High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	
805310N	12.0	4.1	310	8.0	
805360N	12.5	4.4	360	8.4	

The Round Collection™					ES 81
					High Profile 85% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
81325	11.4	4.5	325	8.0	37M325



SMOOTH



SMOOTH
OPAQUE



MICRO
TEXTURED



ANATOMICAL

Corresponding Breast Implant*

Luna XT™					MIA-HR1	
					High Height, Low Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-HR1/315	12.5	13.0	4.2	6.2	315	9.3

Luna XT™					MIA-MR1	
					Moderate Height, Low Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-MR1/280	12.5	12.0	4.2	6.0	280	8.9

Luna XT™					MIA-LR1	
					Low Height, Low Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-LR1/265	12.5	11.0	4.2	5.8	265	8.5

Luna XT™					MIA-MR2		
					Moderate Height, Moderate Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-MR2/290	12.0	11.5	4.7	6.3	290	9.2	GS-XM2 290

Luna XT™					MIA-LR2	
					Low Height, Moderate Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-LR2/280	12.0	10.5	4.7	6.2	280	8.9

Size 2

Width (cm)	Height (cm)	Projection (cm)	Volume (g/cc)
12.8	17.0	4.8	320



ROUND

Corresponding Breast Implant*

The Round Collection™ Smooth					ES 28
					High Profile 85% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
28325	11.4	4.5	325	8.0	37M325
28350	11.7	4.7	350	8.3	37M350
28375	11.7	4.9	375	8.5	37M375

PERLE™					SOR-MR
					Moderate Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
SOR-MR 390	12.6	4.5	390	8.5	GS-MR 390

PERLE™					SOR-HR
					High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
SOR-HR 335	11.1	5.0	335	8.3	GS-HR 330
SOR-HR 360	11.4	5.1	360	8.5	GS-HR 360

The Round Collection™					ES 805N
					Medium - Extra High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	
805335N	12.1	4.3	335	8.2	
805390N	12.8	4.6	390	8.7	

The Round Collection™					ES 810N
					High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	
810320N	11.2	4.8	320	8.2	

The Round Collection™					ES 812N
					Extra High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
812325N	11.0	5.1	325	8.3	37R325

The Round Collection™					ES 801
					Medium Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
801325	12.7	4.0	325	8.2	35R325



ROUND

Corresponding Breast Implant*

The Round Collection™					ES 81
					High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
81350	11.7	4.7	350	8.3	37M350
81375	11.7	4.9	375	8.5	37M375



ANATOMICAL

Corresponding Breast Implant

Luna XT™ 				MIA-MR1 Moderate Height, Low Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-MR1/320	13.0	12.5	4.4	6.2	320	9.3

Luna XT™ 					MIA-LR1 Low Height, Low Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-LR1/345	13.5	12.0	4.6	6.3	345	9.3

Luna XT™					MIA-HR2		
					High Height, Moderate Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-HR2/345	12.5	13.0	4.9	6.8	345	10.0	GS-XF2 350

Luna XT™					MIA-MR2		
					Moderate Height, Moderate Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-MR2/330	12.5	12.0	4.9	6.6	330	9.6	GS-XM2 320

Luna XT™ 				MIA-LR2 Low Height, Moderate Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-LR2/355	13.0	11.5	5.1	6.7	355	9.7

Size 3

Width (cm)	Height (cm)	Projection (cm)	Volume (g/cc)
13.0	17.3	4.9	340



ROUND

Corresponding Breast Implant*

The Round Collection™ Smooth					ES 28
					High Profile 85% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
28350	11.7	4.7	350	8.3	37M350
28375	11.7	4.9	375	8.5	37M375
28400	12.2	4.9	400	8.7	37M400

PERLE™					SOR-MR
					Moderate Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
SOR-MR 430	12.9	4.7	430	8.8	GS-HR 420

PERLE™					SOR-HR
					High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
SOR-HR 360	11.4	5.1	360	8.5	GS-HR 360

The Round Collection™					ES 805N
					Medium - Extra High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	
805360N	12.5	4.4	360	8.4	
805420N	13.0	4.7	420	8.9	

The Round Collection™					ES 810N
					High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	
810350N	11.5	5.0	350	8.5	

The Round Collection™					ES 812N
					Extra High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
812350N	11.3	5.1	350	8.4	35R350
812375N	11.6	5.2	375	8.6	35R375

The Round Collection™					ES 81
					High Profile 85% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
81350	11.7	4.7	350	8.3	37M350
81375	11.7	4.9	375	8.5	37M375
81400	12.2	4.9	400	8.7	37M400



SMOOTH



SMOOTH
OPAQUE



MICRO
TEXTURED



ANATOMICAL

Corresponding Breast Implant*

Luna XT™					MIA-HR1	
					High Height, Low Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-HR1/400	13.5	14.0	4.6	6.7	400	10.1

Luna XT™					MIA-MR1	
					Moderate Height, Low Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-MR1/360	13.5	13.0	4.6	6.5	360	9.7

Luna XT™					MIA-LR1	
					Low Height, Low Projection	
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-LR1/345	13.5	12.0	4.6	6.3	345	9.3
MIA-LR1/385	14.0	12.5	4.8	6.6	385	9.7

Luna XT™					MIA-HR2		
					High Height, Moderate Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-HR2/345	12.5	13.0	4.9	6.8	345	10.0	GS-XF2 350
MIA-HR2/390	13.0	13.5	5.1	7.1	390	10.4	GS-XF2 400

Luna XT™					MIA-MR2		
					Moderate Height, Moderate Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-MR2/375	13.0	12.5	5.1	6.9	375	10.0	GS-XM2 370

Size 4

Width (cm)	Height (cm)	Projection (cm)	Volume (g/cc)
13.4	17.9	5.3	400



ROUND

Corresponding Breast Implant*

The Round Collection™ Smooth					ES 28
					High Profile 85% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
28450	12.5	5.3	450	9.1	37M450

PERLE™					SOR-HR
					High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
SOR-HR 425	12.1	5.3	425	8.9	GS-HR 420
SOR-HR 460	12.4	5.4	460	9.2	GS-HR 460

The Round Collection™					ES 805N
					Medium - Extra High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	
805480N	13.9	4.9	480	9.4	

The Round Collection™					ES 810N
					High Profile
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	
810390N	12.0	5.2	390	8.8	
810435N	12.5	5.4	435	9.2	

The Round Collection™					ES 801
					Medium Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
801400	13.7	4.2	400	8.8	35R400

The Round Collection™					ES 81
					High Profile 85% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
81450	12.5	5.3	450	9.1	37R450

The Round Collection™					ES 811
					Extra High Profile 100% Fill
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
81400	11.8	5.5	400	9.0	37R400



SMOOTH



SMOOTH
OPAQUE



MICRO
TEXTURED



ANATOMICAL

Corresponding Breast Implant*

Luna XT™					MIA-HR2		
					High Height, Moderate Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-HR2/440	13.5	14.0	5.3	7.3	440	10.8	GS-XF2 450

Luna XT™					MIA-MR2		
					Moderate Height, Moderate Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-MR2/415	13.5	13.0	5.3	7.1	415	10.4	GS-XM2 415

Luna XT™					MIA-LR2		
					Low Height, Moderate Projection		
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	
MIA-LR2/405	13.5	12.0	5.3	7.0	405	10.1	



How does the external stackable sizer work?

The external stackable sizer is used in combination with the external sizer in order to simulate a bigger projection. Usually the stackable is placed behind the sizer however, some surgeons put the stackable sizer on top and have good results as well.

External Stackable Sizer			
Width (cm)	Height (cm)	Projection (cm)	Volume (g/cc)
9.90	13.2	1.0	110
EXTERNAL STACKABLE SIZER	EXTERNAL SIZER	COMBINATION	

Example:

SIZE	Corresponding Implant*
SIZE S	SOR-HR 270
SIZE S + External Stackable sizer	SOR-EHR 290

PERLE™				SOR-HR/SOR-EHR	
				High to Extra High Profile	
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
SOR-HR 270	10.0	4.9	270	7.8	GS-HR 270
SOR-EHR 290	10.1	5.1	290	8.1	GS-EHR 300

*These are not medical devices and are intended for an approximate estimation of sizing only. These are not intended to replace a review with the physician on the most appropriate size implant to be used for each individual patient as this can vary due to many factors such as patient anatomy and implant location. They are not intended to replace the use of internal sizers which provide a more accurate representation of the size of the implant.

SURGERY

Aesthetic, Revision and
Total Breast Reconstruction™ Surgery





SURGERY

The **Quality** you Expect.
The **Options** you Need.
The **Peace of Mind** you Demand.

Backed by years of continuous research and development, GC Aesthetics® offers a comprehensive range of silicone implants and innovative solutions for reconstructive and aesthetic surgery.


GCA® has collaborated with the medical community to deliver safe options that establish a newer, higher standard of care for women who seek to recover their breast and nipple shape, look and feel.

GCA® has a broad portfolio:

- Implant Insertion Device
- Round and Anatomical Breast Implants
- Smooth, Smooth Opaque and Micro-Textured Breast Implants
- Breast Sizers
- Nipple-Areola Complex Implant
- Fat Grafting Closed Device
- Tissue Expanders
- Face and Body Implants


A complete range of Breast Implants

By Shape and Surface Roughness

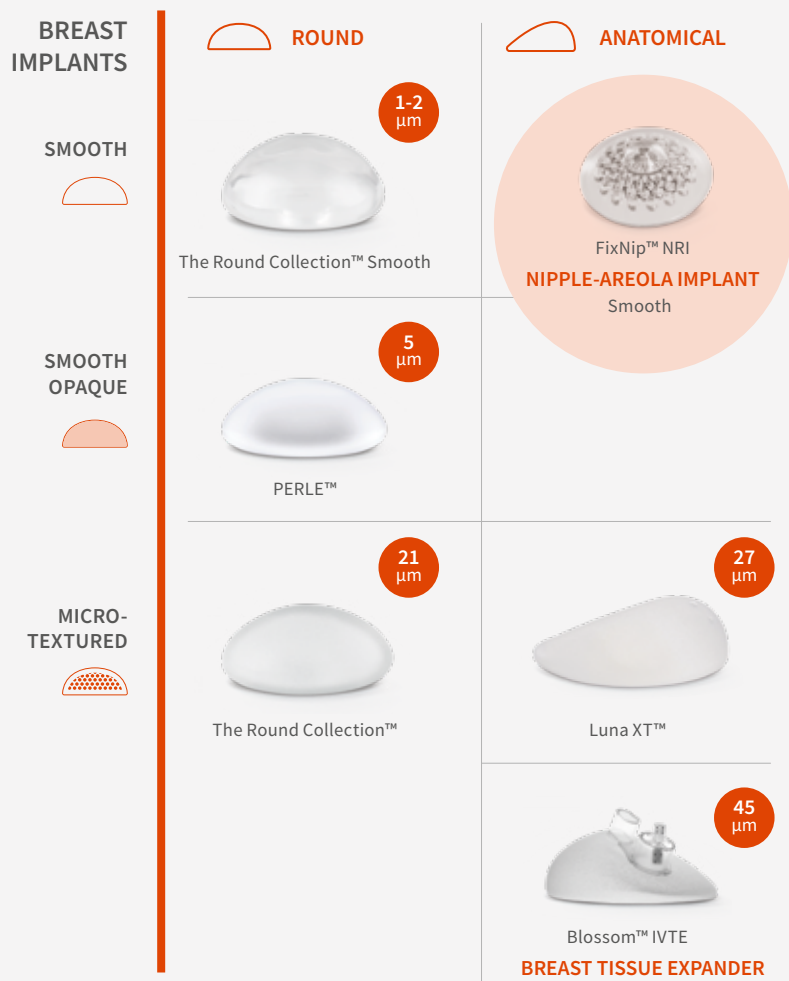

Average Surface roughness of GCA® Round and Anatomical Implants: ISO 14607:2018 (S_a).

By Firmness and Cohesivity

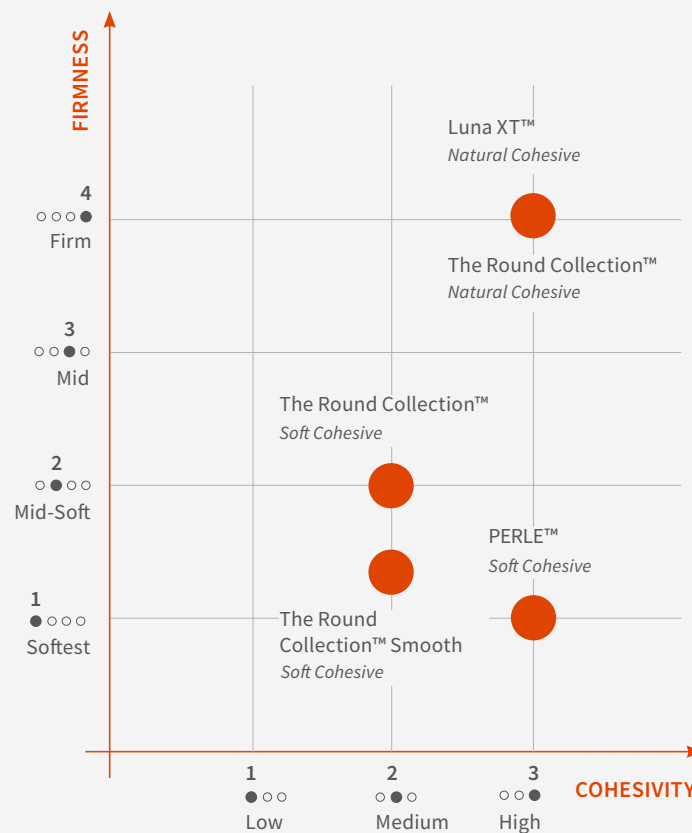

Please note that this scale is subjective, based on the feel of the product, and is not based on actual test results.

Global Product Certifications

Demonstrating global compliance

	MFDS	TFDA	NMPA
APAC	South Korea	Taiwan	China

We have passed the rigorous certifications processes underpinned by proven safety and performance.

EUROPE	CE
--------	----

For over 40 years we have maintained continuous CE marking in Europe. We became the first company to place MDR breast implants on the market under the latest EU Medical Device Regulation (MDR).

	Inmetro	Cofepris	ANMAT	INVIMA	Dlgimed
LATAM	Brazil	Mexico	Argentina	Colombia	Peru

We have product certifications issued to us by Health authorities across the region.

	FDA
US	United States

GCA® supplies products which have passed through the US FDA 510(k) clearance process

GC Aesthetics®

Certifications

With GC Aesthetics® as your partner you can have confidence, we utilize high purity medical silicones and other qualified materials throughout our validated manufacturing processes. All our systems from design and manufacture through to distribution are maintained to the highest global standards by our rigorous certified quality management systems. The pre-clinical and clinical safety and performance data is audited to meet state-of-the-art standards and global regulatory certification requirements.

Quality Systems & Good Manufacturing Practices control our processes

ISO 13485	GMP	ANVISA	KGMP	QSD
International	Brazil	South Korea	Taiwan	



HydroCone™

Surgical support

HydroCone™, the GCA® implant insertion device that will help you to achieve a better outcome, by optimising your time during surgery and reducing risk of complications^{7,8,9,10}.

Product features:

- Sterile, single patient use instrument for the delivery of gel implants made from a clear polymeric material that is soft and flexible¹¹.
- Designed to make inserting all breast gel implants easier¹⁰.
- Inner surface has a coating made of hyaluronic acid, extremely slippery when hydrated with sterile saline solution¹¹.
- Packaged in a sterile double pouch system¹¹.

Key Benefits on using HydroCone™ in surgery:

Safety: Employs minimal-touch delivery technique which may reduce post-operative complications^{7,8,9}.

Convenience: Helps to distribute the force applied to the implant. Ease of insertion which may reduce surgery time or make the procedure more efficient¹⁰.

Broad usage: Designed to work with silicone implants and with different surgical approaches¹¹.



Product Code: REF_01KR





The Round Collection™ Smooth

The Round Collection™ Smooth Breast Implants are engineered and designed with performance in mind. This collection offers a wide variety of sizes and profiles and is designed to meet the needs of surgeons and patients.



Main features:

- Round Smooth Silicone Breast Implant².
- Filled with silicone gel, ParaGel™, which is cohesive and soft².
- Shell is made of multiple cross-linked layers².
- The internal layers of the shell are formulated with ParaGelGard™ 360° barrier layer to minimise gel bleeding¹.
- Wide variety of sizes and profiles².


 The Round Collection™ Smooth
By GC Aesthetics®

The Round Collection™ Smooth

Round Smooth | Soft Cohesive Gel

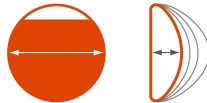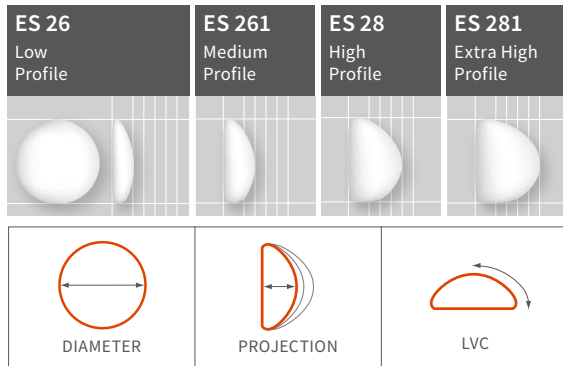

ES 26 Low Profile 85% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
26100	9.6	2.0	100	5.6	35M100
26120	9.7	2.4	120	5.9	35M120
26140	10.2	2.5	140	6.1	35M140
26160	10.4	2.8	160	6.4	35M160
26180	10.9	2.9	180	6.7	35M180
26200	11.6	2.8	200	7.0	35M200
26220	11.8	3.0	220	7.2	35M220
26240	12.0	3.1	240	7.3	35M240
26260	12.2	3.3	260	7.5	35M260
26280	12.5	3.3	280	7.7	35M280
26300	12.9	3.4	300	7.9	35M300
26325	13.2	3.5	325	8.1	35M325
26350	13.5	3.6	350	8.3	35M350
26375	13.8	3.7	375	8.5	35M375
26400	14.2	3.7	400	8.7	35M400
26450	14.7	3.8	450	9.0	35M450
26500	15.3	4.0	500	9.4	

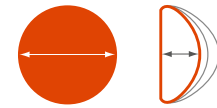


ES 261 Medium Profile 100% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
261100	8.5	2.8	100	5.6	35R100
261120	9.0	2.9	120	5.9	35R120
261140	9.3	3.2	140	6.2	35R140
261160	10.0	3.1	160	6.4	35R160
261180	10.6	3.2	180	6.8	35R180
261200	11.1	3.3	200	7.1	35R200
261220	11.5	3.2	220	7.2	35R220
261240	11.7	3.5	240	7.5	35R240
261260	11.9	3.6	260	7.6	35R260
261280	12.3	3.7	280	7.9	35R280
261300	12.4	3.9	300	8.0	35R300
261325	12.7	4.0	325	8.2	35R325
261350	13.0	4.0	350	8.4	35R350
261375	13.3	4.2	375	8.6	35R375
261400	13.7	4.2	400	8.8	35R400
261450	14.1	4.4	450	9.1	35R450
261500	14.5	4.6	500	9.4	
261550	15.0	4.8	550	9.8	
261600	15.8	4.5	600	9.9	
261650	15.8	4.9	650	10.2	
261700	15.8	5.4	700	10.5	
261750	15.9	5.8	750	10.9	
261800	15.8	6.0	800	11.0	

Depending on the product specifications, 1g is approximately equivalent to 1cc.

The Round Collection™ Smooth

Round Smooth | Soft Cohesive Gel

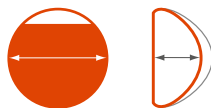

ES 28 High Profile 85% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
28100	7.9	2.9	100	5.4	37M100
28120	8.9	2.8	120	5.8	37M120
28140	8.8	3.2	140	6.0	37M140
28160	9.3	3.3	160	6.3	37M160
28180	9.5	3.6	180	6.6	37M180
28200	9.9	3.7	200	6.8	37M200
28220	10.3	3.7	220	7.0	37M220
28240	10.5	3.9	240	7.2	37M240
28260	10.7	4.1	260	7.5	37M260
28280	10.9	4.3	280	7.7	37M280
28300	11.1	4.4	300	7.8	37M300
28325	11.4	4.5	325	8.0	37M325
28350	11.7	4.7	350	8.3	37M350
28375	11.7	4.9	375	8.5	37M375
28400	12.2	4.9	400	8.7	37M400
28450	12.5	5.3	450	9.1	37M450
28500	12.6	5.8	500	9.5	

ES 281 Extra High Profile 100% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
281100	7.3	3.4	100	5.5	37R100
281120	8.0	3.5	120	5.9	37R120
281140	8.5	3.8	140	6.3	37R140
281160	8.8	3.9	160	6.5	37R160
281180	9.1	4.3	180	7.0	37R180
281200	9.4	4.3	200	7.1	37R200
281220	9.7	4.4	220	7.3	37R220
281240	10.2	4.4	240	7.5	37R240
281260	10.4	4.5	260	7.6	37R260
281280	10.6	4.7	280	7.9	37R280
281300	10.8	4.9	300	8.1	37R300
281325	11.0	5.1	325	8.3	37R325
281350	11.3	5.1	350	8.4	37R350
281375	11.6	5.2	375	8.5	37R375
281400	11.8	5.5	400	9.0	37R400
281450	12.1	5.7	450	9.2	37R450
281500	12.3	6.1	500	9.6	
281550	12.5	6.4	550	9.9	
281600	13.1	6.4	600	10.2	
281650	13.6	6.4	650	10.4	
281700	13.4	7.1	700	10.8	
281750	13.5	7.4	750	11.1	
281800	13.6	7.7	800	11.4	

Depending on the product specifications, 1g is approximately equivalent to 1cc.
The dimensions have a tolerance +/- 4mm on diameter and +/- 3mm for projection.



PERLE™

PERLE™, by GC Aesthetics®, is the novel round breast implant with (R)evolutionary BioQ™-surface technology and highly cohesive Emunomic™ Breast Tissue Dynamic Gel.


PERLE™ is the evolution of smooth opaque breast implants to give patients and surgeons a new alternative for aesthetic and reconstruction procedures.

Main features:

- Round Smooth Opaque Silicone Breast Implant².
- (R)evolutionary BioQ™ -Surface designed with reverse technology for optimal harmonisation with the tissue².
- Dinamic Emunomic™ Gel: 6th Generation Gel, highly cohesive it provides projection, fullness and firmness².
- The implant is 100% gel-filled².
- Enhanced GCA® Ultra-link shell: Optimum cross-linking in our proprietary formula for a resistant yet soft shell².
- Advanced RRE™ Patch: Resistant, Reinforced and Elastic patch/shell interface properties².
- Versatile range of sizes and profiles².



Round Smooth Opaque | Soft Cohesive Gel

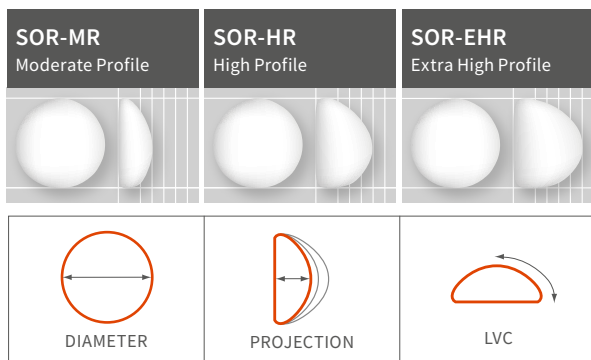

SOR-MR Moderate Profile 100% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
SOR-MR 165	9.2	3.4	165	6.5	
SOR-MR 190	9.9	3.6	190	6.8	
SOR-MR 220	10.3	3.8	220	7.1	GS-MR 210
SOR-MR 255	10.8	4.0	255	7.4	GS-MR 240
SOR-MR 280	11.0	4.2	280	7.7	GS-MR 270
SOR-MR 300	11.4	4.3	300	7.9	GS-MR 300
SOR-MR 340	11.9	4.4	340	8.1	GS-MR 330
SOR-MR 365	12.2	4.5	364	8.3	GS-MR 360
SOR-MR 390	12.6	4.5	390	8.5	GS-MR 390
SOR-MR 430	12.9	4.7	430	8.8	GS-MR 420
SOR-MR 480	13.4	4.9	480	9.1	GS-MR 460

SOR-HR High Profile 100% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
SOR-HR 150	8.3	3.9	150	6.5	
SOR-HR 175	8.6	4.2	175	6.8	
SOR-HR 210	9.3	4.3	210	7.1	GS-HR 210
SOR-HR 235	9.7	4.6	235	7.4	GS-HR 240
SOR-HR 270	10.0	4.9	270	7.8	GS-HR 270
SOR-HR 300	10.6	4.9	300	8.1	GS-HR 300
SOR-HR 335	11.1	5.0	335	8.3	GS-HR 330
SOR-HR 360	11.4	5.1	360	8.5	GS-HR 360
SOR-HR 380	11.6	5.2	380	8.6	GS-HR 390
SOR-HR 425	12.1	5.3	425	8.9	GS-HR 420
SOR-HR 460	12.4	5.4	460	9.2	GS-HR 460
SOR-HR 505	12.6	5.7	505	9.4	
SOR-HR 550	13.4	5.7	550	9.7	

SOR-EHR Extra High Profile 100% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
SOR-EHR 160	8.1	4.3	160	6.7	
SOR-EHR 180	8.4	4.7	180	7.0	
SOR-EHR 210	9.0	4.8	210	7.3	GS-EHR 210
SOR-EHR 230	9.2	5.0	230	7.5	GS-EHR 240
SOR-EHR 260	9.8	5.0	260	7.7	GS-EHR 270
SOR-EHR 290	10.1	5.1	290	8.1	GS-EHR 300
SOR-EHR 330	10.7	5.2	330	8.4	GS-EHR 330
SOR-EHR 360	11.0	5.3	360	8.7	GS-EHR 360
SOR-EHR 400	11.4	5.5	400	8.9	GS-EHR 390
SOR-EHR 435	11.8	5.7	435	9.1	GS-EHR 420
SOR-EHR 460	12.0	5.8	460	9.3	GS-EHR 460
SOR-EHR 520	12.5	6.1	520	9.7	
SOR-EHR 575	12.8	6.4	575	10.2	
SOR-EHR 620	13.1	6.6	620	10.4	

Depending on the product specifications, 1g is approximately equivalent to 1cc.



The Round Collection™

The Round Collection™ breast implants unique micro-textured surface and gel², are clinically proven to be safe and effective and have been shown to have long-term patient satisfaction⁶ and low rates of capsular contracture, rupture and reoperation⁵.



Main features:

- Round Cristalline™ Micro-Texture Silicone Breast Implant².
- High performance shell: Maintains consistent strength and flexibility².
- ParaGel™ Gel: Provides controlled gel distribution and high gel cohesivity. It is available with two types of gels classified as soft or natural².
- ParaGelGard™ 360° barrier layer is designed to minimise silicone leakage while also enhancing shell integrity².
- Wide variety of sizes and profiles².



The Round Collection™

Round Micro-Textured | Natural Cohesive Gel

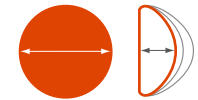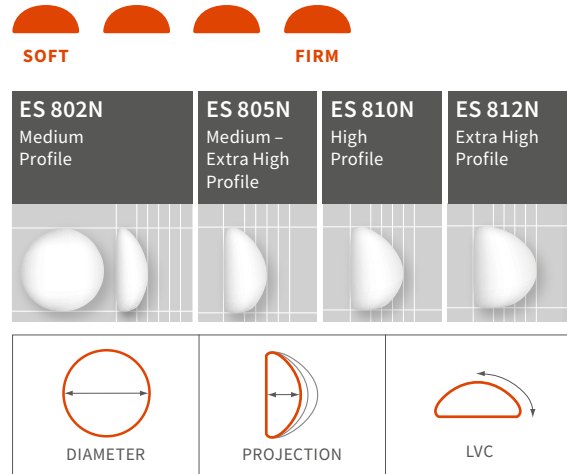

ES 802N Medium Profile 100% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
802100N	8.5	2.8	100	5.6	35R100
802120N	9.0	2.9	120	5.9	35R120
802140N	9.3	3.2	140	6.2	35R140
802160N	10.0	3.1	160	6.4	35R160
802180N	10.6	3.2	180	6.8	35R180
802200N	11.1	3.3	200	7.1	35R200
802220N	11.5	3.2	220	7.2	35R220
802240N	11.7	3.5	240	7.5	35R240
802260N	11.9	3.6	260	7.6	35R260
802280N	12.3	3.7	280	7.9	35R280
802300N	12.4	3.9	300	8.0	35R300
802325N	12.7	4.0	325	8.2	35R325
802350N	13.0	4.0	350	8.4	35R350
802375N	13.3	4.2	375	8.6	35R375
802400N	13.7	4.2	400	8.8	35R400
802450N	14.1	4.4	450	9.1	35R450
802500N	14.5	4.6	500	9.4	
802550N	15.0	4.8	550	9.8	
802600N	15.8	4.5	600	9.9	
802650N	15.8	4.9	650	10.2	
802700N	15.8	5.4	700	10.5	
802750N	15.9	5.8	750	10.9	
802800N	15.8	6.0	800	11.0	




ES 805N Medium - Extra High Profile 100% Fill				
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)
805170N	9.9	3.4	170	6.6
805195N	10.5	3.5	195	6.9
805215N	10.6	3.7	215	7.1
805240N	11.4	3.7	240	7.5
805265N	11.4	3.9	265	7.6
805290N	11.8	4.0	290	7.8
805310N	12.0	4.1	310	8.0
805335N	12.1	4.3	335	8.2
805360N	12.5	4.4	360	8.4
805390N	12.8	4.6	390	8.7
805420N	13.0	4.7	420	8.9
805450N	13.2	4.8	450	9.0
805480N	13.9	4.9	480	9.4


Depending on the product specifications, 1g is approximately equivalent to 1cc.
The dimensions have a tolerance +/- 4mm on diameter and +/- 3mm for projection.

The Round Collection™

Round Micro-Textured | Natural Cohesive Gel

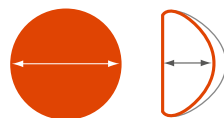

ES 810N High Profile 100% Fill				
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)
810155N	9.0	3.6	155	6.4
810180N	9.5	3.8	180	6.8
810215N	10.0	4.0	215	7.1
810240N	10.3	4.2	240	7.4
810265N	10.6	4.4	265	7.7
810295N	10.9	4.6	295	7.9
810320N	11.2	4.8	320	8.2
810350N	11.5	5.0	350	8.5
810390N	12.0	5.2	390	8.8
810435N	12.5	5.4	435	9.2
810500N	13.0	5.8	500	9.7
810575N	13.5	6.2	575	10.2
810665N	14.0	6.6	665	10.7
810765N	14.5	7.0	765	11.2


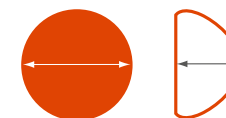


ES 812N Extra High Profile 100% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
812100N	7.3	3.4	100	5.5	37R100
812120N	8.0	3.5	120	5.9	37R120
812140N	8.5	3.8	140	6.3	37R140
812160N	8.8	3.9	160	6.5	37R160
812180N	9.1	4.3	180	7.0	37R180
812200N	9.4	4.3	200	7.1	37R200
812220N	9.7	4.4	220	7.3	37R220
812240N	10.2	4.4	240	7.5	37R240
812260N	10.4	4.5	260	7.6	37R260
812280N	10.6	4.7	280	7.9	37R280
812300N	10.8	4.9	300	8.1	37R300
812325N	11.0	5.1	325	8.3	37R325
812350N	11.3	5.1	350	8.4	37R350
812375N	11.6	5.2	375	8.6	37R375
812400N	11.8	5.5	400	9.0	37R400
812450N	12.1	5.7	450	9.2	37R450
812500N	12.3	6.1	500	9.6	
812550N	12.5	6.4	550	9.9	
812600N	13.1	6.4	600	10.2	
812650N	13.6	6.4	650	10.4	
812700N	13.4	7.1	700	10.8	
812750N	13.5	7.4	750	11.1	
812800N	13.6	7.7	800	11.4	


Depending on the product specifications, 1g is approximately equivalent to 1cc.
The dimensions have a tolerance +/- 4mm on diameter and +/- 3mm for projection.

The Round Collection™

Round Micro-Textured | Soft Cohesive Gel

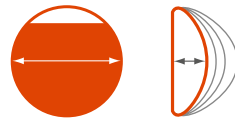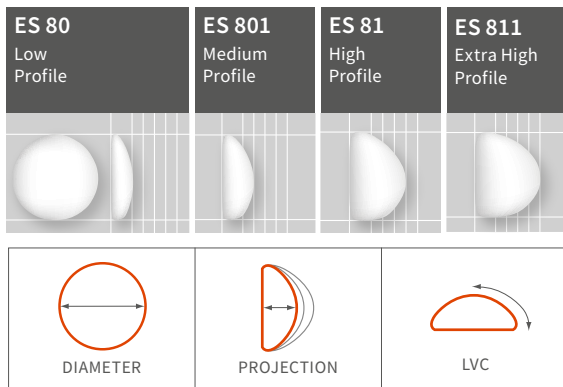

ES 80 Low Profile 85% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
80100	9.6	2.0	100	5.6	35M100
80120	9.7	2.4	120	5.9	35M120
80140	10.2	2.5	140	6.1	35M140
80160	10.4	2.8	160	6.4	35M160
80180	10.9	2.9	180	6.7	35M180
80200	11.6	2.8	200	7.0	35M200
80220	11.8	3.0	220	7.2	35M220
80240	12.0	3.1	240	7.3	35M240
80260	12.2	3.3	260	7.5	35M260
80280	12.5	3.3	280	7.7	35M280
80300	12.9	3.4	300	7.9	35M300
80325	13.2	3.5	325	8.1	35M325
80350	13.5	3.6	350	8.3	35M350
80375	13.8	3.7	375	8.5	35M375
80400	14.2	3.7	400	8.7	35M400
80450	14.7	3.8	450	9.0	35M450
80500	15.3	4.0	500	9.4	

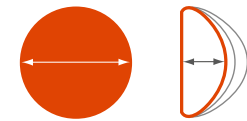


ES 801 Medium Profile 100% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
801100	8.5	2.8	100	5.6	35R100
801120	9.0	2.9	120	5.9	35R120
801140	9.3	3.2	140	6.2	35R140
801160	10.0	3.1	160	6.4	35R160
801180	10.6	3.2	180	6.8	35R180
801200	11.1	3.3	200	7.1	35R200
801220	11.5	3.2	220	7.2	35R220
801240	11.7	3.5	240	7.5	35R240
801260	11.9	3.6	260	7.6	35R260
801280	12.3	3.7	280	7.9	35R280
801300	12.4	3.9	300	8.0	35R300
801325	12.7	4.0	325	8.2	35R325
801350	13.0	4.0	350	8.4	35R350
801375	13.3	4.2	375	8.6	35R375
801400	13.7	4.2	400	8.8	35R400
801450	14.1	4.4	450	9.1	35R450
801500	14.5	4.6	500	9.4	
801550	15.0	4.8	550	9.8	
801600	15.8	4.5	600	9.9	
801650	15.8	4.9	650	10.2	
801700	15.8	5.4	700	10.5	
801750	15.9	5.8	750	10.9	
801800	15.8	6.0	800	11.0	

Depending on the product specifications, 1g is approximately equivalent to 1cc.
The dimensions have a tolerance +/- 4mm on diameter and +/- 3mm for projection.

The Round Collection™

Round Micro-Textured | Soft Cohesive Gel

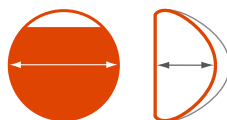

ES 81 High Profile 85% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
81100	7.9	2.9	100	5.4	37M100
81120	8.9	2.8	120	5.8	37M120
81140	8.8	3.2	140	6.0	37M140
81160	9.3	3.3	160	6.3	37M160
81180	9.5	3.6	180	6.6	37M180
81200	9.9	3.7	200	6.8	37M200
81220	10.3	3.7	220	7.0	37M220
81240	10.5	3.9	240	7.2	37M240
81260	10.7	4.1	260	7.5	37M260
81280	10.9	4.3	280	7.7	37M280
81300	11.1	4.4	300	7.8	37M300
81325	11.4	4.5	325	8.0	37M325
81350	11.7	4.7	350	8.3	37M350
81375	11.7	4.9	375	8.5	37M375
81400	12.2	4.9	400	8.7	37M400
81450	12.5	5.3	450	9.1	37M450
81500	12.6	5.8	500	9.5	



ES 811 Extra High Profile 100% Fill					
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	LVC (cm)	Size (Gel-filled)
811100	7.3	3.4	100	5.5	37R100
811120	8.0	3.5	120	5.9	37R120
811140	8.5	3.8	140	6.3	37R140
811160	8.8	3.9	160	6.5	37R160
811180	9.1	4.3	180	7.0	37R180
811200	9.4	4.3	200	7.1	37R200
811220	9.7	4.4	220	7.3	37R220
811240	10.2	4.4	240	7.5	37R240
811260	10.4	4.5	260	7.6	37R260
811280	10.6	4.7	280	7.9	37R280
811300	10.8	4.9	300	8.1	37R300
811325	11.0	5.1	325	8.3	37R325
811350	11.3	5.1	350	8.4	37R350
811375	11.6	5.2	375	8.6	37R375
811400	11.8	5.5	400	9.0	37R400
811450	12.1	5.7	450	9.2	37R450
811500	12.3	6.1	500	9.6	
811550	12.5	6.4	550	9.9	
811600	13.1	6.4	600	10.2	
811650	13.6	6.4	650	10.4	
811700	13.4	7.1	700	10.8	
811750	13.5	7.4	750	11.1	
811800	13.6	7.7	800	11.4	

Depending on the product specifications, 1g is approximately equivalent to 1cc.
The dimensions have a tolerance +/- 4mm on diameter and +/- 3mm for projection.



Luna XT™

The newly CE approved breast implant, under the latest Medical Device Regulation (MDR)⁴ by GC Aesthetics®. Luna XT™ is an industry-leading innovation to meet the needs of both surgeons and patients in breast reconstruction procedures.

The implant technology fulfils all the rigorous standards required for medical devices and has reliable long-term performance.

Main features:

Placement & Insertion:

- Micro-Textured implants have been shown to provide tissue adhesion to decrease rotation or displacement¹².
- Designed with 7 orientation dots for surgeon location and placement guidance¹³.

Designed for safety:



- Shell performance: exceeds the ISO 14607 device standard specifications¹⁴.
- The implant is 100% gel filled with highly cohesive gel¹⁷.
- Patch has high strength which exceeds ISO 14607 mechanical integrity requirements¹⁴.
- A variety of projections, heights, and volumes to suit different patient needs⁴.





Luna XT™

Anatomical Micro-Textured | Natural Cohesive Gel

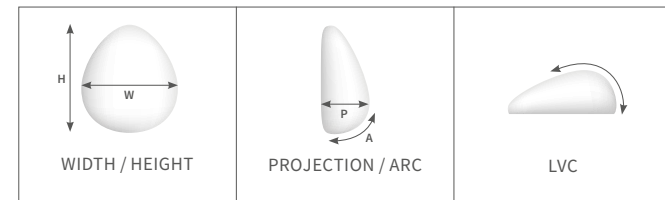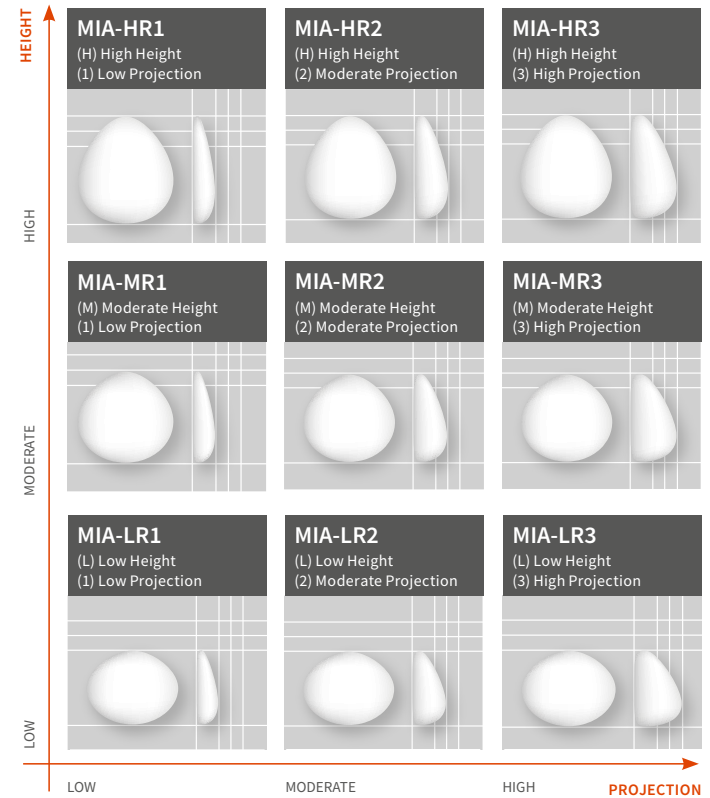

MIA-HR1		High Height, Low Projection				
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-HR1/135	9.5	10.0	3.0	4.6	135	6.9
MIA-HR1/155	10.0	10.5	3.2	4.8	155	7.3
MIA-HR1/185	10.5	11.0	3.4	5.1	185	7.7
MIA-HR1/215	11.0	11.5	3.6	5.4	215	8.1
MIA-HR1/245	11.5	12.0	3.8	5.6	245	8.5
MIA-HR1/275	12.0	12.5	4.0	5.9	275	8.9
MIA-HR1/315	12.5	13.0	4.2	6.2	315	9.3
MIA-HR1/345	13.0	13.5	4.4	6.4	345	9.7
MIA-HR1/400	13.5	14.0	4.6	6.7	400	10.1
MIA-HR1/450	14.0	14.5	4.8	7.0	450	10.5

MIA-MR1		Moderate Height, Low Projection				
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-MR1/125	9.5	9.0	3.0	4.4	125	6.5
MIA-MR1/140	10.0	9.5	3.2	4.6	140	6.9
MIA-MR1/160	10.5	10.0	3.4	4.9	160	7.3
MIA-MR1/195	11.0	10.5	3.6	5.2	195	7.7
MIA-MR1/215	11.5	11.0	3.8	5.4	215	8.1
MIA-MR1/245	12.0	11.5	4.0	5.7	245	8.5
MIA-MR1/280	12.5	12.0	4.2	6.0	280	8.9
MIA-MR1/320	13.0	12.5	4.4	6.2	320	9.3
MIA-MR1/360	13.5	13.0	4.6	6.5	360	9.7
MIA-MR1/405	14.0	13.5	4.8	6.8	405	10.1

MIA-LR1		Low Height, Low Projection				
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-LR1/110	9.5	8.0	3.0	4.2	110	6.1
MIA-LR1/130	10.0	8.5	3.2	4.4	130	6.5
MIA-LR1/145	10.5	9.0	3.4	4.7	145	6.9
MIA-LR1/175	11.0	9.5	3.6	5.0	175	7.3
MIA-LR1/200	11.5	10.0	3.8	5.2	200	7.7
MIA-LR1/230	12.0	10.5	4.0	5.5	230	8.1
MIA-LR1/265	12.5	11.0	4.2	5.8	265	8.5
MIA-LR1/305	13.0	11.5	4.4	6.1	305	8.9
MIA-LR1/345	13.5	12.0	4.6	6.3	345	9.3
MIA-LR1/385	14.0	12.5	4.8	6.6	385	9.7

MIA-HR2		High Height, Moderate Projection					
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-HR2/165	9.5	10.0	3.7	5.2	165	7.6	GS-XF2 145
MIA-HR2/190	10.0	10.5	3.9	5.4	190	8.0	GS-XF2 170
MIA-HR2/215	10.5	11.0	4.1	5.7	215	8.4	GS-XF2 200
MIA-HR2/245	11.0	11.5	4.3	6.0	245	8.8	GS-XF2 230
MIA-HR2/275	11.5	12.0	4.5	6.2	275	9.2	GS-XF2 270
MIA-HR2/310	12.0	12.5	4.7	6.5	310	9.6	GS-XF2 315
MIA-HR2/345	12.5	13.0	4.9	6.8	345	10.0	GS-XF2 350
MIA-HR2/390	13.0	13.5	5.1	7.1	390	10.4	GS-XF2 400
MIA-HR2/440	13.5	14.0	5.3	7.3	440	10.8	GS-XF2 450
MIA-HR2/500	14.0	14.5	5.5	7.6	500	11.2	GS-XF2 505
MIA-HR2/575	14.5	15.0	5.7	7.9	575	11.6	
MIA-HR2/655	15.0	15.5	5.9	8.1	655	12.0	
MIA-HR2/710	15.5	16.0	6.1	8.4	710	12.4	

MIA-MR2		Moderate Height, Moderate Projection					
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-MR2/155	9.5	9.0	3.7	5.0	155	7.3	GS-XM2 130
MIA-MR2/175	10.0	9.5	3.9	5.2	175	7.6	GS-XM2 165
MIA-MR2/195	10.5	10.0	4.1	5.5	195	8.0	GS-XM2 185
MIA-MR2/230	11.0	10.5	4.3	5.8	230	8.4	GS-XM2 215
MIA-MR2/255	11.5	11.0	4.5	6.1	255	8.8	GS-XM2 250
MIA-MR2/290	12.0	11.5	4.7	6.3	290	9.2	GS-XM2 290
MIA-MR2/330	12.5	12.0	4.9	6.6	330	9.6	GS-XM2 320
MIA-MR2/375	13.0	12.5	5.1	6.9	375	10.0	GS-XM2 370
MIA-MR2/415	13.5	13.0	5.3	7.1	415	10.4	GS-XM2 415
MIA-MR2/465	14.0	13.5	5.5	7.4	465	10.8	GS-XM2 445
MIA-MR2/530	14.5	14.0	5.7	7.7	530	11.2	
MIA-MR2/595	15.0	14.5	5.9	8.0	595	11.6	
MIA-MR2/655	15.5	15.0	6.1	8.2	655	12.0	

MIA-LR2		Low Height, Moderate Projection				
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-LR2/135	9.5	8.0	3.7	4.8	135	6.9
MIA-LR2/160	10.0	8.5	3.9	5.1	160	7.3
MIA-LR2/190	10.5	9.0	4.1	5.3	190	7.7
MIA-LR2/215	11.0	9.5	4.3	5.6	215	8.1
MIA-LR2/245	11.5	10.0	4.5	5.9	245	8.5
MIA-LR2/280	12.0	10.5	4.7	6.2	280	8.9
MIA-LR2/320	12.5	11.0	4.9	6.4	320	9.3
MIA-LR2/355	13.0	11.5	5.1	6.7	355	9.7
MIA-LR2/405	13.5	12.0	5.3	7.0	405	10.1
MIA-LR2/440	14.0	12.5	5.5	7.2	440	10.5

MIA-HR3		High Height, High Projection					
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-HR3/190	9.5	10.0	4.5	5.9	190	8.5	GS-XF3 170
MIA-HR3/220	10.0	10.5	4.7	6.2	220	8.9	GS-XF3 200
MIA-HR3/250	10.5	11.0	4.9	6.4	250	9.3	GS-XF3 235
MIA-HR3/295	11.0	11.5	5.2	6.8	295	9.8	GS-XF3 280
MIA-HR3/330	11.5	12.0	5.4	7.1	330	10.2	GS-XF3 320
MIA-HR3/370	12.0	12.5	5.6	7.3	370	10.6	GS-XF3 365
MIA-HR3/415	12.5	13.0	5.9	7.7	415	11.1	GS-XF3 405
MIA-HR3/460	13.0	13.5	6.2	8.1	460	11.6	GS-XF3 465
MIA-HR3/505	13.5	14.0	6.3	8.2	505	11.9	GS-XF3 515
MIA-HR3/580	14.0	14.5	6.4	8.4	580	12.2	GS-XF3 540
MIA-HR3/645	14.5	15.0	6.6	8.7	645	12.6	
MIA-HR3/720	15.0	15.5	6.8	9.0	720	13.0	
MIA-HR3/775	15.5	16.0	7.0	9.2	775	13.3	

MIA-MR3		Moderate Height, High Projection					
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)	Sizer (Gel-filled)
MIA-MR3/175	9.5	9.0	4.5	5.7	175	8.1	GS-XM3 155
MIA-MR3/205	10.0	9.5	4.7	6.0	205	8.5	GS-XM3 190
MIA-MR3/235	10.5	10.0	4.9	6.3	235	8.9	GS-XM3 215
MIA-MR3/275	11.0	10.5	5.2	6.6	275	9.4	GS-XM3 260
MIA-MR3/310	11.5	11.0	5.4	6.9	310	9.8	GS-XM3 300
MIA-MR3/345	12.0	11.5	5.6	7.2	345	10.2	GS-XM3 335
MIA-MR3/395	12.5	12.0	5.9	7.5	395	10.7	GS-XM3 375
MIA-MR3/440	13.0	12.5	6.2	7.9	440	11.2	GS-XM3 440
MIA-MR3/475	13.5	13.0	6.3	8.1	475	11.5	GS-XM3 485
MIA-MR3/545	14.0	13.5	6.4	8.3	545	11.8	GS-XM3 520
MIA-MR3/610	14.5	14.0	6.6	8.5	610	12.2	
MIA-MR3/685	15.0	14.5	6.8	8.8	685	12.6	
MIA-MR3/735	15.5	15.0	7.0	9.1	735	13.0	

MIA-LR3		Low Height, High Projection				
Product Code	Width (cm)	Height (cm)	Proj. (cm)	Arc Lgth. (cm)	Volume (cc)	LVC (cm)
MIA-LR3/160	9.5	8.0	4.5	5.6	160	7.8
MIA-LR3/190	10.0	8.5	4.7	5.9	190	8.2
MIA-LR3/215	10.5	9.0	4.9	6.1	215	8.6
MIA-LR3/260	11.0	9.5	5.2	6.5	260	9.1
MIA-LR3/290	11.5	10.0	5.4	6.8	290	9.5
MIA-LR3/330	12.0	10.5	5.6	7.0	330	9.9
MIA-LR3/385	12.5	11.0	5.9	7.4	385	10.4
MIA-LR3/415	13.0	11.5	6.2	7.8	415	10.9
MIA-LR3/480	13.5	12.0	6.3	7.9	480	11.2
MIA-LR3/525	14.0	12.5	6.4	8.1	525	11.5

Depending on the product specifications, 1g is approximately equivalent to 1cc. The dimensions have a tolerance +/- 2mm for the height and width and ± 3mm for the projection.



CE
0483

FixNip™ NRI

FixNip™ NRI is the innovative biocompatible¹⁵ smooth silicone implant specially designed for a long-lasting projection¹⁶ of the Nipple-Areola Complex (NAC).

Main features:

- Minimally invasive, an office-based procedure¹⁷.
- Subcutaneous pocket, under local anaesthesia¹⁷.
- Easy to use¹⁷.
- Performed on average in 15 minutes or less¹⁷.
- MRI compatible^{16,18}.
- Very low rate of device related adverse events¹⁷.
- Patients are relieved to see their breast fully reconstructed¹⁹.
- Minimum healing time: Procedure can be performed under local anesthesia supporting return to normal activities¹⁷.
- High patient satisfaction levels in terms of aspect and feel/touch¹⁹.
- Tested successfully in different patient profiles: patients with breast implants, expanders, flaps (PAP, DIEP and latissimus dorsi), lipofilling and a combination of those, and in both non-irradiated and irradiated patients¹⁹.





Blossom™ IVTE

GC Aesthetics® Integrated Valve Tissue expanders (Blossom™ IVTE), must be filled with a sterile isotonic saline solution only. A 21 gauge standard needle must be used in order to inflate the device through the injection port.

Main features:

- Anatomical and Micro-Textured surface roughness².
- Re-usable and ergonomic External Locating Device (ELD)².
- Extra strong magnet for ease of location².
- Large integrated injection port with self-sealing technology, which helps to facilitate the transfer of saline solution from the needle to the product².
- Reinforced base layer to ensure unidirectional expansion².
- Available in full, moderate and low height projections².
- Available in a range of sizes to complement GCA® breast implants².

External Location Device:


Supplied sterile with the Integrated-Valve Tissue Expander.





Blossom™ IVTE


Anatomical Micro-Textured

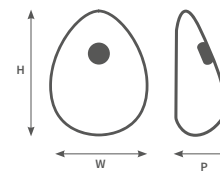

TE-XF3

(F) Full Height (3) High Projection


Product Code	Width (cm)	Height (cm)	Projection (cm)	Volume (cc)
TE-XF3 300	11.0	11.5	5.2	300
TE-XF3 400	12.0	12.5	5.6	400
TE-XF3 500	13.0	13.5	6.2	500
TE-XF3 600	14.0	14.5	6.4	600
TE-XF3 700	15.0	15.5	6.8	700
TE-XF3 800	15.5	16.0	7.0	800

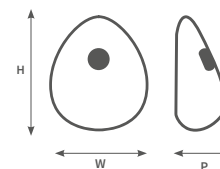

TE-XM3

(M) Moderate Height (3) High Projection

Product Code	Width (cm)	Height (cm)	Projection (cm)	Volume (cc)
TE-XM3 200	10.0	9.5	4.7	200
TE-XM3 300	11.0	10.5	5.2	300
TE-XM3 400	12.0	11.5	5.6	400
TE-XM3 500	13.0	12.5	6.2	500
TE-XM3 600	14.0	13.5	6.4	600
TE-XM3 700	15.0	14.5	6.8	700

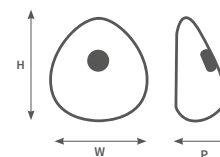

TE-XL3

(L) Low Height (3) High Projection

Product Code	Width (cm)	Height (cm)	Projection (cm)	Volume (cc)
TE-XL3 200	10.0	8.5	4.7	200
TE-XL3 300	11.5	10.5	5.4	300
TE-XL3 400	12.5	11.0	5.9	400
TE-XL3 500	13.5	12.0	6.3	500

All device dimensions correspond to the suggested volume.

POST-SURGERY

Patient After Care Solutions






POST-SURGERY

GCA® expands its catalogue of safe and innovative solutions to contribute even more to women's healthcare and well-being, offering post-surgery support.

Patients must minimize implant movement during the early post-operative period to maintain the dissected pocket boundaries. Surgeons should recommend support garments and educate patient on the importance of post-operative care.

When planning for a breast augmentation, many factors must be taken into account. If we talk about the post-surgery care, another important factor is the warranty on the implants. GCA® offers the most comprehensive warranty in the market. It contains a full coverage plan to provide confidence to all patients.

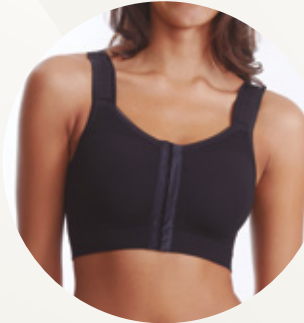

NOA X0



NOA X1



NOA X2



NOA B1

NoaTM

GCA® compression garments designed for comfort to support optimum results after cosmetic or reconstruction surgery.

The appropriate care after surgery may be as important as the surgery itself²⁰. Clinical studies show that compression aids the body's healing process, reducing post-surgical complications, and helping to optimise results after surgery²⁰.

Main features:

- Designed for use immediately following breast surgery, when the skin is still tender and sensitive².
- Seamless design for superior comfort².
- Four-way stretch material that conforms to the breast contour².
- Targeted built in compression zones².
- No underwire or moulded cups².
- Breathable and moisture-permeable material².
- Tested for harmful substances according to Oeko-Tex 100 standard².
- A flexible sizing concept with less references².

Noa
By GC Aesthetics®



Easy-to-use²



Discreet²¹



UV blocking²²



Water resistant²²

Healena™

Medical Grade Scar Treatment

Healena™ reduces and improves the appearance of surgical scars. The product advances the benefits of traditional silicone sheeting with medical grade compression²¹. Designed by a board certified plastic surgeon, Healena™ is a multi-layer adhesive single-use dressing composed of 2 layers: a) bottom layer; perforated silicone hydrocolloid and b) top layer; polyurethane of 25 microns film with a strong acrylate medical glue, which ensure uniform adhesion and continuous contraction and compression simultaneously on the scar, improving cosmetic outcomes²¹.

Indications:

Healena™ is indicated for use on surgical or post traumatic linear scars after primary healing has been reached (after removal of the suture threads). The dressing should only be changed if it comes off spontaneously; it can stay in place for a month or more.

Healena™

By GC Aesthetics®



GCA Comfort Plus™ Warranty

All GCA® Breast Implants are covered by the GCA Comfort Plus™ Warranty, a comprehensive lifetime warranty²³.

Key features:

- Lifetime implants replacement in case of rupture and capsular contracture (grade III & IV on the Baker Scale)²³.
- Reoperation reimbursement in the rare event of BIA-ALCL and free of charge implant replacement²³.
- No time limitation, it is available for the lifetime of the patient²³.



*This warranty is dated 01 July 2020 and replaces and supersedes any prior GC Aesthetics®, Nagor® or Eurosilicone® guarantee or warranty terms.

To learn more about the GCA Comfort Plus™ Warranty visit:
gcaesthetics.com/comfort-plus-warranty



Face and Body Implants





Face and Body Implants

GC Aesthetics® completes its broad implant portfolio with face and body implants and remote valve tissue expanders to offer surgeons a quality surgical option.

The face defines our appearance and how we are perceived. Facial implants are used to bring balance to the appearance of the face.

Body Implants can be intended to replace missing body parts or provide support to tissues. Patients can look and feel better with a perfectly proportioned body.

The Renuvé™ RVTE series is the remote valve tissue expanders by GC Aesthetics®. The device is mainly used in cosmetic or reconstructive plastic surgery to expand an area of healthy, soft tissue, primarily skin, to gain new tissue that can be used to repair a specific area of the body.

Face and Body Implants & Tissue Expanders

CHIN IMPLANTS

SMOOTH



Thin half-moon



Wide half-moon



GLUTEAL IMPLANTS*

SMOOTH



GLUTEAL



TESTICULAR IMPLANTS

SMOOTH



Ovoid



CALF IMPLANTS

SMOOTH



Symmetric



Asymmetric



*Product available in limited markets.
Not Available in Europe.

REMOTE VALVE TISSUE EXPANDERS

Renuv ™ RVTE



Round



Rectangular



Cylindrical straight



Cylindrical curved



Croissant shape






Main features:


- Semi-curved Half-moon Smooth Opaque Implant².
- The shell consists of a very resistant silicone elastomer which is highly resistance to mechanical stress².
- Available in both thin semi-curve and wide semi-lune shapes².
- Wide range of volumes and sizes².

Chin Implants

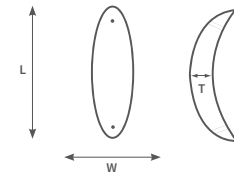

ES 32 / 34 / 36 / 38

Product Code	Length (cm)	Width (cm)	Thickness (cm)
32000	3.4	1.3	0.9
34000	3.9	1.3	0.9
36000	4.2	1.5	0.9
38000	4.5	1.3	0.7

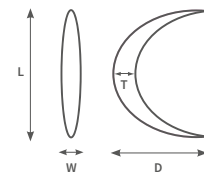

ES 325


Product Code	Length (cm)	Depth (cm)	Width (cm)	Thickness (cm)
32501	6.2	3.0	0.6	0.6
32502	6.5	3.4	0.8	0.8
32503	7.0	3.8	0.9	0.8




Main features:

- Round and anatomical gluteal silicone implants².
- Smooth surface².
- The shell consists of a very resistant silicone elastomer and highly elastic equipped with a barrier layer².
- Implant filled with highly cohesive silicone gel².
- Wide range of volumes and sizes².

Gluteal Implants

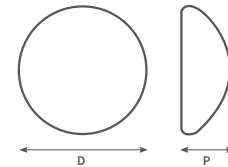

ES 62		Smooth Surface, Round		
Product Code	Diameter (cm)	Projection (cm)	Volume (cc)	
62140	8.4	4.1	140	
62180	9.1	4.3	180	
62200	9.4	4.4	200	
62220	9.7	4.5	220	
62240	10.3	4.7	240	
62260	10.4	4.7	260	
62280	10.6	4.9	280	
62300	10.8	5.0	300	

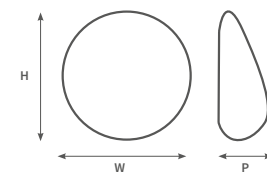


ES 60		Smooth Surface, Anatomical		
Product Code	Height (cm)	Width (cm)	Projection (cm)	Volume (cc)
60150	13.1	9.6	2.3	150
60210	14.6	10.7	2.6	210
60230	12.0	11.2	3.5	230
60275	12.4	11.4	3.8	275
60300	13.9	10.9	3.6	300
60325	14.7	11.3	3.7	325
60375	14.7	11.7	4.1	375
60425	15.0	12.4	4.5	425
60480	15.5	13.2	4.8	480

Product available in limited markets.
Not Available in Europe.



Testicular Implants

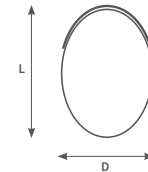

ES 520		Smooth Surface	
Product Code	Diameter (cm)	Length (cm)	Volume (cc)
52006	2.2	2.8	6.0
52012	2.8	3.3	12.0
52018	3.0	4.2	18.0
52022	3.3	4.7	22.0
52026	3.3	5.2	26.0

Main features:

- Ovoid silicone testicular implant².
- Smooth Surface².
- High performance silicone elastomer shell².
- Shell provides high tear resistance with low permeability due to the incorporation of a barrier layer within the shell².
- The shell is filled with highly firm silicone gel².
- The implant is sealed with a disc-shaped closure patch².
- The upper pole of the device is equipped with a fixation tab, which allows the user to attach the implant to the surrounding tissue².
- GCA® Testicular Implants are available in five sizes².




Main features:

- Smooth Surface².
- Filled with highly firm silicone gel².
- A barrier layer has been formulated to minimise permeability of the shell to the silicone gel².
- Disc shaped patch, facilitating the insertion of the implant².
- Available in both symmetrical or asymmetrical shapes².
- Wide range of volumes and sizes².

Calf Implants



ES 54 Symmetric Shape				
Product Code	Length (cm)	Width (cm)	Projection (cm)	Volume (cc)
54035	16.6	5.0	1.0	35
54045	16.3	5.0	1.5	45
54055	17.7	5.2	1.5	55
54060	16.8	3.3	2.0	60
54075	17.3	5.1	2.0	75
54090	17.0	3.7	2.8	90
54110	17.7	5.8	2.1	110
54120	18.9	5.2	2.6	120
54160	19.5	6.9	2.3	160
54180	19.0	5.9	3.1	180
54220	18.7	7.1	3.2	220

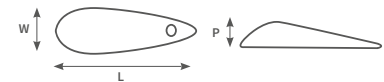

ES 54 Asymmetric Shape				
Product Code	Length (cm)	Width (cm)	Projection (cm)	Volume (cc)
54135	22.0	4.9	2.3	135
54140	18.6	5.3	2.6	140
54150	13.9	6.4	2.9	150
54165	23.9	5.6	2.3	165



Renuvé™ RVTE

The Renuvé™ RVTE series is the remote valve tissue expanders by GC Aesthetics® available for body contouring (not applicable to the breast).

The product is indicated in cosmetic or reconstructive surgery, to expand an area of healthy, soft tissue, primarily skin, and in order to create a flap that can be used to repair a damaged area of the body².

The Renuvé™ RVTE series by GC Aesthetics® may be used in a variety of surgical techniques².





Renuvé™ RVTE

Renuvé™ RVTE collection has significant advantages²:

- A unique port base diameter of 28mm for all shapes.
- Memory-shape injection port prevents escape of sterile physiological saline solution.
- Metal safety needle stop, built-in to the base of the valve, to avoid puncture until the base of the injection port.
- The injection port is well away from the expander avoiding any risks of puncture while inflating the expander.
- Reinforced base of the shell to keep the shape of the expander.
- Medical stainless steel connector helps the tube length customization for every patient.
- Radiopaque tubing and injection port for ease of detection on X-ray.
- The expander exceeds the mechanical ISO standards, which assures its shell integrity resistance.

Main features:

- Smooth Surface Shell¹.
- Available in 5 different configurations²: Round, Rectangular, Cylindrical Straight, Cylindrical Curved, Croissant.
- Volumes from 200cc to 1,100cc².
- Supplied Sterile².

¹All device dimensions correspond to the suggested volume.


²Please consult availability of this product in your country.



Renuvé™ RVTE Round

ES 12


Product Code	Diameter (cm)	Height (cm)	Valve (cm)	Volume (cc)
12200	7.9	5.1	2.8	200
12250	9.0	5.7	2.8	250
12300	10.0	6.3	2.8	300
12375	13.0	5.0	2.8	375
12400	11.0	6.0	2.8	400
12500	12.8	6.3	2.8	500
12600	13.8	6.7	2.8	600
12700	14.8	7.5	2.8	700
12800	15.7	6.7	2.8	800
12900	16.6	6.5	2.8	900
121100	17.4	8.2	2.8	1100

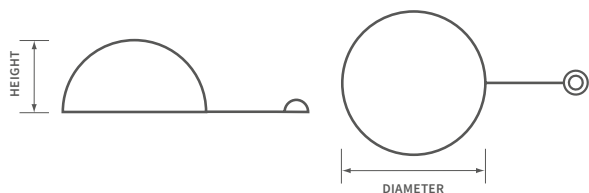

Renuvé™ RVTE Rectangular

ES 14

Product Code	Lenght (cm)	Width (cm)	Projection (cm)	Valve (cm)	Volume (cc)
14200	8.5	5.0	4.9	2.8	200
14300	10.0	6.0	6.3	2.8	300
14400	10.8	5.8	5.9	2.8	400
14520	11.7	7.8	5.4	2.8	520
14650	11.7	9.0	5.8	2.8	650
14720	14.5	7.8	5.8	2.8	720
14750	13.5	8.7	5.9	2.8	750
14800	15.5	8.2	5.8	2.8	800
14950	17.5	8.8	5.5	2.8	950

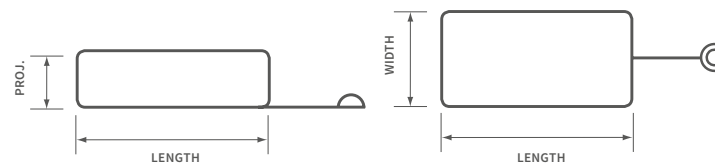





Renuvé™ RVTE Cylindrical Straight

ES 16


Product Code	Lenght (cm)	Width (cm)	Projection (cm)	Valve (cm)	Volume (cc)
16200	14.0	4.5	4.5	2.8	200
16250	10.2	7.2	5.4	2.8	250
16300	16.0	5.0	5.0	2.8	300
16400	17.5	6.0	6.0	2.8	400

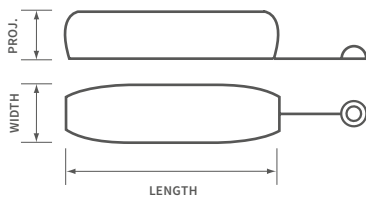

Cylindrical Curved

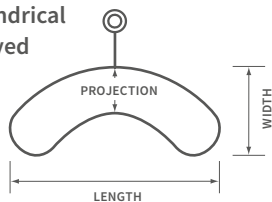

Croissant

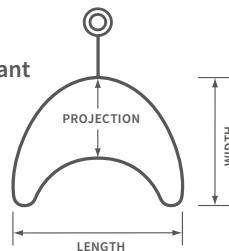

Renuvé™ RVTE Cylindrical Curved

ES 16C

Product Code	Lenght (cm)	Width (cm)	Projection (cm)	Valve (cm)	Volume (cc)
16C250	14.5	7.7	4.5	2.8	250



Renuvé™ RVTE Croissant

ES 19

Product Code	Lenght (cm)	Width (cm)	Projection (cm)	Valve (cm)	Volume (cc)
19250	13.0	10.0	5.7	2.8	250
19500	17.0	12.0	7.8	2.8	500

References

1. GC Aesthetics®. GCA® catalogue, 2023.
2. GC Aesthetics®. Internal data, 2023.
3. Data on file, GC Aesthetics®, Countries Sales Report, 2021.
4. BSI, EU Technical Documentation Assessment Certificate, MDR 746961 R000. 2022-09-24
5. Duteille F, Perrot P, Bacheley MH, Bell E, Stewart S. Ten-Year Safety Data for Eurosilicone's Round and Anatomical Silicone Gel Breast Implants. *Aesthet Surg J Open Forum*. 2019;1(2):ojz012. Published 2019 Apr 27.
6. Kooiman L, Torensma B, Stevens H, van der Lei B. Single Center and Surgeon's Long-Term (15-19 Years) Patient Satisfaction and Revision Rate of Round Textured Eurosilicone Breast Implants. *Aesthet Surg J*. 2022;42(5):NP282-NP292.
7. Newman & Davison (2018). Effect of Keller Funnel on the Rate of Capsular Contracture in Periareolar Breast Augmentation. *Plastic and Reconstructive Surgery*. *Global open*, 6(6), e1834. <https://doi.org/10.1097/GOX.0000000000001834>.
8. Flugstad, N.A. et al (2016). Does Implant Insertion with a Funnel Decrease Capsular Contracture? A Preliminary Report Nicholas. *Aesthetic Surgery Journal*. 26 (5). 550-556.
9. Montemurro, P. et al (2020). Controllable factors to reduce the Rate of Complications in Primary Breast Augmentation: Review of the Literature. *Aesthetic Plastic Surgery*. 2021;45(2):498-505. doi:10.1007/s00266-020-01726-x.
10. Montemurro, P. et al (2019). Implant insertion time and incision length in breast augmentation surgery with the Keller Funnel: Results from a comparative study. *Aesthetic Plastic Surgery*. 43; 881 – 889.
11. HydroCone™. Instructions for use. 2020.
12. Julia A. Chiemi, BSc; and S. Sean Kelishadi, MD, FACS. A Rationale for Micro-textured Breast Implant Augmentation. *Aesthetic Surgery Journal Open Forum* 2022, 1–10
13. Data on file, GC Aesthetics®, Luna XT™ Instructions for Use
14. Data on file, GC Aesthetics®, Luna XT™ Design Output summary report
15. FixNip LTD. Internal data, 2021. Biocompatible test.
16. FixNip LTD. Internal data, 2022.
17. GC Aesthetics®, Internal Data. Product Trial, 2022.
18. FixNip™ NRI. Instructions for use. 2022.
19. GC Aesthetics® Patient Video Testimonials, April 2022. Internal Data.
20. Gladfelter, J (2007). Compression Garments 101. *Plastic Surgical Nursing: Vol. 27, Issue 2*.
21. Healena™. Instructions for use. 2024.
22. Determination of moisture vapour transmission and measurement of UV transmission resistance. Test Report, National Metrology and Testing Laboratory. 2023.
23. GCA Comfort Plus™ Warranty brochure. A50-239 Rev 2_Nov 2023



GC Aesthetics®

A Confident Choice for Life™

GC Aesthetics®, registered in Ireland with number 466501. Registered office: Suite 601, Q House, Furze Road, Sandyford Industrial Estate, Dublin 18, Ireland.
© 2025 GC Aesthetics®. A50-272_Rev03_Jan 2025.

GC Aesthetics® is a registered trademark of Global Consolidated Aesthetics Limited.

For more information about our products please visit:
gcaesthetics.com
or email: **info@gcaesthetics.com**

