


Analysis of herbicides

BPX35

Column part number	054711		
Phase	BPX35	Final temperature	300°C, 5 min
Column	25 m x 0.22 mm x 0.25 µm	Detector	FID, 380°C
Initial temperature	80°C	Injection mode	Split (20:1)
Rate	10°C/min	Carrier gas	He, 100 kPa

Components

1. Eptam®	7. Balan®	13. Simazine	19. Bromacil
2. Sutan®	8. Ro-Neet®	14. Terbacil	20. Oxadiazon
3. Vernam®	9. Propachlor	15. Sencor®	21. GOAL®
4. Tillam®	10. Tolban®	16. Dual®	22. Hexazinone
5. Ordram®	11. Propazine	17. Paarlant®	
6. Treflan®	12. Atrazine	18. Prowl®	

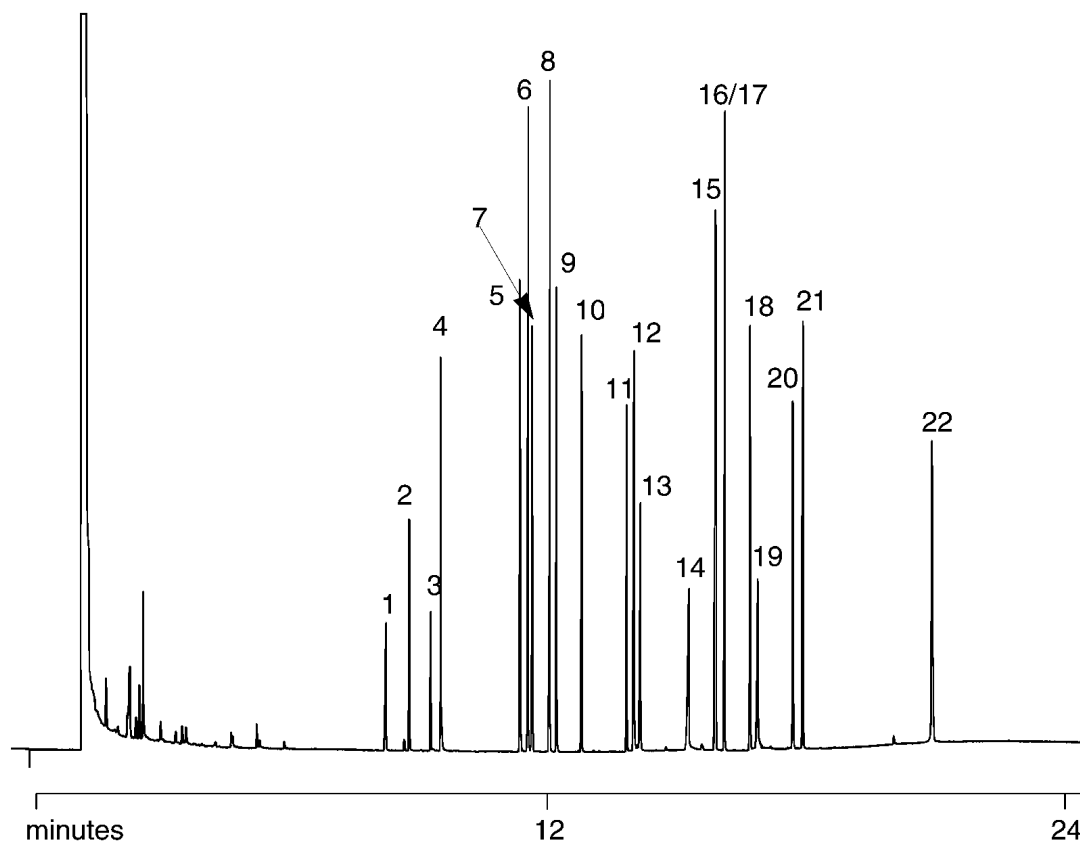

For more information visit www.trajanscimed.com or contact techsupport@trajanscimed.com