



Extended simulated distillation of motor oil

Application Note

Energy & Fuels

Authors

Agilent Technologies, Inc.

Introduction


Extended simulated distillation of a motor oil according to D2887 with an Agilent HT-SimDist CB UltiMetal column.



Agilent Technologies

Conditions

Technique : GC
Column : Agilent HT-SimDist CB UltiMetal, 0.53 mm x 5 m
(df = 0.09 μ m) (Part No. CP7569)
Temperature : -20 °C \rightarrow 430 °C, 10 °C/min, 5 min
Carrier Gas : He, 25 kPa
Injector : PTV
Detector : FID

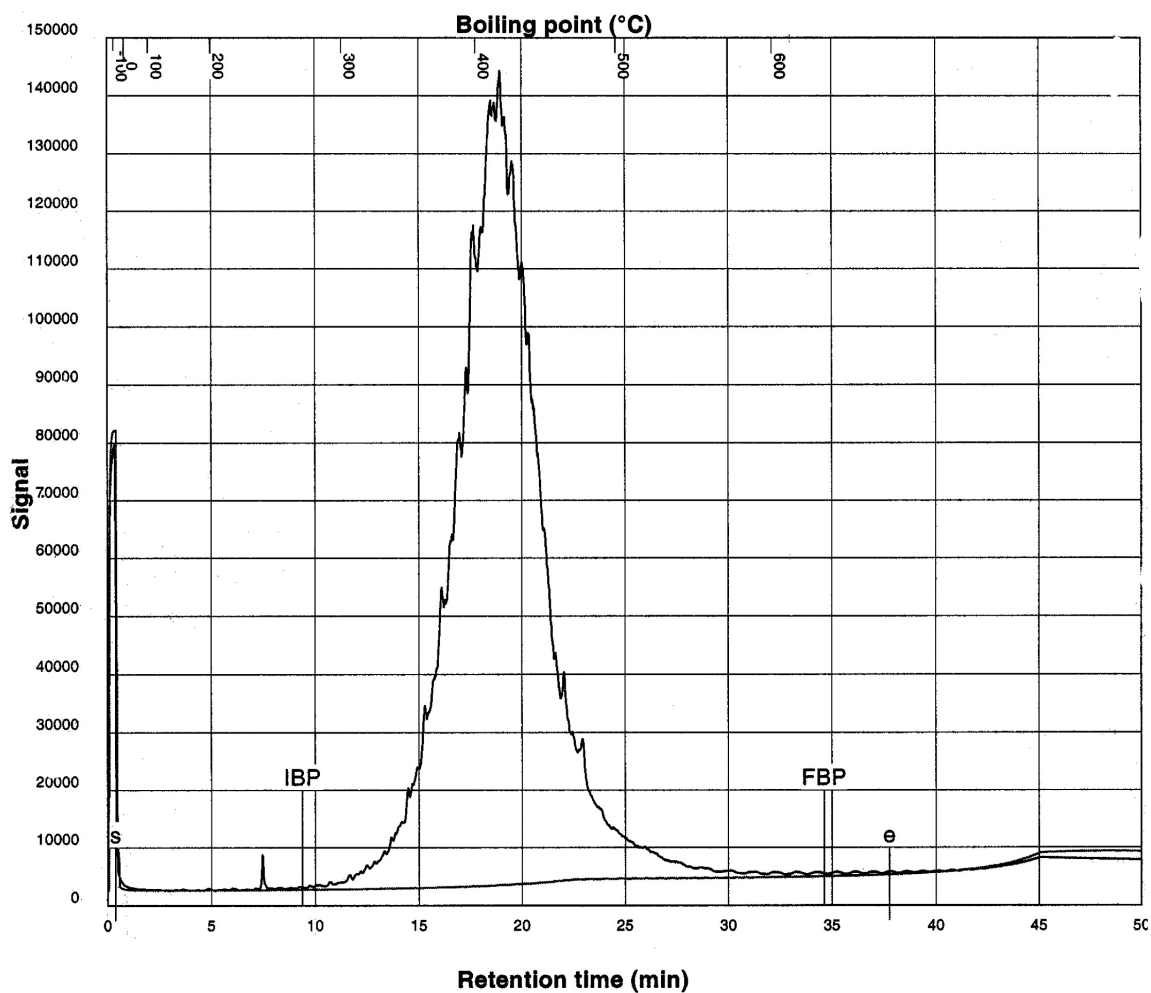

www.agilent.com/chem

This information is subject to change without notice.

© Agilent Technologies, Inc. 2011

Printed in the USA

31 October, 2011

First published prior to 11 May, 2010

A02262



Agilent Technologies