


## “Phthalate-free Quartz Wool”

Phthalate-free Quartz Wool helps prevent scattering of solids from sample cup during the sample cup recovery when Auto-Shot Sampler is used. Also, this product is characterized by having extremely low level of seven phthalates including DIBP, DBP, BBP, DEHP, DNOP, DINP and DIDP, and proven by test certificate. There is no concern about phthalate contaminations before analysis.

### Features

- Ultra clean quartz wool and ready to use without any further treatment.
- Extremely low level of seven phthalates guaranteed by test certificate.

### < Using Phthalate-free Quartz Wool >

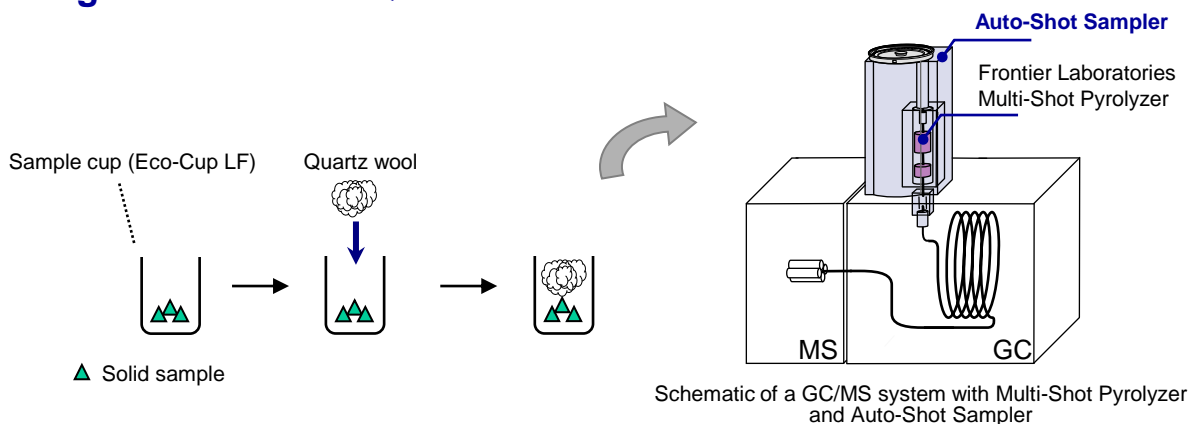

- 1) Place a solid sample in a sample cup.
- 2) Put approx. 2 mg of quartz wool over the sample.
- 3) Load the sample cups to the carousel of Auto-Shot Sampler and start measurements.

### < Specifications >

- Material: Quartz
- Max use temperature: 1050°C
- Level of seven phthalates: Test certificate included.  
(DIBP, DBP, BBP, DEHP, DNOP - less than 1 ppm each)  
(DINP, DIDP - less than 5 ppm each)



| Product name               | Product No. | Quantity                       |
|----------------------------|-------------|--------------------------------|
| Phthalate-free Quartz Wool | PY1-5121    | 200 mg each in 5 glass bottles |

 **FRONTIER LABORATORIES LTD.**

4-16-20 Saikon, Koriyama, Japan, 963-8862  
Phone: 81(24) 935-5100 FAX: 81(24) 935-5102  
<http://www.frontier-lab.com>