


# Use of Surface Plots and No Surface Plots with Peak Markers for GCxGC Data

Life Science and Chemical Analysis Centre  
LECO Corporation, Saint Joseph, Michigan USA

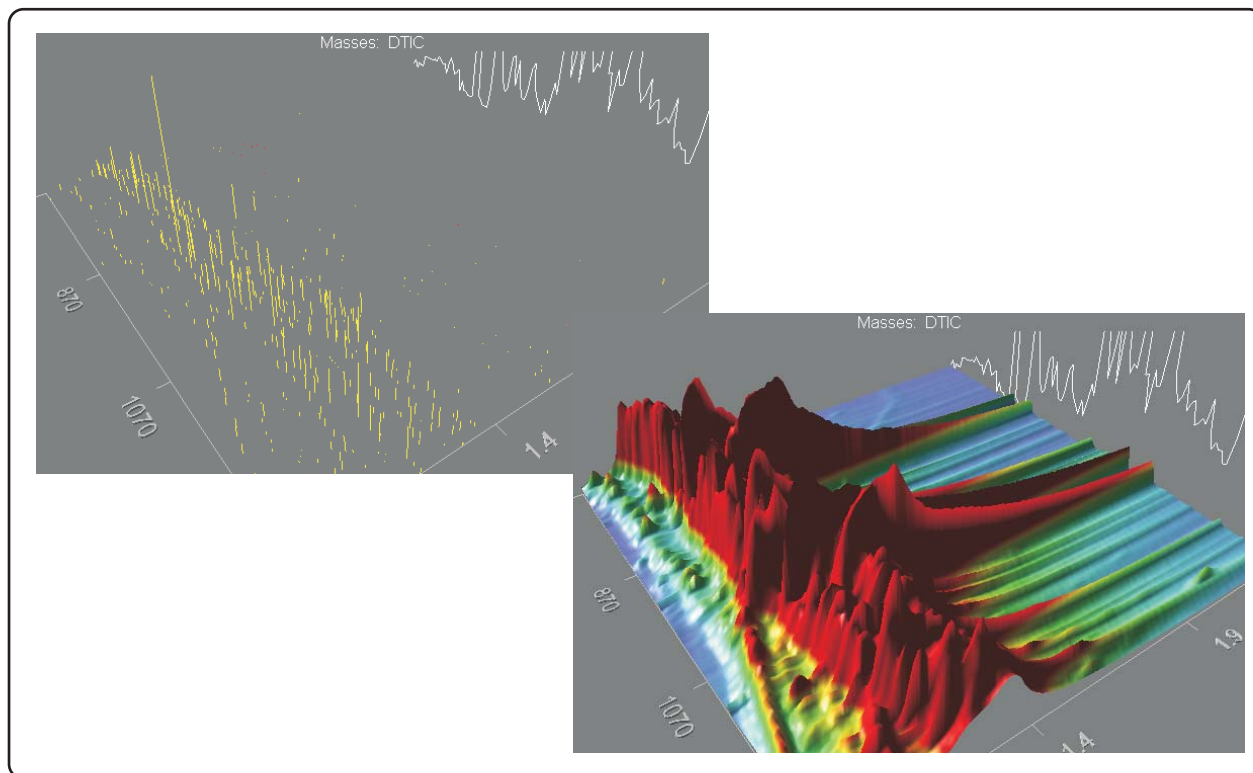

**By using the “No Surface Plot” and displaying Peak Markers, the identified peaks and their intensities can be more easily viewed in a sample with significant overlapping of peaks due to a complex matrix.**

## Columns:

30 m x 0.25 mm x 0.25  $\mu$ m Rtx-1ms  
1.25 m x 0.10 mm x 0.10  $\mu$ m BPX-50

## TOFMS:

45 to 760 m/z @ 100 spectra/s

For questions on this analysis  
e-mail us at: [life\\_science@leco.com](mailto:life_science@leco.com)

For a complete listing of snapshots  
and application notes visit us on the  
web at [www.leco.com](http://www.leco.com)

