



Vector Rehab Power Chair

Power tilt / Power Recline



Vector Rehab Power Chair

Power tilt / Power Recline

Features

- Centre of gravity changes when operating the power tilt/recline seat function (for safety reasons the centre of gravity of the user will always stay above the drive wheels)
- Shear reduction for comfort and prevention of injury
- LED lighting front & rear
- Adjustable seat width 18"/20"/22" and seat depth 16"/18"/20"
- Removable battery case for replacing and maintaining battery easily
- Dual in-line motors for enhanced efficiency, torque, range and performance
- Full front and rear adjustable suspension which provides a comfortable ride
- Power tilt operation in 5°~55° for pressure relief
- Power recline operation 90°~150° for comfort
- Optional seat elevator, power centre foot plate or power elevating leg rests, lateral & abductor supports
- PG R-NET controller standard
- Colours available Red or Ocean Blue



Seat elevator travel range:
30cm/12inch




Tilt



Recline

Dimensions cm /in

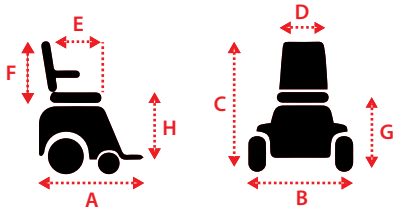

Overall Dimension : A x B x C	98 x 64x 107cm / 39in x 25in x 42in
Seat Dimension : D x E x F	46-56 x 40-51 x 55cm / 18-22 x 16-20 x 22in
Seat Height from Ground : G	47-52cm / 18.5-20.5in
Seat Height from Deck : H	35.5-40.5cm / 14-16in

Weight kg /lb

base	battery/ea	seat	total with battery	capacity
66/145	17/37	44/97	145 / 319	182 / 400
Seat Elevator	Power Footplate			
15/33	10/22			

Technical Specs

Max Speed	8kph/ 5mph
Range	20km/ 12.5mi.
Turn Radius	68cm/ 27"
Ground clearance	8cm/ 2"
Motor	2 Pole DC24V, 200W(Rating)/ 575W(Max)
Controller	PG R-Net 60A
Battery	12v/55-62a/H x 2pcs
Charger	6A off-board
Gradient	8°
Caster	F: 8" / R: 8"
Drive Wheel	14 PU tire
Brake	Intelligent, regenerative, electromagnetic brakes