

voyanceX: Acquisition Software

voyance.health



OUR PROMISE Better Health For All

AI ENABLED END TO END MEDICAL IMAGING SOLUTION




Navigating Your Health Path, Every Step Together.



PATIENT

GP PRACTICE

Powered by Voyance

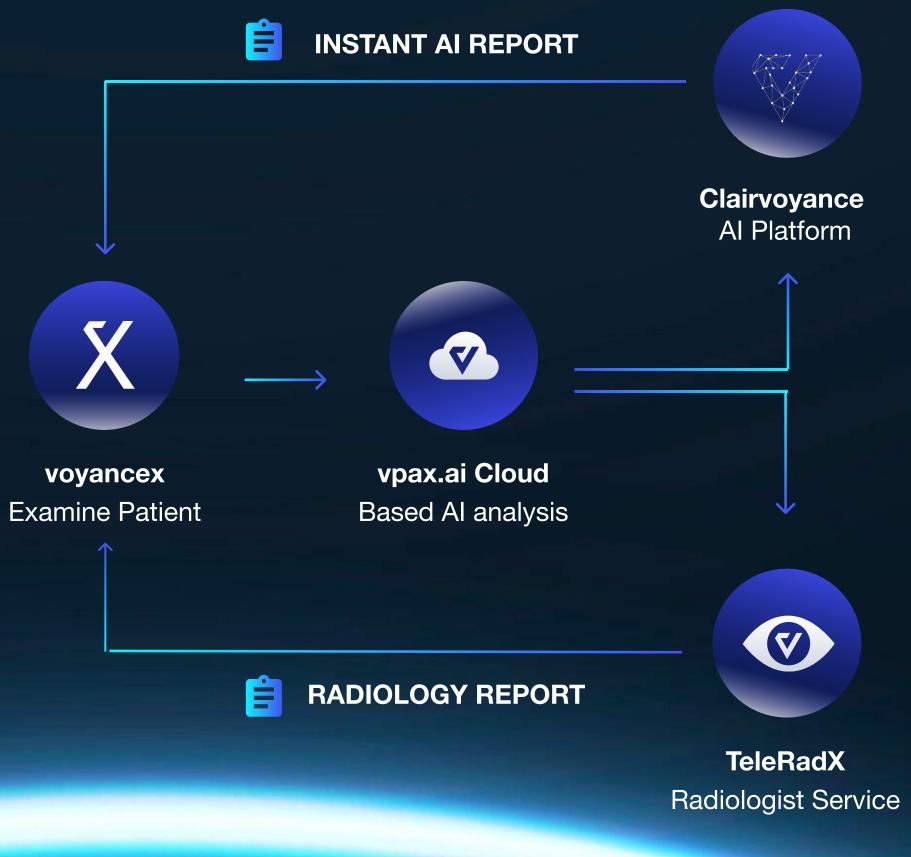

voyancex Acquisition Software

Advanced Acquisition Software for Digital Radiography: X-ray and Fluoroscopy





Intelligent Workflow

					Peterit UID				1/17/2023 🗇 1/1	17/2023 🖽	Default	Ckiar	D M Y	SEARCH +	Doleto
Schedded 91/2023 92.9 PM 91/2023 92.9 PM 91/2023	State	Created	Scheduled	Accession 1	Physician	Patient ID	Pallent nerr	Birth dwie	Age (years and mo	nths) Sex	Count	Exposure date	Description	Last Sending	Status
Schwaled 97/2023 (2.22 PM 97/20 PM 97/20 PM 97/20 P	Gpan i								Thyon						
Passed VIV/2023 S2/S/PM VIV/2023 S2/S/S/M VIV/2023 S2/S/S/S/M S2/S/S/S/S/S/M S2/S/S/S/S/S/S/S/S/S/S/S/S/S/S/S/S/S/S/	Scheduled	917/2023 12:29 PM	99/2023			20230045343	John Wick	\$101968	Sty Dro	Mile	ò				
Schedded VV/2023 G236 PM VV/2023 G236 PM </td <td>Scheduled</td> <td>917/2023 12:22 PM</td> <td>417/2020</td> <td></td> <td></td> <td>202300016543</td> <td>Reatin Tarag</td> <td>\$17,1990</td> <td>12y (m</td> <td>Male</td> <td>0</td> <td></td> <td></td> <td></td> <td></td>	Scheduled	917/2023 12:22 PM	417/2020			202300016543	Reatin Tarag	\$17,1990	12y (m	Male	0				
Deal VV/2022332/VM VV/20223/VM	Dane -	WINDOW COSTIN	Vinison				Dept Clair	Citation:	(and a	Mah		VOIDED 12/16 PM			
Parent V100023 0210 PM V100021 20235500 M MeX Statutore N101014 300 On Male 1 Schedded V100023 5216 PM V100023 20225585523 Kester Roy V101078 45y On Male 0	Scheduled	V1/2023 02:36 PM	917/2023			20232765363	Trella Claivaon	Section 1	35y Ora	Female	9				
Scheduled 1/17/2023 1/2:10 PM 1/17/2023 20235587523 Kastar Day 1/17/978 45y Dra Mala 0	Bien	WWWWWWWWWWW	(Physica)			-202345452824	Taran Witson	11121800	12,00	Permit					
	Panel	M10023 12:00 PM	1910/2020			2023456702M	data Spanne	dannings.	09/09	MAR					
Schedded V172023 5235 FM V172023 2022105423 DahlarRec 11722056 Wr Panado 0	Scheduled	MA2021 12:10 PM	977/2022			< 20236567623 S	Kester Tray	\$12(978	< dig Dre	Mala					

A modern single-page design with Advanced Intelligent automation capabilities that address any manual task or workflow steps to streamline processes, including worklist import, automatic image processing, and auto PACS archiving & sending, increase technologist productivity 1.5x times, and enhance patient care.





Image Enhancement

Delivers superb image quality through a set of intelligent image processing filters and specialized contrast enhancement techniques, tailored to the selected imaging protocol. Key features include:

- Contrast enhancement.
- Smart noise cancellation
- Automatic cropping
- Grid lines suppression











Auto Stitching

Auto and manual stitching for whole body, full leg, and whole spine imaging for orthopedic applications.



Al Bone Suppression

Our algorithm, ClairVoyance, utilizes state-of-the-art deep learning techniques to suppress bones in typical medical X-rays, eliminating the need for additional exposure doses required by the older dual-energy technique. This leads to faster diagnoses and improved patient care.



Stay Secure. Free Protection Against Ransom Attacks.

Ransomware Attacks Cost US Healthcare \$77.5B in Downtime. By 2031, global ransomware damages are forecast to exceed \$265 billion, with a new attack every two seconds. To combat this, Voyance is now offering its auto cloud backup service for free to all its customers and all healthcare facilities in the US. This initiative ensures enhanced security against ransom attacks and easy recovery from natural disasters.



- Automated Data and Settings
 Backup.
- Intuitive Touch-Enabled
 Interface.
- Comprehensive Keyboard Shortcuts.
- Multilingual Support: English, Spanish, Korean, etc.
- Multiple Users and Roles access.
- Storage Space Monitoring.
- Advanced Rejection Analysis: Exposure Trends, Projection Insights, and Repeat Analysis.
- Integrated with Vpax.Al Cloud PACS.



HW Integration

- Simultaneous Multi-Detector connection.
- Multi-Vendor Panel Integration.
- Enhanced Sleep and Wakeup
 Mode Support.
- Real-Time Panel Status Monitoring (Battery, Temperature, etc.).
- Fully Integrated with CPI, DRGEM, and IQflex Generators.

DICOM 3.0

DICOM Communication

- DICOM Store, Storage Commitment.
- DICOM Query/Retrieve.
- RIS Communication.
- Modality Worklist (MWL).
- Modality Performed Procedure Step (MPPS).
- DICOM Printing.
- Flexible DICOM Export Capabilities



- DICOM Printing Compatibility.
- Inkjet Printing Support.
- Free Layout Printing Option.
- True Size DICOM Printing.
- Configurable for Multiple
 Printers.

Image Export

- Multi-Format Image Export (DICOM, JPEG, BMP, TIFF).
- Mini DICOM Viewer Included with Image Export.
- Image Burning to CD/DVD/USB.
- Built-In Email Functionality.
- Auto-send Feature.



Exam Managment

- One-Click Emergency Mode for Immediate Shooting.
- Technique Chart Guidance.
- Auto-Positioning System.
- Exposure Index Tracking.

\mathbf{O}

Viewing Capabilities

- Customizable Layout for Intra-Study Viewing.
- Cross-Patient History Comparison.
- Advanced Image Editing.
- 15 Basic Imaging Tools.
- 32 Precision Measurement Tools.
- 15 Chiropractic Toolset.



Study Management

- Integrated MWL and Manual Patient Registration.
- PDF Report Generation: Includes Images, Study Notes, and Patient Information.
- Seamless Image Copying Between Studies.
- Direct Hard Drive Image Import via Copy-Paste.
- Sending Queue Management.
- Auto-Send Study to PACS.



IMAGE MODIFICATION	Edit Body partSoft and Hard Tissue	ContrastBrightness
TOOLS	 Sharpness 	 Grid Suppres (Vertical & Ho
IMAGE BASIC TOOLS	 Full Screen Mode Selection Mode Screen Layout Crop 	 Rotate Clock Counterclock Angle Measu Zoom Reset L and R Labe
<section-header></section-header>	 Point Line Arrow Free Distance Circle Three Point Circle Circle Angle Vertical Discrepancy Meta-Diaphyseal Angle 	 Angle Cross Angle Middle Angle V-Axis H-Axis Free Axis Horizontal Leve
CHIROPRACTIC TOOLSET	 Spine Label Spine Analysis Cobb Angle Vertebral Line George's Line 	 Cervical Curve Atlas-Axis Arve Atlas (C1) Plate Cervical Curve

ession Iorizontal)	 Grid lines selection Spline vs. Linear Contrast Curve Graph Invert 	
kwise/ ckwise urement t bel	 Two Point Measurement Flip Horizontal Flip Vertical Window Level 	 Free Label Reset Image Delete Image
e evel el	 Free Axis Level Point Ratio Three Point Ratio Cardio-Thoracic Ratio Leg Length Discrepancy Femoral Symmetry Acetabular Angle Iliac Angle 	 Double Line Axis Center Edge Angle PMAA / DMAA Axis Angle IMA / HVA Measurement Calibrati Free Rotate
rve Ingle Iane Angle rve C2-C7	 Cervical Curve C1-T2 Spinolaminar Line Vertebral Angulations Lumbar Curve 	 Spondylolisthesis Pelvic Parameters Pelvis Analysis



ion

Multi Vendor Integration

X ray integration





Flat Panel Detectors



DRTECH





















We'd love to hear from you!

Address

50 Maxwell, Irvine, 92618 CA, USA

Phone Number +1 (949) 284- 5179

Mail info@voyance.health



