

IC222 Specifications

Adjustable Operations



Head Lift
Angle of 62°



Knee Break
Angle of 42°



Height Adjustable
30cm to 61cm from
floor to top of base

Clearance from
floor to rail - 21cm

Overall Bed and Mattresses Sizes

Measured in millimetres	Long Single	King Single
Overall Base Measurement	2150 x 960	2150 x 1120
Mattress Surface Area Measurement	2030 x 900	2030 x 1070

Wheels

The IC222 features 4 x hospital grade swivel castors for ease of moving
All castors lockable.



Mattress Surface

Strong metal mesh under mattress allowing increased airflow



Actuators

Long Single and King Single size
1x 6000 newton - 24 watt height adjust actuators
1x 6000 newton - 24 watt head adjust actuator
1x 6000 newton - 24 watt knee adjust actuator

Hand Piece

- Backlit for night time use
- Easy press buttons
- Wired stretch cord
- Hook attachment for hanging to side of bed or convenient place



Fabric

- Hard wearing commercial grade Warwick fabric
- Ardo Onyx range (BLACK) 82% Polyester, 18% Olefin
- Ada Stone range (GREY) 70% Polyester, 30% Acrylic



Warranty

- Steel Frame of Bed - 5 Year Limited Warranty
- Actuators - 5 Year Limited Warranty
- Hand Piece - 2 Year Limited Warranty
- Control Box - 5 Year Limited Warranty
- Upholstered Surround and Fabric - 3 Year Limited Warranty
- Castor wheels - 5 Year Limited Warranty
- For Terms & Conditions - please refer to the User Manual and Warranty Information booklet

Safe Working Load

- Long singles and king singles - 180kg
- The safe working load includes any model of i-care mattress and 2 side rails

Mattresses Retention


- All beds come with a mattress retainer which slots into the foot end of the bed
- There is also a 25mm lip around the mattress from the frame surround which holds the mattress in place

Headboards

- Upholstered Headboards attaches on to the back of the IC222 with 4 hand bolts
- Available in matching fabric
- Easy to remove

Companion Bases

- Companion bases are a standard non mechanical base that can be pushed up against the IC222. These do not move or adjust. They cannot be attached or joined to the IC222 because of the high low function. They must have a separate headboard to the IC222. There is a small gap between the mattresses when this is pushed against the IC222. However this is required for bedding and to allow for the IC222 to adjust height.

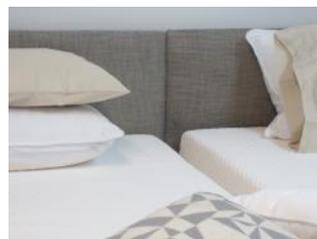

Side rails and poles

High Side Rails

Height 72cm from top of base

Width: 36cm

2 brackets can be fitted anywhere along the side of the bed

These hand clamp on the bed frame surround and can also be screwed for extra strength

The rail then can slide in or out of brackets

Easy to remove rail or brackets

Low Side Rails

Height 42cm from top of base

Width: 68cm

2 brackets can be fitted anywhere along the side of the bed

These hand clamp on the bed frame surround and can also be screwed for extra strength

The rail then can slide in or out of brackets

Easy to remove rail or brackets

Bedstick

Height 75cm from top of base

1 bracket can be fitted anywhere along the side of the bed.

These hand clamp on the bed frame surround and can also be screwed for extra strength

The rail then can slide in or out of bracket

Easy to remove rail or bracket

Handgrip at the top of bedstick

Full Length Fold Down Rail

Height 34cm from top of base

Width: 138cm

Brackets can be fitted anywhere along the side of the bed

These hand clamp on the bed frame surround and can also be screwed for extra strength

The rail then can be folded down below mattress with a simple pull catch

Easy to remove rail

3 horizontal rails

Overbed pole

Free Standing

Overall Height - 175cm

Adjustable handle height

Overbed table

Height adjustable

Castor wheels