


VFC400 Software Setup



Open Control Solutions VTMC

by clicking on the icon on your desktop

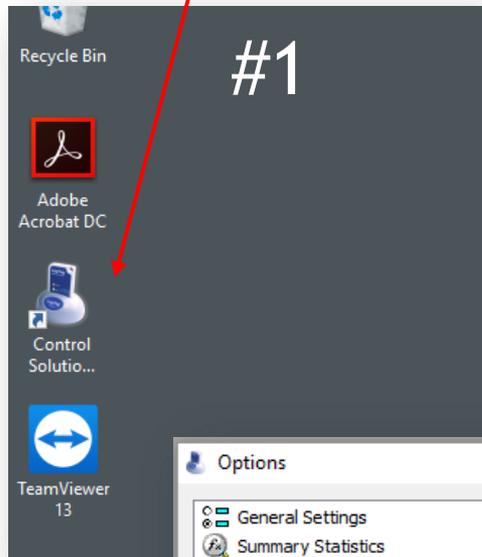


#1

#2


Next, open the Edit menu and choose Options

#3

From the Options menu choose General Settings and select your preferred temperature scale: F° or C°*

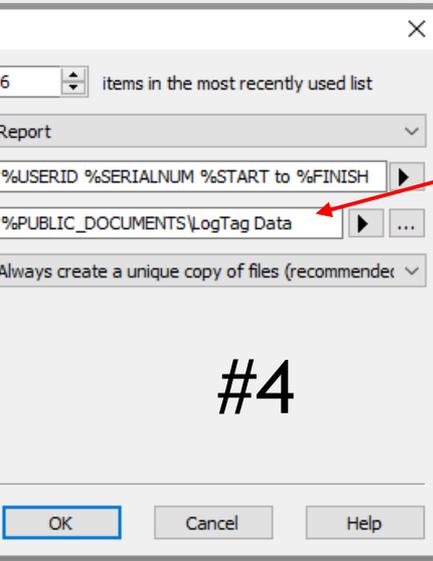

#4

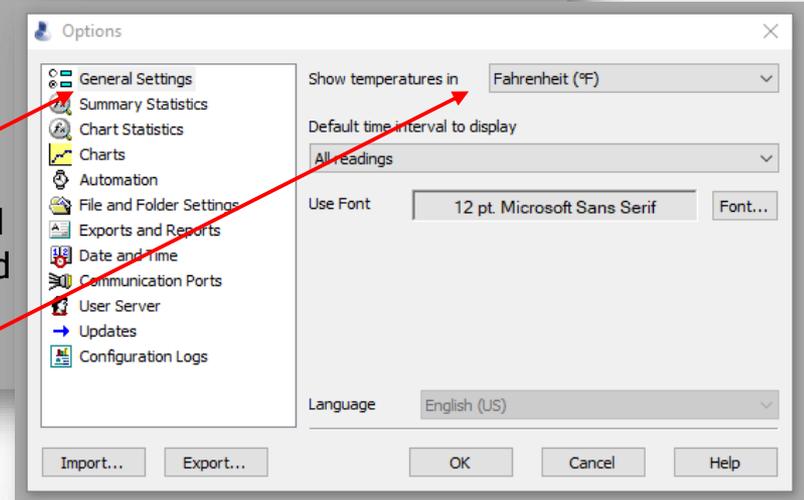

Data will automatically be saved to the newly created folder Public>Documents>LogTag Data every time your logger is placed in the docking station. Anyone logged into the PC can access the data.

If you prefer to save the data to a network or shared drive, click the button with 3 dots(...). A box will appear saying "Locate Folder". Simply go to the folder you want to save the data to and click "OK".




Software Configuration Continued.

VTMC customers will need to follow these steps.

The automation settings will verify cloud connectivity.

1. Click "Automation"

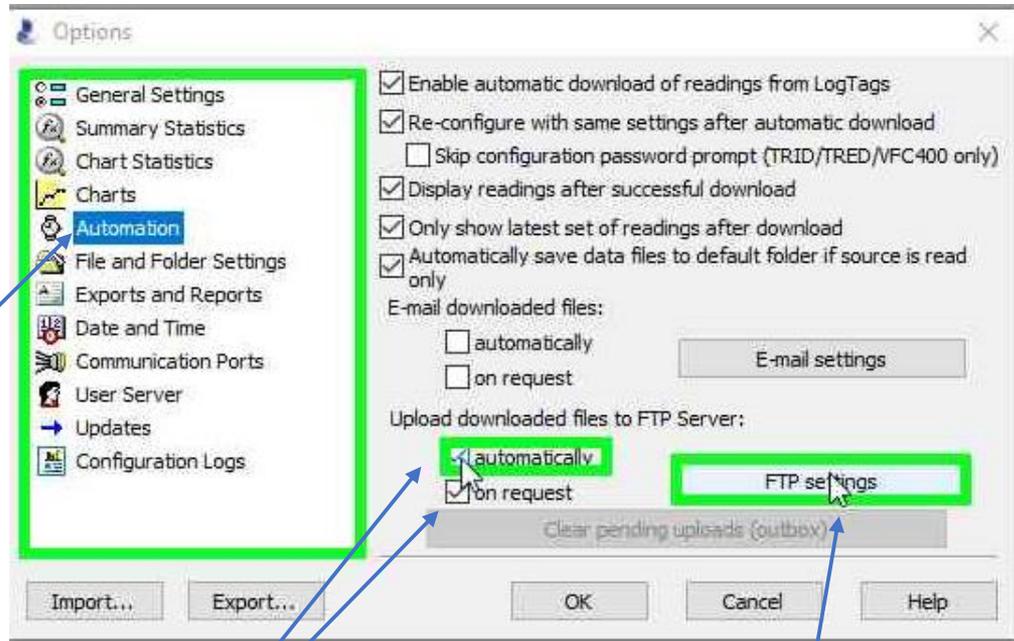

2. Click Check both of these boxes.

4. Click FTP settings



Software Configuration Continued.

VTMC customers will need to follow these steps.


The automation settings will verify cloud connectivity.

1. These three fields should match your software.

2. Click only "CSV Comma delimited"

3. Click "Test Connection"

4. You will receive this message if you have successfully connected to the cloud.

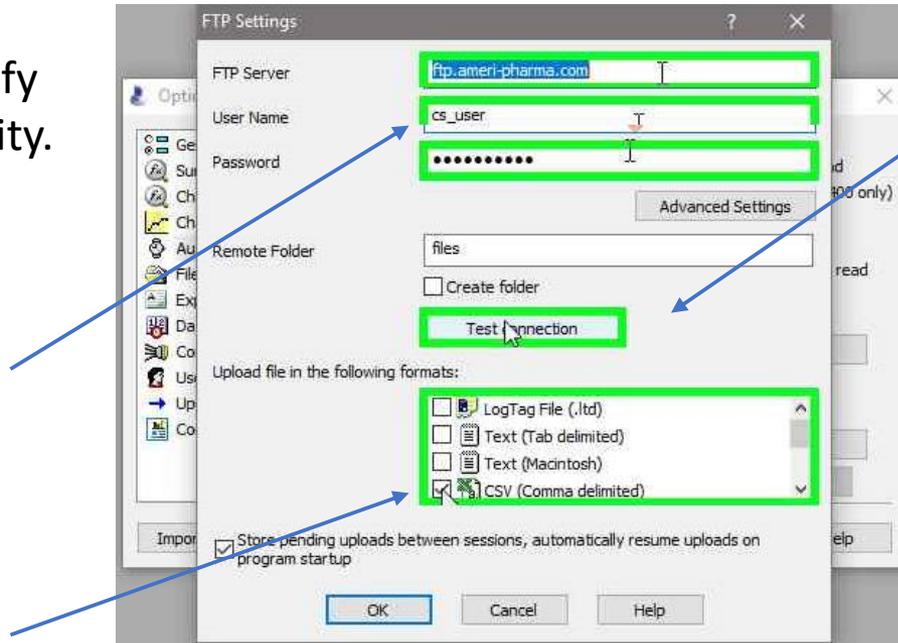

Choose OK

to complete software configurations

the remaining options do not require changes

Other options are configurable but not necessary

